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
INSTRUCTIONS BEAUTY MOLE REMOVAL SWEEP SPOT PEN

OPERATING GUIDE

1. Power-on key (M): long press (M) key for 2 seconds to turn on and off the machine, and then the gear position will be displayed (default: 1).
2. After power on, press the power button (M) to turn on the led light. Press the light again to turn it off.
3. Short press plus or minus to switch the output power, with 9 gears it can be adjusted, 1 to 9 gradually increase.
4. Work button (OUT): press the button (OUT) to start working, loosen it to stop working.
5. Protection function: long press (OUT) for 2 minutes to stop the device from losing out, if you need to continue working, release the (OUT) key after again press once.
6. Electric quantity display and charging instructions: the electric quantity is displayed in four grids. When charging, the indicator will beat in turn. When fully charged, the indicator light will stop flashing.

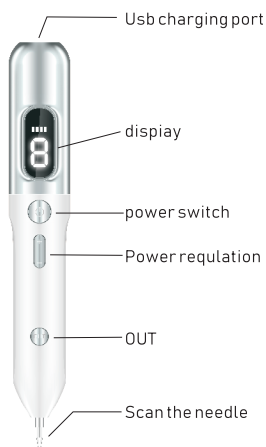function Introduction

To use high and new technology, ion carbonization function, instant sweep freckles, spots, flat mole And so on, can achieve fast, safe without side effect effect under the circumstance of not bleeding.

Usage

In addition to the spot:

Including the contact spot location moves back and forth slightly, can be clearly found that the melanin in electric acupuncture cue after Saul zoom black fall off, if precipitation shallow spots, area is lesser, cana removal. If a larger area of pigmentation is deeper, it must be 2 or 3 times to deal with more appropriate, should also pay attention to the hairdresser when using electric acupuncture treatment, besides keeping customers pose, according to the hairdresser with the hand of electric acupuncture should balance the jitter, lest cause burns when processing, shades, influence wound healing, leave scar phenomenon.

Dot mole

When dot mole, pay special attention to distinguish the mole extent, some mole grows on the skin surface layer, and root growth of moles, deep subcutaneous tissue, including using the needle at a time when the end of the cut surface of moles, the best way to burn in small squares, if growth relatively deep mole, besides should avoid inflammation after processed had better be in 2 to 3 days before the wound healing nourishing cream coated with a small living cells, can prevent skin concave hole, when wiping nutrition frost, should pay attention to the wound with inflammation of the situation, if is not to wipe.

In addition to granulation:

The electric acupuncture mouth slightly contact granulation, derivative granulation root cells to radium radio needle burned to death, about a week after granulation can fall off on their own, such as granulation is small, can use laser radio disposable needle burning breeding ground for granulation, such as granulation is big, the duty to burn for many times, lest affect wound healing.

In addition to tattoo and tattoo:

Generally burn tattoos and tattoo, its area is larger, appropriate USES large needle, the method of burning tattoo tattoo, roughly the same as the scanning spot, operating time may be slightly longer, the area is too large (more than 3cm) should deal with many times, everyone time.

Otherwise, accidentally bacteria infection, consequence is unimaginable, beautician and customer should be especially careful if you have severe inflammation of the situation, as soon as possible to the hospital for treatment.

Note

1. It is recommended to use fresh fruit test equipment when using this equipment for the first time, when fully cooked. Practice equipment performance and operation procedures before use.
2. Due to the effect of this equipment can produce high temperature in an instant, when not in use, please properly keep away from children and do not know the equipment use and to play, to avoid danger happening.
3. The use of this equipment had better not too clean, affected within 30 days (simple) under the care not scrub, not artificially make hard skin fall off.
4. when the scab place itchy (tender meat phase), don't by hand grasp, for 20 to 30 days right, scab will fall off on their own.
5. Recovery period is 3 months or so, it is suggested that during the recovery period is unfavorable eat ginger (affect the skin smoothness beef (produces light red) soy sauce (produces light black) in convalescence, the good things above, the skin can be restored to its original 95% of skin burn pigment directly by machine. Regardless of any color is OK. Then thin scab skin, the skin pigment will all concentration below the scabs, scab off the pigment in the skin such as follows scabs all fall out, absolutely not infection.
6. If long time not to use this equipment should every three to four months to charge, keep the battery using life.
7. The device failure occurs, please don't open trying to repair in order to avoid risk, suggest to contact the original factory maintenance.
8. It is strictly prohibited to water or clean this equipment directly, like liquid to avoid entering the device internal cause damage or danger happening.

Product parameters

Product name: beauty mole removal sweep spot pen
Product model: dz-01 充电输入: DC5V=1A
Battery capacity: 1200mA 输出功率: 4W
Product size: 190*28*29 净重: 85g