Accutrend Plus

User's Manual

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Manual do Utilizador

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Accutrend[®] Plus

User's Manual

©2011, Roche Diagnostics. All rights reserved

The contents of this document, including all graphics, are the property of Roche Diagnostics. Information in this document is subject to change without notice. Roche Diagnostics shall not be liable for technical or editorial errors or omissions contained herein. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Roche Diagnostics.

Please send questions or comments about this manual to your local Roche representative.

ACCU-CHEK, ACCU-CHEK MULTICLIX, ACCUTREND, COBAS and SAFE-T-PRO are trademarks of Roche.

All other trademarks are the property of their respective owners.

This instrument has been constructed and tested in accordance with EN 61010-1 ("Safety requirements for electrical equipment for measurement, control and laboratory use; Part 1: General requirements") and left our factory in a safe condition.

To ensure that this remains so and to ensure safe operation the user must observe the instructions and warnings given in this manual.

Installation, use and maintenance of the Accutrend Plus meter lie in the full responsibility of the user.

On the packaging and on the identification plate of the instrument you may encounter the following symbols, shown here with their meaning:

Caution, consult accompanying documents. Refer to safety-related notes in the manual accompanying this instrument.

Temperature limitation (Store at)

Use by

Manufacturer

Batch code/ Lot number

Catalogue number

In vitro diagnostic medical device

0123

This product fulfills the requirements of the European Directive 98/79/EC on *in vitro* diagnostic medical devices.

Consult instructions for use

c US LISTED

The system fulfills the Canadian and U.S. safety requirements (UL LISTED, in accordance with UL 61010A-1:02 and CAN/CSA-C22.2 No. 61010-1-04)

1	Introduction	7
	Before you start.....	7
	Intended use.....	7
	Important information regarding use.....	7
	If you need help.....	7
	Test principle.....	8
	Checking the contents.....	8
	Important safety instructions and additional information.....	9
	Safety information.....	10
	Disposal of the system.....	10
	General care.....	11
	Electromagnetic interference.....	11
	Operating conditions.....	12
	Quality control.....	13
2	The Accutrend Plus Meter	14
	Overview of the meter elements.....	15
	Display information and icons.....	16
	Power supply.....	18
3	Starting Up	19
	Inserting the batteries.....	20
4	Meter Setup	23
	Brief overview of the meter settings.....	23
	General procedure for setting up the meter (set mode).....	24
	Setting the date format.....	26
	Setting the date.....	26
	Setting the time format.....	28
	Setting the time.....	28
	Setting the beeper.....	29
	Setting the lactate display.....	30
5	Performing a Test	31
	Brief overview of the steps to be carried out.....	33
	Code strips.....	35
	Powering on the meter.....	36
	Inserting the code strip.....	37
	Checking the stored code numbers.....	38
	Sample material.....	39
	Testing by healthcare professionals.....	40
	Preparing to perform a test.....	41
	Performing a test.....	43
	Recommendations for the collection and measurement of capillary blood.....	45
	Blood collection.....	46
	Blood application inside the meter.....	47

Blood application outside the meter.....	48
Starting the measurement	49
Display of results	50
Disposal of used items after testing with blood	51
Flagging test results.....	52
6 Control Testing	55
Preparing to perform a control test.....	56
Performing a control test.....	56
Applying the control solution	59
Starting the measurement	60
Display of results	61
Disposal of used items after control testing	61
7 Memory	63
Reviewing results in memory.....	63
Deleting results from memory	66
Deleting the last result.....	67
Deleting all results of a test parameter.....	68
Deleting all results.....	70
8 Cleaning	71
Recommended cleaning/disinfecting solutions	71
Cleaning the exterior (meter housing).....	71
Cleaning the measurement chamber cover/test strip guide	72
Cleaning the optical measuring system.....	73
9 Troubleshooting	75
10 Product Specifications	79
Operating conditions and technical data.....	79
Sample material	79
Storage and transport conditions	80
Ordering	80
Product limitations.....	80
Contact Roche	81
Index	83

1 Introduction

Before you start

Intended use

The Accutrend Plus meter is used for the quantitative determination of four blood parameters: **glucose, cholesterol, triglycerides, and lactate**. The reflectance photometric measurement is performed using test strips specific for each of these blood parameters. For detailed information about each test, see the package inserts of the respective test strips.

The meter is suitable for professional use as well as for self-testing.

Important information regarding use

This manual contains all the information you need to operate and care for the Accutrend Plus meter. **Read this entire manual carefully** before using the system for the first time.

Pay particular attention to the *Important safety instructions and additional information* section in this chapter before you start testing. Also read the package inserts of the test strips and the lancing device to be used.

Healthcare professionals: In addition, observe carefully the precautions and procedures applying to professional use provided in the relevant sections of this manual.

Note: Before you use the meter for the first time (after you have first inserted the batteries), you must set the date and time correctly to allow you to perform measurements properly. Each time you replace the batteries you need to check (and, if necessary adjust) the date and time.

If you need help

If you have any questions about the Accutrend Plus meter, please contact your local customer support and service centre. You will find contact details on page 81.

Test principle

The meter reads the lot-specific characteristics of the test strips currently in use by means of a code strip. This information is stored and therefore needs to be read only once per test strip container. To run a test, an unused test strip is taken from the container and inserted into the meter. When a strip is inserted, the application area of the test strip is illuminated by an LED (light-emitting diode) from below. Before sample application, the reflection behaviour of the test strip is determined by means of the light which is reflected from the application area (blank value).

The blood sample is then applied to the application area and the measurement chamber flap is closed. The constituent to be determined in the applied sample undergoes an enzymatic reaction with formation of a dye. The amount of dye formed increases with the concentration of the substance to be determined.

After a certain period of time (dependent on the test parameter), the colour intensity is measured by illuminating the application area again from below using the LED. The intensity of the reflected light is measured with a detector (reflectance photometry). The measured value is determined from the signal strength of the reflected light, using the previously measured blank value and the lot-specific information from the code strip. Finally, the result is displayed and simultaneously saved to memory.

Checking the contents

- Accutrend Plus meter
- User's manual
- Four batteries (1.5 V, AAA)

Important safety instructions and additional information

This section explains how safety-related messages and information related to the proper handling of the system are presented in the Accutrend Plus User's Manual. Read these passages carefully.

The safety alert symbol by itself (without a signal word) is used to promote awareness to hazards which are generic or to direct the reader to related safety information

WARNING

Indicates a hazardous situation which, if not avoided, could result in death or serious injury.

CAUTION

Indicates a hazardous situation which, if not avoided, could result in minor or moderate injury.

NOTICE

Indicates a hazardous situation which, if not avoided, may result in damage to the system.

Important information that is not safety relevant is presented against a grey background (without a symbol). Here you will find additional information on correct use of the meter or useful tips.

Illustrations in this manual show two different kinds of hands:

Hand without glove

Hand with glove

Safety information

Protection against infection

Any object coming into contact with human blood is a potential source of infection.

Patient Self-Testers: There is a potential risk of infection if you let other people use your meter or lancing device. Do not share your meter or lancing device with any other person.

Protection against infection in a healthcare environment

There is a potential risk of infection. Healthcare professionals using the Accutrend Plus meter must be aware that any object coming into contact with human blood is a potential source of infection. Healthcare professionals must also be aware that any cross-contamination is a potential source of infection for patients.

- Use gloves.
 - Use an auto-disabling single-use lancing device for each patient.
 - Dispose of used lancets in a sturdy sharps container with lid.
 - Dispose of used test strips according to your institution's infection control policy.
 - Follow all health and safety regulations in force locally.
-

Disposal of the system

Risk of infection

During testing the meter may come into contact with blood. Used meters therefore carry a risk of infection. Dispose of your used meter - after removing the batteries and cleaning the meter - according to the regulations applicable in your country. For information about correct disposal, contact your local council. The meter falls outside the scope of the European Directive 2002/96/ EC (Directive on waste electrical and electronic equipment (WEEE)).

For healthcare professionals: Infection by a potentially biohazardous instrument

Dispose of the meter in line with your facility's procedures for dealing with potentially biohazardous waste.

Exploding batteries

Do not throw used batteries into an open fire. They may explode.

Disposal of used batteries

Do not dispose of the batteries with normal domestic waste. Dispose of used batteries in an environmentally responsible manner and in accordance with applicable local regulations and directives. Contact your local council/local authority or the manufacturer of the used batteries for advice on correct disposal.

General care**NOTICE**

Clean the meter only with the solutions recommended (see page 71). Using other solutions may result in incorrect operation and possible system failure. Do not let cleaning solution enter the instrument. Make sure that the meter is thoroughly dried after cleaning or disinfecting.

Electromagnetic interference

Electromagnetic interference

Strong electromagnetic fields may interfere with the proper operation of the meter. Do not use the meter near strong electromagnetic fields.

Operating conditions

To ensure that the meter functions properly, please observe the following guidelines:

- Operate the meter only within the acceptable temperature range*. This range is test-dependent:
 - For cholesterol 18-35 °C
 - For glucose 18-35 °C
 - For triglycerides: 18-30 °C
 - For lactate 15-35 °C
- Only use the meter at a relative humidity between 10% and 85% (no condensation).
- When testing, place the meter on a level, vibration-free surface, or hold it so it is roughly horizontal.

* **Please note:** Different temperature ranges apply for Accutrend Controls. Refer to the package insert for details.

Quality control

The Accutrend Plus meter comprises numerous integrated or available control functions such as the following:

- An automatic check of the electronic components and functions when the meter is powered on.
- An automatic check of the ambient temperature before and during the measurement.
- An automatic check of the test strip to make sure that the code information necessary for measurement exists.
- A check of the optical system and the overall function by means of control solutions.

2 The Accutrend Plus Meter

Overview of the meter elements

A **Display**

Shows results, information, icons and results recalled from memory.

B **M (memory) button**

Press this button to recall results from memory and (together with the **Set** button) to change the meter settings.

C **On/Off button** ①

Press this button to power the meter on or off.

D **Measurement chamber flap**

To apply the sample, open this flap. To start the measurement, close the flap.

E **Test strip guide**

Insert the test strip here.

F **Infrared (IR) window**

This infrared interface is for Roche internal use only. You cannot use it to download results from the meter to a computer.

G **Set button**

By pressing this button, you access the various meter settings which are changed using the **M** button. You also use this button to switch between the different test parameters to show the currently stored code numbers (before measurement), or to review results (when in memory mode).

H **Battery compartment cover**

Provides access to the battery compartment (4 AAA 1.5 V alkaline manganese batteries).

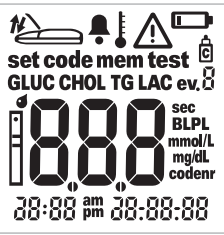
I **Serial number**

on the identification plate.

J **Measurement chamber cover (with test strip guide)**

You can remove this cover to clean the test strip guide.

Display information and icons

Icon	Meaning
	<p>Every time the meter is powered on, you can check the display. The meter temporarily shows all elements that may appear in the display.</p> <p>Regularly check that all of the display elements are displayed properly. Results may be misread if a segment is missing.</p>
	Close measurement chamber flap
	Open measurement chamber flap
	Beeper activated
	Temperature warning
	Error (see: <i>Troubleshooting</i>)
	Battery warning (batteries almost empty)
	Flagged as a control test with control solution
	Flag for specific events (event 0–9)
	Set mode
	Memory mode
	Code display
	Test mode (measurement)
	Test strip flashes: insert test strip static: test strip is inserted
	Test strip and blood drop: apply blood
	Test parameter: glucose
	Test parameter: cholesterol

Icon	Meaning
LAC	Test parameter: lactate
TG	Test parameter: triglycerides
BL	Lactate displayed as a whole blood value
PL	Lactate displayed as a plasma value
mmol/L	Standard unit for lactate and (in some countries) for glucose, cholesterol and triglycerides
mg/dL	Unit for glucose, cholesterol and triglycerides (in some countries)
sec	Measurement time in seconds
codenr	Display of the code number
am	In the morning (if using the 12h time format with a.m./p.m.)
pm	In the afternoon (if using the 12h time format with a.m./p.m.)

Power supply

To save power, the Accutrend Plus meter automatically powers off after 2 minutes unless a button has been pressed or a new test strip has been inserted. When the meter powers off, all results obtained up to that point remain in memory.

With a set of fresh batteries, you will normally be able to perform at least 1,000 measurements. When the battery warning is displayed for the first time, approximately 50 measurements can still be performed. In this case, replace the batteries as soon as possible.

When replacing the batteries you must insert the new batteries within one minute of removing the old ones, to keep the date and time settings. If you take longer than this, you must re-enter the date and time. Use only alkali-manganese batteries type AAA.

Results are retained in memory together with the date and time even when no batteries are inserted. All other settings are retained, too.

Think of the environment. Dispose of used batteries sensibly.

WARNING

Exploding batteries

Do not throw used batteries into an open fire. They may explode.

3 Starting Up

Before using the meter for the first time, perform the following steps:

- 1 Insert batteries.
- 2 Set date, time, and beeper.
- 3 Select how lactate test results are to be displayed (blood or plasma value).
- 4 Insert code strip (can also be done directly before performing a test).

Inserting the batteries

- 1 With the meter powered off, turn it over.
- 2 Press the battery compartment cover release tab.

- 3 Lift off the cover.

- 4 Insert the four batteries in the battery compartment as indicated. Use only alkaline manganese batteries (1.5 V, AAA). Always replace all four batteries at the same time because batteries with different capacities may impair the function of the meter.

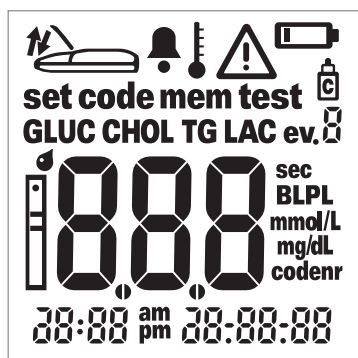
Do not use rechargeable batteries.

Pay attention to the positions of “+” (top of battery) and “-” (flat end).

- 5 Close the battery compartment cover.
- 6 Power the meter on to test the function of the new batteries.

- 7 Check that all of the display elements are displayed properly (see illustration). Results may be misread if a segment is missing.

If the display check when you power on the meter is too quick, you can put it “on hold” by pressing and holding down the **On/Off** button ① the next time you power on the meter. For as long as you keep the button pressed, the display remains frozen.

This page intentionally left blank.

4 Meter Setup

Brief overview of the meter settings

The following table provides an overview of the available settings.

Setting	Options	Default setting *
Date format	Day.Month.Year (31.12.00) Month-Day-Year (12-31-00)	Day.Month.Year
Date		31.12.00
Time format	24-hour time format (24h) 12-hour time format (12h), with am/pm	24h
Time		0:00
Beeper	On Off	On
LAC result display	Blood (BL) Plasma (PL)	BL

* "Default setting" means the setting that is in the meter when it leaves the factory.

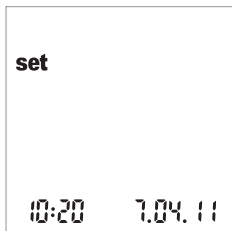
General procedure for setting up the meter (set mode)

You make all settings using the **Set** and **M** buttons as described below. Please note that the meter must be powered off before you can activate the set mode.

- If you power on the meter without setting the date and time, dashes will appear where the date/time information should be. If no date/time setting is made, no date/time information can be saved with the test results in memory.
- If the meter has lost previously entered date/time settings (e.g., after a prolonged period without batteries), you need to enter date/time settings again.

- 1 Press the **Set** button (on the left side of the meter) to power the meter on in set mode.

- 2 The date and time as well as the icon *set* now appear on the display. To make or change settings, press the **Set** button again.

To leave set mode at any time, press the **On/Off** button **ⓘ**.

3 If the displayed setting is correct (e.g. the date is correct and you want to change the time only), you can continue directly to the next setting by pressing the **Set** button **or**:

4 Press the **M** button to change the currently flashing setting. Press the **M** button as often as needed (or keep it pressed) until the desired setting (number) is reached.

Settings with only two options (date/time format, beeper, LAC display and unit) can be toggled on and off with the **M** button.

5 Press the **Set** button again to confirm (save) the current setting and go to the next setting.

You can only move forward through the settings. Moving backwards is not possible. Corrections can only be made by repeating the settings.

The setting procedure can be terminated at any time by pressing the **On/Off** button **ⓘ**. All settings made up to that point are saved.

Setting the date format

In the first step, you set the date format (the entire date flashes).

When elements in the display are surrounded by a halo, it indicates that they are flashing in the meter's display.

The following date formats are available:

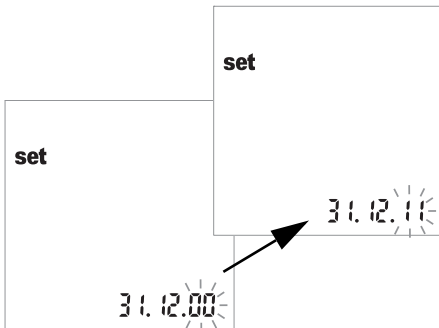
- **31.12.00 (= default setting) Day.Month.Year**
- 12-31-00 Month-Day-Year

- 1 Press the **M** button to select the date format. Each time you press the button, the (flashing) formats *31.12.00* and *12-31-00* appear in alternation.
- 2 Press the **Set** button to save your chosen date format. The display then automatically moves on to the date setting.


Setting the date

The next three settings involve entering first the **year**, then the **month** and finally the **day**.

- 3 Press the **M** button to change the number that is currently flashing and set the current year.
- 4 Press the **Set** button to save the year. The display then automatically moves on to the month setting.

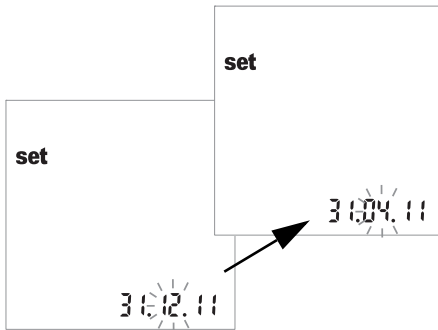If you use the meter without setting the date, no date information will be saved with the test results.


- 5 The default month flashes. Press the **M** button until the correct month is displayed.
- 6 Press the **Set** button to save the setting. The display then automatically moves on to the day setting.

- 7 The default day flashes. Press the **M** button until the correct day is displayed.
- 8 Press the **Set** button to save the setting. The display then automatically moves on to the time format setting.

Setting the time format

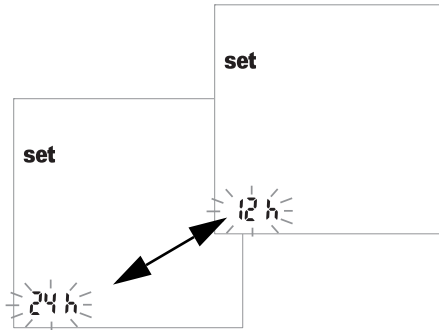
Now select the time format. You can choose between the 24-hour time format (default setting) and the 12-hour time format with “a.m.” or “p.m.”.

- 9 Press the **M** button toggle between 24-hour and 12-hour format.
- 10 Press the **Set** button to save the chosen time format and continue with the time setting. The display automatically moves on to the time setting.

Setting the time

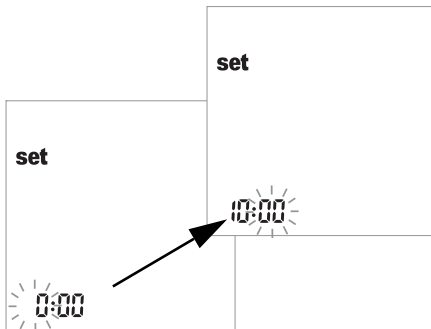First you set the hour, then the minutes.

- 11 Press the **M** button to change the number that is currently flashing. Press the **Set** button to confirm. You can now set the minutes (again using the **M** button). If you have selected the **12-hour** time format, *am* changes to *pm* and vice versa when you reach “12:xx”.
- 12 Press the **Set** button to save the correct setting and to continue with setting the beeper.

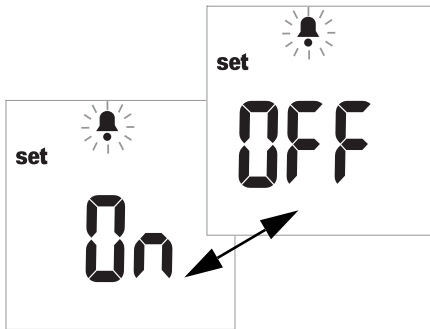
If you use the meter without setting the time, no time information will be saved with the test results.

Setting the beeper

After setting the time, you can set the beeper by choosing *On* or *OFF*.

If the beeper is activated, you will hear a beep in the following situations:

- When the meter detects that the test strip is inserted.
- When the result appears.
- When an error occurs.

13 Press the **M** button to toggle between *OFF* and *On* (*On* is the default setting).

14 Press the **Set** button to save the selected setting. The display automatically moves on to the next setting option.

We recommend that you keep the beeper activated at all times.

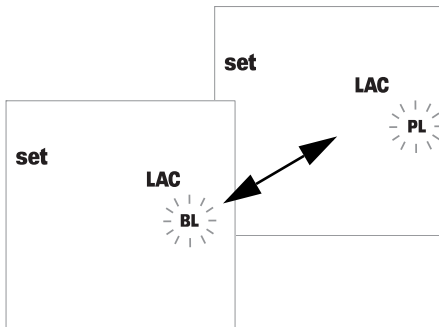
Setting the lactate display

In the default setting, lactate values are represented based on a whole-blood measurement (*BL* is displayed). If you want to display the result as plasma referenced (*PL*), you can switch the representation. This setting affects **only** the (internally converted) representation of the value.

Risk of incorrect results with plasma samples

When performing a **lactate test**: Do not use plasma for testing. Using plasma will lead to incorrect results. Although the test result can be displayed as plasma-referenced values, only fresh or heparinized capillary blood is permitted as sample material.

15 Press the **M** button to toggle between *BL* (blood) and *PL* (plasma).

16 Press the **Set** button to save the setting. The setup procedure is now complete. The meter automatically moves to test mode.

5 Performing a Test

What you need:

- Accutrend Plus meter
- Test strip for selected test parameter* with matching code strip:
 - Accutrend Glucose
 - Accutrend Triglycerides
 - Accutrend Cholesterol
 - BM-Lactate
- Lancing device:
 - **Patient Self-Testers:** Use a lancing device appropriate for personal use (for example, the Accu-Chek Multiclix lancing device from Roche).
 - **Healthcare Professionals:** Use an auto-disabling single-use lancing device for each patient. The lancing device must be intended for use by healthcare professionals in a multiple patient setting (for example, the Accu-Chek Safe-T-Pro Plus lancing device from Roche). Follow the manufacturer's instructions for use.
- Alcohol or cellulose pad, if required

* Note: Not all tests are available in all countries.

Protection against infection

Any object coming into contact with human blood is a potential source of infection.

Patient Self-Testers: There is a potential risk of infection if you let other people use your meter or lancing device. Do not share your meter or lancing device with any other person.

Healthcare Professionals:

- See safety message *Protection against infection* on page 10.
 - See safety message *Potential risk of infection* on page 40
-

Always ...

- ... operate the meter only within the acceptable temperature range. This range is test dependent (refer also to the insert of the respective test strip):
 - For cholesterol 18-35 °C
 - For glucose 18-35 °C
 - For triglycerides: 18-30 °C
 - For lactate 15-35 °C
- ... place the meter on a level surface or hold it steady in your hand.
- ... make sure that all display elements are displayed during the display test.
- ... read the test strip inserts.
- ... keep the test strip guide and housing clean (see description starting page 71).

Never ...

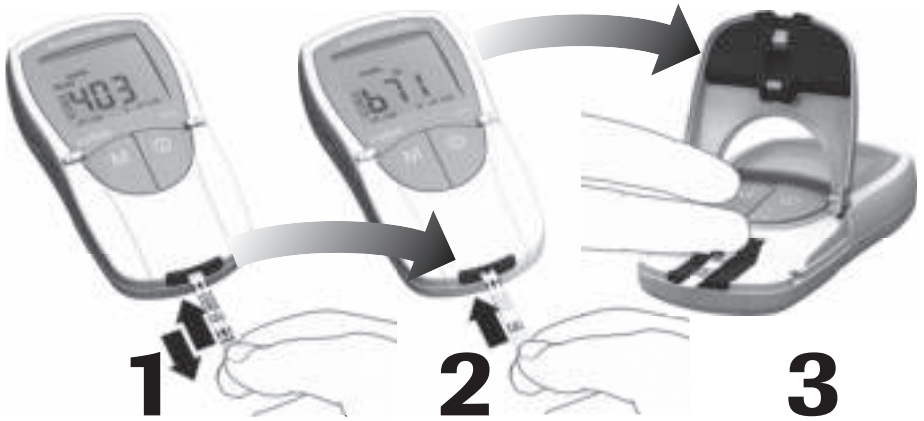
- ... touch or remove the test strip during actual measurement (although this is possible before starting the measurement when blood is applied outside the meter).
- ... delay starting the measurement after blood application.
- ... subject the meter to sudden movements during a measurement.
- ... store the meter and strips at extreme temperatures (see *Product Specifications* on page 79 and test strip inserts).
- ... store the meter and test strips under humid or damp conditions without suitable protection (see *Product Specifications* on page 79 and test strip inserts).

Accuracy/precision of measured results

Failure to comply with the above may lead to inaccurate results. An incorrect result may lead to an error in diagnosis, therefore posing danger to the patient.


Brief overview of the steps to be carried out

Code the meter (once per test strip container)

Insert the test strip

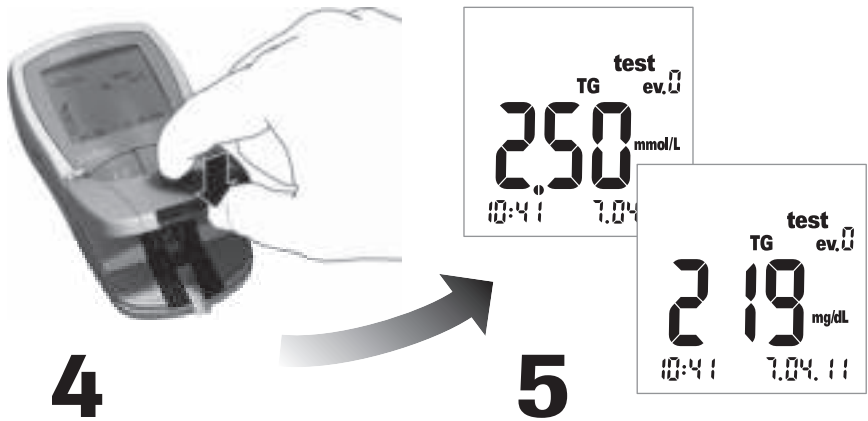Apply blood


Start measurement by closing the measurement chamber flap

Read the result

Test strips

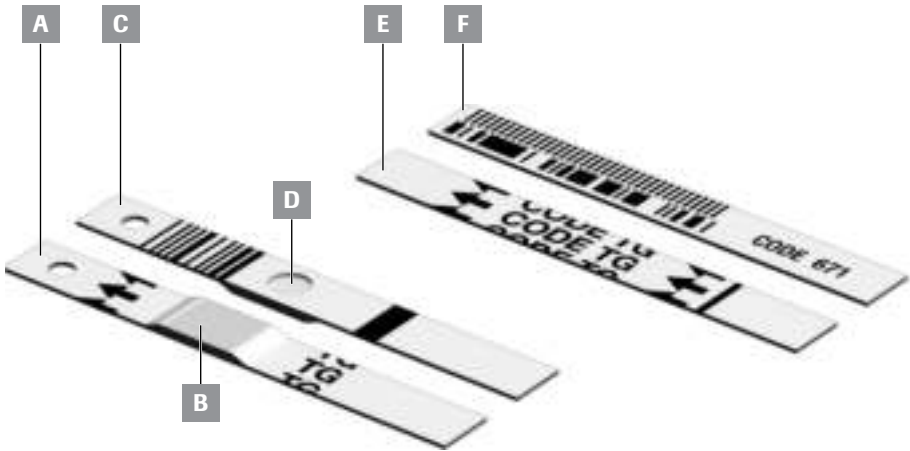Code strips

- A Test strip** (top of strip, TG in the example)
Contains the application area.
- B Sample application area**
Apply the sample here.
- C Test strip** (back of strip)
The imprinted barcode is used to identify the type and lot of the test strip.
- D Round window**
After performing a blood test, turn the strip over and visually check that the round window is evenly coloured.

- E Code strip** (top of strip, TG in the example)
Supplied with each test strip container.
- F Code strip** (back of strip)
The imprinted barcode contains lot-specific information which is read and stored in the meter.

Code strips

The code strip provides the meter with important information on the production-specific properties of the respective test strip lot. The code strip is always required when a new test strip container is opened and before the strips are used. The properties of this lot of test strips are then stored in the meter. The meter always stores the data of only **one** code strip per test parameter (i.e. a total of four codes at a time).

- There is a new code strip in each test strip container. **Code the meter with this code strip before performing the first measurement with a new lot of strips.** Once the code strip data is stored in the meter, you normally do not need the code strip again.
- However, we recommend keeping the code strip. Store the test strip container together with the code strip in the external packaging. That way you always have the code strip at hand if you should need to recode the meter (e.g., after a prolonged period without batteries).

Do not store the code strip in the test strip container

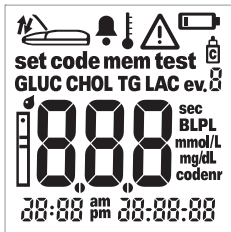
Store the code strip in the external packaging and **not** in the test strip container. If code strip and test strips are stored together, this may adversely affect the quality of both and may lead to error messages or even incorrect test results.

Powering on the meter

- 1 Place the meter on a level surface or hold it in your hand. Power on the meter by pressing the **On/Off** button ①.


Note: To power the instrument off again, press and hold the **On/Off** button ① until the meter powers off.

- 2 Check that all of the display elements are displayed properly. Results may be mis-read if a segment is missing.
- 3 Check the battery level. If the battery icon appears outside the display test, you will only be able to perform a few more tests.

When the display test is complete, the meter displays the most recent code stored. If no code has been stored in the meter so far, you will see the display shown on the left.

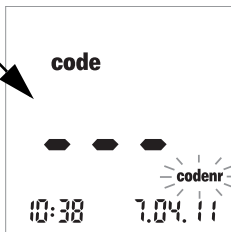The flashing icon *codenr* instructs you to insert a code strip.

Inserting the code strip

- 1 Hold the code strip with thumb and index finger in the white area located at the end of the strip. Do not touch the printed area (behind the black bar).
- 2 Insert the code strip smoothly into the test strip guide up to the stop in the direction of the printed arrows.

Withdraw it immediately. Leave the measurement chamber flap closed during this procedure.

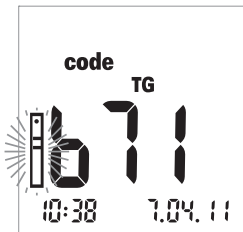
If the meter reads the barcode information correctly, successful coding is confirmed with a short beep (if beeper is activated).

The three-digit code number (which is also printed on the reverse side of the code strip and on the test strip container) appears in the display.

If the code number has a leading zero (e.g. 039), the meter might not display the leading zero.


If any problems occur during reading, an error message will be displayed (see *Troubleshooting*). In this case repeat reading the code strip after a few seconds.

The flashing test strip icon prompts you to insert the test strip.

Checking the stored code numbers

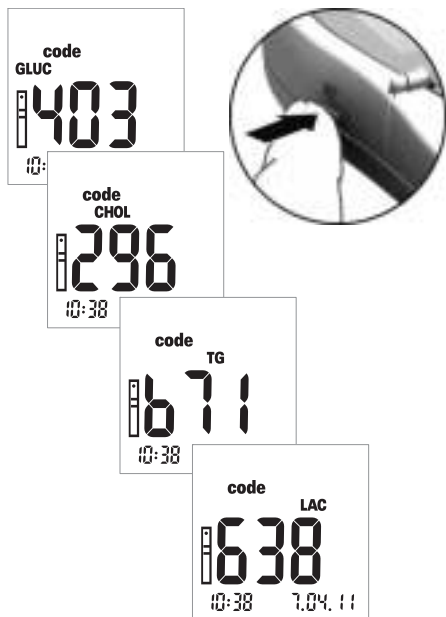You can display all of the stored code numbers and move back and forth between them.

- 1 If the last stored or used code number is displayed after powering on the meter, press the **Set** button to switch to the next stored code number. Each time you press the button, the display shows the code numbers of the test strips in succession GLUC > CHOL > TG > LAC (if already stored).

This display is for information purposes only. You do not need to select the code of the test parameter first before performing a test.

Note: You only need to code the meter the first time you use a new lot of test strips. If you insert the code strip every time you test, you may damage the barcode on the code strip. Then it may not work when you need it, for instance, after changing the batteries. See *Code strips* on page 35.

Sample material

Fresh capillary blood is used as sample material. You need a hanging blood drop to perform a test. For more information on sample material and use of heparinized blood, refer to the appropriate test strip package insert.

Protection against infection

Any object coming into contact with human blood is a potential source of infection.

Patient Self-Testers: There is a potential risk of infection if you let other people use your meter or lancing device. Do not share your meter or lancing device with any other person.

Healthcare Professionals:

- See safety message “Protection against infection” on page 10
 - See safety message “Potential risk of infection” on page 40
-

Testing by healthcare professionals

Potential risk of infection

Healthcare professionals using the Accutrend Plus meter to perform tests in more than one patient must be aware that any object coming into contact with human blood is a potential source of infection. Healthcare professionals must also be aware that any cross-contamination is a potential source of infection for patients. (See: Clinical and Laboratory Standards Institute: Protection of Laboratory Workers from Occupationally Acquired Infections; Approved Guideline - Third Edition; CLSI document M29-A3, 2005).

Any patient with an infection or suffering from an infectious disease and any patient who is a carrier of a multiresistant organism must be assigned his/her own meter. This applies also if a patient is suspected of being one of the above and, for as long as the suspicion exists, the meter must not be used to test any other patient.

- Follow the infection control procedures of your institution when handling blood testing equipment.
- Wear gloves. Change gloves between patients, even if patient dedicated testing devices and single-use, self-disabling lancing devices are used.
- If used for testing multiple patients, the Accutrend Plus meter must be properly cleaned and disinfected after every patient following the guidelines provided in this manual (see section *Cleaning* starting on page 71).
- The gloves worn during cleaning /disinfection should be removed and hands washed thoroughly with soap and water before performing the next patient test.
- Use an auto-disabling single-use lancing device for each patient. The lancing device must be intended for use by healthcare professionals in a multiple patient setting (for example, the Accu-Chek Safe-T-Pro Plus lancing device from Roche). Follow the manufacturer's instructions for use.
- Use capillary pipettes to apply blood outside the meter when performing a test on multiple patients if the test is suitable for the use of capillary pipettes (see test specific package insert).
- If using capillary pipettes, avoid damaging the sample application area with the edge of the pipette during sample application.

- When performing several measurements on the same patient, apply blood outside the meter (see page 48).
- Dispose of used lancets, capillary pipettes, and test strips in line with the disposal policy of your hospital, institute or medical practice.
- Follow all other locally applicable health and safety regulations.

Preparing to perform a test

- 1 Get the test strip container for the required test (e.g. triglycerides).
- 2 Check the expiration date of the test strip. Always use test strips **before** their expiration date has passed.
- 3 Code the meter with the code strip belonging to these test strips (unless the meter has already been coded with this code strip).

Damage to strips by environmental influences

Environmental factors (e.g. humidity and light) may damage the test strips and lead to incorrect results or error messages. Do not remove the test strips from the packaging until immediately before performing a test.

- 4 Prepare the lancing device.

Self-testers: Do not lance yourself before being instructed to do so in the course of this description.

Patient Self-Testers

We recommend the Accu-Chek Multiclix lancing device for self-testing.

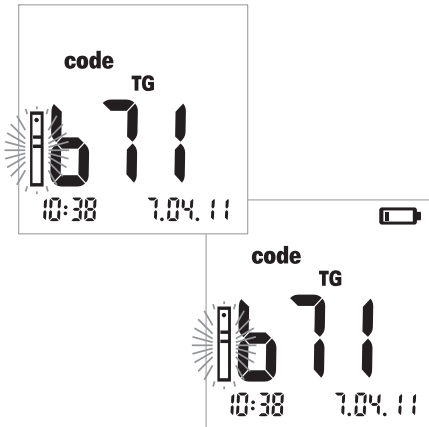
Testing by healthcare professionals

Use an auto-disabling single-use lancing device for each patient. The lancing device must be intended for use by healthcare professionals in a multiple patient setting (for example, the Accu-Chek Safe-T-Pro Plus lancing device from Roche). Follow the manufacturer's instructions for use.

Performing a test


After it is powered on and coded, the meter expects you to insert a test strip. Using the barcode on the back of the strip, the meter detects which test parameter is to be measured and which code is needed for this test strip. If the code strip has not been read yet, an error message is displayed after inserting the test strip.

- 1 Check the following before performing the test:
 - Check that all of the display elements are displayed properly. Results may be mis-read if a segment is missing.
 - Are date and time correct? If you want the test result to be saved with date/time information, enter the correct settings (see description starting on page 26).
 - Check the battery level. When the battery icon appears for the first time (outside the display test), you will only be able to perform a few more tests. Replace the batteries as soon as possible (see description starting on page 20).
- 2 Take a test strip from the test strip container.

NOTICE

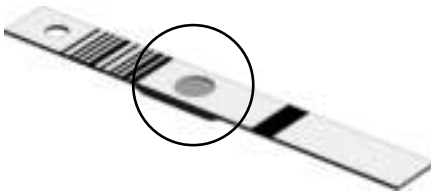Close the container immediately after opening to protect the desiccant, otherwise the test strips may become unusable before their expiration date. Prevent liquids from entering the test strip container. Do not mix up the caps of different test strip containers!


For cholesterol and glucose test strips:

Before performing a test, please check the colour in the round window on the back of the test strip for discoloration. If you detect discoloration, this test strip is unusable. Refer to the package insert for details.

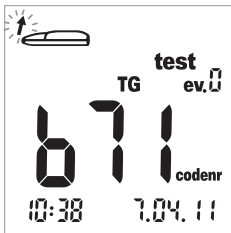
- 3 Hold the test strip with thumb and index finger so that the the lettering (for the test parameter) and the sample application area are facing upwards.
- 4 Insert the test strip into the test strip guide up to the stop. When the test strip reaches the correct position, you will hear two beeps (short – long; if beeper is activated).

A flashing arrow now instructs you to open the measurement chamber flap to apply the blood. (For a description of how you can add additional information to test results, see *Flagging test results* on page 52.)

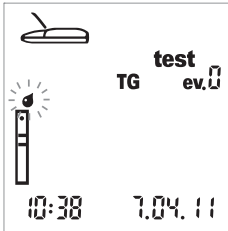- 5 Open the measurement chamber flap. The flap firmly locks into place when it reaches a vertical position.

The flashing drop icon (above the strip icon) now instructs you to apply the blood.

Blood can be applied to the sample application area of the test strip either:

- with the test strip inside the meter (when used by a single user) or
- with the test strip outside the meter (when tests are performed in the professional sector with, for example, heparinized capillary pipettes).

See description on pages 47 and 48.

Recommendations for the collection and measurement of capillary blood

To obtain a suitable blood drop:

- Wash your hands with warm water.
- Ensure hands are warm and dry before lancing. Massage your fingertip.
- After lancing the finger, discard the first drop of blood. Then try to obtain a sufficiently large hanging drop without excessive pressing or squeezing.
- We recommend obtaining the capillary blood from the side of the fingertip as this part is the least sensitive to pain.

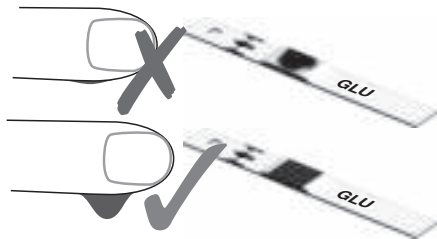
Potential risk of incorrect results due to fatty residues

When performing a **triglycerides** test: Nearly all creams and many soap products (e.g. shower gel and shampoo) contain fat. If residues of these fatty substances on the skin come into contact with the test strips, the measurement will be incorrect. It is therefore particularly important to wash your hands very thoroughly and rinse them with plenty of clear water when performing this test.


Blood collection

- 6 Lance the outer side of the fingertip using the lancing device to obtain a **large hanging drop of blood**.

Blood application inside the meter

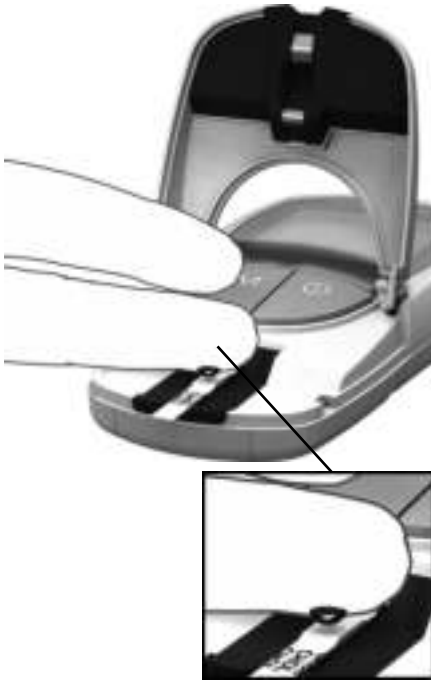- 7 Apply a large hanging drop of blood directly from the finger to the yellow sample application area of the test strip.

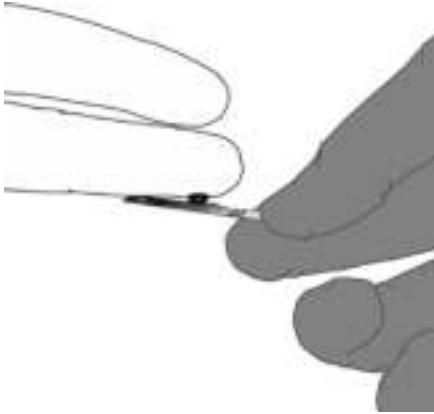
Do not touch the sample application area with the finger to avoid damaging the application area.

Potential risk of incorrect result

Apply the drop of blood to the test strip **immediately** after lancing the fingertip. Ensure that the blood drop is sufficiently large and hanging. Re-dosing is not possible. Blood which is applied later may lead to an incorrect result.

Blood application outside the meter

- 8 Withdraw the test strip after opening the flap (and leave it open).
- 9 Apply a large hanging drop of blood directly from the finger to the yellow sample application area of the test strip. Heparinized capillary pipettes can be used to apply blood. Refer to the insert of the respective test strip.

Do not touch the sample application area with the finger or the pipette to avoid damaging the application area.

- 10 With the measurement chamber flap open, insert the test strip back into the meter.

Checking the applied blood:

A

Make sure the application area is completely covered with blood (example A). Otherwise incorrect results may be obtained.

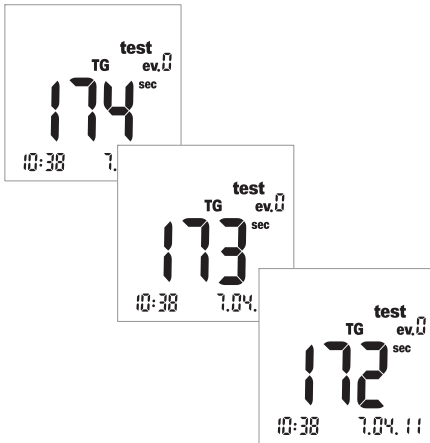
B

If too little blood is applied (example B), **do not try to spread it or to apply more blood** as this may lead to an incorrect result. Repeat the test with a new test strip.

Starting the measurement

- 11** Close the measurement chamber flap. This automatically starts the test.

Measurement starts. The time needed to run a test depends on the test parameter. This time is shown in the display and counted down to "0". The measurement times of the test parameters are:

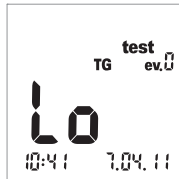
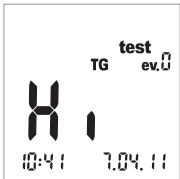
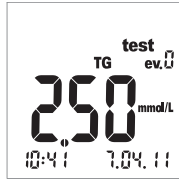
- Glucose: 12 seconds
- Cholesterol: 180 seconds
- Triglycerides: max. 174 seconds
- Lactate: 60 seconds

For Glucose, Cholesterol, Lactate only:

The meter signals that it has reached the last four seconds of measurement time by emitting one short beep per second (if beeper is activated). Completion of the test and the subsequent display of results is indicated by a longer beep.

Note: If you have set the beeper to *OFF*, the meter will not beep when the result is displayed.



Display of results

When measurement is finished, the result is displayed. If you are measuring glucose, please carry out the plausibility check described below.

Test results which fall outside the measurement range are displayed as *Hi* (above the measurement range) or *Lo* (below the measurement range).

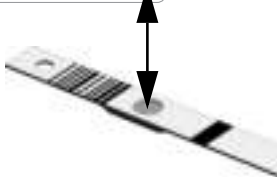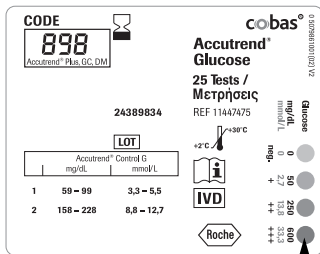When the result is displayed, the *ev. 0* event (no event) is displayed. For an explanation on how to attach additional information to a result, see section *Flagging test results* on page 52.

Always carry out a plausibility check after performing a glucose test.

- 12 Remove the test strip and turn it so that the back is facing upwards.
- 13 Compare the colour of the round window on the back of the test strip with the colour field scale on the label of the test strip container.

The colour of the round window must approximately match the colour assigned to your measurement result. If the colours do not match very well, perform a control test. Refer to the Accutrend Glucose test strip package insert for more information.

If the displayed result (particularly when blood glucose is measured) does not match the way you feel or seems to be unusually high or low, perform a control test using control solution and a new test strip (see description starting on page 55). If this check confirms proper functioning of the meter, please read again the instructions on how to perform a test (starting on page 31). Perform another test using a new test strip. If this result does not seem plausible either, please consult your doctor.

If you do not want to flag the measurement result with an event or as a control test, the measurement is now complete.

- 14** Open the measurement chamber flap and remove the test strip.
- 15** Keep the **On/Off** button  pressed until the meter powers off.

Clean the meter if necessary (see description starting on page 71).

Disposal of used items after testing with blood

Patient Self-Testers

Discard the used test strip with your regular household waste.

Dispose of used lancets carefully (e.g. use a sturdy sharps container with lid) so that the needles cannot cause injury to yourself or to others.

Healthcare Professionals

-
- Dispose of used consumables in line with the disposal policy of your hospital, institute or medical practice.
 - See precautions in section *Testing by healthcare professionals* starting on page 40.
 - See safety message *Protection against infection in a healthcare environment* on page 10.
-

Flagging test results

Flags can provide additional information about the conditions when you performed the test (e.g., meals, sport, illness, etc). In addition, you can flag tests carried out with control solutions as control tests. You can flag a test result at different points in time:

- At the start of the test after inserting the test strip.
- When the test result is displayed.

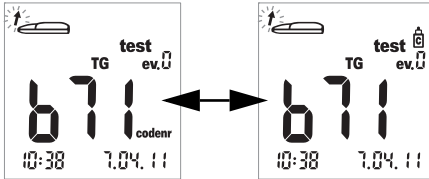
You **cannot** flag a result while a test is running.

Flagging the measurement as a control test:

You flag measurements as control tests when they have been performed using a control solution (as described in the following chapter).



- 16 Press the **M** button (after inserting the test strip or when the result is displayed) to flag the measurement as a control test (bottle icon). Pressing the **M** button again turns the flag off.

Assigning additional information to a measurement:

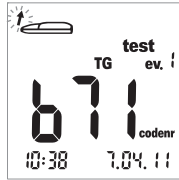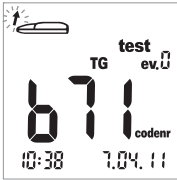Optionally, you can assign one of nine different events to a test. The event “0” represents “no flagging”. If used, you must define the events and their respective event numbers yourself. Make sure that your assignments are unique and reproducible.

- 17** Press the **Set** button (after inserting the test strip or when the result is displayed) to flag the measurement with an event. Pressing the **Set** button increases the displayed event number in increments of one. After event “9”, the event display is reset to “0”.

This page intentionally left blank.

6 Control Testing

Perform a control test regularly to check that the meter and test strips are working properly and that you are performing the test correctly. For each test parameter a separate control solution is available. Always perform a control test in the following situations:

- When you open a new test strip container.
- After replacing the batteries.
- After cleaning the meter.
- When you doubt that the measured values are correct.

A control test is performed in the same way as a regular measurement with the difference that control solutions are used instead of blood.

When performing a control test, operate the meter only within the acceptable temperature ranges for the control solutions. This range is test dependent (refer to the package insert of the respective control solution):

- For Accutrend Control CH1 (Cholesterol): 18-30 °C
- For Accutrend Control G or G2 (Glucose)*: 18-32 °C
- For Accutrend Control TG1 (Triglyceride): 18-30 °C
- For BM-Control Lactate: 15-35 °C

Healthcare Professionals: Refer to the package insert. Follow the applicable government regulations and local guidelines for quality control.

What you need:

- Accutrend Plus meter
- Test strip for selected test parameter with matching code strip:
 - Accutrend Glucose
 - Accutrend Triglycerides
 - Accutrend Cholesterol
 - BM-Lactate
- Control solution for selected test:
 - Accutrend Control G or G2 (Glucose)*
 - Accutrend Control TG1 (Triglyceride)
 - Accutrend Control CH1 (Cholesterol)
 - BM-Control Lactate

* Note: Availability of G or G2 may vary from country to country.

Preparing to perform a control test

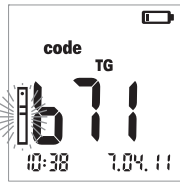
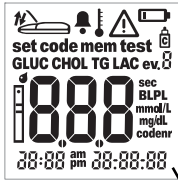
- 1 Get the test strip container for the required test (e.g. triglycerides).
- 2 Code the meter with the code strip belonging to these test strips (unless the meter has already been coded with this code strip).
- 3 Get the correct control solution for the test strips.

The following description assumes that the meter has already been coded for this lot of test strips. If this is not the case, see *Inserting the code strip* on page 37.

Performing a control test

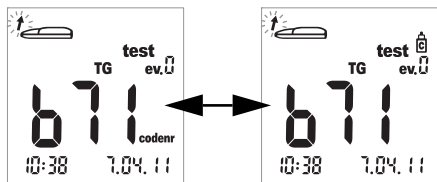
- 1 Place the meter on a level surface or hold it in your hand. Power on the meter by pressing the **On/Off** button .

- 2 Check the following before performing the control test:
 - Check that all of the display elements are displayed properly. Results may be mis-read if a segment is missing.
 - Are date and time correct? If you want the test result to be saved with date/time information, enter the correct settings (see description starting on page 26).
 - Check the battery level. When the battery icon appears for the first time (outside the display test), you will only be able to perform a few more tests. Replace the batteries as soon as possible (see description starting on page 20).
- 3 Take a test strip from the test strip container. Close the container again immediately after removing the strip, otherwise the test strips may become unusable before their expiration date.
- 4 Hold the test strip with thumb and index finger so that the printed measurement parameter is facing upwards.
- 5 Insert the test strip into the test strip guide up to the stop. When the test strip reaches the correct position, you will hear two beeps (short – long; if beeper is activated).

- 6 Press the **M** button to flag the test result as a control test (bottle icon).

- 7 Open the measurement chamber flap. The flap firmly locks into place when it reaches a vertical position.

Applying the control solution

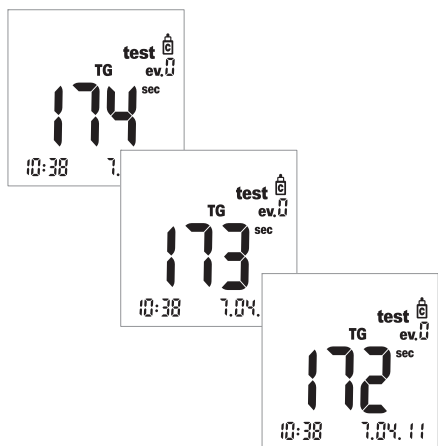
- 8 Apply a large hanging drop of control solution directly from the bottle to the test strip. Ensure that neither the bottle nor the finger touches the application area. The application area must be completely covered.

You can also apply the control solution outside the meter. The procedure is the same as when you apply a blood sample. Refer to the description starting on page 48.

Starting the measurement

- 9 Close the measurement chamber flap. This automatically starts the test.

Measurement starts. The time needed to run a test depends on the test parameter. This time is shown in the display and counted down to “0”. The measurement times of the test parameters are:

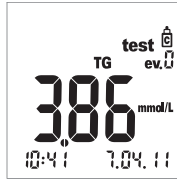
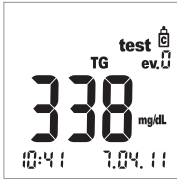
- Glucose: 12 seconds
- Cholesterol: 180 seconds
- Triglycerides: max. 174 seconds
- Lactate: 60 seconds

For Glucose, Cholesterol, Lactate only:

The meter signals that it has reached the last four seconds of measurement time by emitting one short beep per second (if beeper is activated). Completion of the test and the subsequent display of results is indicated by a longer beep.

Note: If you have set the beeper to *OFF*, the meter will not beep when the result is displayed.

Display of results

When measurement is finished, the result is displayed.

Check whether the displayed result falls within the acceptable range.

Target values which should be obtained with control solutions can either be found on the labels or on the package inserts for the test strips or control solutions. If the value is outside the range, repeat the control test. If the second result is again outside this range, please contact your local customer support and service centre (see page 81).

Disposal of used items after control testing

Patient Self-Testers

Discard the used test strip with your regular household waste. For more details, refer to the package insert.

Healthcare Professionals

- Refer to the package insert.
- Dispose of used consumables in line with the disposal policy of your hospital, institute or medical practice.

This page intentionally left blank.

7 Memory

The Accutrend Plus meter has four memory areas, each of which can be used to save up to 100 test results together with date, time and flags.

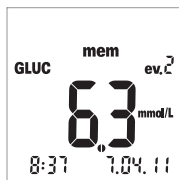
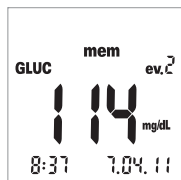
If you use the meter without setting the date or time, no time/date information will be saved with the test results.

Reviewing results in memory


- 1 Power on the meter directly in memory mode by pressing the **M** button, or by pressing this button when the meter is in coding mode.

You can exit memory mode by pressing the **On/Off** button .

After the usual display test has been performed when powering on the meter, the last saved test result is displayed.

The icon *mem* appears in the display and indicates that the meter is in memory mode. The date and time on the display match the time at which the measurement was performed, and not the current time. In memory mode, the colon between the hour and the minutes of the time display does not flash.



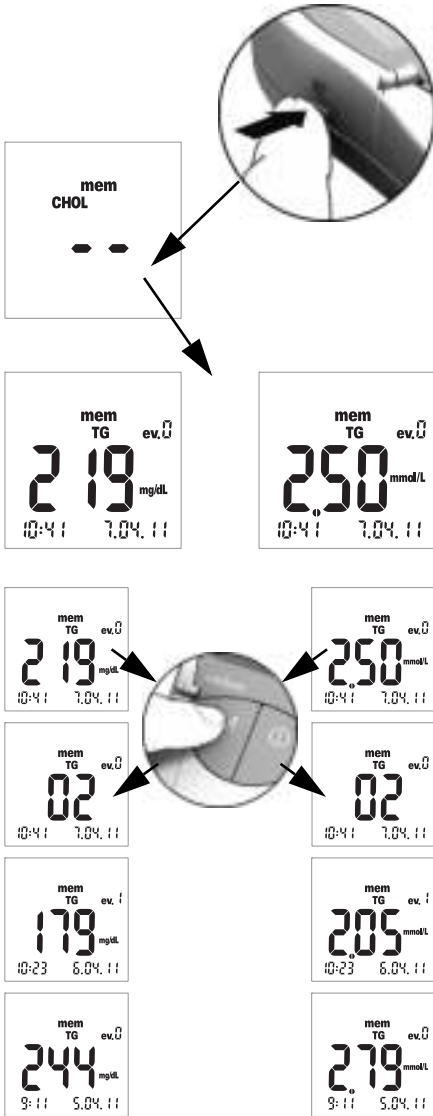
- 2 Press the **Set** button to switch between the memory areas of the four test parameters. Each time you press the **Set** button, the last stored test result of each memory area is displayed (GLUC > CHOL > TG > LAC > GLUC >...). If the area contains no stored results, two dashes are displayed.

To navigate within one memory area, use the **M** button.

- 3 Press the **M** button to show the next oldest result within one memory area. When you press the button, the place of the result within the memory is displayed (here: 02); the actual value is not shown until you release the button. When you hold down the **M** button, older memory places (02 ... 03 ... 04 ... and so on) continually appear in the display until you release the button. Then the result in the last displayed memory place is shown.

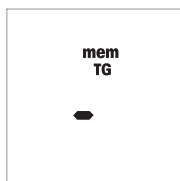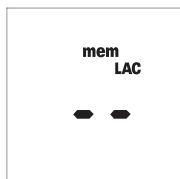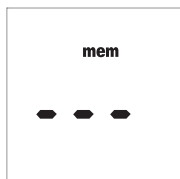
Special cases:

- If memory is empty because the meter has not performed any tests so far, three dashes (– – –) will appear in the display. If memory is empty because of a manual delete process, two dashes (– –) will appear in the display.
- If the entire memory area of a test parameter is empty, two dashes (– –) will appear in the display.
- If the meter is not able to save a result to memory correctly*, a single dash (–) is displayed in memory mode.
- If you delete the last result manually (see page 67), this entry is completely cleared from memory. Nothing, not even a dash, will be displayed.

* This is very unlikely to happen. For example, it could occur if the meter suddenly loses power (due to a drained battery) while the result is being saved to memory.

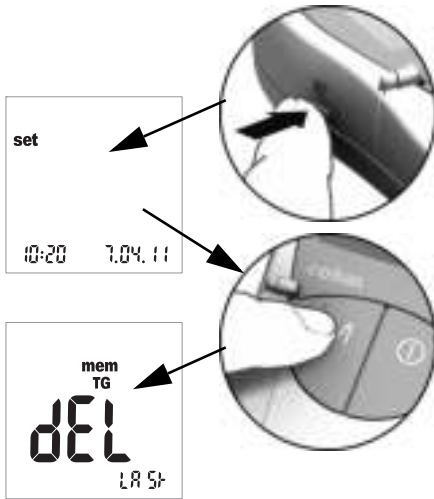Deleting results from memory


- 1** Power on the meter by pressing the **Set** button (on the left side of the meter).
The date and time as well as the icon *set* now appear on the display.
- 2** Now press the **M** button to display the options for deleting results from memory. You can exit delete mode at any time by pressing the **On/Off** button **ⓘ**.
All the following descriptions start with this display (*dEL LAST*).

You can choose between three different deletion options:

- You can delete the last result.
- You can delete all results of a single test parameter at once.
- You can delete all results of all test parameters at once.

It is not possible to delete a single result (unless it is the last result).

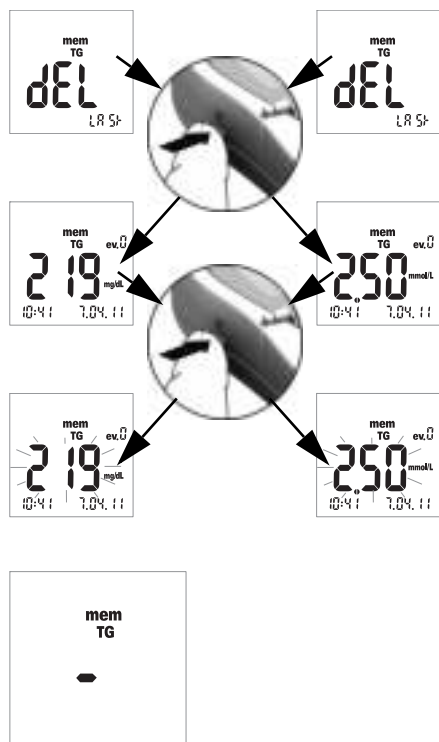
Deleting the last result

- To delete the **last result** (*dEL LAST* display), press the **Set** button. The related test parameter is displayed.
- To delete the displayed result from memory, **press and hold down** the **Set** button **for more than 3 seconds**. While you hold down the button, the selected result flashes. At the same time you will hear a short beep once a second.

After 3 seconds the result is deleted from memory and a dash appears briefly to confirm deletion (see illustration on the left). Once you release the **Set** button, you automatically exit the delete mode.

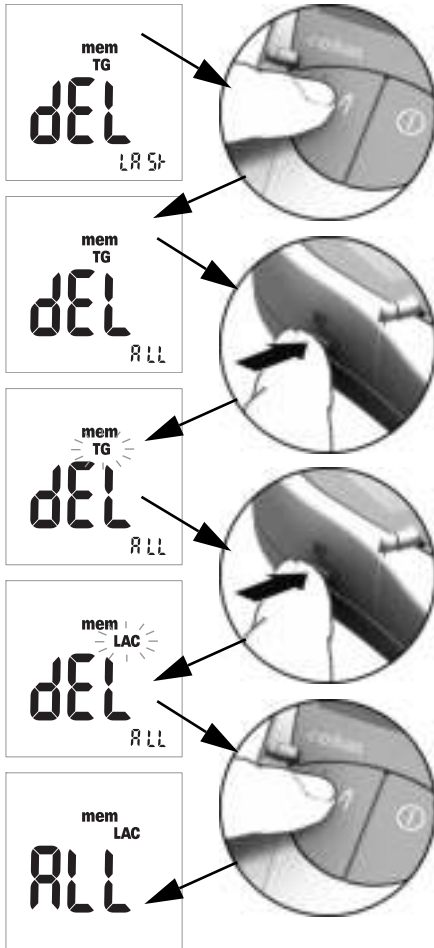Deleting all results of a test parameter

5 If you want to delete **all results** of a test parameter, press the **M** button. The display shows *dEL - ALL*, and the currently selected test parameter (here: TG) is displayed.

6 If the displayed test parameter is not the parameter for which you want to delete all results from memory, press the **Set** button.

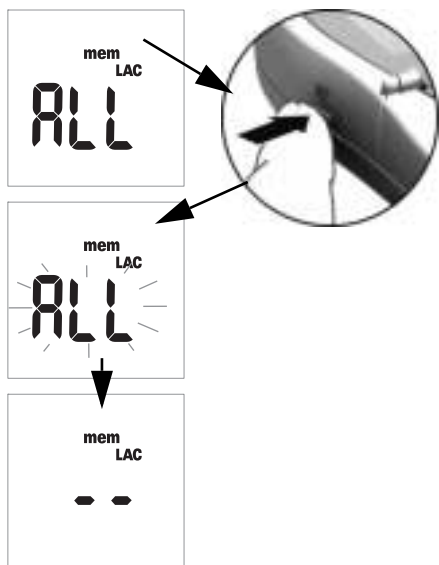
The currently selected test parameter flashes.

7 Press the **Set** button again to select a different test parameter.

Each time you press the button, the next test parameter appears in the display (GLUC > CHOL > TG > LAC ...).

8 To select the displayed test parameter for deletion, press the **M** button.

The display switches to *ALL*.


- 9 To delete the results of the currently displayed test parameter from memory, **press and hold down** the **Set** button for **more than 3 seconds**.

While you hold down the button, *ALL* flashes in the display. At the same time you will hear a short beep once a second.

After 3 seconds all results for the selected test parameter are deleted; this is confirmed by the display as shown on the left (- -). Once you release the **Set** button, you automatically exit the delete mode.

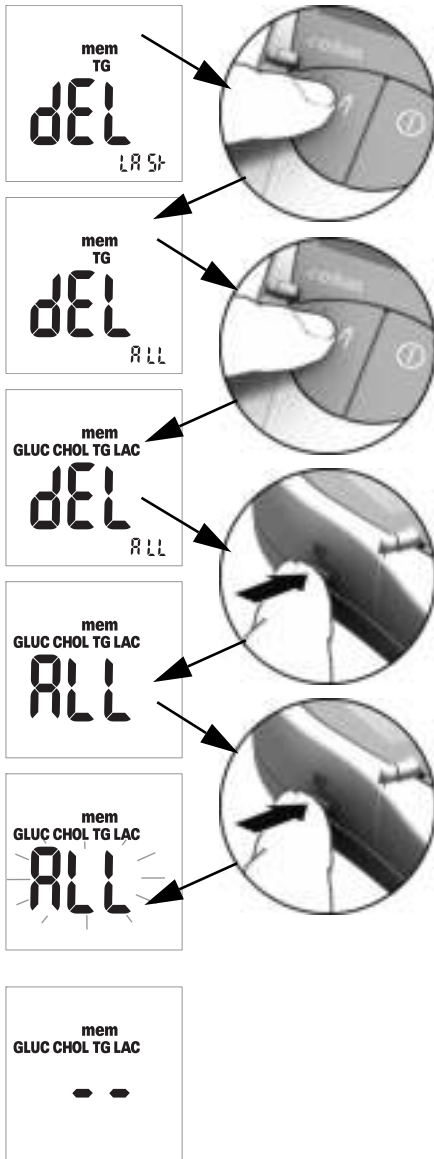Deleting all results

10 If you want to delete **all saved results for all test parameters**, press the **M** button twice.
The display shows *dEL - ALL*, and all test parameters are shown.

11 Press the **Set** button to select this delete mode.

The display switches to *ALL*.

12 To delete all results from memory, **press and hold down** the **Set** button for **more than 3 seconds**.

While you hold down the button, *ALL* flashes in the display. At the same time you will hear a short beep once a second.

After 3 seconds all results are deleted; this is confirmed by the display as shown on the left (- -). Once you release the **Set** button, you automatically exit the delete mode.

8 Cleaning

A clean optical measuring system is a basic prerequisite for determining precise results. Therefore clean the meter regularly and immediately after it becomes dirty. Always power off the meter before cleaning it.

Use only the following items for cleaning:

- Ordinary lint-free cotton buds
- Ordinary lint-free tissues
- Ordinary disinfecting tissues

Recommended cleaning/disinfecting solutions

Mild soap suds as well as 70 % ethanol or isopropyl alcohol are suitable for cleaning.

Healthcare Professionals: When used in the professional sector (e.g. doctor's office), a mixture consisting of 1-propanol (400 mg/g), 2-propanol (200 mg/g) and glutaraldehyde (1.0 mg/g)¹ is recommended.

NOTICE

Damage to the instrument through moisture

- Ensure that no liquid enters the meter. If moisture enters the meter, it may cause malfunction of the instrument.
 - Do not use any disinfectant sprays or tissues or cotton buds which are dripping wet as the liquid may enter the meter and damage it.
-

Cleaning the exterior (meter housing)

- Ensure the meter is powered off.
- Wipe the outside of the meter with a lightly moistened, lint-free cotton cloth.

1. Sold in some countries under the name "Bacillol Plus."

Cleaning the measurement chamber cover/test strip guide

- 1 Open the measurement chamber flap.
- 2 Remove the measurement chamber cover (including the test strip guide) by slightly pushing it to the centre of the meter and then pulling it upwards.

- 3 If necessary, you can rinse the measurement chamber cover/test strip guide (separately from the meter) under warm running water. Dry the measurement chamber cover/test strip guide with a clean cloth.

Cleaning the optical measuring system

- 4** Clean the easily accessible areas of the optical measuring system with a lint-free pad or a moistened cotton bud. Make sure that no liquid enters the meter. Do not insert any objects into the meter.
- 5** Allow the meter to dry thoroughly.

- 6** Do not put the measurement chamber cover/test strip guide back into the meter until it is completely dry. Press it down gently at the front until you feel it click into place.
- 7** Close the measurement chamber flap.
The meter is now ready for operation again. Perform a control test (see description starting on page 55).

This page intentionally left blank.

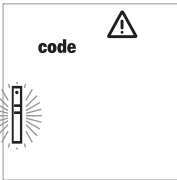
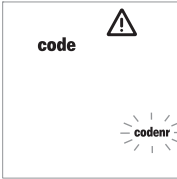
9 Troubleshooting

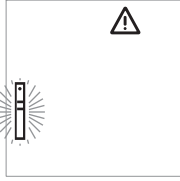
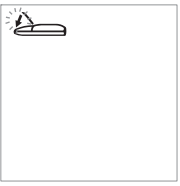
In certain circumstances, the meter may display an error message. A troubleshooting table follows that will help you when the meter is not performing as expected. Most concerns can be resolved quickly by referring to this table for help.

Take the following steps when the meter displays an error message :

- Find the displayed message or condition in the Troubleshooting Table.
- Take the action suggested under the column headed *Description/Solution*.

If the problem persists, please contact your local customer support and service centre (see page 81).

Error/Display	Description/Solution
<p>Error: Unknown barcode</p> 	<p>The barcode of test strips or code strips could not be read.</p> <p>Solution</p> <ul style="list-style-type: none"> ■ All strips: Remove the strip and check the barcode for soiling. ■ Code strips: Repeat the procedure with the code strip. Remember to insert the strip and remove it immediately. ■ Test strips: Repeat the procedure with a new strip. ■ Check whether there may be interference caused by electromagnetic fields in the vicinity of the meter. Move meter away from such sources (e.g., X-ray equipment).
<p>Error: Strip code does not match stored code</p> 	<p>The test strip belongs to a different test strip lot than the lot which was last coded for this parameter.</p> <p>Solution</p> <ul style="list-style-type: none"> ■ Remove the strip. Repeat the test using a test strip from a lot which matches the last lot coded into the meter for this parameter. ■ Code the meter with the code strip corresponding to the test strips being used.

Error/Display	Description/Solution
<p>Error: Test strip is unusable</p> 	<p>The test strip has already been used or is dirty.</p> <p>Solution</p> <ul style="list-style-type: none"> ■ Remove and discard the strip. Repeat the test using a new test strip.
<p>Error: Temperature</p> 	<p>The meter is too hot or cold, or the ambient temperature does not fall within the acceptable range for the selected test</p> <p>Solution</p> <ul style="list-style-type: none"> ■ Move the meter to an environment which has an appropriate temperature and repeat the test after several minutes. Do not artificially heat or cool the meter by any means.
<p>Error: Flap open (after meter has been powered on)</p> 	<p>The measurement chamber flap is not closed properly.</p> <p>Solution</p> <ul style="list-style-type: none"> ■ Close the measurement chamber flap.
<p>Error: Flap open (during measurement)</p> 	<p>The measurement chamber flap is not closed properly.</p> <p>Solution</p> <ul style="list-style-type: none"> ■ Close the measurement chamber flap.

Error/Display	Description/Solution
Error: Flap not closed after blood application	The measurement chamber flap was not closed properly after blood application; the test did not start in time.

Solution

- Remove the strip and repeat the test using a new test strip.

Error: Internal error (example 142)	The meter has detected an internal error.
--	---

Solution

- Power the meter off and on again. If the error still persists, the meter is defective. Please call your local customer support and service centre (see page 81).

This page intentionally left blank.

10 Product Specifications

Operating conditions and technical data

Temperature range for measurements with patient samples	Dependent on test parameter: <ul style="list-style-type: none"> ▪ For cholesterol: 18-35 °C ▪ For glucose: 18-35 °C ▪ For triglycerides: 18-30 °C ▪ For lactate: 15-35 °C
Temperature range for measurements with control solutions	Dependent on test parameter: <ul style="list-style-type: none"> ▪ For cholesterol: 18-30 °C ▪ For glucose: 18-32 °C ▪ For triglycerides: 18-30 °C ▪ For lactate: 15-35 °C
Relative humidity	10 to 85%
Measurement range	Blood glucose: 20–600 mg/dL (1.1–33.3 mmol/L) Cholesterol: 150–300 mg/dL (3.88–7.76 mmol/L) Triglycerides: 70–600 mg/dL (0.80–6.86 mmol/L) Lactate: 0.8–21.7 mmol/L (blood value), 0.7–26 mmol/L (plasma value)
Memory	100 test results, optionally with date, time, and additional information per test parameter
Battery operation	4 x AAA 1.5 V alkaline manganese batteries
Number of tests per set of batteries	More than 1,000 tests (with new batteries)
Safety class	III
Dimensions	154 x 81 x 30 mm
Weight	Approx. 140 g

Sample material

Sample type	Fresh capillary blood (for detailed information see package insert)
Sample size	A hanging drop of blood
Interactions	Refer to test strip package insert

Storage and transport conditions

Temperature range	-25 °C to +70 °C
Relative humidity	10 to 85 % (no condensation)

Ordering

Please contact your specialist supplier.

Item	Description
Accutrend Glucose 25	25 test strips for determining blood glucose
Accutrend Control G	Control set to be used with Accutrend Glucose test strips
Accutrend Control G2	Control set to be used with Accutrend Glucose test strips
Accutrend Cholesterol 25	25 test strips for determining cholesterol in blood
Accutrend Cholesterol 5	5 test strips for determining cholesterol in blood
Accutrend Control CH1	Control solution to be used with Accutrend Cholesterol test strips
Accutrend Triglycerides 25	25 test strips for determining triglycerides in blood
Accutrend Control TG1	Control solution to be used with Accutrend Triglycerides test strips
BM-Lactate 25	25 test strips for determining lactate in blood
BM-Control Lactate	Control solution to be used with BM-Lactate test strips
Accutrend Plus meter	

Note: Not all items are available in all countries.

Product limitations

Refer to the package inserts accompanying the test strips or control solutions for detailed product data and limitations.

Contact Roche

For all questions about the Accutrend Plus system that are not answered in this manual, contact your Roche representative. If you do not already have contact details, a list of Roche Diagnostics offices is provided below or visit our website at www.Roche.com. Select “Roche Worldwide” at the top of the page and then select your country to find the appropriate local office contact information.

Australia	Roche Diagnostics Australia Pty Ltd. ABN 29 003 001 205 31 Victoria Avenue Castle Hill, NSW, 2154 Telephone 02-9860 2222 or 1800 645 619
Canada	Roche Diagnostics 201 Boulevard Armand-Frappier Laval, Québec H7V 4A2 Telephone 1-450-686-7050
New Zealand	Roche Diagnostics N.Z. Ltd, 15 Rakino Way, Box 62-089 Mt. Wellington, Auckland Telephone +64 9 276 4157 Toll Free 0800 652 634
South Africa	Roche Products (Pty) Ltd. South Africa Diagnostics Division 9, Will Scarlet Road / Ferndale P.O. Box 1927, Randburg 2125 Telephone +27 -11 504-4600
United Kingdom	Roche Diagnostics Ltd. Charles Avenue Burgess Hill, RH15 9RY UK Freephone: 0808 100 99 98 Rep. of Ireland Freephone: 1800 509 586
USA	Roche Diagnostics 9115 Hague Road Indianapolis, IN 46256 Service information For technical assistance or information, contact the Accu-Chek® Customer Care Service Center at 1-800-440-3638, 24 hours a day, 365 days a year.

This page intentionally left blank.

Index

A

Addresses (Roche) 81

B

Batteries

Inserting 20–21

Battery compartment cover 15

Beeper 29

Blood application 47–48

checking 48

Blood collection 46

Blood parameters 7

Buttons

M button 15

On/Off button 15

Set button 15

C

Capillary blood 39

Capillary blood

recommendations for collection 45

Cleaning 71–73

cleaning agents 71

meter components 71–72

optical measuring system 73

Code

switching the display 38

Code strips 35

inserting 37

overview 34

Contents (packaging) 8

Control solution 55

applying 59

Control testing 55–61

display of results 61

flagging 52

performing 56–61

preparing 56

starting the measurement 60

what you need 55

D

Date 26–27

Date format 26

Display 15

icons 16–17

Display of results 50–51

Display test 21, 36

Disposal 10–11, 51

E

Error messages 75–77

Event 53

I

Icons 16–17

Infection (protection against) 10, 40

Infrared window 15

L

Lactate display 30

Lancing device 31, 42

LED (light-emitting diode) 8

M

- M button 15
- Measurement chamber cover 15
- Measurement chamber flap 15
- Memory 63–70
 - deleting test results 66–70
 - displaying test results 63–65
- Meter
 - Overview 14–18
 - powering on 36

O

- On/Off button 15
- Operating conditions 12, 79
- Ordering information 80
- Overview
 - Meter elements 14–18

P

- Powering on 36
 - in memory mode 63
 - in set mode 24
- Putting the meter into operation 19–21

Q

- Quality control
 - functions 13

R

- Reflectance photometry 8
- Relative humidity 12
- Roche (Address) 81

S

- Safety
 - Protection against infection 10, 40
- Sample material 39, 79
- Set button 15
- Settings
 - date 26–27
 - date format 26
 - general procedure 24–25
 - lactate display 30
 - overview 23
 - time 28
 - time format 28
- Storage conditions 80
- Symbols
 - on packaging and identification plate 4

T

- Technical data 79
- Temperature range 12
- Test 31–53
 - brief overview 33
 - display of results 50–51
 - flagging 52–53
 - performing 43–53
 - plausibility check 50
 - preparing 41–42
 - starting 49
 - what you need 31
- Test principle 8
- Test strip guide 15
- Test strips
 - inserting 44, 57
 - overview 34
- Time 28
- Time format 28
- Transport conditions 80
- Troubleshooting 75–77

ACCU-CHEK, ACCU-CHEK MULTICLIX, ACCUTREND, COBAS,
and SAFE-T-PRO are trademarks of Roche.

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Roche Diagnostics GmbH