
# Good4You<sup>®</sup>

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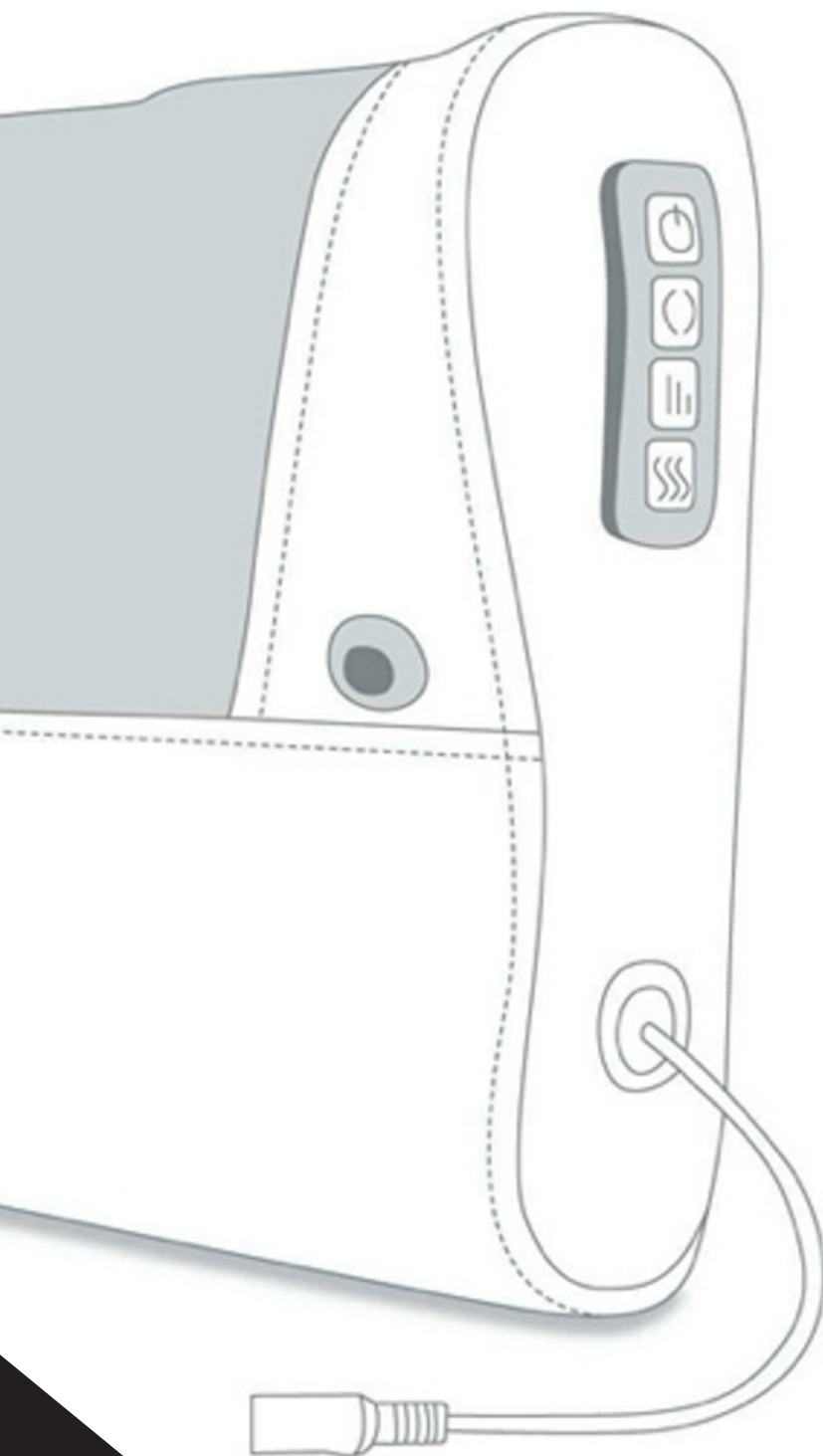## Introduction

Thank you for choosing our product! Please always follow the safety instructions for safe use. Should there be any questions about the product please feel free to contact us via email: [info@good-4you.nl](mailto:info@good-4you.nl)

## Warnings and Safety Instructions



Please make sure the packaging material does not get into the hands of children. They may suffocate! Read the operating instructions carefully before use and keep them in a safe place!


- This device is intended for private use only. Commercial or professional use voids the warranty. The supplier cannot be held liable if the device is used for anything other than what is indicated in these operating instructions.
- If the device is not working properly, has been dropped or damaged, do not use it any longer.
- Do not disassemble, repair, or restore the device.
- Never leave the device unattended when it is connected to power. When not in use, turn the device off and unplug the power cord.
- Do not pull the cord to unplug the device.
- Do not move the appliance by pulling the cord or twist the cord.
- Periodically check both power cord and plug for signs of damage or wear. If there is any, have repairs done only by qualified professionals.
- Stay awake while using the unit.
- This product is not suitable for children and should be kept out of reach of children.
- Do not use the appliance if you are pregnant.
- Do not use the device if you suffer from the following health problems: sensitivity to heat, blood circulation problems, varicose veins, fever, bruises or bleeding, wounds or other skin damage, skin irritation, inflammation of muscles or tendons or any injury.
- Consult a physician before using the device if you have a pacemaker, artificial joints, or electronic implants, have unexplained pain or are undergoing medical treatment.
- If you experience pain or discomfort while using this device, stop using it immediately and consult a physician.
- To avoid excessive strain on muscles and nerves, do not exceed the recommended massage time of 15 minutes.
- Cleaning the device: Before cleaning, make sure the device is not connected to electricity. Clean it with a slightly damp cloth and then dry thoroughly with a tea towel. Under no circumstances should the device be connected to power before it is completely dry.
- Do not cover the device when it is on. Never place a blanket or pillow against it. This can cause a fire hazard, electric shock, personal injury, or damage to the appliance.
- Do not use oil or lotion on the skin before or while using the machine.
- Do not allow the device to get wet and use it only when it is completely dry.
- Do not allow the device to come into contact with any liquid.
- Keep the device away from moisture, heat and areas where the temperature is high.
- Do not touch the device with wet or damp hands.
- Do not use the device just before going to sleep because the stimulation may make it difficult to fall asleep.

## Operating Instructions

1. Connect the adapter cable to the device plug.
2. Connect the adapter to an electrical outlet.
3. Press the “⏻” button to turn on the massage cushion.
4. Press the “↻” button to adjust the massage direction.
5. Press the “▮▮▮” button to adjust the intensity. The massage speed goes from low to medium and then to high.
6. Press the “🌀” button to turn the heat function on/off.
7. Press the “⏻” button to turn off the massage cushion. Unplug from the electrical outlet after use.
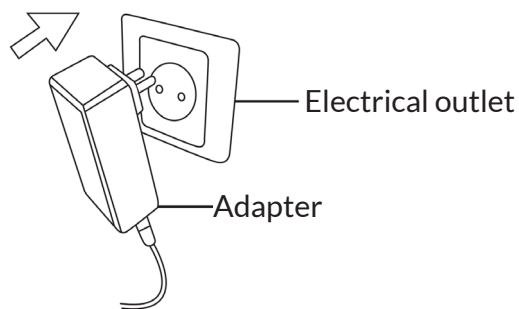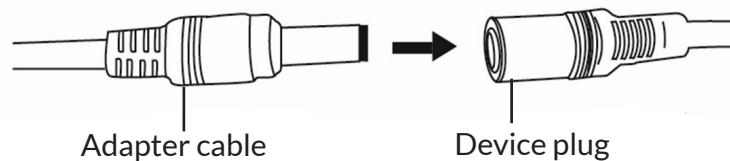
## Bedieningspaneel
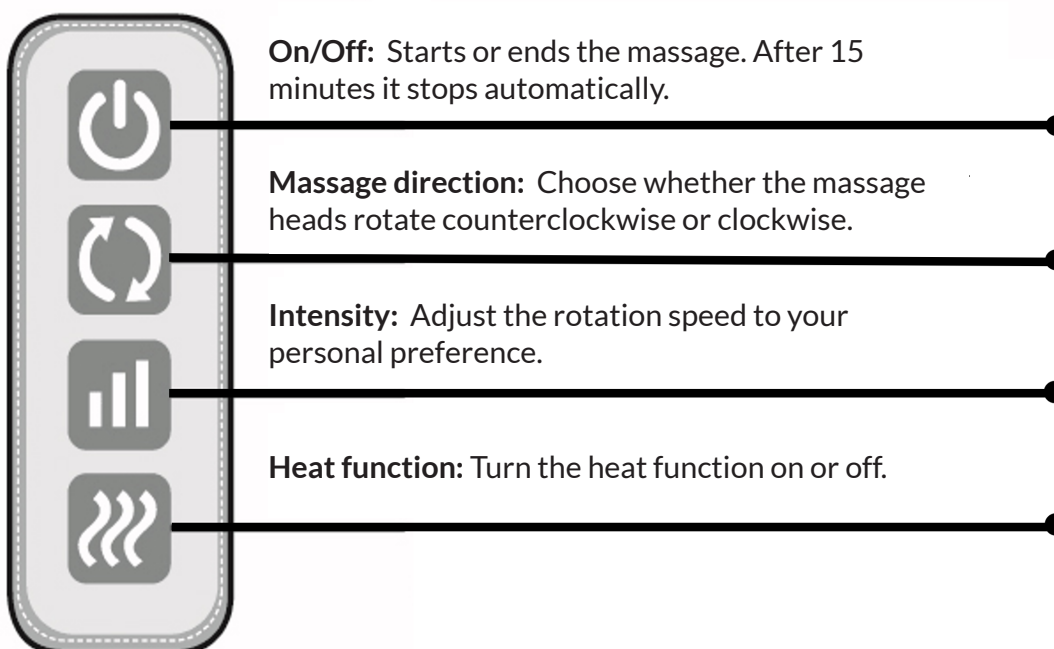
## Technical Data

Size	30 x 22 x 9 cm
Carrying bag size	43 x 13 x 18 cm
Weight	1,3 kg
Material	PU leather
Input - Frequency - Output	Input: 200-240V - 50HZ Output: 12V - 1600 mA
Power	24W
Suitable for (body part)	Upper back, lower back, neck, shoulders, abdomen, legs, calves, feet, hands, arms, waist
Suitable for (functions)	Muscle aches, stiff muscles, stimulation blood circulation, making the muscles supple and flexible
Automatic shut-down	After 15 min.
Overheat protection	Yes
Massage Type	Shiatsu massage
Number of Massage Nodes	16
Variable Intensity	Ja
Number of Intensity Modes	3
Heat Function	Yes

## Contents of Package

- Good4You Massage Pillow - Model RP3

### *Attached Accessories*

- User Manual
- Car Adapter
- Power Adapter

## Cleaning and Maintenance

- Remove the device from the charger before cleaning.
- Use a soft, damp cloth and non-abrasive cleaning agents.
- Do not immerse the device in water or other liquids and keep it away from solvents and harsh detergents.
- Store the device in a dry, clean place. Do not store the device in extremely high or low temperatures.

## Disposal and Recycling

- Think of the environment and contribute to a cleaner living environment! At the end of its life, do not dispose of the device with normal household waste, but take it to a government designated collection point for recycling.

## Service and Warranty

- If you need service or information regarding your product, do not hesitate to contact us.
- We provide a one-year warranty on our products. The customer must be able to show proof of purchase and the warranty must be valid within a period of one year from the date of purchase. Please note that the warranty does not apply: for damage caused by improper use or repair; for parts subject to wear and tear; for defects of which the customer was aware at the time of purchase; for defects caused by the customer's own fault; for damage caused by third parties.

## Disclaimer

- Changes subject to change; specifications subject to change without notice.

## Problems

<b>Device does not work</b>	Check that all cables are connected properly and that the plug is in the wall socket.
<b>The device switches off during operation</b>	<p>For safety reasons, the device is set to stop automatically after 15 minutes. Check that the plug is still in the socket and the cables are properly connected.</p> <p>If the above does not work: switch off the device and wait one hour for it to cool down. The device will automatically turn off if it overheats. If the device continues to turn off, do not use it and contact the seller for advice.</p>

# Good4You<sup>®</sup>

Made in PRC