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Nederland

Product information

Let your face shine with our skin rejuvenation device. With the latest skin care technologies and proven results, this device offers you the key to timeless beauty.

Experience a transformative experience as you tackle fine lines, wrinkles, acne and pigmentation.

Important :

This product uses microcurrent ion technology, so light current (EMS) sensation is a normal occurrence.

Usage & Features

Press and hold the power button for more than 1.5 seconds to turn on the device.

Important :

- The product cannot be used while charging!
- The product itself is not waterproof, avoid excessive contact with water!

Clean Mode (Green)

To improve the function of cells or oxygen utilization in the skin, facilitate skin microcirculation or accelerate lymph circulation.

1. Apply makeup remover or other cleansing lotion evenly to the skin.
2. Press and hold the button for about 1 second to enter "CLEAN" mode (indicated by the green LED color).
3. Follow the steps explained in the illustrated image below*.

EMS Mode (Red)

This mode is used to stimulate and relax the skin, and helps hydrate the deeper layers of skin. Combined with red light, it can penetrate deep into the skin, revitalize collagen, reduce dark spots and reduce fine lines and wrinkles.


1. Apply makeup remover or other cleansing lotion evenly to the skin.
2. Hold the power button for about 1 second to start and press the power button again to enter "EMS" mode (indicated by red light).
3. Follow the steps explained in the illustrated figure below*.

Warm Mode (Blue)

This mode can gently activate and relax the skin, restoring skin radiance and creativity. Combined with blue light, it facilitates collagen synthesis, revitalizes the skin, firms up sagging skin, and skin pores are improved.

1. Apply makeup remover or other cleansing lotion evenly to the skin.
2. Hold the power button for about 1 second to start and press the power button twice to enter "HOT" mode (indicated by blue light).
3. Follow the steps explained in the illustrated image below*.

Illustrated image*

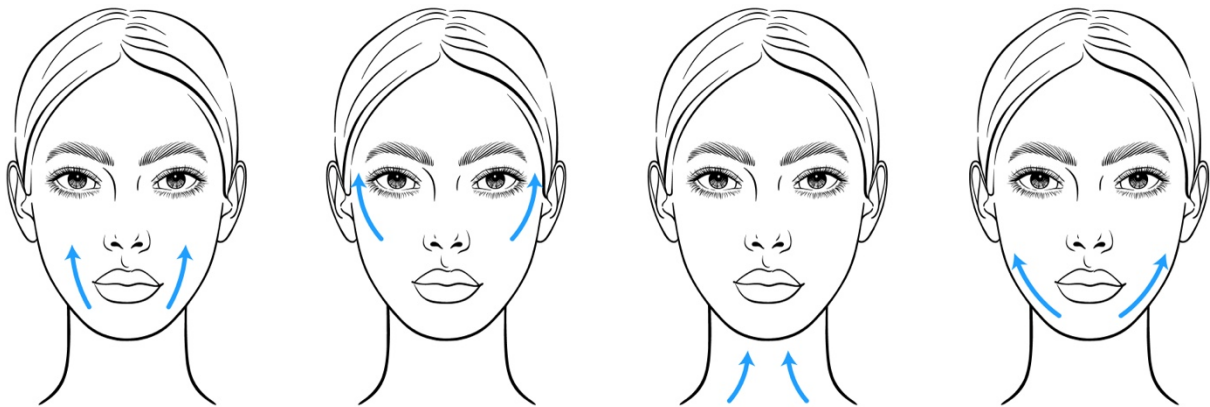1. Rub upward, from chin to cheek.
2. Rub upward from cheek, outward.
3. Rub up, from neck to chin.
4. Rub up, from chin to upper jaw

Safety measures

1. Do not use with plastic surgery, implants or pregnant individuals.
2. Do not store or use in damp areas.
3. Each mode can only be used once a day and it is recommended to use three modes in succession each day.
4. Do not wash the device with water, but wipe with a damp cloth.

Disclaimer

By using this manual, you agree to this disclaimer. TLS Commerce B.V. makes every effort to keep the information in this manual and related web pages as complete, accurate and current as possible. If unclear or incorrect information, please feel free to contact the customer service of TLS Commerce B.V. All information on this website and manual is intended for personal use. To reserves, TLS Commerce B.V. accepts no responsibility for damages of any kind arising from use, incompleteness or inaccuracy of the information provided on the manual or information provided on the website when used carelessly. The information is subject to change without prior warning or notification.

Company details

TLS Commerce B.V. | Nuvance

Chamber of Commerce: 85524271

VAT no: NL863652591B01

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6004 HT Weert

The Netherlands