
NORTHWALL

NORTHGUN FLEX

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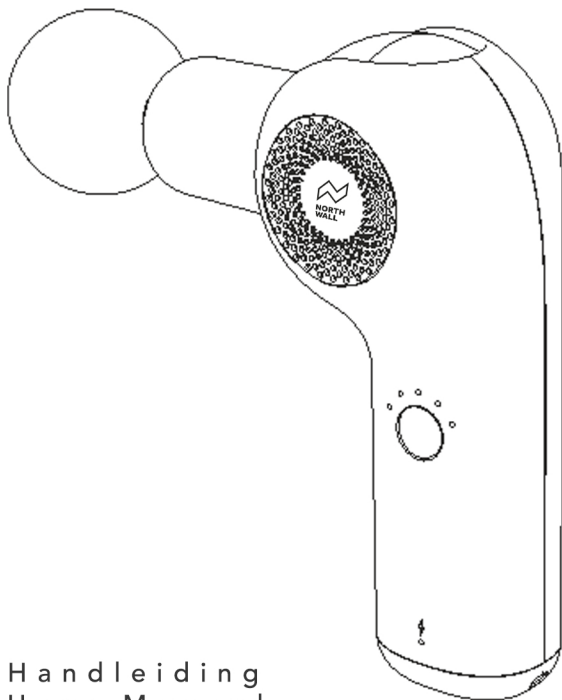READ THIS BEFORE USE

Warning

- Never use this device near water.
- Keep the ventilation grille free from dust and dirt.
- Do not remove any screws and never try to open the device yourself.
- Do not use this appliance for more than one hour at a time, rest 30min. before using again.
- Never charge the device overnight and only charge under supervision.
- Disconnect the charger from the battery before use.
- Suitable for adults only. Do not use this device on injuries. Consult a doctor before using this device.

Warning

DANGER: To reduce the risk of electric shock, injury to limit fire, this appliance may only be used in accordance with instructions below:

- Only suitable for adults.
- Never leave the Northwall massage gun unattended while charging.
- Never use this massage gun in any other way than indicated.
- Check the appliance carefully for defects before each use.
- Charge only with the supplied Northwall USB-C cable.
- Never use water near this massage gun.
- Keep hair, fingers and other body parts away from the shaft of the massage heads because these can get pinched.
- Bruises may occur when using a massage gun.
Check treated parts of the body regularly and stop immediately in case of pain or discomfort.
- Use the Northwall massage gun only on soft body tissues, without discomfort or pain. Don't ever use these on head and bone.
- Use the Northwall massage gun only on dry skin or clothes in an up and downward movement for about 60 to 90 seconds per muscle.

DISCLAIMER

READ THIS BEFORE USE

Never use the Northwall massage gun without approval from your doctor in the cases below:

Pregnancy, diabetes, use of pacemakers, recent surgeries, epilepsy, migraines, hernias, spondylolisthesis, spondylolysis, spondylosis, joint replacements, metal plates in the body or other concerned present inc. your physical health. These cases do not exclude the use of the Northwall massage, but advice from a medical expert is needed.

GUARANTEE

If a manufacturing error or defect occurs within two year of purchase, then Northwall will repair or replace the massage gun free of charge. In case of defect please contact us by mail: info@northwallproducts.com

CONTACT

For questions about the Northwall massage gun or any other Northwall product you can always contact Northwall customer service. We are happy to help and we will do our absolute best to come to a solution together.

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
Specifications

IMPORTANT: If the massage gun is not charging than try to use a USB port from a laptop or pc. Not all chargers are suited. Apple adapters, chargers with an output lower than 2A and fastchargers are not suited.

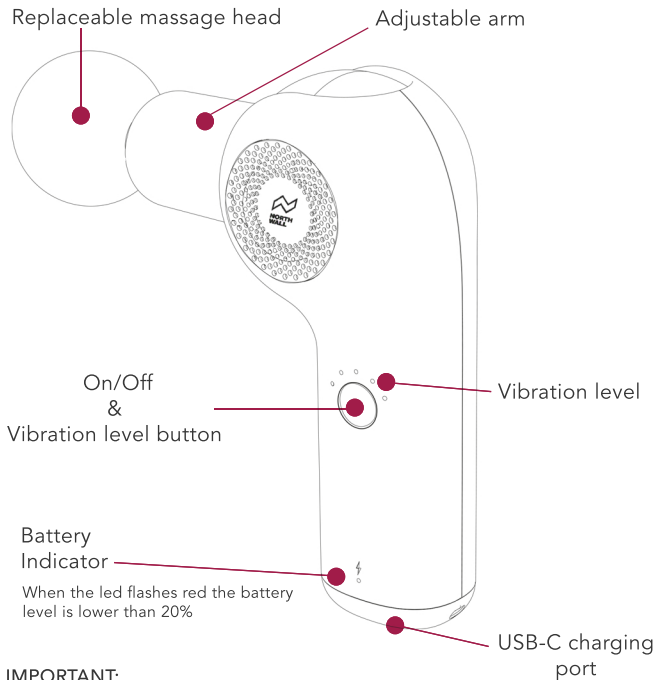IMPORTANT:

The massage gun shuts down after 10m. of use as a safety measure.

**Bal**

Gentle massage for large muscles such as arms, legs and back

**Bullet**

For a deep massage round muscle knots and joints

**Fork**

For neck massage, achilles and back

**Flat head**

Firm massage for large muscles such as arms, legs and back

Manual

- PAY ATTENTION! The Northgun switches itself off after 10m. for security.
- Charge the battery for up to 4 hours before first use.
- The average working time ranges between 1 and 2 hours depending on the vibration level and the applied pressure.
- Attach the massage head before turning on the device. You can push in the head and pull it out after use.
- Switch the device on and off by pressing the On/Off button for 2 sec.
- Short press once to go to a higher vibration mode.

How to use the Northgun Flex:

www.northwall.be/pages/how-to

Technical Data

Brand:	Northwall
Model:	Flex
Battery Capacity:	2600mAh
Power:	25W
Voltage:	7.3V
Weight:	0.4Kg
Charging time:	3 hours
Vibration levels:	5

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