**WELLCARE®**

# Mobile Heating Belt

3-in-1 | Cordless

**Instruction Manual**

## **MHB10**

**Mobiler Wärmegürtel**

3-in-1 | Drahtlos

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**Mobiele Verwarmingsgordel**

3-in-1 | Draadloos

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**Ceinture Chauffante Mobile**

3 en 1 | Sans fil

**Mode d'emploi**

## MODEL : MHB10

<b>EN</b>	Mobile Heating Belt 3-in-1   Cordless Instruction Manual	..... 03
<b>DE</b>	Mobiler Wärmegürtel 3-in-1   Drahtlos Gebrauchsanleitung	..... 19
<b>NL</b>	Mobiele Verwarmingsgordel 3-in-1   Draadloos Gebruiksaanwijzing	..... 37
<b>FR</b>	Ceinture chauffante mobile 3 en 1   Sans fil Mode d'emploi	..... 55

B820.101\_ G1020

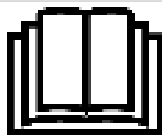
**WELLCARE CO., LTD.**

16F, No. 738, Zhong Zheng Rd., Zhong He District, New Taipei City 235, Taiwan

Tel. 886 2 8226 2899 Fax. 886 2 8226 2889

[www.wellcare-global.com](http://www.wellcare-global.com)

## Explanation of symbols

Read instructions

Do not use when folded or creased

Do not pierce with needles

Not to be used by very young children  
(0-3 years)

Do not dry clean

30°C extra gentle washing program

Do not bleach

Do not tumble dry

Do not iron

Thank you for purchasing the Mobile Heating Belt, MHB10, from Wellcare.

To enjoy the most comfortable experience and the greatest satisfaction, please take a few minutes to read this instruction carefully before using this product.




---

## **IMPORTANT SAFETY INSTRUCTIONS RETAIN FOR FUTURE USE**
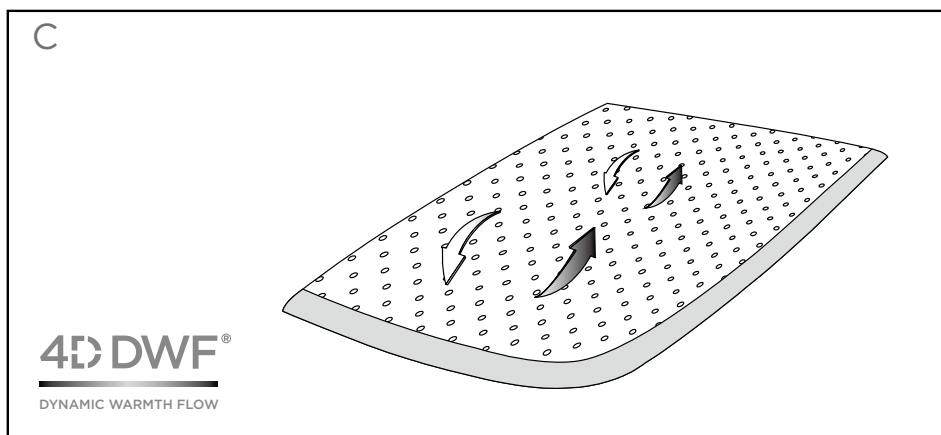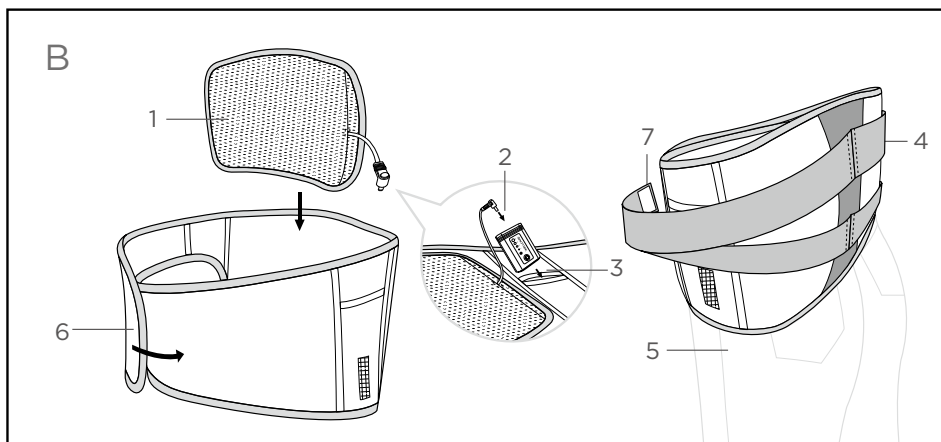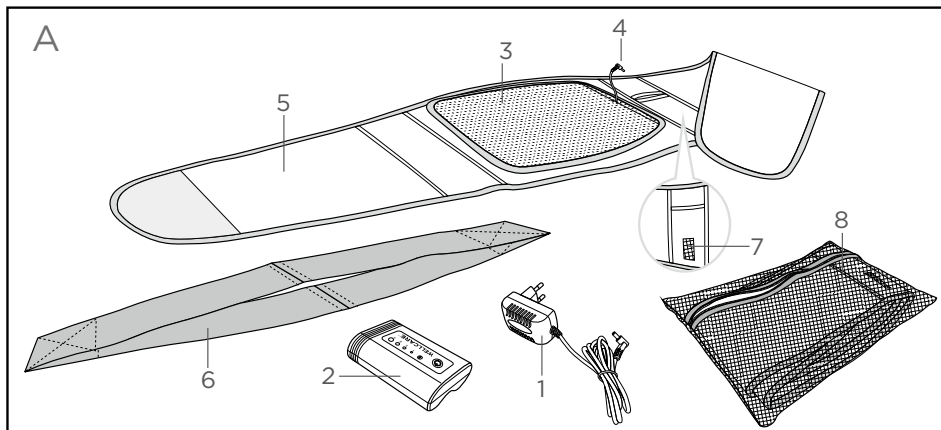
### **KEEP CAREFULLY IN CASE IT IS NEEDED!**

The mobile heating belt is dedicated to protecting the waist area and providing warmth which helps stimulate the circulation in areas of tense muscle on the lower back. It can be used indoor and outdoor and is designed to be used for many years if it is properly used. Read the instruction carefully to prevent damage or danger caused by incorrect use. Keep this instruction manual at an easy-to-access location so you can refer back when necessary.

If you let somebody else use this product, please pass on this instruction manual as well and make sure the user has read it before using.

#### **4D DWF® Technology** (please see fig. C)

A specially-designed heating module is used in all of Wellcare's mobile heating products. This heating module utilizes WELLCARE's recent and the most innovative technology-4D DWF®. While this technology pioneers in the flexible heating industry, its production process is also the state of the art.

Our physical living environment consists of three dimensions, width, length, and depth; hence, the term "3D". As Wellcare introduces to you the 4D DWF® technology, the "temperature" brought to you by dynamic flow of warmth is the definition of the fourth dimension. 4D DWF® employs concepts of breathing system which creates full-body heat circulation for heat to disperse evenly, thoroughly, and quickly. The heat is generated by heating elements, builds up and gathers around several hundred uniquely designed circular holes. Through these holes, the heat flows bi-directionally to create a gentle heat flow and the thermal energy can flow rapidly to warm up human bodies.

There are many benefits of this innovative technology :

- Better breathability
- Rapid heating
- Higher temperature
- Uniquely designed holes to provide 360° warmth circulation
- Greater heat distributing area
- Effective heat transmission
- Strands of warmth flow retain air moisture
- Warmth kept at a constant temperature
- Energy saving

#### **This appliance is not intended for medical use in hospitals**

- Please read carefully and follow the instructions.
- Do not use this appliance if there is visible damage.
- Examine the appliance frequently for signs of wear or damage. If there are such signs or if the

heating belt has been misused or does not work, return it to the supplier before switching it on again.

- Use the appliance only as described in this instruction manual.
- Before connecting the charger to the wall socket, please make sure the wall socket is with the same mains voltage rating as stated on the rating label on the charger.
- The appliance must not be used when folded or crumpled.
- This appliance must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
- Do not use the appliance for persons who cannot move without help, small children, persons who cannot feel heat, or animals. Persons who can feel heat to some extent may only use the appliance with supervision.
- Do not operate this mobile heating belt unsupervised.
- Caution! Never fall asleep while the appliance is in operation.
- Prolonged use can cause burns to the skin.
- When the appliance will be used for an extended period (falling asleep), the temperature setting should be adjusted to the lowest level (setting 1) for continuous use to avoid overheating on the same body part, which may lead to burns to the skin.
- Do not pierce the appliance with needles, pins or other sharp objects.

- Do not tug, twist or make sharp bends to the cables.
- Do not expose the charger nor rechargeable battery pack to moisture or liquids.
- Do not use the appliance when it is damp or wet.
- If the appliance is damaged, it must be replaced by the manufacturer or qualified persons to prevent hazards.
- The rechargeable battery pack should not be placed on or beneath or covered by any sources of heat when the mobile heating belt is in use or when charging.
- Avoid sharp bending or bundling of the appliance.
- While the appliance is switched on,
  - Do not place any large or heavy objects on it such as suitcases, etc.,
  - Do not place any sources of heat on it such as hot water bottles, heating pads or suchlike.
- Do not fold the appliance when it is hot.
- Children shall not play with this appliance.
- Children below the age of three (0-3) are not allowed to use this device due to their inability to react to overheating.
- This device may not be used by young children over the age of three(3) unless the control unit has been pre-set by a parent or guardian or the child has been adequately instructed on how to operate the control unit safely.

- This appliance may be used by children over the age of eight(8), and persons with reduced physical, sensory or mental capabilities or lack of experience or knowledge if they have been supervised or instructed concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision.
- If you have any concerns about health issues, please consult with your doctor before using the appliance.
- If you use a pacemaker, please consult your doctor and the manufacturer of the pacemaker before using this product due to the EMF (electromagnetic fields) emanated from this electric product may interfere with functions of your pacemaker under certain circumstances.
- Do not use this appliance on pets or other animals.
- Do not try to open the rechargeable battery or the charger nor to drop them to the floor.
- If you do not use the heating belt daily, disconnect the DC power plug from the battery pack.
- When the product will not be used for a period of time, we recommend that you store it in the original packaging.
  - Allow it to cool down before the storage.
  - Do not crease the heating pad by placing items on top of it during storage.
  - Store it in a dry place away from any liquids

or sources of heat.

- The rechargeable battery pack must be fully charged before being used for the first time (Please see : Introduction of rechargeable battery pack).
- When the charging is finished, always disconnect the battery pack and the charger from the wall socket.
- The device can only be used with rechargeable temperature control (Model no. : WE-221a) and charger (Model no. : RSS1002-101084-W2E).
- Use only the original charger (Model no. : RSS1002-101084-W2E) to charge the rechargeable battery. The charger should not be used for any other batteries, especially non-rechargeable batteries.
- Do not use the heating belt directly on your skin. Only use it on top of clothing.
- Do not use the heating belt while lying down.

## **ITEMS INCLUDED IN THE PACKAGE**

1 Belt with Velcro fastener

1 Detachable heating pad

1 Elastic compression belt

1 Rechargeable battery pack/Temperature control

1 Charger

1 Storage bag

1 Instruction manual

1 Product catalogue

## 1 Service instruction

### **DESCRIPTION OF THE DEVICE** (please see fig.A)

1. Charger
2. Rechargeable battery pack with temperature control
3. Detachable heating pad with DC power plug and fastening ring
4. DC power plug
5. Belt with Velcro fastener
6. Elastic compression belt
7. Battery pocket with window
8. Storage bag

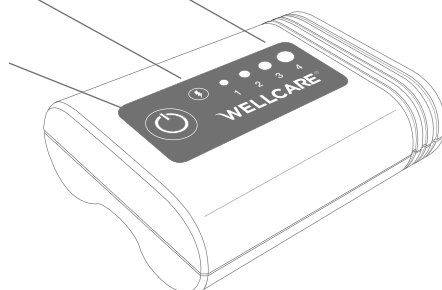
### **INTRODUCTION OF RECHARGEABLE BATTERY PACK**

The rechargeable battery pack(2) works as a temperature control unit with power switch button. Blue LEDs display as 4 temperature settings and Green LED shows the state of charge.

Blue LED (1-4) : temperature level

Green LED (  ): state of charge

Power switch button

1. The rechargeable battery pack must be fully charged (for at least 5 hours) before being used for the first time or not having been used for a long time. One single fully charge allows you to use this heating belt 5 hours at the highest temperature setting.
2. Always charge the battery pack in a dry environment at normal room temperature and do not let the battery pack be covered by any objects.
3. In order to charge the battery pack, connect the charger(1) to the battery pack(2) first, then connect the charger(1) to the wall socket. During the charging process, green LED will flash.
4. When the battery pack(2) is fully charged, the green LED will keep lighting up. Remove the charger(1) from the wall socket and disconnect the battery pack (2) from the charger(1).
5. If you want to check the state of charge of the battery pack(2), you may press the power switch button. Please make sure the battery pack(2) is not connected to the heating pad(3) nor the charger(1) :
  - If the green LED lights up constantly for 3 seconds, the state of charge is more than 25%
  - If the green LED flashes 3 times in 1 second, the state of charge is below 25%. Please recharge the battery pack before using it.
6. Please note that the battery pack must be protected against extreme weather (such as placing it under clothing in very cold weather) to maintain its performance.
7. The battery pack should be charged at least one

time per year to main its capacity.

8. We highly recommend to use and store the battery pack :
  - Operating temperature : 0 ~ 45 °C
  - Storage temperature : -10 ~ 45 °C, best recommended : -10 ~ 25 °C

## **PREPARING FOR USE** (please see fig.B)

Always make sure the detachable heating pad is not folded during use.

---

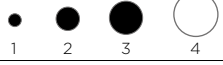
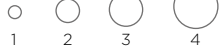
1. Attach the detachable heating pad to the middle area of the inner side of the belt. Make sure the heating pad is without being folded and the DC power plug is on the side of the battery pocket.
2. Connect the DC power plug to the battery pack. Insert the rechargeable battery pack into the battery pocket. Make sure the connection socket on the battery pack is put upward and the LEDs can be seen through the window of the battery pocket.
3. Set the desired temperature level (1-4) by pressing the power switch button repeatedly and the device is now ready for use.
4. Elastic compression belt is optional. If you want a stronger support to the body, attach it to the middle of outer side of the belt.
5. Place the heating pad area around your lower back or your abdomen.
6. Secure the belt by fastening the Velcro.

7. Adjust the elastic compression belt to fit different body contours.

This device can be also used without the detachable heating pad to enjoy the nature warmth from the body temperature, or to keep the body contour for daily use. Please put the heating pad flat in the storage bag; keep it in a dry place and without any contact to liquids, heat or any heavy objects on the top of it.


## TEMPERATURE CONTROL

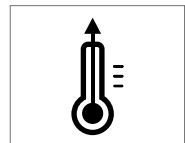When the mobile heating belt is in use, you can increase the temperature level by pressing repeatedly the power switch button on the rechargeable battery pack(2). The blue LEDs will display as :

Press the 1 <sup>st</sup> time		1 <sup>st</sup> temperature level
Press the 2 <sup>nd</sup> time		2 <sup>nd</sup> temperature level
Press the 3 <sup>rd</sup> time		3 <sup>rd</sup> temperature level
Press the 4 <sup>th</sup> time		4 <sup>th</sup> temperature level
Press the 5 <sup>th</sup> time		Switched Off

## RAPID HEATING

The mobile heating belt is with the rapid heating function for user to feel warmth within 7-10 minutes after it is switched on at temperature setting 4.

## TIMER FUNCTION

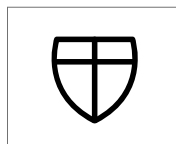
Now you can enjoy the pleasant warmth. After

approx. 90 minutes, the mobile heating belt will be switched off automatically. If you want to continue using the appliance, press the power switch button and select your desired temperature level.

## OVERHEAT PROTECTION

This mobile heating belt is equipped with a safety temperature control that regulates the temperature automatically and prevents the heating pad from becoming too hot.

## CAUTION

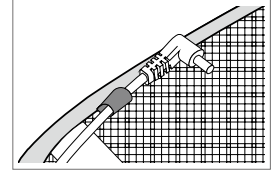
1. Do not use the mobile heating belt when the heating pad is folded, creased or drawn together.
2. The rechargeable battery pack may be slightly warm during the operation. Therefore, the battery pack must not be covered or placed on top of any heating sources.
3. You should not gauge the heating belt's temperature by hands because the skin on your body is more sensitive than hands. The temperature may feel lukewarm by hand but could be unpleasantly hot to your body.
4. The repairs to the mobile heating belt may only be undertaken by qualified persons by using original replacement parts. Repairs carried out by unqualified persons may cause substantial risk to the user.

## CLEANING AND CARE

This mobile heating belt is machine washable

- The charger(1) and the rechargeable battery pack(2) CANNOT be washed.
- Before cleaning, always disconnect the battery pack(2) from the heating pad(3) and remove it from the battery pocket(7).
- Use a dry cloth to clean the battery pack(2) only instead of using any chemical or corrosive cleaning agents.
- Minor stains can be removed with a cloth or damp sponge with a little liquid mild detergent.
- Please note that the mobile heating belt should not be dry-cleaned, wrung out, machine-dried, pressed or ironed.

- Set washing machine to a gentle washing cycle at 30 °C and wash the mobile heating belt together with other textiles. Use a mild washing detergent for delicate fabrics and dose it according to the manufacturer's instructions.
- Always use a laundry bag for better protection of the product.
  - Put the DC power plug through the fastening ring on the detachable heating pad.
  - Close the belt by fastening the Velcro.
- For dark color products, always wash separately.
- In order to dry it, lay the appliance flat on a moisture-absorbent towel, mat or on a clothes horse.
- Please note that excessive washing will wear out the heating pad. For this reason, the appliance should only be washed a maximum of 15 times throughout its entire life.
- Do not use pegs or similar items to attach the heating pad to the clothes horse, otherwise the heating pad may be damaged.
- Only use the mobile heating belt when it is complete dry.
- Never switch on the appliance to dry it.

## STORAGE

If this appliance will not be used for a long time, we highly recommend :

- Disconnect the heating pad and the rechargeable battery pack.
- Always let the appliance to cool down first before storage.
- Store the appliance in the storage bag; put it in a dry and cool place.
- Best recommended temperature: -10 ~ 25 °C, relative humidity : 60 ± 25%.
- Do not place any objects on the top of the mobile heating belt.

## BATTERY REMOVAL AND DISPOSAL

The battery cells must be removed from the battery pack before being scrapped. Disconnect the battery pack, open the covers of battery

pack, cut off the PCB of battery cells and remove them out. The battery pack can be disposed safely. Never put the battery pack into the fire or close to heat to avoid explosion. Dispose the battery cells to the special recycle location only.

Be careful of leaking liquid from the battery. Do not let it touch your hand or any other part of the body. Use any other tool to remove the battery if the leakage occurs.

This appliance must not be disposed together with domestic waste. Please dispose of this appliance in accordance with Directive 2002/96/EC on Waste Electrical and Electronic Equipment (WEEE). If there is any question, please contact the municipal authority for more information about disposal in your area.

## WARRANTY

Please note Wellcare provides two-year limited warranty for this product. Wellcare warrants that this product will be free from defects in materials and workmanship for a period of two years beginning on the date of purchase.

During the two-year period beginning on the date of purchase, Wellcare will replace defective products covered under this limited warranty. Please note a clear copy of the purchase receipt for the product is required. To request limited warranty service, please return the defective product to the store where you purchased it with the clear copy of purchase receipt within the limited warranty period. The replacement product will assume the remaining limited warranty of the original product.

We advise you to call our helpdesk for any problems regarding the product before returning it to the store where you purchased it.

However, please note that our limited warranty does not cover any damage due to shipment of the products to you and external causes, including accident, abuse, misuse, problems with electrical power, usage not in accordance with product instructions, failure to perform reasonable maintenance, problems caused by use of parts and components not supplied by Wellcare, and problems caused by unauthorized repair of Wellcare's product. Furthermore, our limited warranty does not cover any damage due to the defects known to exist

by you when you purchase the product.

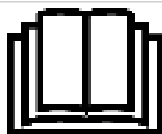
This limited warranty gives you specific legal rights, and you may also have other rights that vary from country to country (or jurisdiction to jurisdiction). Wellcare's responsibility for malfunctions and defects in product is limited to replacement as set forth in this limited warranty statement. All express and implied warranties for the product, including but not limited to any implied warranties of an conditions of merchantability and fitness for a particular purpose, are limited in duration to the limited warranty period set forth above and no warranties, whether express or implied, will apply after such period. Wellcare does not accept liability beyond the remedies set forth in this limited warranty statement or liability for indirect, incidental or consequential damages.

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