



FASCIAL GUN

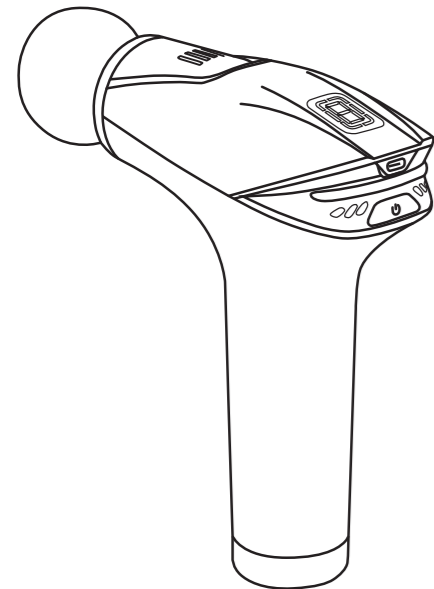


an instruction manual

Product overview

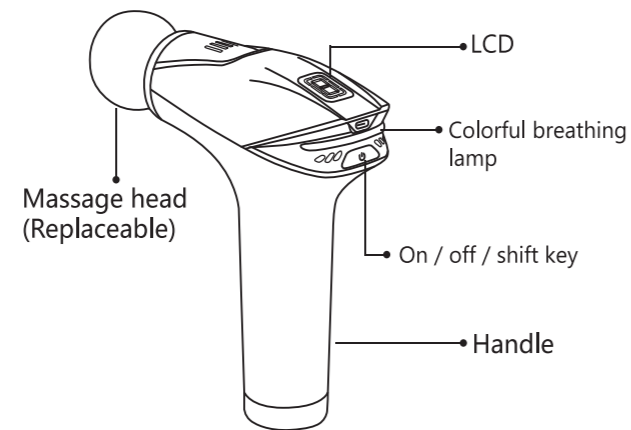

Massage head

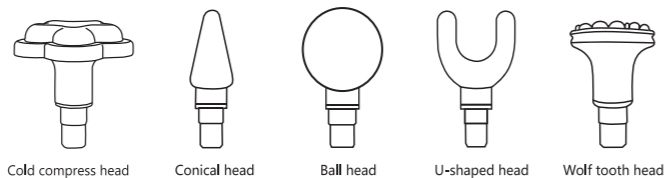

It is suitable for massaging sore muscle parts to quickly relieve pain after exercise and prevent inflammation.

It is suitable for impacting deep tissues, plantar acupoints, such as meridians and joints.

It is suitable for massaging various muscle parts, such as palm, sole, etc.

It is suitable for massaging the neck and Achilles tendon of spine.

It is suitable for massaging large muscle groups, legs, hips and other deep stimulation.

Product use

On / Off (key)

Press and hold the "⏻" key for two seconds to start the machine. Long press the "⏻" on key for two seconds to shut down.

Strength adjustment


In the power on state, click "⏻" once to add one file, The maximum gear is 8th gear. When you reach the eighth gear,click "⏻" again to turn to the first gear.

Application method of cold compress head


Put the cold compress head into the refrigerator for one hour or half an hour, and then cold compress, the effect will be better.

Precautions: in order to protect bones and muscles and avoid damage caused by long-term fixed-point massage, open the intelligent protection mode that stops automatically for 15 minutes.

Replace massage head

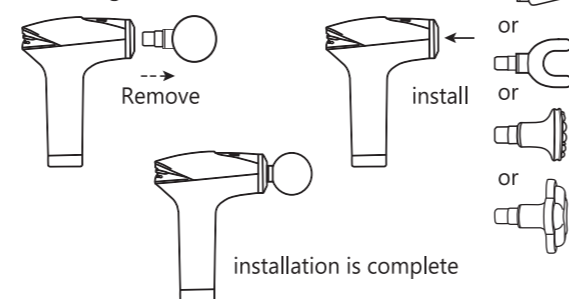

About charging

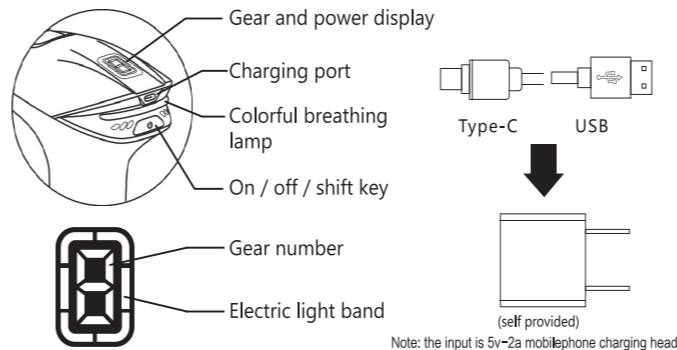

1. Power on / off button (long press to start / long press to shut down)
2. The massage head can be removed and replaced. When removing the massage head, just pull it out with force.
3. Please use 5V charging equipment for charging
4. For the first use, please charge for three hours. Please connect the 5V power supply line with the data line for charging
5. The charging screen shows that 1-6 light strips are composed of "⏻" and all 6 are on, indicating full charge


♥ indication

- Pain and spasm due to muscle contusion, sprain and strain
- Help the edema fluid flow in the swollen area
- Loose thickened connective tissue and fascia
- Reduce the accumulation of lactic acid in muscle
- Increase joint range of motion
- Eliminate muscle fatigue

taboo

- Aneurysms, bleeding, use of blood thinners
- Heart disease, pacemaker or defibrillator
- Pregnancy, cancer
- The device cannot be used within 3 inches of internal fixation within 90 days
- Sensitive areas: head, face, cervical spine, vertebral body
- Tibia, dorsum of foot, dorsum of hand, etc

detailed list

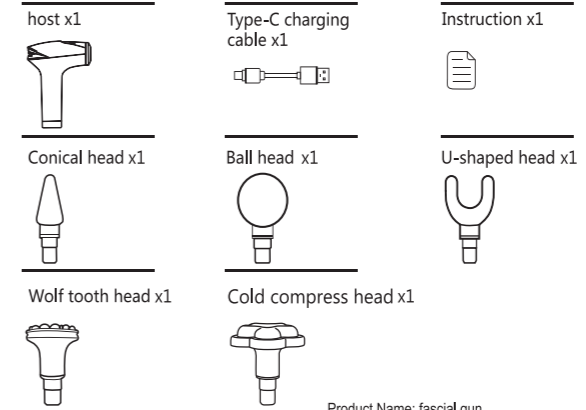

Product Name: fascial gun
 Rated voltage: 7.4v
 Motor power: 30W
 Interface type: type-C interface
 Revolutions: 1800-3200
 Gear: 8-gear force adjustment
 Number of massage heads: 5
 Material: electronic components + ABS plastic
 Product weight: 620g
 Product size: 155x66x198mm

Warranty Card

Product Model _____

Serial Number _____

Purchase Date _____

Customer Name _____

Customer Tel _____

Customer Email _____

Warranty service terms:

- (1) The product is guaranteed for one year and the built-in battery is guaranteed for 6 months
- (2) The warranty starts from the first purchase date of the product.
- (3) In case of quality problems requiring warranty, please this warranty card. The purchase receipt is sent back to the dealer or engineering maintenance service together with the product.
 The following conditions are not covered by the warranty:
 (1) Unable to provide formal purchase certificate and warranty card with detailed information.
 (2) Damage caused by man-made or irresistible factors.
 (3) Irreversible damage caused by man-made machine dismantling gun.

The company reserves the right to make final interpretation of the above after-sales terms and restrictions

FC CE RoHS Appearance patent counterfeiting must be investigated Made in China