Shyineyou

USER MANUAL
RED LIGHT THERAPY LAMP


PACKAGE INCLUDING

Protective Goggles×1

Power Adapter×1

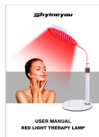User Manual×1

Controller×1

Red Light Lamp×1

Thank you for purchasing Shyineyou red light lamp. You are going to experience this red light therapy device designed to improve skin and relieve the pain. Please read the user manual carefully before use. Check whether the product is damaged or not when you received the product.

PRODUCT SPECIFICATIONS

Item Number	HGD-LAMP-140LED
Color	White
Wavelength	660nm
LED QTY	140PCS
Timer	15/30/45mins
Brightness level	20%-100%
Power	12.5W
Input voltage	AC100-240V
Output voltage	5V3A

LIGHT SENSITIVITY TEST

Please do a skin sensitivity test first if you are the first time using red light therapy. Turn on the light first and expose the skin of your hand to the light, then turn off the light in two minutes. Observe the areas of skin where the light hits in the next five minutes. Please stop using the product if any redness or rash reaction is found.

TIMER FUNCTION

BRIGHTNESS“-”100%-20%

BRIGHTNESS“+”20%-100%

3 TIMING SETTINGS

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POWER BUTTON

Turn on/off

HOW TO INSTALL

WARM TIPS

- 1.Sensitive skin or allergic reaction can cause redness and itching.If these occur, seek medical care from a physician.
- 2.Red light therapy presents as a safe therapy, it is not a heat lamp so they do not cause thermal injuries, only feels a little warm during use.
- 3.Please wear eye mask or goggles to avoid eyes discomfort.
- 4.Indoor use only not waterproof.

DO NOT USE IN THE FOLLOWING SITUATIONS

1. Open wounds and damaged skin.
2. Over the gonads or anesthetic areas, Over fractured bones, the eye directly.
3. Parts of the body with metal implants, like in a total knee.
4. Near or over an implanted electrical stimulation device.
5. Tissues previously treated with deep Xray or radiation.
6. Near or over a pacemaker.
7. During the use of photosensitizing medications (i.e. lithium, melatonin, phenothiazine antipsychotics, and certain antibiotics) or your skin sensitive to sun or light specific wavelength of light (660nm/850nm)
8. People with a history of skin cancer, thyroid conditions and Malignancies and precancerous lesions.
9. For hemophiliacs not covered by factor replacement.
10. Diseases that involve the retina of the eye or systemic Lupus erythematosus.
11. On patients with vascular abnormalities, i.e. deep vein thrombosis, emboli, severe atherosclerosis.
12. Women who are pregnant or clients with epilepsy should consult their physician before beginning IR light therapy treatments.
13. You must wait five days after Botox or cosmetic fillers.

RECOMMENDATIONS FOR USE

Discover the ultimate way to maximize your experience by following our expert recommendations. We suggest starting from the lowest gear to ensure a smooth initiation into this exciting journey. Remember, timing is everything! Begin with 10 minutes per session, allowing your body to adapt and embrace the transformative power. After just three days, you'll be ready to ascend to the next gear, unlocking new levels of intensity and satisfaction. Trust your instincts and personal progress, as you have the freedom to tailor both gear selection and duration according to your unique journey. Embrace the adventure and experience the incredible possibilities that await you.

WHAT RED LIGHT THERAPY CAN DO

- 1.Improve wound healing.
- 2.Reduce stretch marks.
- 3.Reduce wrinkles, fine lines and age spots.
- 4.Improve facial texture.
- 5.Improve psoriasis, rosacea and eczema.
- 6.Improve scars.
- 7.Improve sun-damaged skin.
- 8.Improve acne.
- 9.Increase blood circulation to the tissue.
- 10.Reduce inflammation in cells.

WARRANTY

We offer a 1-year warranty from the date of purchase for all our products. We warrant each device to be free from defects in materials and workmanship.

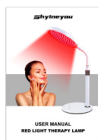
If you have any questions about our products and/or warranty, please feel free to contact us via the email below. We typically respond to all emails within 24 hours.

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