# LUVION®

SONIC BABY TOOTHBRUSH

HANDLEIDING

INSTRUCTION MANUAL

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NEDERLANDS.....4

ENGLISH.....8

FRANÇAIS.....12

DEUTSCH.....16

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#### PACKAGE CONTENTS:

- Electric sonic baby toothbrush
- Replaceable baby brush head (3-18 months)
- Replaceable toddler brush head (18-36 months)
- AAA (R03) battery (non-rechargeable)
- This instruction manual

#### PRODUCT OVERVIEW

##### *Electric sonic baby toothbrush (Fig. 1)*

1. Bristle
2. Replaceable brush head
3. Lamp
4. Power button
5. AAA battery
6. Battery cover
7. Battery cap
8. To detach: Turn brush head anti-clockwise and pull upwards.
9. Turn battery cap anti-clockwise to remove. Insert a new AAA battery.

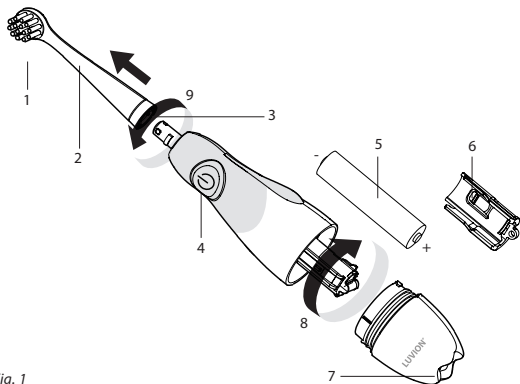
Fig. 1

#### INTRODUCTION

Thank you for purchasing this Luvion sonic baby toothbrush. Make sure you always keep your original proof of purchase. For warranty cases you will always be asked to show a copy of your proof of purchase. For warranty claims you can contact the Luvion dealer or the local Luvion representative. More information can be found at [www.luvion.com](http://www.luvion.com)

This manual contains all the information you need to be able to use this Luvion sonic baby toothbrush safely and with pleasure. Before using the toothbrush, we recommend that you read this manual carefully. We ask for your extra attention for the safety instructions.

#### PRECAUTIONS

- Please read the instruction carefully before use;
- If you do not understand the instructions for use or precautions listed in this manual, do not use the device. Incorrect use of the device may damage it;
- Adult supervision is required at all times;
- Children should be supervised while using the device;
- Store out of reach of children;
- This device is not a toy. Children shall not play with appliance;
- Cleaning and users maintenance shall not be made by children;
- Adult assembly is required;
- Keep small parts out of reach of children;
- Children should not play with, clean or repair the product;
- The toothbrush can be used by children and people with little experience as long they are supervised by a person responsible for their safety;
- Toothbrushes can be used by children and persons with reduced physical sensory or mental capabilities or lack experience or knowledge if they are supervised or have been given instructions concerning the use of the appliance in a safe way and understand the hazards involved;
- If using the device causes your child discomfort or bleeding from the gums, consult your dentist before continuing to use it.


#### IMPORTANT SAFETY PRECAUTIONS

Please note the following precautions before usage. Failure to comply may be dangerous.

- Do not soak the product in water for a long period of time;
- Do not open, repair, or alter the product without authorization;
- Do not expose the product to high temperatures or direct sunlight;
- Do not alter the product in any way other than a battery or brush head change;
- Non-rechargeable batteries are not to be recharged;
- Remove rechargeable batteries before charging;
- Stop using the device if it's clearly damaged;
- Do not use other accessories other than those recommended by the manufacturer;
- Immediately remove exhausted batteries from the appliance and dispose of properly;
- Do not try to recharge the included non-rechargeable battery;
- Do not replace the battery of the device when your hands are wet;
- Please insert the battery as indicated by the polarity markings.



## HOW TO USE

1. First use: Remove the plastic insulator from the battery by holding the toothbrush upside down, thereafter turn the cap (7) counter-clockwise to open and peel off the plastic insulator from the battery compartment. Place battery cap back in place by turning clockwise;
2. Wet the brush and apply appropriate amount of toothpaste;
3. Press on the power button (4). You will see the toothbrush light up;
4. Divide the teeth into up, down, left and right sections. See picture below (fig. 2). This toothbrush is equipped with a 2-minute-total timer and 30-second pulse reminders that signals when it's time to move on to the next quadrant of the mouth;

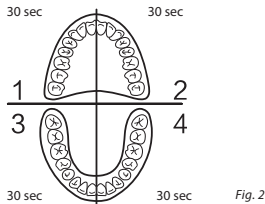5. Place toothbrush against teeth at a slight angle towards the gum line;
6. Press the power button (4) again to turn on the toothbrush. Gently move the brush head across the teeth and gums with small up and down movements;
7. Move to the next section when the brush signals (pulses), (every 30 seconds) until all four sections of the mouth have been brushed;
8. The toothbrush will power off automatically after 2 minutes. If you wish to proceed, simply press the power button (4) again to continue.

## CHANGING THE BRUSH HEAD

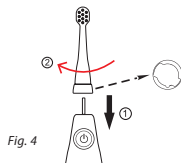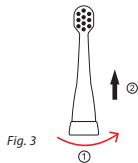
**Step 1:** Detach the brush head counter-clockwise (fig. 3) and pull upward.

**Step 2:** To attach: align the new brush head and turn it clockwise (fig. 4).

**Please note:** The notch on the neck should align with the inner notches of the brush head in order to lock in place.

## REPLACING THE BATTERY

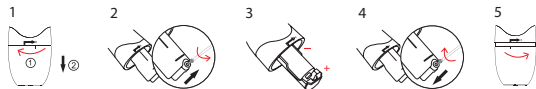**Step 1:** Hold the toothbrush upside down, turn the battery cap counter-clockwise to open;

**Step 2:** Unscrew the battery compartment with a small philips screwdriver;

**Step 3:** Insert a new AAA battery, pay attention to plus and minus poles;

**Step 4:** Place the battery cover door and tighten the screw. Do not over-tighten

**Step 5:** Align line marker on battery cap with the line marker on brush handle. Place the battery cap back in place by turning clockwise.

## SPECIFICATIONS

Material:	DuPont Nylon / ABS / TPE
Power supply:	1 x AAA (R03) battery (included)
Waterproof:	IPX7
Frequency:	16.000 movements per minute

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