COMEDO SUCTION BEAUTY DEVICE

Operation Manual

Please read this manual carefully before using and keep this in a safe place for future reference

Function

Microcrystalline skin-smooth

The probe with microcrystalline particles has a large surface area with high adsorption. Without hurting the skin, this probe can massage and exfoliate.

Comedo suction, skin tightening

The large circular probe has strong suction. It is best used for oily spots and blackheads. It has a unique mechanism that uses vacuum pressure to massage the skin by micro-lift which deeply effects the underlying skin. Stimulates micro-circulation operation, promotes blood and lymphatic circulation, let the skin rejuvenate again.

Cleansing pores

The small circular probe's suction is weaker. Place the round hole on the blackhead acne to remove. If the skin gets hurt or swollen, it is not recommended to use it. If the blackhead is deep in the pores, use the larger size probe. This small size probe is more suitable for sensitive skin.

Smooth fine lines

The special oval probe is designed to have a micro puffing massage for the canthus, ala nasi, the corner of the mouth and other areas that have fine lines, areas difficult to care for, or sensitive skin. It stimulates the micro circulation, preventing the appearance of wrinkles.

Tips

- 1 Tightening function shapes and pulls the skin, which may cause injury or scarring. Thus, it is not recommended to keep the probe in the same place for too long, generally 3-5 seconds is appropriate.
- 2 It is recommended that you use the microcrystalline particle probe depending on your skin type: once a week for dry skin, twice a week for oily skin, and no more than 5 minutes each time per type.
- 3 Please don't press the probe firmly on your skin; the machine has strong suction.
- 4 After skin care, the skin may feel a little tingling and flushed because of blood circulation. Depending on your skin type, the redness phenomenon will automatically fade away between 0.5-24hours.
- 5 Please pay more attention to moisturizing and sun care after skin care.
- 6 Please do not use scrub cream or exfoliator before or after skin care.

Working Principle

Skin-smooth and tighten function

After cleansing the face, open the pores with a hot towel or water vapor and keep your face dry before use;

- 1 Change to the microcrystalline probe or large probe;
- 2 Start with your chin and move the device from bottom to top;
- 3 Use the device on your cheeks and move it from the inside to the outside;
- 4 Use the device in the forehead area and move it from the center to the sides;
- 5 Use the device in the T-zone and move it from bottom to top.

Clean the pore

- 1 Change to large or small size circular probe (depending on your skin type);
- 2 Place the round hole on the blackhead acne to remove it. The levels can be changed according to your need;

Method: when working on a blackhead, please move slowly in the same direction. Repeat several times to avoid redness of the skin. Please don't stay in one location for too long.

Smooth fine lines


Change to the oval probe, use on canthus, alanas, the corner of the mouth and other areas where it is easy to have fine lines. The levels can be changed according to your need;

Method: Suction on the areas of need to smooth fine lines, pull away quickly and gently again and again. It stimulates the micro circulation, preventing the appearance of wrinkles.

(Attention: Please don't suction on the same area all the time, it will make the skin swell. If the skin swells, stop using immediately.

Product Feature

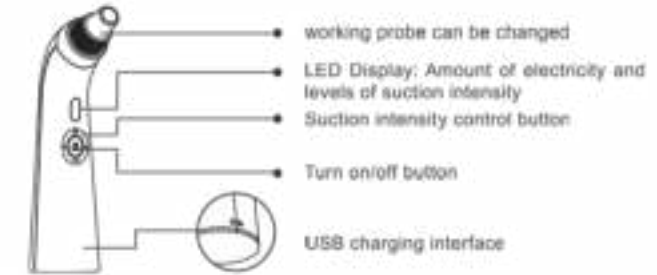Using area

Microcrystalline Probe:

Any area on the face except the area around the eyes. Also suitable for elbows, knees and ankles.

Big Circular Probe:

Any area on the face except the area around the eyes

Small Circular Probe:

Sensitive areas like the area around the eyes or any other areas of the face.

Oval Probe:

Sensitive areas like the area around the eyes or any other areas of the face.

Attention:

- 1 Please dry your hands before using the device. The water may get into the device from its air exit hole, which will damage the device;
- 2 Before using it, please make sure the sponge has been put into the filter cup. The sponge must be dry, or the liquid will damage the device;
- 3 Please don't wash the device;
- 4 This product is not for people who are physically or mentally disabled or lack experience and knowledge (including children) without a responsible parent or guardian present.
- 5 This product is not suitable for body parts that have had cosmetic surgery or rhinoplasty. Pregnant women and heart patients should not use this product.

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