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
MEDISANA FTF THERMOMETER WITH FLEXIBLE SENSOR TIP

Operating Instructions. Please read carefully!
Thank you for your confidence!

To ensure the best results and long-term satisfaction with your MEDISANA FTF Thermometer, we recommend that you read the following operating and maintenance instructions carefully.

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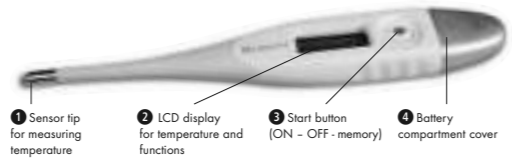SAFETY INFORMATION
1. The thermometer contains small parts that could be swallowed by children. Do not let children handle the thermometer without adult supervision.
2. Only use the thermometer for measuring body temperature.
3. The unit is not designed for commercial use.
4. If a fault occurs, please do not try to repair the unit yourself. Discontinue further use of the unit. Contact your dealer and only have repairs carried out by MEDISANA service personnel.
5. Please take note of the operating and storage conditions to ensure that your MEDISANA FTF thermometer will continue to function correctly for a long time.

USEFUL INFORMATION ON TAKING TEMPERATURE WITH THE MEDISANA FTF THERMOMETER

What is fever?
Our normal body temperature is between 36°C and 37,5°C. Fluctuations are quite normal, depending on time of day and the body's condition, for example in the female menstruation cycle. However, when body temperature rises above 38°C, the phenomenon is known as fever. If the fever should exceed 39°C and still persist after 1 to 2 days, you should definitely contact your doctor.

Taking a temperature reading - quickly and accurately!
A fever is often associated with fatigue, pains in the limbs, shivering and chills. The course of an illness can be best monitored by taking a temperature reading 2 to 3 times a day at 4 to 6 hour intervals. Readings should always be taken at the same part of the body and approximately the same time so that the results can be compared.

The MEDISANA FTF Thermometer enables you to quickly and simply determine the current body temperature. Typical errors, such as the duration of measurement being too short can be almost completely eliminated. A precise result is displayed after just 10 - 20 seconds, and this result is stored until the next reading is taken. Temperature readings can be taken at the following parts of the body:
a) Under the tongue (oral)
b) In the armpit (axillary)
c) In the rectum (rectal)

GETTING STARTED

Inserting / removing battery:
The MEDISANA FTF high-speed thermometer is ready to use right away, as the battery is already installed in the unit. However, if the „▼“ symbol should appear on the right of the display after the unit has been in operation for some time, you need to replace the battery. To do this, remove the battery compartment cover (1) and lift the compartment opening. 0.5cm out of the thermometer. Use a suitable tool, e.g. the tip of a ballpoint pen, to carefully lever out the battery. Insert the new battery with the „+“ sign facing upwards. Push the battery compartment back into the thermometer and replace the compartment cover.

OPERATION

Taking a reading
Press the start button (1) once to switch the unit on. All display symbols will appear for approx. two to three seconds and you will hear a short tone. If the temperature is below 32°C, you will see an „L“ and a flashing „°C“ on the LCD display (2). The MEDISANA FTF thermometer is now ready to take a temperature reading. Insert the tip of the sensor (1) either orally, into the armpit or into the rectum. The thermometer will now automatically begin to take a reading. During the reading, the current temperature will be shown continuously. As soon as the variation in temperature is below 0,1 °C in 8 seconds, you will hear a tone and the reading taken will be shown on the display (2). The thermometer will switch off automatically ten minutes after taking the reading. To prolong the life of the battery, however, you should switch off the thermometer manually by pressing start button (1) once you have taken the reading.

Taking temperature readings at various parts of the body
The MEDISANA FTF thermometer is suitable for oral (in the mouth), axillary (in the armpit) and rectal (in the rectum) temperature readings. **Oral reading takes approx. 60 seconds.** To ensure safe and accurate measurement, you should take account of the following information when positioning the thermometer:
ORAL READINGS:
Place the sensor tip (1) of the thermometer under the tongue and close your mouth. Breathe through your nose while taking the reading in order to avoid inaccurate results.

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RECTAL READINGS:

Many patients find rectal measurement uncomfortable, but it delivers precise results. Thanks to the sensor's flexible tip, the MEDISANA FTF thermometer makes measuring temperature in this part of the body much easier. To quickly and reliably determine body temperature, carefully insert the sensor tip (1) approx. 3,5 cm into the rectum. After approx. 60 seconds, you can read the result on the display.

USING THE PROTECTIVE PROBE COVERS
The MEDISANA FTF Thermometer package includes 3 protective probe covers intended for single use (e.g. for rectal use). If you wish to use a protective cover, pull it gently apart on the side marked and insert the tip of the thermometer into the opening. Now hold the upper end of the protective cover tightly against the thermometer and slowly pull the paper strip and the projecting material downwards. The thermometer is now ready to use with the protective cover. After taking the reading, pull off the protective cover and dispose of it.

CLEANING AND DISINFECTION
The MEDISANA FTF thermometer is waterproof, making it easy to clean and disinfect. Clean the thermometer before and after each use with a cloth moistened in a mild solution containing alcohol. You may also soak the thermometer in a suitable disinfectant solution. Your dealer will be happy to recommend suitable cleaning products and give you advice on their use.

USEFUL INFORMATION ON TAKING TEMPERATURE WITH THE MEDISANA FTF THERMOMETER

What is fever?
Our normal body temperature is between 36°C and 37,5°C. Fluctuations are quite normal, depending on time of day and the body's condition, for example in the female menstruation cycle. However, when body temperature rises above 38°C, the phenomenon is known as fever. If the fever should exceed 39°C and still persist after 1 to 2 days, you should definitely contact your doctor.

Taking a temperature reading - quickly and accurately!
A fever is often associated with fatigue, pains in the limbs, shivering and chills. The course of an illness can be best monitored by taking a temperature reading 2 to 3 times a day at 4 to 6 hour intervals. Readings should always be taken at the same part of the body and approximately the same time so that the results can be compared.

The MEDISANA FTF Thermometer enables you to quickly and simply determine the current body temperature. Typical errors, such as the duration of measurement being too short can be almost completely eliminated. A precise result is displayed after just 10 - 20 seconds, and this result is stored until the next reading is taken. Temperature readings can be taken at the following parts of the body:
a) Under the tongue (oral)
b) In the armpit (axillary)
c) In the rectum (rectal)

GETTING STARTED

Inserting / removing battery:
The MEDISANA FTF high-speed thermometer is ready to use right away, as the battery is already installed in the unit. However, if the „▼“ symbol should appear on the right of the display after the unit has been in operation for some time, you need to replace the battery. To do this, remove the battery compartment cover (1) and lift the compartment opening. 0.5cm out of the thermometer. Use a suitable tool, e.g. the tip of a ballpoint pen, to carefully lever out the battery. Insert the new battery with the „+“ sign facing upwards. Push the battery compartment back into the thermometer and replace the compartment cover.

OPERATION

Taking a reading
Press the start button (1) once to switch the unit on. All display symbols will appear for approx. two to three seconds and you will hear a short tone. If the temperature is below 32°C, you will see an „L“ and a flashing „°C“ on the LCD display (2). The MEDISANA FTF thermometer is now ready to take a temperature reading. Insert the tip of the sensor (1) either orally, into the armpit or into the rectum. The thermometer will now automatically begin to take a reading. During the reading, the current temperature will be shown continuously. As soon as the variation in temperature is below 0,1 °C in 8 seconds, you will hear a tone and the reading taken will be shown on the display (2). The thermometer will switch off automatically ten minutes after taking the reading. To prolong the life of the battery, however, you should switch off the thermometer manually by pressing start button (1) once you have taken the reading.

Taking temperature readings at various parts of the body
The MEDISANA FTF thermometer is suitable for oral (in the mouth), axillary (in the armpit) and rectal (in the rectum) temperature readings. **Oral reading takes approx. 60 seconds.** To ensure safe and accurate measurement, you should take account of the following information when positioning the thermometer:
ORAL READINGS:
Place the sensor tip (1) of the thermometer under the tongue and close your mouth. Breathe through your nose while taking the reading in order to avoid inaccurate results.

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