


From Flat2Fit™ in Minutes.



Scan to watch fitting video

SOVA Aero Night Guard gives a true custom fit to help relieve the pain from bruxism. SOVA allows you to remold until you have a perfect fit.

Starts flat

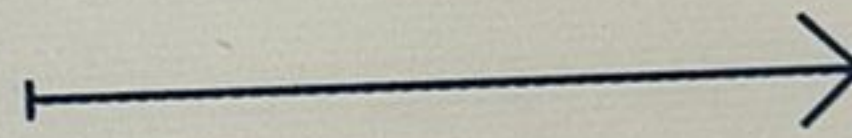

Perfect fit



DO NOT FIT ON YOUR OWN, BUT CONSULT A DENTIST BEFORE USE IF YOU HAVE: • Orthodontics — you must have the SOVA Night Guard fitted by your orthodontist or dentist. • Missing teeth. • Clicking of your jaw. • Sores in your mouth or bleeding of the gums. • If you have been diagnosed with temporomandibular joint (TMJ) problems. WARNING — DO NOT USE: • With artificial dentures. • If you have loose teeth, loose fillings, loose caps or cavities with no fillings. • An improperly fitted night guard. • As an athletic mouthguard. Instead please go to www.sisuguard.com to see our athletic line of SISU® Mouthguards. CONSULT A DENTIST IF: • The same symptoms continue after several months of use. • If you experience pain while using the SOVA Night Guard. First, try remolding using less suction. DISCONTINUE USE IF: • Soreness, bleeding, or other reaction inside of your mouth. • You notice new symptoms such as pain in your jaw, teeth, ear, neck stiffness, headache or joint clicking after use. • You have loose teeth or a change in your bite that lasts more than a few minutes after removing. • The SOVA Night Guard has been worn or bitten through or show signs of excessive wear.



MADE IN THE USA

Akervall Technologies Inc.

1512 Woodland Dr. Saline, MI 48176 USA
1.800.444.0570 | sovanightguard.com
customerservice@sovanightguard.com
sovanightguard.com/patents



FDA Cleared