


英语

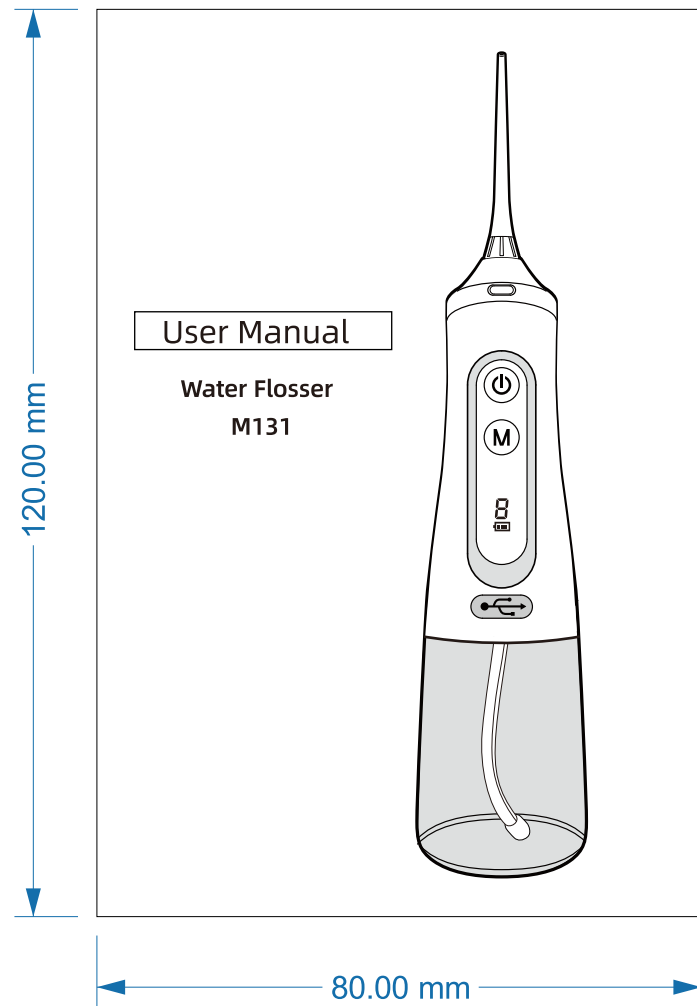

User Manual

Water Flosser
M131

Please keep this manual in a safe place and read it carefully before use.

Safety Notices

When using electrical appliances, especially children users, should pay attention to the basic safety precautions, read the operation guide before use, and the following obvious marks and warnings have been made for important information to reduce the risk of burning, electric shock, fire or injury to human body.

WARNING Incorrect handling may result in death or serious injury

- When going out, the power plug of the product should be unplugged.
- Children and disabled persons should use it under the supervision of adults.
- Do not store liquid in the water tank when carrying or storing outside.
- This product is only used for the purpose indicated in this manual. Do not use accessories not recommended by the manufacturer.
- This product cannot be used under the following circumstances:
 - The power cord or plug is damaged.
 - Not working properly.
 - In case of falling or damage.
 - Keep the power cord away from the overheated enclosure.
 - Do not insert any sundries into the tube or product.
 - Do not use power supply and frequency other than those indicated.
 - When the product is used close to other electrical products, electrical appliances may interfere with each other. Please use them separately.
 - This product can only be used by one person. Cross infection or aggravation may occur when multiple people use it. If multiple people use it, please replace the nozzle assembly.


2

Basic precautions:

- Read the relevant materials in the operation guide before use.
- Do not use accessories not officially designated and recommended.
- Before connecting the power supply, please check whether the marked charging voltage is consistent with the supply voltage.
- If the product is damaged in any form (nozzle/host machine), please stop using it.
- Please keep the oral irrigator away from high temperature, high heat, fire source, heat source, strong light source and direct sunlight.
- Please ensure that the safe socket is used for charging, and do not force the plug into the socket.
- In order to ensure safety, people (including children) who are physically or mentally impaired or lack relevant experience and knowledge need to be supervised or guided for safe use.
- This product contains detachable accessories such as nozzle and water tank. Please take care of children to prevent accidental ingestion of nozzle accessories.
- The oral irrigator is a personal item and cannot be used by dental clinics or institutions for multiple patients.
- Please cut off the power when the mouth rinser is not charged or before cleaning, and do not clean this product in the dishwasher.
- Do not let hair care products (such as gel, oil) or hand sanitizer contact this product to avoid damaging the surface.

3

Enjoy brighter, cleaner teeth and healthier gums in just 1 minute a day

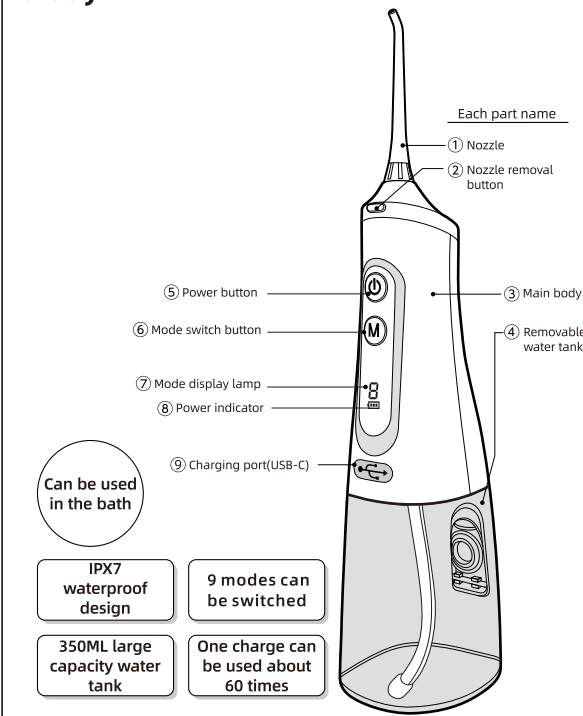

Can be used in the bath

- IPX7 waterproof design
- 350ML Large capacity water tank
- 9 modes can be switched
- One charge can be used about 60 times

4

About charging

If the machine suddenly stops working or the battery icon flashes at the same time when you press the power button during use, it indicates that the remaining battery is low and needs to be charged.

How to charge

- Open the silicone cover of the charging port of the machine, connect the Type-C charging cable provided by the machine, insert it into the Type-C port, and connect the other end of the charging cable to the charging adapter with 5V/1A output, and then connect it to the power supply.
 - *If there is water in the charging port of the host, please wipe it with a dry cloth, then connect the charging cable, and start charging after it is dry.

- When the charging power supply is connected, the battery icon indicator lights flash in turn, indicating that charging is in progress.
- When all battery icon indicators are on, it indicates that the battery is fully charged.

How to use

- Nozzle mounting**
Insert the nozzle into the nozzle hole on the top of the host and hear it click, indicating that it is installed.

5



Company website

We provide lifetime customer service for this product, if there is any problem with the product, please feel free to contact us. We will do our best to serve you.

- Add cleaning solution:**
Open the water tank cover, place the body horizontally, fill in an appropriate amount of water, and cover the water tank cover.

Cleaning fluid available:

- On the right is a list of liquids that can be put into the tank and used for cleaning. Please do not use any other cleaning fluid, which may cause failure.
- Tap water
 - Warm water (40°C or lower)
 - Mouthwash

- Mode setting:**
Tap the mode button to select the required gear; Each time you press it to cycle from gear 1 to gear 9, the display will display the corresponding gear.
1 → 2 → 3 → 4 → 5 → 6 → 7 → 8 → 9

- Power on**
Put the head of the nozzle into the mouth, and press the power button to start the machine.

7

If no water sprays out for a long time, please unplug the nozzle and press the water valve in the water outlet intermittently with a tool with a diameter of less than 6mm until water sprays out.

- Cleaning**
During the cleaning process, keep the body forward and slightly open the mouth during the cleaning process, so that the food residues in the mouth can be flushed out.
*Please note that if you open your mouth vigorously, the washing water will splash everywhere.
● The main body is equipped with a waterproof mechanism, so you can use it with confidence even if it gets wet with washing water.
● The direction of the nozzle can be freely rotated 360°. Please adjust it to an easy to use direction when using. Tilt for ease of use.

Point of care

The nozzle is aimed between the teeth and gums, and then slowly spray along the gums from the molar (the back teeth) to the front teeth, as well as the surface and back of the teeth.

8

- Stop**
When you want to pause use, tap the power button to shut down; The machine will automatically shut down after each gear is set to run for 2 minutes.

- How to remove the nozzle**
While pressing the nozzle release button, pull the nozzle up.

About maintenance

- After use, open the water tank cover and remove the accumulated water in the water tank;
- If there is dirt inside the water tank, you can pour a small amount or less of 40°C warm water into the water tank, cleaned it slightly, and then the above operations are carried out.
- To keep it clean, rinse the nozzle with clean water and dry it after each use.
- When cleaning this product, please gently wash away the dirt with clean tap water (flowing water).
- After cleaning, wipe off the water with a dry cloth and dry it. If the nozzle is damaged, replace it with the attached spare nozzle.

9

Troubleshooting

Problem	Possible Causes	Solution
Can not be used	Power button not pressed.	Please press the power button.
	Power LEDs blink.	Please charge the battery.
Water pressure is very low	The nozzle is worn or deformed. Nozzle is clogged with dust.	Replace with a new nozzle.
	There is no water in the tank.	Fill the tank with water.
Water does not come out	The tip of the tube inside the tank is above the water surface.	Adjust the direction of the main unit so that the tip of the tube is below the surface of the water.
	Clogged water tank air inlet hole.	Please clean the air inlet hole foreign matter.
Short usage time after charging	The charging cable is not properly connected.	Make sure that the charging cable is securely connected to the main unit and the USB output charging adapter, and that the charging adapter is connected to power outlet.
	Battery life	Even if charging is completed normally, the usage time is short, it may be the end of the battery life. Contact service center for support.

10