CARMEN

ILL 11 / CIL1001W • INFRARED LAMP

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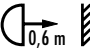
Explanation of symbols

	Read the instructions!		Manufacturer
	Caution: hot surface		This product satisfies the requirements of the applicable European and national directives.
	The device has double protective insulation and therefore complies with protection class 2.		Dispose of packaging in an environmentally friendly manner.
	WARNING; Warning of risk of injury or danger to your health.		Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).
	CAUTION; Safety information about possible damage to the device/ accessories.		NOTE; Refers to important information.
	Keep away from rain		Serial number
	Distance between infrared heat lamp and irradiated body part		

Important safety instructions - read carefully and keep for later reference!

WARNING

- Non-observance of the following notes may result in personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information is not only intended to protect your health and the health of others, it should also protect the product. For this reason, pay attention to these safety notes and include these instructions when handing over the product to others.
- This infrared lamp is only intended for radiating the human body.
- Always use the protective goggles supplied during usage.
- Always use the infrared lamp with the protective cover in place.
- Never remove the protective cover during use. Only remove the protective cover when changing the infrared lamp.
- Before use, ensure that all packaging materials are removed and that there is no visible damage to the unit or accessories. When in doubt, do not use the unit and contact our Customer Services department.
- Always ensure that the appliance is on a secure and even surface.
- The lamp housing and the infrared lamp become very hot during operation. There is a risk of burning if these are touched! Always let the lamp cool down first before touching it.
- Inflammable objects must be kept at a distance of at least 1 meter from the infrared lamp. Do not hang the appliance on the wall or from the ceiling.
- Connect the unit only to the mains voltage listed on the type plate.
- **Do not immerse the appliance into water and do not use it in wet rooms.**
- This device is not intended for use by persons (including children) with restricted physical, sensory or mental capabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety or they are given instructions on how to use the device.
- Do not use on babies and small children.

- Children should be supervised to ensure that they do not use the device as a toy.
- This device is not intended for use on changing units, cots or similar.
- Do not use it on persons sensitive to heat.
The sensitivity to heat can be reduced or increased in the following cases:
 - diabetic patients;
 - persons with sleepiness, dementia or difficulties concentrating;
 - persons with illness-related skin changes or persons with scarred skin areas in the application area;
 - persons with allergies;
 - children and older persons;
 - after consuming medications or alcohol.
- Do not operate the device if you are wearing an implant, particularly if it is located in the area of the body to be treated. If you have any doubts, consult your doctor.
- If you are pregnant, do not operate the device near your abdomen, in order to prevent this area from overheating. If you have any doubts, consult your doctor.
- In case of acutely inflammatory processes, radiation should only be carried out after consulting a doctor.
- Always limit the period of use and check the reaction of the skin.
- Medications, cosmetics or food can lead to a over-sensitive or allergic reaction of the skin under certain conditions. Radiation must be ended immediately in this case.
- An appropriate distance between the infrared lamp and the desired treatment area (not less than 60 cm / 24 inch) must be maintained during treatment so that no skin irritation, such as tingling, itching, burning or too high heat sensitivity occurs. Stop the treatment if this occurs. Alternatively, increase the radiation distance and/or reduce the radiation duration.
- This device is not intended for radiation from above.
- The minimum distances described in these instructions for use must be adhered.
- If the appliance is operated continuously, caution and attention are advisable - never fall asleep during radiation!
- Prolonged use can cause burns to the skin.
- Do not touch the infrared lamp or unscrew it when it is connected to the mains.

- Children do not understand the risks associated with electrical appliances. Make sure the appliance is not used by children without appropriate supervision.
- The appliance may be used only under supervision.
- Keep children away from packaging materials (risk of suffocation).
- Do not cover up or pack away the unit while it is warm.
- Always pull out the plug from the socket and let the appliance cool off before you touch it.
- Do not touch the appliance with wet hands when it is connected to the mains; water should not be sprayed on the appliance. Operate it only when it is completely dry.
- Protect the unit against impacts.
- Do not pull the mains plug from the socket by the cord.
- In the event of a defect or malfunction, switch the device off immediately and disconnect it from the power supply. If the power cord of the device is damaged, it must be replaced by the manufacturer or an authorised customer service technician in order to avoid hazard.
- Disconnection from the mains supply is only guaranteed if the mains plug has been removed from the socket.
- There is a risk of strangulation for babies and children who become trapped in the mains cable.
- This product is not suitable for lighting a room.
- This unit is not intended for commercial or clinical use, but only for individual, private household use!
- If you have health concerns of any kind, contact your general practitioner!
- Lighting is excluded from the warranty.
- Do not open the device. Please do not attempt to repair the unit yourself. This could result in serious injuries. Failure to observe this will invalidate the warranty.
- In case the cable or the housing is damaged, please contact our Customer Services department, as special tools are required for repairing it.
- If you should have any questions about using our devices, please contact our Customer Services department.

Appliance description

1. Infrared lamp
2. Protective cover
3. Adjustable lampshade (5 positions)
4. Base
5. Protective goggles

Power consumption: 100W

Mains connection: 220-230V ~ 50-60 Hz

Max. Radiation intensity: 680W/m² at 60 cm

Operating conditions:

Temperature: +5°C until +35°C

Relative humidity: 15%-93%

Expected operating life of the device: 7 years

Subject to technical modifications. This device meets the requirements of the European Medical Devices Directive 93/42/EEC with amendments as well as the German Medical Devices Act.

The serial number is located on the device or in the battery compartment.

Instructions for use

This infrared lamp is only intended for radiating the human body.

Infrared radiation causes warmth to be transported to the human body. The blood supply in the radiated skin is increased and the metabolic turnover in the field of temperature is enhanced. The mode of action of infrared light encourages the body to heal; healing processes are thus assisted in a specific area. Infrared light can be used as a parallel treatment in diseases of the ear, nose and throat and to assist in facial and beauty care, especially in case of impure skin.

However, first ask your family doctor whether application in the individual case is advisable from a medical standpoint.

Initial use

Before you use the appliance for the first time, please do as follows: Carefully unpack the appliance and remove all the packaging material. Keep the material (plastic bags and cardboard) out of reach of children. Before use, check whether the infrared lamp is properly screwed in. Also check that the protective cover is securely in place. Always ensure that the appliance is on a secure and even surface.

The appliance has an adjustable lampshade. The 5 positions of inclination make it possible for you to adjust the infrared lamp according to your needs. Set the angle of the fluorescent screen in accordance with your individual requirements. The distance between the infrared lamp and the radiated part of the body **should not be less than 60 cm (24 inch)**. Always use the protective goggles supplied during usage. See chapter protective goggles for more information. The protective goggles corresponds to the requirements of EN 171:2002 and serves to protect the eyes.

CAUTION: Never look directly into the infrared light without protective goggles.

The recommended duration of treatment is **maximum 12 minutes**.

We recommend selecting a shorter radiation period to start with. However, the skin of some people may show an over-sensitive (e.g. with extreme redness, blistering, itching, heavy sweating) or allergic reaction appear. Examine the radiated section of the body regularly. Discontinue radiation immediately if any signs of an over-sensitive or allergic reaction appear and consult a physician.

Switching on

Push the plug completely into the mains socket. Connect the unit only to the mains voltage listed on the type plate.

NOTE

- Ensure that there is a socket near the place where the lamp will stand.
- Lay the mains cable in such a way that no one can trip over it.

Switching off

Pull the mains plug from the socket.

NOTE

The lamp is hot after being used. Let the lamp cool off first long enough before you put it away and/or pack it up! Please note that the mains cable must not be wound round the infrared lamp either for storage or when in use.

Protective goggles

The product complies with the PPE regulation (EU) 2016/425 (reference: OJ L 81, 31/03/2016, p. 51-98), the harmonised standard EN 166:2001 (reference: OJ C 209, 15/06/2018, p. 18) and the standard EN 171:2002.

1.1 Warnings and safety notes

Read these instructions for use carefully! Non-observance of the following information may result in personal injury or material damage. Store these instructions for use and make them accessible to other users. Make sure you include these instructions for use when handing over the product to third parties.

WARNING

- Keep packaging material away from children. There is the risk of suffocation.
- Children must not play with the protective goggles.
- Before using the protective goggles, check it for signs of external damage (e.g. cracks or holes). If the protective goggles is damaged, you must NOT continue using it.
- When wearing the protective goggles, people with sensitive skin may have an allergic reaction to the materials of the goggles that may come into contact with the skin.
- Make sure that the protective goggles does not come into contact with grease, oils or alcohol (e.g. cleaning products or facial toners).
- Please note that the protective goggles is opaque, which means it will restrict your spacial awareness.
- Do NOT operate the infrared lamp whilst wearing the protective goggles.
- Do NOT move about whilst wearing the protective goggles.
- Make sure that you do NOT reduce the distance between you and the infrared lamp whilst wearing the protective goggles.
- Use the protective goggles only in quiet environments, so that you are not disturbed during use.

1.2 Intended use

The protective goggles is exclusively intended for protecting your eyes against infrared radiation while using an infrared lamp. The protective goggles provides protection against infrared radiation only. The protective goggles does NOT provide protection against other types of radiation. Use the protective goggles only with this type of infrared lamps. Never use the protective goggles with infrared lamps from other manufacturers.

Always use the protective goggles directly over your eyes, not over another pair of glasses. The protective goggles is intended for domestic/private use only, not for commercial use.

The protective goggles is not suitable for use as eye protection against mechanical hazards (e.g. impacts).

It is free from pollutants as per the legal requirements. The material of the protective goggles may initially smell. This smell will dissipate after a few uses.

1.3 Usage

- Position the infrared lamp as described in the instructions for use.
- Prepare the desired treatment distance.
- Switch on the infrared lamp as described in the instructions for use for the infrared lamp.
- Sit down in preparation for the application.
- Hold the protective goggles to your eyes and pull the rubber band over your head. To make the rubber band longer/ shorter, pull on the strap at the side.
- To end the application, first take off the protective goggles. Then switch off the infrared lamp.

Replace the product after a period of use of 3 years.

1.4 Storage

If you do not plan to use the protective goggles for an extended period of time, we recommend storing it in its original packaging in a dry environment, ensuring that it is not weighed down by other items.

1.5 Cleaning

Clean the protective goggles from time to time. Use a slightly damp, soft cloth to clean it. Do not use any sharp or abrasive cleaning products.

1.6 Technical data

Model	EP 10
Protection level	4 -10 in accordance with EN 171:2002
Protective eyewear designation	4-10 BEU CE The code for the lenses is comprised of the following: <ul style="list-style-type: none">- Initial number for IR protection filter (1st part of code)- Protection level number (2nd part of code)- Manufacturer's code (3rd part of code)- CE mark (4th part of code)
Light transmission level	Max. 0,023%, min. 0,0085%
Maximum average spectral transmittance value with infrared	780 - 1400 nm: 0,050% 780 - 2000 nm: 1,3%

Changing the lamp

The device must always be switched off, disconnected from the mains and allowed to cool down before changing the infrared lamp.

The procedure for replacing the infrared lamp is as follows:

1. To remove the protective cover, unscrew the screw on the side using a screwdriver.
2. Carefully pull the protective cover off the device.
3. Unscrew the old infrared lamp from the device.
4. Fully screw the new infrared lamp into the device.
5. Place the protective cover back on the device.
6. Tighten the screw on the side of the protective cover using a screwdriver.
7. Check that the protective cover is firmly attached to the device.

Please note that only infrared lamps of the original type, R95E (100W), may be used as replacements.

WARNING

- Do not touch the infrared lamp or unscrew it when it is connected to the mains.
- Unplug the device, to avoid the risk of electrocution.
- In the interest of environmental protection, do not dispose of the device in household waste.

Cleaning & maintenance

To avoid odours caused by fluff and dust, we recommend regularly cleaning the device as follows:

- Before every cleaning, the appliance must be switched off, disconnected from the mains and must be cool.
- Clean using a damp cloth, to which you can apply a little detergent if necessary.
- Ensure that no water penetrates the appliance. Use the appliance again only when it is completely dry.

NOTE

- Do not use abrasive detergents and never immerse the unit in water.
- Do not wash the unit in a washing machine!
- Do not touch the unit with wet hands while it is plugged in; do not allow any water to be sprayed onto the unit. The unit must be operated only when it is completely dry.

Disposal

Legislation requires that all electrical and electronic equipment must be collected for reuse and recycling. Electrical and electronic equipment marked with the symbol indicating separate collection of such equipment must be returned to a municipal waste collection point.

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Ordering Accessories - EN

Go to www.carmen.nl/support to purchase accessories or spare parts.

Warranty conditions - EN

A 24 month warranty is provided for this product. Your warranty is valid if the product has been used in accordance with the user instructions and for the purpose for which it was made. Also the original proof of purchase must be submitted, stating the date of purchase, the name of the retailer and the item number of the product. For extensive warranty information you are referred to our website: www.carmen.nl/support

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Carmen
Postbus 159
6920 AD DUIVEN
The Netherlands

support@carmen.nl

Download the most recent
manual from www.carmen.nl