Comforheat Model 99360

Handleiding 2-in-1 warmte- en massagekussen

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www.trebs.nl

IMPORTANT SAFETY INSTRUCTIONS

WARNING: this product contains magnets. It can be harmful to those who have a pacemaker, defibrillator or other medical implants. Failure to follow this warning can result in serious injury.

CAUTION: handle with care. Always lay flat or fold carefully. Crushing or crimping of the heating wires can lead to premature failure

DANGER: to reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions.

This pad may not be used by children under the age of 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been supervised/instructed and understand the hazards involved, they may use. Children shall not play with the pad. This is not a toy. Cleaning and user maintenance shall not be done by children unless they are older than 8 and supervised. Children under the age of 3 may not use this pad due to their inability to react overheating.

1. Read all instructions
2. Control that your mains voltage corresponds to that stated on Ensure the rating plate of the pad.
3. Do not use while sleeping.
4. Do not use on children under 3 years.
5. This may not be used on or by an invalid, sleeping or unconscious person, or a person with poor blood circulation unless carefully attended.
6. Turn off the power when not in use and unplug the supply cord from the main outlet.
7. Keep dogs, cats and other pets away from the pad.
8. It's for indoor use only
9. The pad has a heated surface so persons insensitive to heat must be careful when using the pad. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
10. Burns can occur regardless of control setting so check the skin under the heating pad frequently.
11. Never use the heating pad without a cover in place.
12. Do not use pins or other metallic means to fasten this heating pad in place.
13. Do not sit on or crush the heating pad - avoid sharp folds.
14. Never pull this heating pad by the supply cord. And do not use the cord as a handle.
15. Carefully examine the inner cover before each use. Discard the heating pad if the inner covering shows any sign of deterioration.
16. This pad must not be used by persons sensitive to heat and other vulnerable persons who are unable to react to overheating.
17. When storing the pad, allow it to cool down before folding.
18. Do not crease the pad by placing items on top of it during storage.
19. This pad is not intended for medical use in hospitals.
20. Do not use when it is wet.
21. Too long use can result in skin burns, so be careful.
22. When continuous use; adjust always a low level.
23. The switch or control unit must not be allowed to get wet during washing or drying and make sure the cord must be positioned to ensure that water does not flow into the switch or control.
24. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
25. In order to avoid a hazard due to inadvertent resetting of the thermal cut-out, this pad must not be

supplied through an external switching device such as a timer, or connected to a circuit that is regularly switched on and off by the utility.


26. When not in use, let the pad cool down before folding. Fold neatly but not tightly (or roll) without sharp bends in the heat element. Do not store other objects upon it.
27. Do not use the pad in areas of the body that are inflamed, injured or swollen. If in doubt, you should seek medical advice before use. This also applies for diabetic persons.
28. This pad can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the pad in a safe way and understand the hazards involved. Children shall not play with the pad. Cleaning and user maintenance shall not be made by children without supervision.
29. Do not leave packaging components (plastic bags, cardboard, polystyrene, etc.) within the reach of children, as it could be dangerous.
30. Do not let the cord hang over the edge of a table or counter, or touch hot surfaces.
31. The cord should not come into contact with water or other liquids.
32. Do not place on or near a hot gas, electric burner or on ovens.
33. Never use the pad outside and always place it in a dry environment.
34. To avoid any risk of electrical shock do not immerse cord, plug, switch or pad in water or other liquid.
35. Never use the pad when it is wet. Do not use the electric pad while standing on wet surface or with wet hands.
36. Unplug the pad before removing additional parts.
37. **Prolonged use can cause burns on the skin. Caution!** Avoid to fall asleep while the warming heating pad is in operation.

Warning: Check the pad frequently for signs of wear or damage. If there are such signs or if the pad has been misused, it must be returned to the manufacturer or dealer for check prior to any further use.

1. If the mains cable is damaged, it should only be replaced at a repair workshop approved by the manufacturer, as special tools are required for the repair.
2. Repairs to the electrical pad must be carried out only by trained personnel. Improper repair could result in considerable risks for the user.
3. When the pad is not in use, the plug should be removed from the mains and the pad.
4. Unplug the cord from the mains:
 - when malfunction occurs during use
 - before cleaning
 - after use and for storing
5. Do not operate the pad with a damaged cord, if it malfunctions or has been damaged in any manner. Ask a (*) competent qualified electrician to check and repair it.
6. Should the cord be damaged, it must be replaced by a competent qualified electrician(*) to avoid any danger. All repairs should be made by a competent qualified electrician(*).

(*)Competent qualified electrician : after sales department of the producer or importer or any person who is qualified, approved and competent to perform this kind of repairs in order to avoid all danger. In any case you should return the pad to such electrician.

USING INSTRUCTIONS

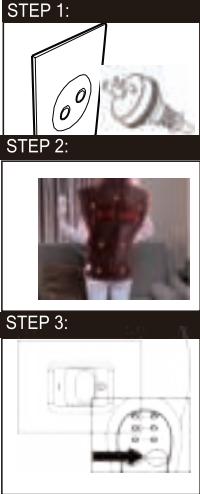
Plug the heat and massage pad into a 220V-240V outlet.

Place the pad on your body as shown on the picture so it will fit your body perfectly and fasten it with the clip closure.

Turn on the switch of the controller's backside and then select the desired heat or massage function on the front of the controller. **IMPORTANT:** if the heat indicators do not light up or stay on, unplug from the electrical outlet.


IMPORTANT: allow several minutes for the pad to preheat. If there is a red blinking "High" light, this means an error has occurred. Unplug the pad from the outlet and inspect the condition of the cord and the function of your main outlet. If blinking continues, stop using the pad and contact the Trebs customer service.

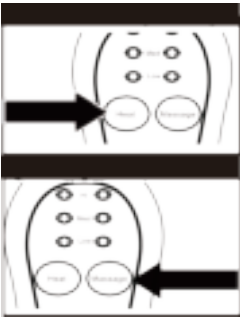To apply heat, turn on the switch at the controller's backside, then press the heat activation button and select your desired temperature: HI, MED or LOW
Press Low for low heat
Press Med for med heat
Press Hi for high heat

To activate the massage feature, turn on the switch at the controller's backside, then press the massage activation button to select your desired speed:
Press 1x - low pulse - Press 4x - med steady
Press 2x - low steady - Press 5x - hi pulse
Press 3x - med pulse - Press 6x - hi steady

Press OFF when stop using.

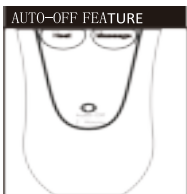The auto-off feature is included in the control on the heating pad. The auto-off allows to conserve energy by turning off after 2 continuous hours of use. The massage feature will automatically shut off after 15 minutes of use. **IMPORTANT:** the auto-off timer will reset after a power outage.

WASHING INSTRUCTIONS

DISCONNECT THE POWER BY UNPLUGGING FROM THE MAIN OUTLET.

NOTE: ONLY THE COVER IS WASHABLE! DO NOT WASH ABOVE 40°C. DO NOT WASH THE HEATING PAD! DO NOT DRY CLEAN! DO NOT USE BLEACH. DO NOT WRING.

1. Unzip the cover on the side and then carefully remove the cover.
2. Wash the cover in warm water (40°C) on gentle cycle, tumble dry on medium heat and remove quickly from dryer.
3. The heating pad can be hand washed in warm, soapy water. Hang to dry. Never fasten the pad with pins.
4. When dry, place the heating pad back inside the cover.

DISPOSAL OLD PADS

Please observe the environment! Dispose the packaging according to the local garbage regulations. Dispose paper and carton separately from plastic bags into the corresponding containers. Dispose the pad at the end of its lifespan according to environmental considerations. (Special waste collection WEEE of your local waste management authorities). Please dispose of the electric pad in accordance with the Ordinance on Waste Electrical and Electronic Equipment (WEEE)2012/19/EU or the applicable local regulations. In case of inquiries, please contact the municipal authorities responsible for waste disposal in your area.

GUARANTEE

The guarantee period for this pad is 24 months. The guarantee starts from the date of purchase (cash receipt). In case of malfunction you may send the pad to your dealer, add information about the malfunction and the original purchase receipt (which is obligatory for a repair free of charge)

Within the guarantee period we repair defects of the pad, which are based on material or manufacturing faults, in our choice by repair or exchange. These obligations do not extend the period of guarantee. No claim on a new period of guarantee results from it.

An exact description of the complaint reduces the processing time. We apologise for any inconvenience, that any claims of good-will are excluded:

- after expiry of the lifetime due to the construction and especially for batteries, accumulators, electric lamps...
- for consumption accessories and wear parts such as driving belts, motor coils, tooth brushes, grinding attachments, dough hooks, beaters, mixer attachments, round cutters etc.
- In case of breakage, as e.g. glass, porcelain or plastic
- in case of damages in transport, misuse, as well as in case of lack of care, cleaning and maintenance and in case of non-observance of operation or mounting instructions
- in case of chemical and/or electrochemical effects
- in case of operation with wrong type of the current voltage and in case of connecting to unsuitable power sources
- in case of abnormal environmental conditions and in case of strange operation conditions
- as well as if repairs or illegal abortions are made by persons, who are not authorised by us or

- if our pads are equipped with spare parts, supplementary parts or accessories, which are not original parts or if contact with unsuitable material have caused the damage

TECHNICAL DATA

Model: 2-in-1 heat and massage pad, 99360

Voltage: 220V- 240V ~ 50/60Hz

Power: 100 W

Product dimension: 66x85 cm

SERVICE

If you have any further questions about the product after reading this manual, please contact the Trebs customer service:

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