LA VAGUE

- FLASH TRAVEL Elektrische Reisezahnbürste mit USB
- FLASH TRAVEL Electric travel toothbrush (USB)
- FLASH TRAVEL Brosse à dents électrique de voyage (USB)
- FLASH TRAVEL Spazzolino elettrico da viaggio (USB)
- FLASH TRAVEL Cepillo de viaje eléctrico (USB)

VEHNSGROUP

Language

Page

Deutsch

3

English

11

Français

18

Italiano

26

Espanol

34

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Version 1.0 – Jan 2018

Table of contents

1	Features.....	12
2	Getting started.....	14
3	How to use	15
4	Cleaning and storage	16
5	Trouble shooting	16

1 Features

- 1 travel toothbrush with sonic technology with rechargeable lithium-ion battery
- 1 USB cable
- 1 Upper casing
- 2 Dupont brush heads

READ ALL INSTRUCTIONS BEFORE THE FIRST USE

Safety guide

Do not place or store FLASH TRAVEL in wet places.

Do not immerse FLASH TRAVEL in water or any other liquid.

Do not use when showering or bathing.

If FLASH TRAVEL has fallen into water or other liquids, unplug the power cord immediately from the wall outlet if this is safe. It is safe to take the toothbrush out of the water only after unplugging it.

Warnings

To reduce the risk of injury including burns, fire, or electrocution:

DO NOT place FLASH TRAVEL components in the dishwasher for cleaning.

DO NOT use any accessories other than those recommended by the manufacturer.

Discontinue use if FLASH TRAVEL appears damaged in any way (brush head, handle).

Keep the cable away from all heat sources.

Do not use a brush head with damaged bristles. Crushed or bent bristles can break off during brushing. Replace the brush head every 3 months or earlier if signs of wear occur. Please note that improper cleaning may damage the bristles. (See Cleaning)

Use FLASH TRAVEL only for its intended use as described in this manual. This product is designed to clean teeth, gums and tongue. If pain or discomfort occurs, stop using FLASH TRAVEL and contact your dentist or physician as this may indicate other health problems.

Persons (including children) with insufficient experience and knowledge or impaired mental, sensory or physical abilities should not use this product unless a person responsible for their safety has given appropriate instructions for the use of FLASH TRAVEL.

Children should be supervised to ensure that the product is used properly.

If your toothpaste contains peroxide, baking powder or other bicarbonate often found in toothpastes, clean your brush head and treat it thoroughly with water after each use. This prevents plastic fractures due to prolonged exposure to these chemicals.

Medical warnings

Consult your doctor before using FLASH TRAVEL if you have any medical problems.

If you have had jaw or gum surgery in the last two months, consult your dentist before using FLASH TRAVEL.

If excessive bleeding occurs after using FLASH TRAVEL or if bleeding persists after one week of use, contact your dentist as this may be a sign of other health problems.

If you have a pacemaker or other implanted device, consult your doctor before use. FLASH TRAVEL has been tested and meets the safety standards for electromagnetic equipment.

Disposal

The battery inside FLASH TRAVEL cannot be replaced, but must be removed from the handle for recycling.

If you have any questions about the proper disposal of electronic waste, please contact your local waste disposal office.

Warranty & service

Defects caused by faulty materials and/or workmanship will be repaired or replaced by LA VAGUE within the warranty period if the proof of purchase is provided during the warranty period. Please contact the customer service of your contractual partner.

Not covered by warranty:

- brush heads
- damage caused by improper use, misuse, neglect, alteration or unauthorized repair.
- Normal wear and tear, including scratches, discoloration or fading.

2 Getting started

Before using the FLASH TRAVEL toothbrush for the first time, it must be charged. It takes about 8 hours to fully charge the toothbrush for the first time or when the battery is fully discharged.

Charging process



1. Open the USB cover on the underside of the toothbrush.
2. Plug the USB end of the USB cable into the end of the toothbrush.
3. Plug the standard USB end of the USB cable into a USB adapter (any USB adapter can be used for this purpose, e. g. from all common smartphones.) with power connection.
4. Plug the adapter into a live power outlet.

Installing the brush head

Put the brush head on the neck of the toothbrush. Align the guide mark on the back of the brush head with the beginning of the guide mark on the back of the brush body.

Rotate the brush head so that the guide marks on the back of the brush head and neck are aligned.

Dentists recommend replacing the brush heads every 3 months to achieve optimum results.

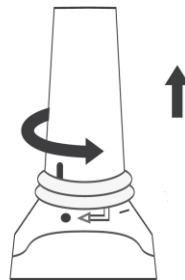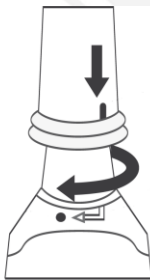3 How to use

Brushing Instructions

1. Put the brush head under water and add toothpaste.
2. Gently place the bristles of the brush on your teeth at an angle of approx. 45° to the gum line.
3. Press the power button to activate the vibration.
4. Light pressure allows FLASH TRAVEL to achieve optimal results.
5. Move the brush head slowly and gently over your teeth and gums in a reciprocating motion throughout the entire brushing time. Make sure that the bristles reach between the teeth and massage the gums.
6. Brush your teeth for 30 seconds on each quadrant of your mouth (4/30 method).
7. After brushing all four quadrants, the vibration stops automatically after 2 minutes of continuous use the cleaning time recommended by the dentist.
8. To stop at any time, press the power button. Note: FLASH TRAVEL is safe to use on braces and dentures (fillings, veneers, crowns, bridges, etc.). When used on braces, the brush heads wear out faster.

Vibration modes

The first time you use FLASH TRAVEL, it starts automatically with the default mode. For later use, the previously selected mode is used. Select the appropriate mode by pressing the Mode Selection button. Together with the corresponding brush heads, FLASH TRAVEL offers three modes for different oral care needs.

Choose from the following options:

- Whiten - Removes stains on the surface
- Clean - For exceptional everyday cleaning
- Massage - Gentle massage of the gums

4 Cleaning and storage

DO NOT CLEAN ANY PARTS OR COMPONENTS IN THE DISHWASHER.

Brush heads

1. Rinse the brush head and bristles after each use.
2. Shake the brush to remove excess water.
3. Rinse the brush head connection weekly with warm water.

Cleaning

1. Remove the brush head and rinse the neck of the toothbrush with warm water.
2. Use a clean, damp cloth to wipe the entire surface of the toothbrush.

Storage

If you do not want to use FLASH TRAVEL for an extended period of time, unplug the charger cable from the wall outlet, clean the toothbrush and store the entire device in a cool, dry place and avoid direct sunlight.

5 Trouble shooting

The toothbrush cannot be switched on

You may need to charge FLASH TRAVEL. Charge the toothbrush for 8 hours. (See Getting Started > Charging)

Tingle sensation when brushing

Due to the vibrations, the first use of FLASH TRAVEL may cause tingling sensation. This feeling will decrease over time as you get used to the vibrations.

The brush head rattles against the teeth while brushing

Keep the brush head at a slight angle to the gum line. Adjust the angle so that the brush head does not touch the surface of your teeth.

The vibrations appear to be weak

You may need to fully charge FLASH TRAVEL. You may need to replace your brush head. Note also that each model has distinct vibration pattern.

Manufacturer:

VEHNSGROUP

VEHNS GROUP GmbH
Oberanger 24
80331 Munich
Germany

www.vehnsgroup.com

service@vehnsgroup.com

www.lavague.de

Subject to alterations and errors.

The manufacturer assumes no liability in case of faulty installation, improper or incorrect use.

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© VEHNS GROUP GmbH