beper

- MASSAGGIATORE ANTICELLULITE - MANUALE DI ISTRUZIONI
- ANTI-CELLULITE MASSAGER - USE INSTRUCTIONS
- MASSEUR ANTICELLULITE - MANUEL D'INSTRUCTIONS
- ANTICELLULITE-MASSAGEGERÄT - BETRIEBSANLEITUNG
- MASSAJADOR ANTICELULÍTICO - MANUAL DE INSTRUCCIONES
- ΣΥΣΚΕΥΗ ΜΑΣΑΖ ΚΑΤΑ ΤΗΣ ΚΥΤΤΑΡΙΤΙΔΑΣ - ΕΓΧΕΙΡΙΔΙΟ ΟΔΗΓΙΩΝ
- SISTEM PROFESIONAL DE LIPOMODELARE - MANUAL DE INSTRUCȚIUNI

Cod.: 40.500

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e-mail: assistenza@beper.com

GENERAL SAFETY ADVICE

READ THE OPERATING INSTRUCTIONS CAREFULLY BEFORE FIRST USE OF THE APPLIANCE.

Read the following instructions carefully, since they give you useful safety information about installation, use and maintenance and help you avoid mishaps and possible accidents.

Remove the packing and make sure that the appliance is intact, with special attention to the supply cable. The elements of the packing (plastic bags, polystyrene, etc) do not have to be accessible to the children, in order to preserve them from danger sources.

It is recommended not to throw these packing material in the domestic waste, but to deliver them to the appropriate station of collection or destruction refusals, asking eventually information to the manager of the Sanitation Department of your own municipality.

Before connecting the appliance make sure that the mains voltage corresponds to the data specified on the label.

If the plug does not suit the socket, let it be replaced with another suitable plug by qualified personnel only, which shall also ensure that the section of the socket cables is suitable for your appliance power absorption.

Never pull the supply cord to unplug the appliance.

Make sure that the supply cord never comes into contact with hot or cutting surfaces

Do not use this appliance if the supply cord is damaged. In case of damage, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.

Only connect this appliance on alternating current

Generally, the use adaptors multiple plugs and/or extensions is not recommended; should they be necessary, please use only simple or multiple adaptors and extensions in compliance with the safety rules in force, making sure that the capacity absorption limit marked on the simple adaptors and the extensions and the maximum power limit marked on the multiple adaptor are not exceeded.

This appliance is to be destined to its special purpose only.

Any other use is to be considered as improper and thus dangerous. The manufacturer cannot be responsible for any damages deriving from improper, wrong or incautious use.

To avoid any heating risk, it is recommended to unroll completely the supply cord and to unplug the appliance, when not used.

Before cleaning, always unplug the appliance.

To protect against electric shock, do not immerse cord, plug or any parts of the appliance into water or other fluids.

Do not hold the appliance with wet hands or feet.

Do not hold the appliance with naked feet.

Do not expose the appliance to atmospheric agents such as rain, wind, snow, hail.

Do not use the appliance near explosives, highly flammable materials, gases, burning flames, heaters.

Do not allow children, disabled people or every people without experience or technical ability to use the appliance if they are not properly watched over. Close supervision is necessary when any appliance used by or near children.

Make sure that children do not play with this appliance.

If you decide not to use the appliance any longer, make it unserviceable by simply cutting the supply cable after disconnecting the plug from the socket.

ALSO MAKE ANY DANGEROUS PARTS UNSERVICEABLE, ESPECIALLY FOR CHILDREN WHO COULD PLAY WITH THE APPLIANCE.

SOME WARNINGS:

ASK YOUR DOCTOR IF YOU SUFFER FROM ANY HEALTH PROBLEM, FOR EXAMPLE, IF YOU USE A PACE MAKER.

DO NOT USE THE MASSAGER IF YOU ARE SUFFERING FROM ANY SKIN RASH OR INFLAMMATION.

IF YOU SUFFER FROM VARICOSE VEINS, DO NOT USE THE INFRARED LIGHT IN THAT AREA, AS THE INFRARED LIGHT EMITS A LOW HEAT WHICH COULD AGGRAVATE YOUR CIRCULATORY PROBLEM.

DO NOT USE THE MASSAGER IN THE ABDOMEN AREA IF YOU ARE PREGNANT OR BELIEVE THAT YOU MIGHT BE PREGNANT.

Description

Main Parts

1. Main body of the appliance
2. Off/massage/massage + IR switch
3. Variable massage intensity control
4. Infrared light emitter
5. Massage head adaptor
6. Cord
7. Finger massage head
8. Spiky roller massage head
9. Ball massage head
10. Wavy surface massage head
11. Protective fabric for wavy surfaced massage head

Selecting the Massage Head:

Head Types and Functions

•"Finger" massage head (7)

For use on painful cellulite and muscular pains.

Do not use on the back.

•"Spiky roller" massage head (8)

This head can be used to eliminate cellulite, promote firmness, stimulate blood flow and return, tone muscles, drain the lymphatic system, reduce size, warm up the muscles, provide relaxation, relieve stress, etc..

Do not use on the back.

•"Massage balls" head (9)

Stimulates and tones the muscles, helping to firm the skin. It improves circulation in the abdominal area, eliminating localised toxins and improving the intestinal tract, preventing constipation.

•"Wavy surface" massage head (10)

Tones and improves blood circulation. For use in the thigh and buttock areas. For greater comfort, cover the wavy surface massage head (10) with the protective fabric (11).

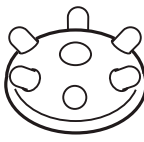
•Use of the "head protector" cloth head (11)

The cloth accessory (11) can only be used with the "Wavy surface" massage head (10).

Types of Cellulite:

Identify the type of cellulite you suffer from, to help you choose the most effective massage head.

- Oedematous Cellulite: Feels spongy, it does not hurt and moves when touched or when you change position.
- Painful Cellulite: Pain is felt when the area is pressed.
- Cystic Cellulite: Grainy and hard to the touch. Usually accompanied by stretch marks. Its temperature feels lower compared to the rest of the body.

ZONE	CONDITION	HEAD	INTENSITY	
LEGS AND BUTTOCKS	Oedematous Cellulite Flaccidity Cystic Cellulite	1st STEP 		
		2nd STEP 		
		3rd STEP 		
LEGS AND BUTTOCKS	Painful Cellulite Tiredness			
		Fluid Retention		
ZONE	CONDITION	HEAD	INTENSITY	
ABDOMEN	Oedematous Cellulite Painful Cellulite Flaccidity Constipation			
		Cystic Cellulite 		
ARMS	Oedematous Cellulite Painful Cellulite Flaccidity			
		Cystic Cellulite 		

General Recommendations for Use

- The intensity control allows you to regulate the intensity of the massage as required. Start at a gentle speed and progressively increase it, depending on your sensitivity to the massage.
- If using special treatment creams (firming and anti-cellulite creams, for example), apply these after the massage, which prepares the skin and its tissues to absorb the ingredients.

We always recommend you use a moisturising lotion after the massage. You can apply the cream using either the “wavy surface” head (10).

- Do not apply creams when using the protective fabric (11) on the “wave- shape” head.

Recommended Treatment Duration

- Apply a 2-3 minute massage on the abdomen and a 5 minute massage on the legs and buttocks. You may establish the duration of the massage depending on your sensitivity.

The total treatment must be approximately 15 minutes.

- Do not use the appliance for more than 15 minutes. Using the appliance continuously for more than 15 minutes may cause it to overheat, reducing its life. Let it cool down before continuing use.
- Use this appliance every day until you notice the results, which may be subsequently maintained by using it three days a week. To maintain the results over the long-term, adapt your lifestyle habits and food.

Use

This appliance has been designed to provide therapeutic massages to help prevent cellulite, fluid retention and flaccidity. Do not use it for any other purpose.

Massaging with this appliance generates a thermal effect in the treated area; this thermoactive effect helps to mobilise the accumulated fats and toxins. To gently intensify the treatment, the appliance has an infrared light that slightly increases the thermal action of the appliance.

Plug in the appliance (fig. 1).

Move the switch (2) to the massage position (fig.2 B) or to the “massage + infrared position (IR) position” (fig.2 C) according to the desired treatment.

Fit the correct head, (7), (8), (9), (10) for the treatment you require to the adaptor (5) (fig. 3).

To remove the head, pull firmly at the edge to release the fastening tabs. (fig. 3).

To choose which massage head to use, follow the instructions in the “Treatment Type” table.

You may alternate the use of different heads during the same massage.

Set the intensity of the massage using the intensity control (3) and place the massager on the area on which you are going to work.

Set the intensity of the massage using the control (3), according to the area to be treated and its sensitivity to massage (see the “Treatment Type” table).

To do so, turn the control to one side or another with your finger, depending on the desired intensity (fig. 4).

Hold the massager firmly with both hands and move it over the skin in the direction of the blood circulation (see “Massage Direction”).

Without pressing too hard, simply glide the massager over the area to be treated.

When you have finished using the massager, move the switch (2) to the “off” position (fig.2 A).

Unplug the appliance.

Massage Direction

Always glide the massager over the skin in the direction of the blood circulation.

• Legs

Place the appliance on the thigh, close to the knee, and glide it gently upwards towards the hip (fig. 5).

• Buttocks

Combine upwards and circular movements. You can end an upward massage on the legs with downwards and circular movements on the buttocks (fig. 6).

• Abdomen

Place the appliance on the navel and move it gently outwards (fig. 7).

Never make circular movements around the navel.

• Arms

Place the appliance close to the elbow and glide it gently towards the armpit. Repeat the massage from the armpit towards the elbow (fig. 8).

NEVER apply a massage on a bone (knee, ankle, elbow). The protruding parts of the massage head could knock the bone, causing an injury.

Only use the massage heads indicated in the “Treatment Type” table.

Cleaning and upkeep

Do not use solvents, detergents or abrasive products.

Unplug the appliance before clean the appliance.

Use a dry cloth to clean the outside.

When the appliance is not being used, unplug it .

To prevent the cable from becoming damaged, do not twist it, pull it or wind it around the appliance. If it is damaged, have it replaced immediately.

TECHNICAL DATA

Power 25W
Voltage 230V~ 50 Hz

For any improvement reasons, Beper reserves the right to modify or improve the product without any notice.

The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment. The crossed out “wheeled bin” symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance’s efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

GUARANTEE CONDITIONS

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- The appliance has not been tampered with, otherwise it is not maintainable.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) the aesthetic parts are excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, carelessness, wrong or improper installation, damage during transport and any other damage not attributable the supplier.

For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge.

In any case, if the part to be replaced for defect, breakage or malfunction is an accessory and/or a detachable part of the product, Beper reserves the right to replace only the very part in question and not the entire product

TECHNICAL ASSISTANCE

Even after the guarantee period we will always pay attention to repairs for defected appliances.

For technical assistance and/or repairs after the guarantee period you can directly contact the below address:

CONTACT YOUR DISTRIBUTOR IN YOUR COUNTRY OR ATER SALES DEPARTMENT BEPER. E-MAIL assistenza@beper.com WHICH WILL FORWARD YOUR ENQUIRIES TO YOUR DISTRIBUTOR.

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