DREAVER

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Hand Massager

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USER MANUAL

Hand Massager

Read these instructions carefully before use

1 About

These are important instructions designed to help you enjoy your hand massager properly and safely. Please read them carefully before use.

2 Specifications and features

Voltage DC3.7 V.

Power: 5W

Weight: 1.15 kg

Program time: 15 min

Charging time: 5 hours

Working time: 3 hours

- Ergonomic design in accordance with (palm) acupuncture points and hand curves
 - Pneumatic segmented massage
 - With heat function for colder days
-

3 Before use

Charge the device before first use to protect the battery life. Fully charging the device takes 5 hours, then it will work for 3 hours straight. You can use a phone charger for charging. While charging, the light will turn on, when the device is fully charged, the light will turn off.


It is important to start at intensity level 1, so you can experience how the device feels. After you are comfortable with this, you can try out the other programs.

4 Target audience



This device is ideal for the following target groups:

- People who use computers for a long time
- People with poor blood circulation in the palm
- People who want to keep their skin young and tight

5 Components

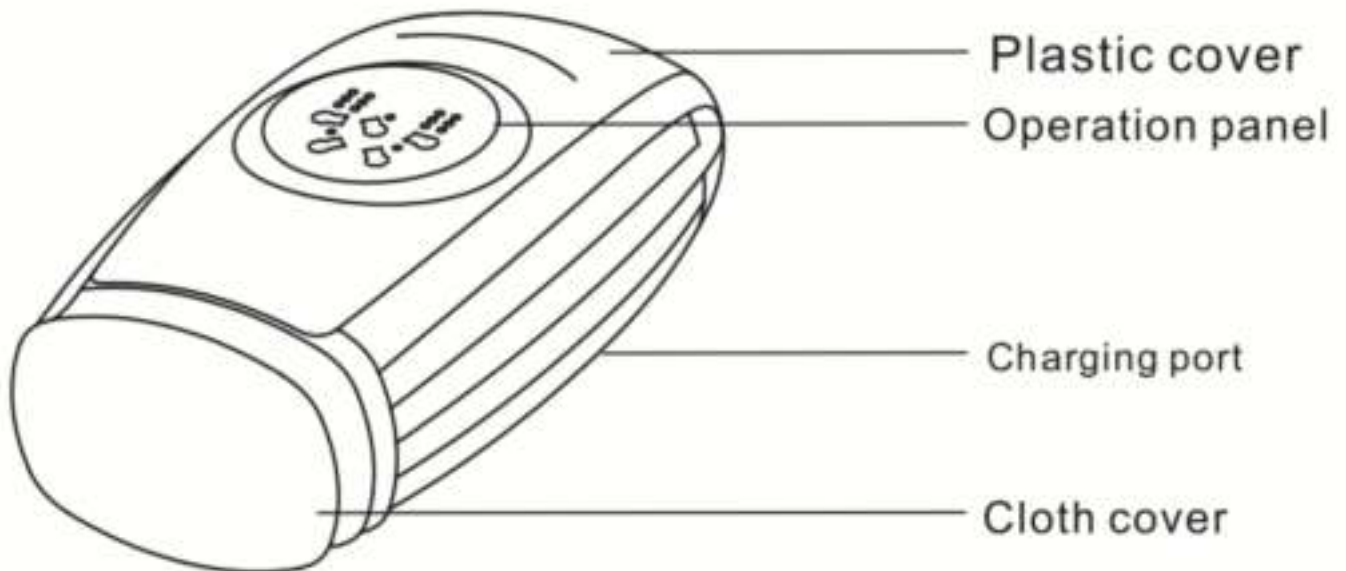6 Features

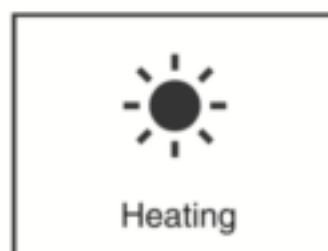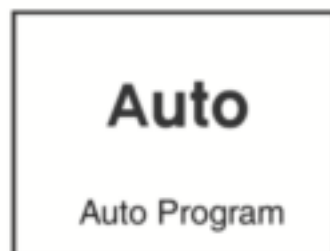7 Instructions

On off switch

Press and hold the power button for two seconds to turn on the device. Then press the intensity button. The machine will now start working automatically. The default massage time is 15 minutes, after these 15 minutes the device will turn off automatically.

Mode button

Now choose a massage mode. The device will start with level P1 by default. In total there are 6 available modes. You can select your favorite massage mode. The order is: P1 - P2 - P3 - P4 - P5 - P6 - P1 etc.

P1: A mix of all different modes (P2 to P6).

P2: The massage area focuses alternately on the hand part and the wrist part.

P3: All areas are working, but the focus is on the full palm area first.

P4: All areas are working, but the focus is on the wrist.

P5: All areas are working, but the focus is on the hands.

P6: The focus is on the hand circulation massage.

Intensitybutton

Then choose an intensity level. The default value is L3, with a total of 6 different intensity levels. You can select the different levels consecutively: L3 - L4 -L5 -L6 - L1 - L2 etc. L1 is the weakest and L6 is the strongest.

Vibrationbutton

Press this key to turn the vibration on or off (optional). There are 2 different vibration modes: red light = partial vibration. Green light = constant vibration.

Heating button

Press to turn the heating on or off (optional). The heater needs 5 minutes to reach its maximum temperature.

8 Problems and solutions

1) Device does not turn on

Possible causes:

- Empty batteries
- The device is not turned on correctly

Solution:

- Recharge the device
- Press the on / off button for 2 seconds

2) The appliance makes a strange noise

Possible noises:

- Buzzing sound
- Ticking sound

Causes:

- The air pump can make a buzzing noise
- The valves can cause a clicking noise

Solution:

These noises are normal. There are no problems with the device and you can continue to use it.

3) The device leaves prints on the hand and arm
(See example in photo below)

This is normal after using the hand massager. The prints will disappear by itself after a few minutes.

4) All lights will illuminate and you will hear a beeping sound

Cause:

The hand massager is (almost) empty.

Solution:

Charge the device.

9 Warnings

These are important instructions designed to help you enjoy your hand massager properly and safely. Please read them carefully before use. Please keep these instructions in a safe place.

This device is designed for private use, not intended for medical or commercial purposes. Children must be supervised to ensure that they do not use this aid as a toy.

Important warnings

- DO NOT immerse the device in water for cleaning.
- DO NOT use the device if you have a wound or after surgery and have not fully recovered from your hands.
- DO NOT use this device by hand for more than 15 minutes continuously. It is recommended to take a 10 to 15 minute break on each hand every 15 minutes of use.
- DO NOT use this product while you are sleeping and avoid falling asleep while the product is in operation.
- Never use the device near a heater, fireplace, stove or other heat source.
- Always supervise and observe while children, patients, the disabled and the elderly use this device. This device is not intended for children under 18 years of age.
- Close supervision is necessary when this device is used by, near, or next to children, invalids, or persons with disabilities.
- Never allow children, the mentally challenged, and those who have difficulties operating this device to access the power source.

- Extended use can cause the product to overheat and shorten its life.
- DO NOT use this product as a substitute for medical attention or treatment.
- This product is designed for home use only.
- Do not wear jewelry while using the hand massager.

PAY ATTENTION

This massage device may not be used in the following cases:

- The massage should be comfortable and pleasant. Should any pain or discomfort occur, discontinue use and consult your physician.
- This product should never be used by a person with a medical condition that would limit the user's ability to operate the controls.
- Never use this product directly on swollen or inflamed areas of the skin or areas where skin eruptions are present.
- DO NOT work after taking painkillers, sedatives, or alcohol. Consult your physician if under these medications.