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
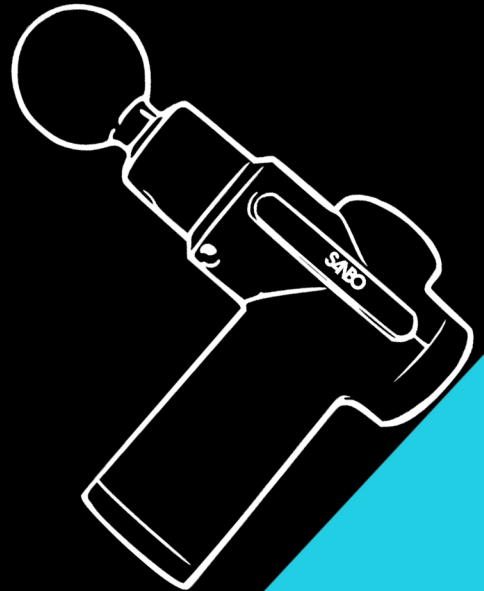
# Congratulations with your new SANBO Massage Gun

You are one step closer to improving your freedom of movement, preventing and curing aches and pains and becoming an improved version of yourself.

Whether you are a top athlete who wants to improve his or her performance, a weekend warrior who wants some extra recovery after his Sunday football game or someone with an office job that wants to improve his or her posture.

The SANBO Massage Gun Mini is designed to help you achieve your personal goals.
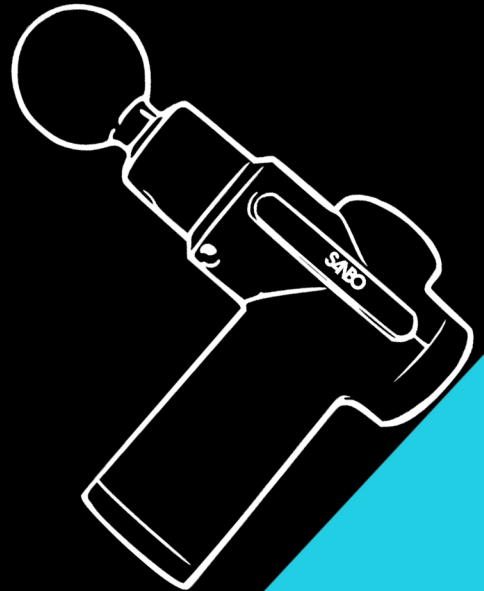
# SANBO

## QUICK START

# Pain-free movement made easy

1


### Turn on the device

Turn on the SANBO Massage Gun pressing and holding the on/off button on the screen.

2

### Find the right speed setting for you


Switch between 6 different speed settings by gently pressing on the on/off button on the screen of the SANBO Massage Gun.

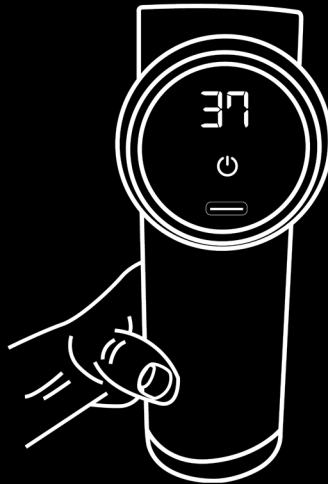**3**

## Choose the right attachment

Choose the right attachment for the type of massage. You “push” the attachment into the massage gun. When replacing the attachment, it is allowed to pull hard.

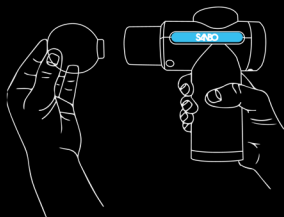**4**

## Adjust the speed


Select the desired speed by pressing the + and - buttons.

**5**

## Charge

Only use the supplied SANBO charger to maintain the warranty.

When the screen indicates 100, the SANBO Massage Gun is fully charged.

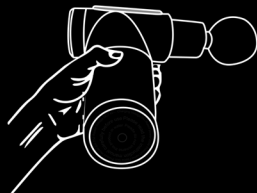**6**

## After care

Use a damp, clean cloth and wipe the Massage Gun and attachments. Make sure the device is turned off and the charger is not connected.

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## Round Ball

Perfect to use for large and small muscle groups. This is an allround attachments.

## Air Plug

Suitable for muscle all muscle groups. Perfect for getting acquainted with massage.


## Y - Fork

Suitable for the back, biceps, calves, achilles tendon and ankles.

## Bullet

The attachment to imitate "deep tissue" massage. This attachment will massage knots out of your muscles in the most precise way.
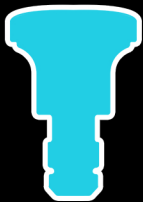
## Dampener

Suitable for applying tapotement massage. Use this attachment especially on the larger muscle groups.

# Package contents

-SANBO Massage Gun Mini- EAN 8720589320208

Attached accessoires

- 5x Attachment
- 1x EU-adapter
- 1x Storage Case
- 1x Manual

# Product Specifications

<b>Size</b>	13,8 x 14,2 x 5,5 CM
<b>Weight</b>	800 gr
<b>Color</b>	Black
<b>Power supply -Power - Voltage</b>	DC 7.4V - Charge 5V - 25W
<b>Type of batterij</b>	Lithium-Ion 1800 mAh
<b>Battery life</b>	2-4 Hours
<b>Speed</b>	0 - 3200 RPM
<b>Noise</b>	35 dB

## Storage

-Store the SANBO Massage Gun in a dry, clean place. do not store the device in extremely high or low temperatures.

## Throw away / recycle

-Please contribute to a cleaner living environment by not throwing your SANBO Massage Gun away with normal household waste at the end of its lifecycle, but by recycling it. This can be done at a collection point for batteries and accumulators designated by the government. This is completely free for now!

## Service and Warranty

- If you need service or information about your SANBO product, do not hesitate and send an e-mail to [customerservice@sanbo.nl](mailto:customerservice@sanbo.nl) or call 0800 3336669.

-There is a standard 1-year manufacturer's warranty on your SANBO product. The customer must be able to show proof of purchase and the warranty must be made within a period of 1 year from the date of purchase. Note the warranty does not apply; in case of damage caused by improper use or repair; wearing parts; for defects of which the customer was aware at the time of purchase; in the event of defects caused by the customer's own fault, in the event of damage caused by third parties.

## Disclaimer

- Specifications are subject to change without notice.

SANBO App

# The massage guide in your pocket

Get guidance with your  
massage.

Read about information about  
percussion therapy and the  
different attachments.

GET IT ON  
Google Play

Download on the  
App Store

SANBO Care

# A Team that is here for you!

Please contact our team if anything is unclear or if you  
have any questions about the product.