
Lifeproducts

Live a Better Life

# Elektrisch Warmtekussen XL

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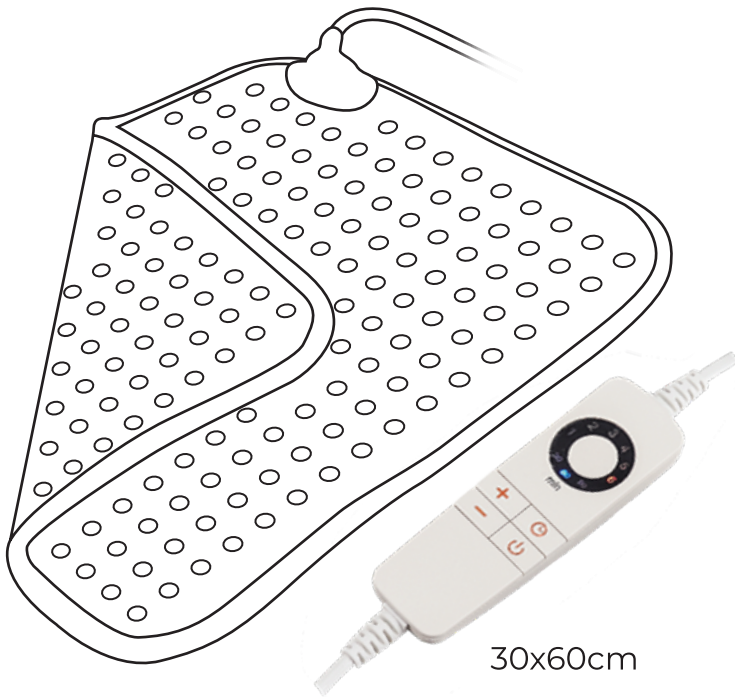
Live a Better Life

# Electric Heating Pad XL

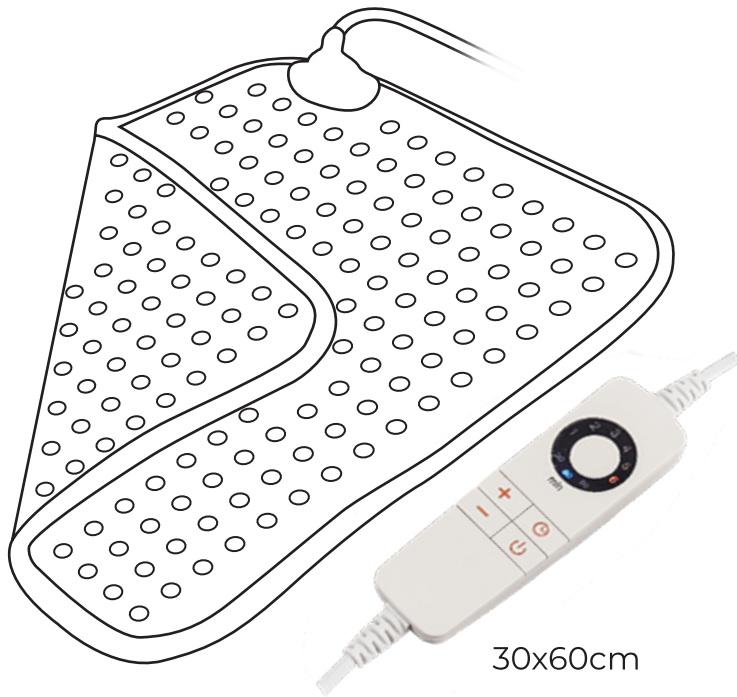Thank you for purchasing your new heating pad. These operating instructions will help you use it properly and safely. We recommend that you spend some time reading this instruction manual so that you can fully understand all the operational features it offers. Read all the safety instructions carefully before use and keep this instruction manual for future reference.

EK-02-NL

## **IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!**

Read the instruction manual carefully before using this product, especially the safety instructions, and keep the instruction manual for future reference. Should you give this device to another person, it's vital that you also give them these instructions.

Don't puncture the product!

Don't use the product when it's folded!

The product shouldn't be used by very young children (0-3 years old)!

Only use the product indoors!

The product can be washed by hand!

Don't bleach the product!

Don't dry the product in a drier!

Don't iron the product!

Don't dry clean the product!


# 1. SAFETY WARNING

- For your own safety and to reduce the risk of injury or electric shock, please follow all the safety precautions listed below.
- Read all the instructions carefully before using this product and keep them for future reference.
- Use the product only for its intended purpose. Improper use may lead to hazards.
- Don't use attachments or accessories which are not recommended by the manufacturer. They may cause malfunction or injury.
- Don't leave the product unattended during use. It must be completely cooled down before storage.
- Don't place heavy objects on top of the product.
- Don't immerse the product or the mains cable or plug in liquids. Never hold the product under water!
- Don't use the product if it's damp or wet.
- Use this product in dry indoor areas only and never outdoors.
- This product is for household use only. Don't use it for nurseries, kindergartens or hospitals. Don't use the product for any purpose not expressly specified in this manual.
- Before using this product ensure that the voltage of your electricity supply is the same as that indicated on the rating plate of the product.
- The product be only completely disconnected from the mains by unplugging from the power socket.
- Unplug the product from the mains supply 1) if you notice any sign of damages, 2) when not in use, and 3) before cleaning any part of the product.
- Only unplug from the power socket by the plug itself. Don't pull on the cable!
- Ensure that the mains cable doesn't touch any hot surface and that it doesn't hang over the edge of your bed.
- Keep the product away from an open fire and spray products. This can cause fires under certain circumstances.
- Examine the appliance frequently to check that there are no signs of wear or damage. If any damage/wear is spotted, don't use the appliance and seek specialist advice.
- Don't operate the product with a damaged mains cable or mains plug. Seek qualified technicians for assistance. If the mains cable is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- Caution! In order to prevent hazards from unintentional resetting of the protective temperature limiter, this product must not be powered by an external switching device such as a timer, or connected to a circuit which is regularly switched on and off by a device.
- People with pacemakers should consult their physicians or the manufacturer of the pacemaker before using this product. The electronic and magnetic emission from this product can affect the pacemaker's function under certain circumstances.

- This product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning the use of the product in a safe way and understand the hazards involved. Children shouldn't play with the appliance!
- The product shouldn't be used by young children over the age of three unless the controls have been pre-set by a parent or guardian, or unless the child has been adequately instructed on how to operate the controls safely.
- The product mustn't be used by people who are sensitive to heat and other very vulnerable persons who are unable to react to overheating.
- Children under the age of three shouldn't use this appliance due to their inability to react to overheating.
- Cleaning and user maintenance of this product shall not be made by children without supervision.
- Don't place items on top of the product during storage. This way you prevent creasing.
- The product should be stored in a cool and dry place.

## 2. PRODUCT OVERVIEW

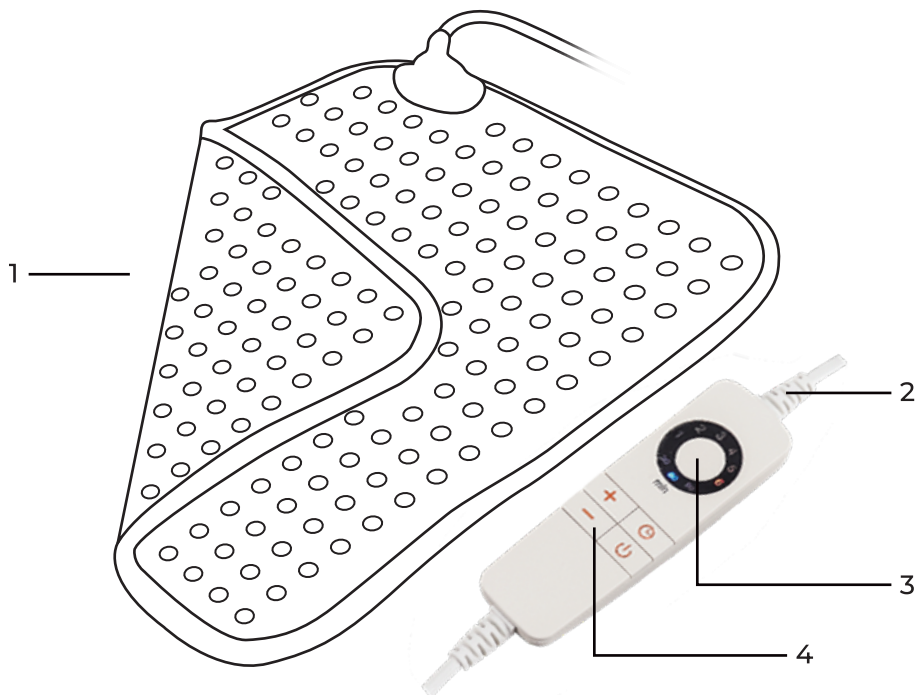1. Heating pad
2. Mains cable
3. Operating indicator
4. Switch 0-9
5. Automatic shutdown 30, 60, 90 min.

### **3. FIRST TIME USE**

Inspect the product for completeness and damages before using it for the first time. Should the delivery contents be incomplete or damaged, please contact us at

klantenservice@lifeproductsstore.nl

### **4. OPERATION**

- 1) Connect the power plug to a suitable power socket.
- 2) Switch on the product.

### **5. TEMPERATURE SETTING**

Easily set the temperature from low heat to maximum heat through 9 different heat settings.

NOTE: Setting 1 is recommended for continuous use. The heating pad isn't designed for use in hospitals or for commercial purposes. It mustn't be used with babies, toddlers, heat sensitive people, disabled people or animals.

### **6. AUTOMATIC SWITCH-OFF FUNCTION**

The heating pad switches off automatically after approx. 90 minutes of continuous operation. Switch off the heating pad (setting 0) or unplug the power plug from the power socket if you're not going to use the heating pad after the automatic switch-off.

### **7. CLEANING AND MAINTENANCE**

- Remove the cable and controller from the heating pad. Clean the heating pad according to the laundry care symbols on the label.
- Don't machine-wash, spin or machine dry, chemically clean, wrung out, press or iron the heating pad.
- Remove smaller stains with a cloth or moist sponge together (with a little liquid detergent for delicate fabrics, if necessary).
- Only use the heating pad again once it's completely dry.
- It's recommended to wash the heating pad only if it's dirty or dusty. If you wash the heating pad too often, it could wear quickly.

## 8. SPECIFICATION

- Lifeproducts Electric Heating Pad XL
- Power: 100 W
- Voltage: 220-240 V~ 50 Hz
- Size: 30x60cm
- Material: super cozy fleece pad
- Automatic switch-off after 90 minutes
- 3 heat settings with LED indicator
- Fast heating and overheating protection
- CE/GS approval
- RoHS/REACH test

This symbol on the product or in the instructions means that your electrical and electronic equipment should be disposed, separately from your household waste, at the end of its lifespan. There are separate collection systems for recycling in the EU.

For more information, please contact us at

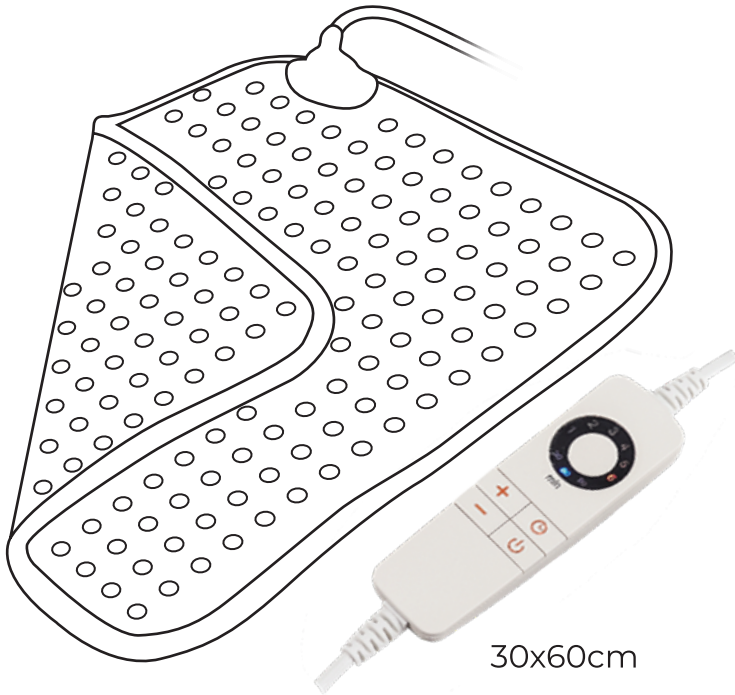
[klantenservice@lifeproductsstore.nl](mailto:klantenservice@lifeproductsstore.nl)

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Lifeproducts B.V.