

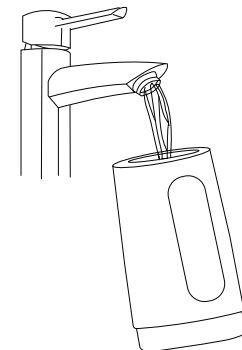
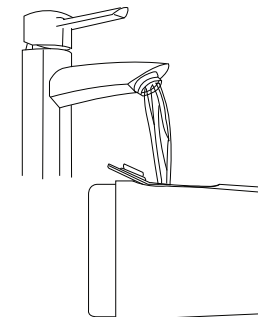
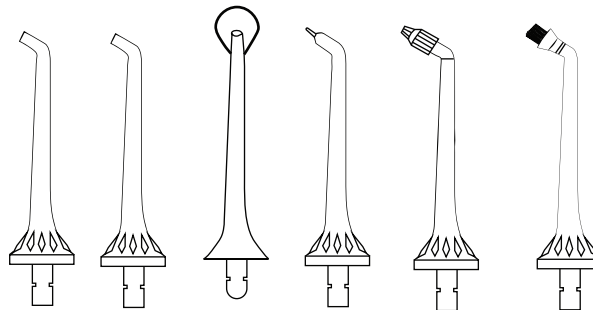


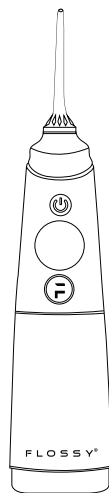

# FLOSSY<sup>®</sup> PORTABLE ORAL IRRIGATOR User Guide

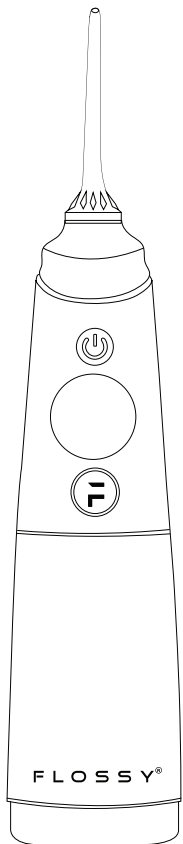


### Step 1

Choose the right nozzle, press the release button, and insert the nozzle into the spray head hold until you hear a "click" ring.


### step 2

Adding water. Open the water tank lid and fill it with the water in it.

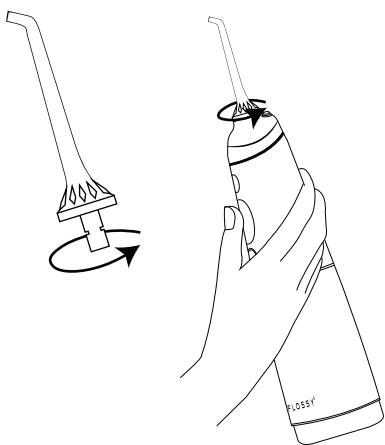


### step 3

Holding the device downward, aim the jet tip to the teeth by rotating the tip base before switching on the device.

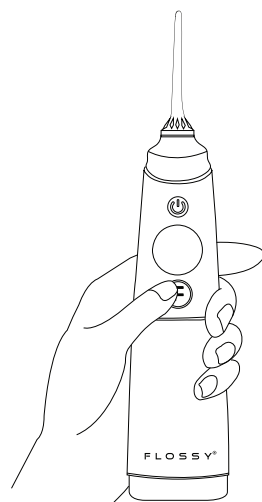


### step 4

Turn on the water flosser, press the mode button continuously to switch among the sequence: SOFT-MEDIUM-STRONG-PULSE-DIY.

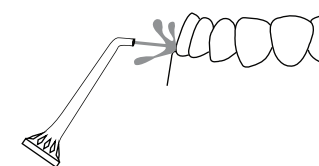




### step 5

When the water flow and the gums are at a 90-degree angle, slightly close your mouth in order to avoid the water spraying out

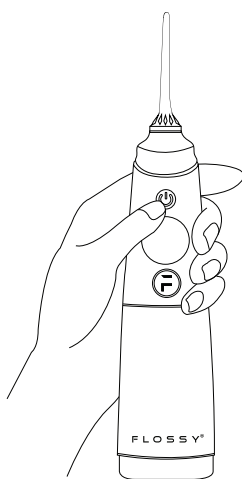

### step 6

After use, please switch off the power, as well as empty the water tank.

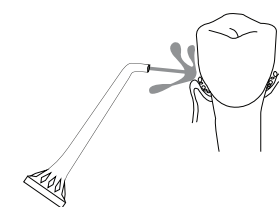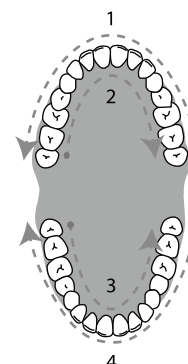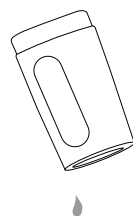

### Tips

In order to achieve the best results. Please start from the molar area to the incisors. Slide the water flow along the gums, as well as stay a while between teeth, until the area around the gums and the spaces between the teeth are fully cleaned.