

BEFORE USE

1. Before using, make sure that the filter sponge is attached correctly and is dry. If water or other liquids is sucked into the air pump of the body this will cause blockages and damage the device.
2. When using this product, all accessories (such as necklaces, earrings, etc.) and contact lenses should be removed to avoid damage to the machine and your body.
3. This product should not be used on eyeballs, eyelids and necks. Otherwise, it may cause damage to skin and body.
4. If you feel pain or abnormality in your skin or body, stop using it immediately.
5. When cleaning and maintenance, please make sure the device is turned off.
6. Do not store this product in the bathroom or other damp rooms.
7. Do not allow children to use this product.
8. Do not use this product near flammable items (alcohol, gasoline, nail polish remover, spray, etc.).
9. This product is not waterproof, so it is strictly forbidden to suck water into the machine or submerge in water or wash it directly under a tap.
10. Do not use the product if it is damaged.
11. Do not disassemble or repair this product by yourself. For maintenance, please contact the shop from where the product was purchased.
12. Do not charge during usage.
13. It is forbidden to use it when bathing.
14. Please use the matching USB cable to charge. To ensure that the charging voltage meets the product requirements.



This manual contains safety advice and instructions for the correct use of this product.

Therefore, please read the instructions carefully to avoid damage and harm caused by improper use.

This product is a personal care product. It is not designed for commercial use.


This product is not suitable for the following people:

1. People with allergic symptoms like skin diseases, dermatitis, skin allergies and severe sunburn.
2. People with skin injury, eczema and swelling.
3. People who are recovering from surgery.
4. Persons with physical, sensory or mental problems, unless directed by a guardian.

Product characteristics

- **Small and easy to carry.**
Easy to Charge, compact, lightweight, easy to carry.
- **Large capacity rechargeable battery**
Built-in 500mAh large capacity rechargeable lithium-ion battery. The USB portable charging connector can be charged through computers, mobile power supply and other devices with USB interface.
- **Replaceable design of attachment head**
Different attachment heads can be used according to the required cosmetic effect. – last entry
- **Easy to clean**
The different attachment heads can be easily removed and cleaned. The filter sponge effectively prevents dirt and aging keratin being sucked into the machine interior.
- **Device usage**
Safe, easy to operate, with no edge effects and side effects.

Component description

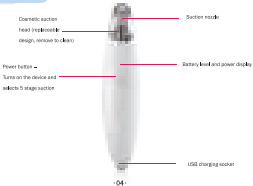

Specification parameter

Product Name	Electric Blackhead Remover
Product Model	X5M
Charging Voltage	5V = 1A
Battery Voltage	3.7V 400mAh
Battery capacity	400 mAh
Standard	GB4706.1 GB4705.15

Introduction of Cosmetic Suction Head



Large Angled Head (High Suction)

Large angled cosmetic head, features: strong suction. Mainly for blackheads, greases, dust, cosmetic residues, etc. Five minutes once a week, depending on the skin condition. It is not recommended to use around the eyes. It can be used in rough and wrinkled skin on the face, as well as on the body and elbow.



Medium Flat Head (Moderate suction)

Medium flat cosmetic head, features: moderate suction. Mainly for blackheads, greases. Once a week, five minutes at a time, depending on the skin condition.



Silicone Heads (Low Suction)

Silicone cosmetic heads have weak suction, and can be used to absorb grease, dust, cosmetic residues and so on. For normal skin.



Silicone Heads (Low Suction)

Silicone cosmetic heads have weak suction, and can be used to absorb grease, dust, cosmetic residues and so on. For normal skin.




Medium Elliptical Head (Low Suction)

Medium Elliptical cosmetic head, can be used for daily skin care, such as wrinkled corners of the eyes and mouth, etc. Pull off skin quickly after suction, then suck and pull it off again, repeat this multiple times. You can also lift your face from the bottom up. Used on cheeks, corners of eyes, side of the nose, etc.

HOW TO USE:

1. Open up the pores by steaming or bathing.
2. Select the desired attachment head based on the area you wish to target. Turn on the device and select the desired power setting. If unsure which power setting to use then please start with the lowest and gradually increase if needed.
3. Place the device on the target skin and move it around slowly and this will extract the blackheads.
4. Please don't stay in one place for more than 2 seconds when using in order to prevent from damaging your skin.
5. After cleaning the pore you can use some skin toning to close the pore. Use something like a facial mask after or cold water.
6. If the skin is red do not worry, it is normal. It will fade away within one or two days

Safety precautions

When you discard this product it should not be disposed of with the garbage used by ordinary households, because it will endanger the environment and human health. Instead, it should be placed in a special environmental recycling bin. 

TIPS

-  Before first use, please charge for 2 hours to prolong the battery life.
-  Long press power button to turn on.
Press the power button once to select 1st suction setting.
Press again for 2nd suction setting.
Press again for 3rd suction setting.
Press again for 4th suction setting.
Press again for 5th suction setting.
Long press power button to turn off the device.

Cleaning and Maintenance

1. Protect the instrument from dampness and liquid.
2. Do not put it in dishwasher or microwave oven.
3. Please turn off the machine before cleaning.
4. Do not use corrosive cleaning solutions.
5. Clean the device with a soft towel
6. Do not rinse the instrument directly under water, because it will cause damage to the machine.
7. Do not leave charging unattended.

Safety advice

1. Do NOT use on scars.
2. Please don't press the cosmetic head onto the skin, the machine itself has a strong suction force.
3. After using the device your skin will appear a little reddish, which is a normal and will disappear.
4. Do not use abrasive plaster or other exfoliating products before using this device.
5. After using this device please moisturize the same areas.

Warranty

Thank you for your interest in the products and services. This Limited Warranty applies to physical goods, and only for physical goods, purchased from (the "Physical Goods").

What does this limited warranty cover?

This Limited Warranty covers any defects in material or workmanship under normal use during the Warranty Period. During the Warranty Period, will repair or replace, at no charge, products or parts of a product that proves defective.

What will we do to correct problems?

will either repair the Product at no charge, using new or refurbished replacement parts or supply a replacement device. The replacement device could be of equal condition.

