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Gezondheid & fitness

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Bluetooth

Lichaamsvetweegschaal

iOS App

Android App

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Apple Health

Google Fit

Mobiel & Apparaat

iOS 7.0

Android 4.3

Bluetooth 4.0

12 parameters in 1 meting

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FitGo - Ideas for health
Health & fitness

User Guide

Bluetooth Body Fat Scale

iOS App

Android App

One-click sync with

Apple Health

Mobile&Device

iOS 7.0

Android 4.3

Bluetooth 4.0

12 parameters in 1 measurement

Multi Users

8 users share one user account in the device.

Target weight

Users can set their own target weight to customize their own health plan.

For First Use

1. Battery Installation

Open the battery door, put into batteries with correct polarity.

2. App Installation


Scan the QR code on the above, or download "FitGo" from App Store or Android Google Play.

3. Registration

Turn on bluetooth on mobile device, run the App"FitGo", follow the App instructions.

APP Instruction


1. Target Weight Edit/Family Member Add(total 8 members)

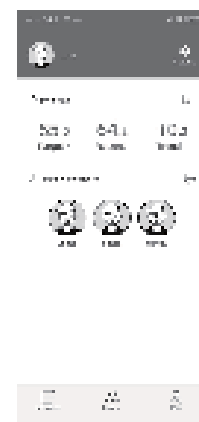2. Click"Me"->"Setting",then set unit/health index/language as you want.


3. Then choose the user and stand on the scale.

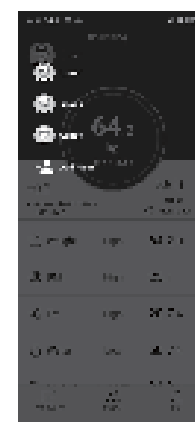4. 12 kinds of measuring data will show as below:

5.Click"Charts",then know history measuring data.

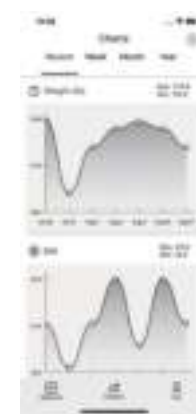Attention

1. Do not use the scale on people who have body implants such as a pacemaker, artificial limbs, contraceptive devices, or metal plates or screws. It may cause the device to malfunction or produce an inaccurate result. When in doubt, consult your physician.
2. Please make sure good contact between bare feet and the metal electrodes, otherwise it will display your weight only.
3. Ensure the scale is on a flat, even and firm surface. Avoid carpet or soft surfaces.
4. The surface of the scale will be slippery if wet. Please make sure both the glass and your feet are dry before use.
5. This scale is not water proof. Do not expose to moisture.
6. Clean the scale with a slightly damp cloth. DO not use any chemical/abrasive cleaning agents to clean the scale.
7. Please remove the batteries, if the scale is not to be used for a long time.
8. The appliance is intended solely for domestic use within the home.
9. Children should be supervised to ensure that they do not play with the appliance.

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