
INFINITY  
G O O D S

# Water Flosser M139

User Manual

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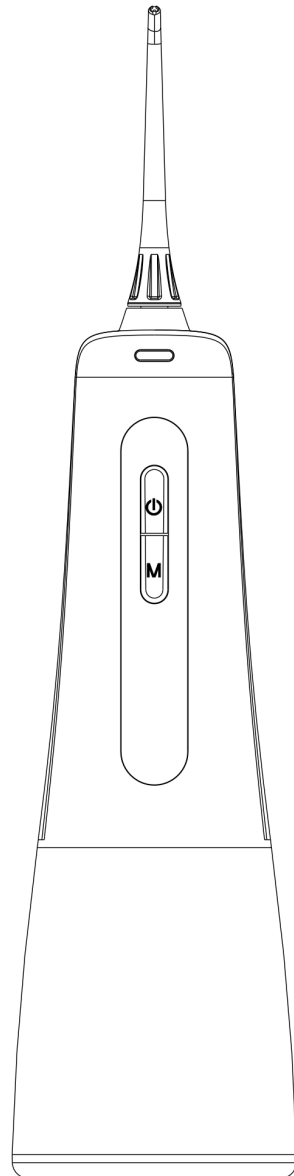Please keep this instruction in a safe place  
and read it carefully before using it.

## Safety Notices

These "Safety Notices" are designed to prevent injury and damages to customers and third parties. Indicates precautions to be followed for safe use of this product.

### WARNING


- Do not disassemble, modify or repair this product except when disposing of it.
- Do not throw this product into fire or heat it. While charging, do not plug or unplug the charging adapter or charging cable with wet hands.
- In the event of an abnormality failure, please stop using it immediately.
- Use and store out of reach infants and pets.

## Basic precautions

1. Read the relevant materials in the operation guide before use.
2. Do not use accessories not officially designated and recommended.
3. Before connecting the power supply, please check whether the marked charging voltage is consistent with the supply voltage.
4. If the product is damaged in any form (nozzle/host machine), please stop using it.
5. Please keep the oral irrigator away from high temperature, high heat, fire source, heat source, strong light source and direct sunlight.
6. Please ensure that the safe socket is used for charging, and do not force the plug into the socket.
7. In order to ensure safety, people (including children) who are physically or mentally impaired or lack relevant experience and knowledge need to be supervised or guided for safe use.
8. This product contains detachable accessories such as nozzle and water tank. Please take care of children to prevent accidental ingestion of nozzle accessories.
9. The oral irrigator is a personal item and cannot be used by dental clinics or institutions for multiple patients.
10. Please cut off the power when the mouth rinser is not charged or before cleaning, and do not clean this product in the dishwasher.
11. Do not let hair care products (such as gel, oil) or hand sanitizer contact this product to avoid damaging the surface.

# Enjoy brighter, cleaner teeth and healthier gums in just 1 minute a day!

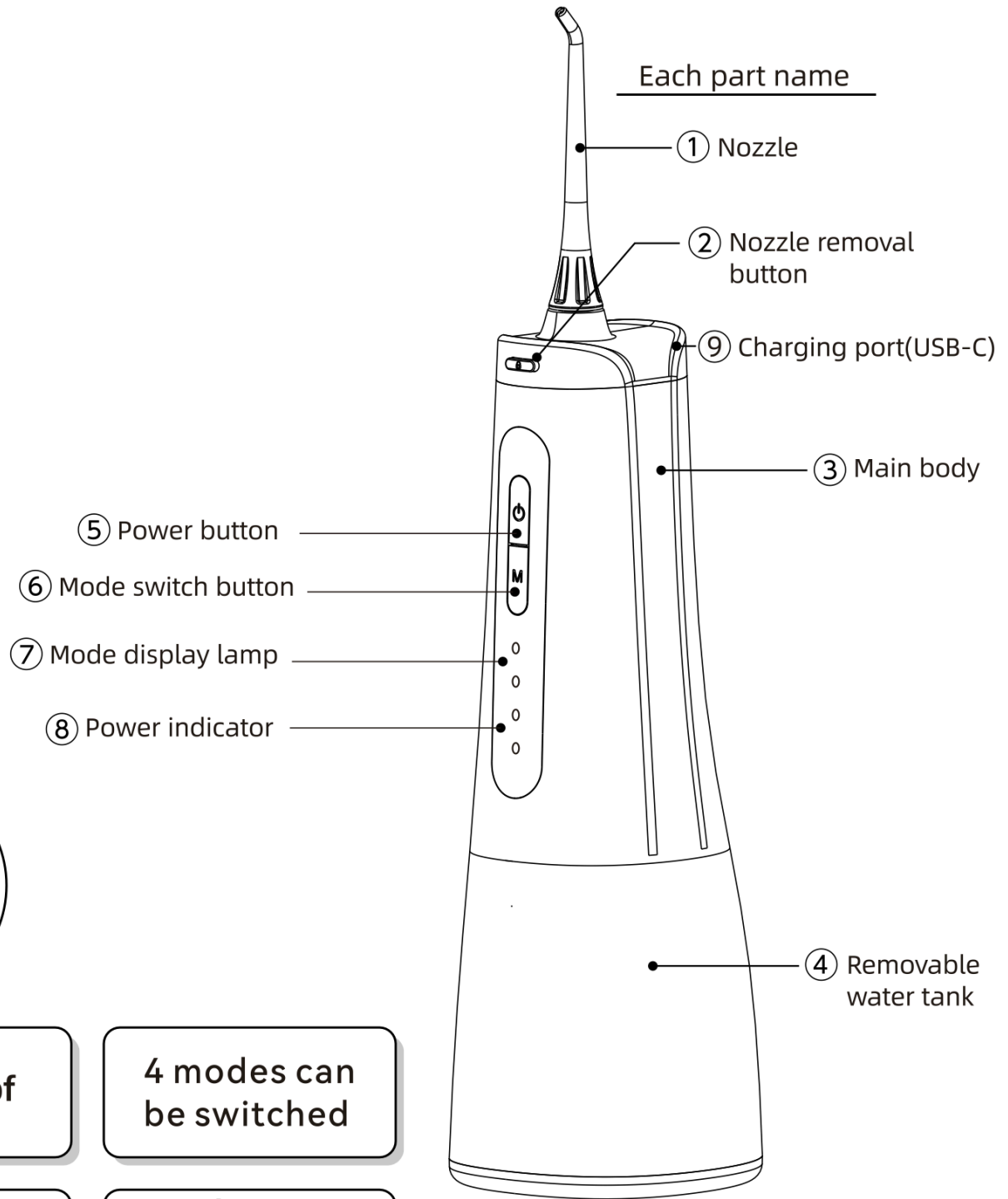Can be used  
in the bath

IPX7  
waterproof  
design

4 modes can  
be switched

300ML large  
capacity water  
tank

One charge can  
be used about  
60 times

This page has been intentionally left blank.


## About charging

If you stop or press the power button during operation and the four mode indicator flash at the same time, the remaining battery power is low and needs to be charged.


## How to charge

- 1 Connect the USB-C port of the supplied charging cable. Please remove the cap of the charging port of the main unit and connect it to the charging port. Connect the other end (USB-A) connector of the charging cable to your USB output charging adapter.

\* If the charging port of the main unit gets wet with water, please wipe it with a dry cloth before connecting the charging cable and start charging after drying.

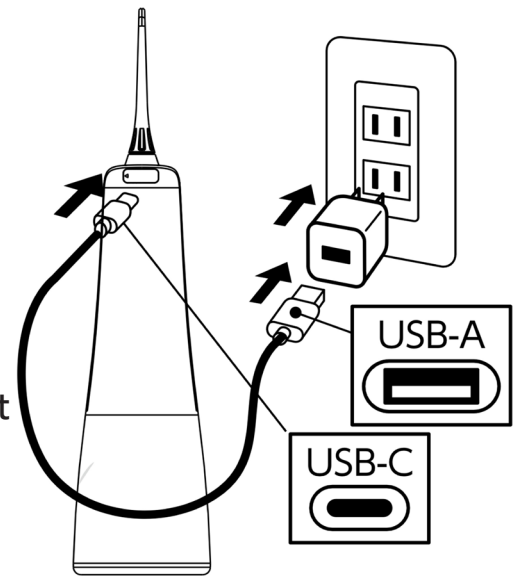
- 2 When you start charging, the four mode indicator lamps light up one after another to indicate that charging is in progress.

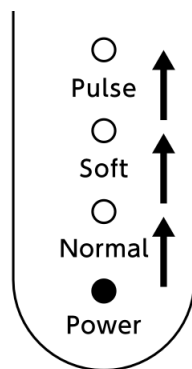Charging

- 3 When charging is completed, the four mode display lamps will light up at the same time.

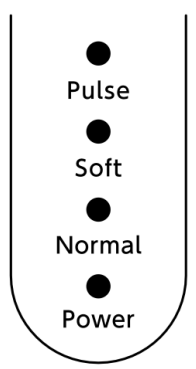
Fully charged

## How to use

- 1 **Nozzle mounting**

Insert the nozzle into the nozzle holder on the top of the main unit until it clicks.

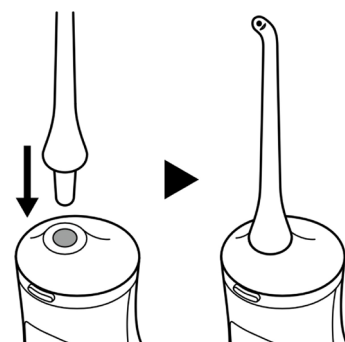
## 2 Injection of cleaning solution

Open the water tank as shown in the method on the right, fill the tank with water, and close the tank after filling with water to start using it.


On the right is a list of liquids that can be put into the can and used for cleaning. DO NOT use any other liquids as this may cause malfunction.

- **Tap water**
- **Warm water (40° or lower)**
- **Mouthwash**

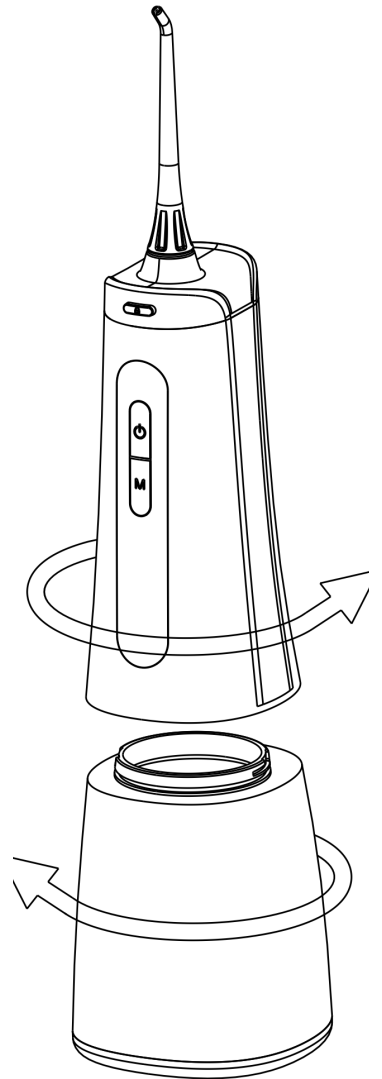
## 3 Mode settings

Press the mode switch button on the front of the main unit and select the desired water pressure. Every time it is pressed, it switches in the order of Pulse, Soft, Normal and Power. Check the set mode with the corresponding mode indicator lamp.

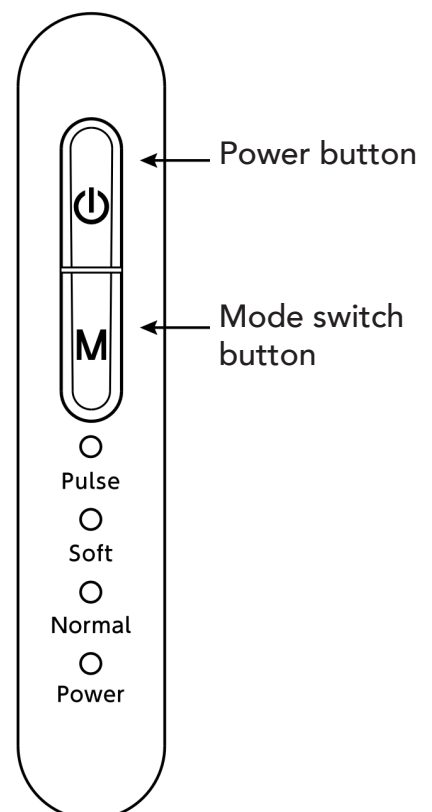
#### 4 Power on

Insert the tip of the nozzle into the mouth and press the power button. The indicator lights up to indicate that cleaning has started.

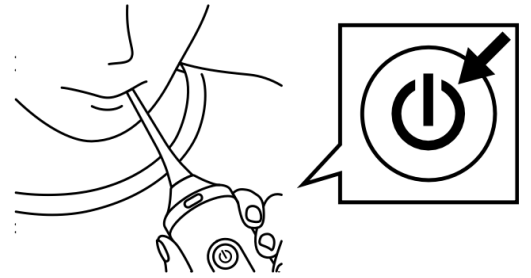*If water does not spray out for an extended period, please disconnect the nozzle and intermittently press the water valve inside the water outlet using a tool with a diameter of less than 6 mm until water begins to spray.*

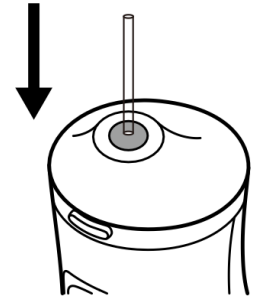
#### 5 Cleaning


During the cleaning process, keep the body forward and slightly open the mouth during the cleaning process, so that the food residues in the mouth can be flushed out.

\*Please note that if you open your mouth vigorously, the washing water will splash everywhere.

- The main body is equipped with a waterproof mechanism, so you can use it with confidence even if it gets wet with washing water.
- The direction of the nozzle can be freely rotated 360°. Please adjust it to an easy to use direction when using. Tilt for ease of use.

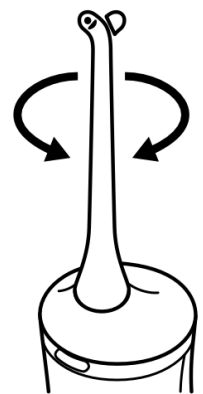
*The nozzle is aimed between the teeth and gums, and then slowly spray along the gums from the molar (the back teeth) to the front teeth, as well as the surface and back of the teeth.*

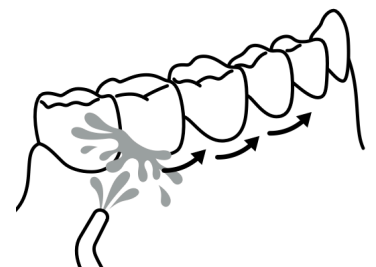
## 6 Stop

When you want to pause use, tap the power button to shut down; The machine will automatically shut down after each gear is set to run for 2 minutes.
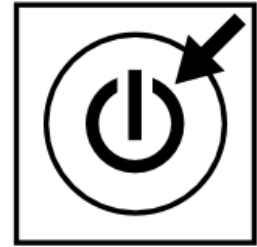
## 7 How to remove the nozzle

While pressing the nozzle release button, pull the nozzle up.
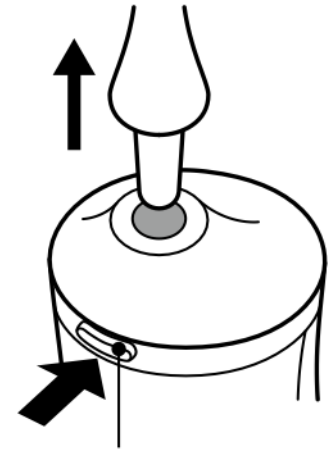
## About maintenance

- After use, open the water tank cover and remove the accumulated water in the water tank;
- If there is dirt inside the water tank, you can pour a small amount or less of 40°C warm water into the water tank, cleaned it slightly, and then the above operations are carried out.
- To keep it clean, rinse the nozzle with clean water and dry it after each use.
- When cleaning this product, please gently wash away the dirt with clean tap water (flowing water).
- After cleaning, wipe off the water with a dry cloth and dry it. If the nozzle is damaged, replace it with the attached spare nozzle.

## Troubleshooting

Problem	Possible Causes	Solution
Can not be used	Power button not pressed	Press the power button. <b>Step 4</b>
	Power LEDs blink	Please charge the battery. <b>About charging</b>
Water pressure is very low	The nozzle is worn or deformed. Nozzle is clogged with dust.	Replace with a new nozzle. <b>Step 1</b>
Water does not come out	There is no water in the tank	Fill the tank with water. <b>Step 2</b>
	The tip of the tube inside the tank is above the water surface.	Adjust the direction of the main unit so that the tip in the tube is below the water surface. <b>Step 2</b>
	Clogged water tank air inlet hole.	Please clean the air inlet hole.
Short usage time after charging	The charging cable is not properly connected.	Make sure that the charging cable is securely connected to the main unit and the USB output charging adapter, and that the charging adapter is connected to power outlet. <b>About charging</b>
	Battery life	Even if charging is completed normally, the usage time is short, it may be the end of the battery life.

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