# KIMODIRECT

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## Instructions


1. Please fully charge it before using this Product for the first time.
2. When the battery power is low, the charging indicator flashes 5 times, and then it will automatically shut down.
3. During the charging operation, insert the charging cable into the corresponding jack on the host panel, and then plug the USB port at the other end of the charging cable into a mobile power supply or a computer USB port (5VDC power supply) to charge, and at the same time the charging indicator is flashing, fully charged. The light is always on.
4. Before each use of this Product, be sure to clean the parts that touch the body.
5. Do not store this product near high temperatures or open flames, or in direct sunlight, otherwise it may cause discoloration and damage.
6. Avoid using corrosive liquids such as benzene, thinner and other chemical solvents on this product and accessories.

## Placement

1. Keep out of reach of children.
2. Do not store in a place exposed to direct sunlight, high temperature and humidity.
3. Store in a dry and ventilated place.
4. Please do not disassemble, repair or modify this product without permission.

## Steps of use

1. Long press the "ϕ" key to turn on the machine. If heating or vibration is required, short press the corresponding function button. The timer time is 30 minutes, and it will automatically stop working after 30 minutes.
2. Short press "ϕ" to select three vibration modes, and the "1, 2, 3" indicator lights will light up accordingly. Cycle in turn: 1-2-3-Vibration stop-1-...
3. If you want to adjust the temperature, short press the "⌘" key to select three heating modes, and the "1, 2, 3" indicator lights will light up accordingly. Cycle in turn: 1-2-3-heating stop-1-2-3-. (1 represents low temperature, 2 represents medium temperature, 3 represents high temperature.  
Note: It is recommended to choose gears 1 and 2 for use in summer; gears 2 and 3 for use in winter.
4. Long press "ϕ" for about 2 seconds, the host enters the shut down mode.
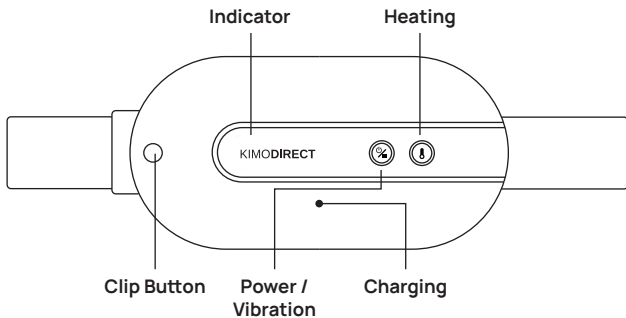
## Cleaning and Maintenance

1. When the front of the main unit is dirty, it can be wiped clean with a clean dry soft cloth or dry towel.
2. When the back of the main unit (the surface covered with hot cloth) is dirty, it can be wiped clean with wet paper towels. Do not add too much water or other liquids to prevent it from penetrating into the main unit. It needs to be dried after wiping. Can be used again.
3. Do not use any steel wire mesh, iron brush, sharp utensils, etc. and contain volatile and corrosive solvents for cleaning.
4. Normally, when not in use, please put this Product in a storage box and store it in a dry and cool place.

## Key Features

1. Designed according to ergonomics, beautiful and comfortable.
2. Using silica gel PVC heating technology, set low (1), medium (2) and high (3) three levels of temperature, which can be selected according to your needs and ambient temperature.
3. With vibration function, 3 kinds of vibration modes can be set, which can be selected according to personal comfort.
4. It adopts LED indicator to display and buzzer Prompt, which is more convenient to use.
5. USB fast charging, built-in large-capacity lithium battery power supply, rechargeable, easy to use.
6. Built-in 30-minute timer.

## Use of rechargeable batteries

1. In the process of using, if the heating is not obvious or the power is turned off after selecting the heating mode, the buzzer will sound multiple times, and the charging indicator will flash quickly, indicating that the built-in lithium battery is low. Connect the charging cable to the host, and then Plug the USB port on the other end of the charging cable into a mobile power supply or a computer USB port (5VDC power supply) to charge.
2. When charging, the host charging indicator light is continuously flashing, and when the battery is fully charged, the charging indicator-light is always on. It takes about 3-4 hours to fully charge, and it can work immediately.


## Environmental protection and resource

This product uses lithium batteries. Please make sure to dispose of used batteries in a designated place in accordance with local regulations.

## How to remove the built-in rechargeable battery

This product contains a rechargeable battery. Before consumers decide to discard this product, be sure to remove the battery from the product.

1. Before disassembling the Product and removing the battery, be sure to turn off the power.
2. Follow the sequence of steps ①—④ shown in the figure below.
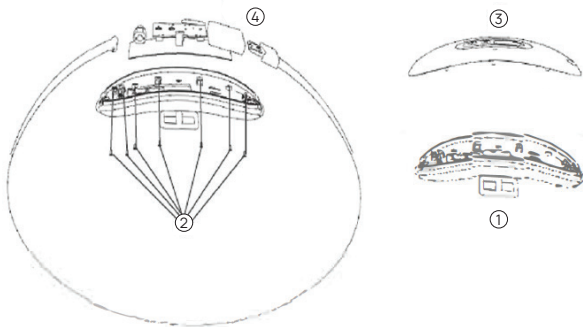
## Precautions

1. It is forbidden to disassemble and assemble this product privately to avoid collision and water ingress to avoid malfunction.
2. It is recommended that each use time not exceed 30 minutes (the product will automatically shut down after 30 minutes of use). If you need to continue using it, please restart the product after 30 minutes.
3. Long continuous use may cause the product to overheat and shorten the service life. If the product is overheated, please stop using it and wait for the product to cool down before using it again.
4. When disposing of this Product, please follow the city's regulations.

## Specifications

|                     |                   |
|---------------------|-------------------|
| Product size        | 160 x 90 x 35mm   |
| Net weight          | About 156g        |
| Power supply        | 5V (DC) 0.8 A     |
| Input power         | 5W                |
| Lithium battery     | 3.7 VDC 1800mAh   |
| Heating temperature | +39°C~+50°C       |
| Charging time       | About 3 - 4 hours |

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