MEDISANA®

home of wellness

DE Massagematte MM 825

GB Massage mat MM 825

FR Tapis de massage MM 825

IT Tappetino massaggiante MM 825

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NL Massagemat MM 825

FI Hierontamatto MM 825

SE Massagematta MM 825

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
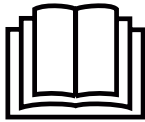
IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

Explanation of symbols

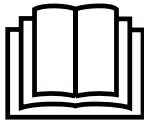This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user.

CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

Lot number

Manufacturer

**about the power supply****Safety Information**

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only insert the mains plug into a socket when the appliance is switched off.
- Keep the cable and the appliance away from heat, hot surfaces, moisture and liquids. Do not touch the mains plug when you are standing in water and always take hold of the plug with dry hands.
- Appliance components carrying voltage may not come into contact with liquid.
- Do not take hold of a device if it falls into water. Pull out the power plug from the socket immediately.
- Connect the unit so that the power plug is freely accessible.
- Always switch off the appliance after use with the buttons on the control unit and unplug the power plug.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull or twist the device on the cable.
- You must not continue to use the unit if the lead or the plug are damaged.
- If the power supply cable is damaged, it may only be replaced by **MEDISANA**, an authorised dealer or by qualified personnel.
- Ensure that the cable will not cause tripping. The cable must not be kinked, trapped or twisted.

specific groups of people

- This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.
- Children must be supervised to ensure that they do not play with the device.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- Do not use the massage mat **MM 825**, or only after consulting your doctor, if:
 - you are pregnant,
 - you have a cardiac pacemaker, artificial joints or electronic implants,
 - you suffer from one or more of the following illnesses and complaints:
Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.

- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Persons, which are insensitive to heat, should be careful using this appliance.
- Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.

before using the unit

- Check the plug, the cable and the massage mat for damage before each use. Do not attempt to use a defective unit.
- Do not use the appliance if it, or the power cord show signs of damage, if it is not in perfect working order, or if it has been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.

operating the unit

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using the device.

Only use the massage mat indoors.

Do not use the massage mat in damp rooms (e.g. when bathing or showering).

- The device must only be used for its intended purpose as described in the instruction manual.
- Using the unit for any other purpose invalidates the warranty.
- Never leave the unit unattended if it is connected to the power supply.
- Do not use the appliance while you are asleep or when you are lying in bed.
- Never use the massage mat if it is folded up or crumpled.
- Do not allow the massage mat to come into contact with pointed or sharp objects.
- Do not stand on the unit.
- The maximum operating time for one session is 15 minutes. Prolonged use reduces the service life of the appliance.
- Very long periods of use can lead to overheating. So after continuously operating for 15 minutes, user should stop it for **at least** 15 minutes before using it again.
- Never put the unit down next to an electric stove or other sources of heat.

- Do not cover up the appliance when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.

maintenance and cleaning

- The unit is maintenance free. If a fault does occur, just check whether the mains power plug is properly connected.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all guarantee rights. Ask your specialist dealer and have repairs carried out by the service only.
- Do not immerse the device in water or other fluids.

Do not wash!

Do not dry clean!

2 Useful Information

Thank you very much

Thank you for your confidence in us and congratulations!

You have gained a high-quality product from **MEDISANA** with the purchase of this massage mat **MM 825**. To derive the most success and to enjoy your **MEDISANA** massage mat **MM 825** for a very long time, we recommend that you read the following information on the use and care of the appliance carefully.

2.1 Items supplied and packaging

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use it and contact your dealer or your service centre.

The following parts are included:

- 1 **MEDISANA** massage mat **MM 825**
- 1 Instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required.

If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

2.2 Range of application

The use of the **MEDISANA massage mat MM 825** will improve your wellbeing after a stressful day. The massage effect is particularly useful for treating tense muscles and tired tissue. The appliance provides pleasant relaxation for overcoming everyday stress. 4 massage zones (Upper back, Lower back, Thighs, Lower legs) can be regulated separately. The combination of massage and additional heating function provides an effective and relaxing treatment. Their settings can be comfortably adjusted by means of the easy-to-handle control unit.

CAUTION

Make sure that the maximum operating time of 15 minutes is not exceeded.

3 Operating

3.1

Operation

1. Connect the appliance to the mains outlet.
2. Press the On-/Off-button "⏻" **1** and all massage zones and the heat function will be activated.
3. Press the zone selection buttons repeatedly (**7** to **10**) according to your personal preference to activate/deactivate the corresponding massage for the selected zone. You may select only one zone or multiple zones for the massage.
4. Press the button multiple times changes the massage intensity or deactivates the massage zone:
 - Weak massage intensity:** The left LED over a zone selection button lights up in green.
 - Strong massage intensity:** The right LED over a zone selection button lights up in green.
 - Deactivated massage zone:** None of the LEDs is alight.
5. **Massage mode selection:**

With the massage mode selector button **6** you can choose one of 5 possible massage modes:

 - 1. Normal:**
Regular vibration massage with constant intensity and vibration frequency.
 - 2. Pulsating:**
Pulse vibration massage, automatically changes the intensity and frequency of the massage.
 - 3. Stream vibration massage:**
Automatically changes intensity and vibration frequency, feels like a massage by "running water".
 - 4. Alternation mode A:**
Automatic alternation of massage zones in short intervals:
Upper & Lower back vibration massage first, then automatically switches to Thighs & Lower legs and so on....
 - 5. Alternation mode B:**
Like alternation mode **A**, but alternates between massage zones Upper back / Lower Legs and Lower back / Thighs.

The respective massage mode is indicated by the green LED lights over the Massage mode selector button **6** (The 1st LED for Mode 1, the 2nd LED for Mode 2 and so on).

Each time you change the massage mode, the massage intensity for all massage zones will automatically be set to intensity 1 (weak intensity).

6. HEAT function button:

Press this button at any massage mode to additionally switch on the heating function. When the heating function is activated, the red LED over the Heat function button **5** will light up. The heating can only be applied to the Lower back area. Press the heat function button again to deactivate the Lower back heating.

7. Press the On-/Off-button "  " **1** to switch off the appliance. The appliance will switch off automatically after approx. 15 minutes of continuous use.

4 Miscellaneous

4.1 Cleaning and maintenance

- Before cleaning the appliance, make sure that the appliance is switched off and the power supply is disconnected. Allow the unit to cool down.
- Clean the massage mat with a slightly damp sponge only. Never use brushes, aggressive cleaning agents, petrol, thinners or alcohol. Rub the unit dry with a clean dry cloth.
- Never immerse the appliance in water to clean it and ensure water cannot get into the appliance.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- Take care to roll up the mains lead to avoid damaging it.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.

4.2 Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Consult your municipal authority or your dealer for information about disposal.

4.3	Name and model:	MEDISANA Massage mat MM 825
Technical	Power supply:	100 - 240 V~ 50/60Hz 1.2A
specifications	Power consumption:	14.4 W
	Automatic switch-off:	after approx. 15 minutes
	Massage zones:	4 massage zones
	Massage intensity levels:	2 massage intensity settings for each massage zone
	Operating conditions:	Only use in dry rooms as described in the operating instructions
	Storage conditions:	In a cool, dry place, storage temperature 0°C - max. 40°C
	Dimensions (L x W x H):	approx. 180 x 62 x 10 cm
	Weight approx.:	1,8 kg
	Article no.	88955
	EAN no.	40 15588 88955 4


In accordance with our policy of continual product improvement, we reserve the right to make technical and optical changes without notice.

Warranty and repair terms Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.
The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

 **MEDISANA AG**
Jagenbergstraße 19
41468 NEUSS
Germany
eMail: info@medisana.de
Internet: www.medisana.com

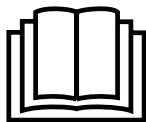The service centre address is shown on the attached leaflet.

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