AVANTA

e.volve[®]
SLIMME WEEGSCHAAL
Model: S5

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Dutch	P1 - P18
German	P19 - P36
English	P37 - P54
French	P55 - P72

PAKKETINHOUD

E.Volve® Slimme Weegschaal

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2 x 1.5V AAA-batterijen

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2 x 1.5V AAA-batteries

TABLE OF CONTENTS

General information	P. 40
Parts / Components	P. 41
Quickstart guide	P. 42
Using the smart features	P. 43
Frequently asked questions & solutions	P. 51

GENERAL INFORMATION

USE THE FITDAYS APP FOR MAXIMUM USABILITY

Scan the QR-code below and download the Fitdays app to enable usage of all 'smart' features of the E.Volve® Smart body weightscale. Alternatively, search for the Fitdays app in your Google play-, Apple- or Android app store.

Thank you for choosing the E.Volve® Smart body weightscale. Paired with the special Fitdays-app this scale allows you to conveniently measure body weight, Body balance, BMI, Body fat rate, Visceral fat, Body water, Skeletal muscle rate, Muscle mass, Bone mass, Protein, BMR, Body age and more.

IMPORTANT: Please read and follow the instructions in this manual before using the device. Failure to comply with the instructions given in the manual, or using this scale in ways other than the ones mentioned in this manual, may result in personal injury or damage to the product itself.

Standard functions (without using the Fitdays app)

- Measuring of general body weight
- Changing the displayed weight units (Kg, Lb, St)

Additional features (while using the Fitdays app)

- Bluetooth pairing
- Combined app-usage up to 24 individuals
- Tracking progress, weight gain and loss
- Standard weight
- BMI (Body Mass Index)
- Body fat ratio
- Fat-free body weight
- Subcutaneous fat ratio
- Visceral fat levels
- Body water ratio
- Skeletal muscle ratio
- Muscle mass
- Bone mass
- Protein
- BMR (Basal Metabolic Rate)
- Body age
- Baby modus

Other (manual) functionalities of the Fitdays-app

- Growth / length tracking
- Body circumference tracking
- Synchronisation with Fitbit and other health apps and health devices other than this scale

GENERAL SAFETY

- Do not use with medical implants such as pacemakers
- Maximum load capacity: 180kg
- Do not drop the scale, or drop products / liquids on the scale
- Store in a dry environment with a stable surface
- Do not place/store other products on the scale, as this may damage the sensors
- Only clean the scale with a slightly damp cloth
- Do not immerse in water and refrain from using chemical cleaning agents as they can damage the sensors
- Not recommended for use during pregnancy
- Not recommended for use by babies and children up to 16 years of age without guidance of an adult
- Obtained measurements are meant for reference only and should not be interpreted as a medical opinion
- Always consult with your primary care doctor or physician prior to changing your diet, exercise regime or other forms of physical activity

- Only use the scale on a hard, dry and flat surface
- Do not step on the edges of the scale but distribute your weight evenly across the sensors
- Do not jump on the scale
- Make sure your feet are dry before stepping on the scale
- Replace batteries as indicated on the safety sticker, located at the bottom of the scale (2x 1.5V AAA)
- Remove the batteries when the scale will be out of use for prolonged periods of time to avoid a possible battery leakage
- This product contains small parts and batteries which could be harmful for children when swallowed

PRODUCT SPECIFICATIONS

Size: 300x260x20mm
LCD-unit: 74x28mm
Available weight units: Kg / Lb / St
Power supply: 2x 1.5V AAA
Maximum number of users: unlimited
Maximum number of users in the Fitdays app: 24 users
In-app Baby modus: allows you to measure the weight of your baby
In-app Height record: allows you to track growth/length for users up to 16 years of age
Min - max operating temperature: 10-40°
Min - max operating humidity: 20% - 90%

SUPPORTED SMARTPHONES/TABLETS/DEVICES


iOS 8.0 or higher
Android 6.0 or higher

Please note: in order to use the scale in combination with the Fitdays app, your device must be equipped with working bluetooth functionality.

PARTS / COMPONENTS

- 1 Weight sensors
- 2 LCD-screen

- 3 Battery protector
- 4 Battery compartment with lid
- 5 Unit selection button (Kg, Lb, St)


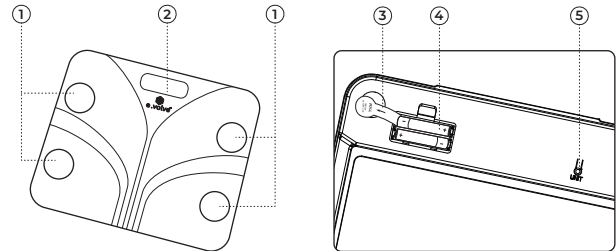
QUICKSTART GUIDE

Please follow the steps on
E.Volve User Guide

POSSIBLE LCD-INDICATORS

Lo	The scale is low on batteries, please replace the batteries (2x 1.5V AAA).
Err	Overload, please allow a maximum load capacity of 180 Kg only, step away from the scale and try again.
C	An error occurred during the analyzing process, step away from the scale and try again.

EXTRA: CHANGE THE DISPLAYED WEIGHT UNIT

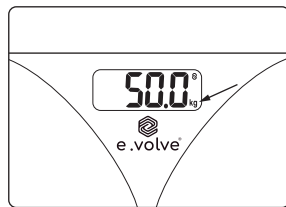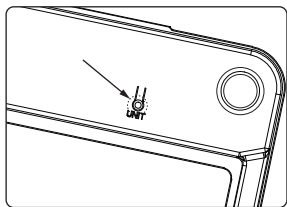STEP 1: Activate the scale by tapping one of the sensors with your hand or foot

STEP 2: Change the weight unit by pressing the dedicated selection button (located in the battery compartment) 2 to 3 times in order for changes to appear in the right corner of the LCD-display

POINTS OF ATTENTION

- Place the scale on a flat and even surface such as tiles and floorboards, do not place on carpets or rugs
- Distribute weight evenly across the scale
- Use with bare feet (in order to activate the sensors and use the smart features)

USING THE SMART FEATURES

STEP 1

Install the Fitdays app on your mobile device or tablet

Scan the QR-code and download the Fitdays app to enable usage of all 'smart' features of the E.volve smart Body Weightscale.

Alternatively, search for the Fitdays app in your Google play-, Apple- or Android app store.


STEP 2
Activate bluetooth-connection

Ensure bluetooth is activated on your mobile device to allow pairing with the scale.

STEP 4
Register & enter your personal data

Follow the registration instructions in the app and add your personal data. Please ensure that you enter these accurately. Failure to comply information regarding gender, height, age and year of birth may result in deviated / incorrect results and measurement data.

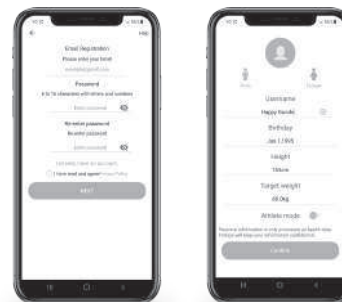STEP 3
Activate the app

Open the Fitdays-app and register as a new user, or log in as an existing user would you already have an active account. Registration and usage of this app are both free of charge.


STEP 5.1
Pair the scale with your app

After registration, click on the 'account' icon at the bottom of the screen.


STEP 5.2

Click on the 'device' icon which appears at the top of your screen after which a new (blank) screen appears which is called 'my device'. Click on the '+' icon at the top of this screen.

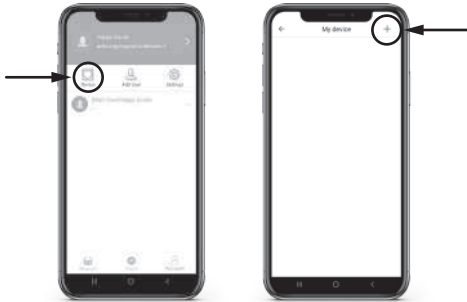STEP 5.4

A new screen will appear, where the smart scale is listed as "MY_SCALE". Click on this item and hit confirm in the pop-up that appears next. Your scale will now be paired with the Fitdays app.

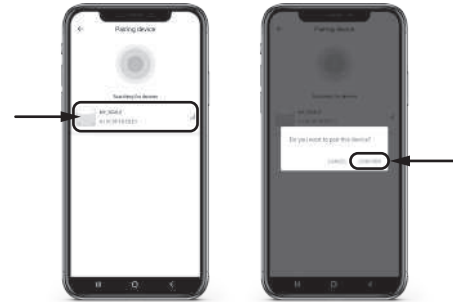STEP 5.3

In the menu which appears next, select the Bluetooth option and directly after, tap on one of the sensors on the scale in order to activate.

STEP 5.5

The Fitdays app will show you a short user guide. You can either watch it, or ignore by clicking on the 'X'.

STEP 5.6

The Fitdays app will ask for access to update your Health data (depending on your type of mobile device / tablet).

When enabled: you can find your results both in the Fitdays app, on the display of the scale AND any other health apps on your device.

When disabled: you can still find your all your measuring data in the Fitdays app and the LCD display of the scale.

STEP 5.7

If you have never used the Fitdays App before the screen will display another quick user guide, similar to the below image.

Click on 'start measuring' to obtain your first smart results!

Reminder:


- Place the scale on a flat surface.
- Do not wear socks or shoes.
- First click on start measuring on your smartphone, then step onto the scale.
- The scale will take approximately 10-30 seconds to obtain all data: do not step off of the scale as you will have to start all over if you step off too quickly.

STEP 6

Finding your full results


The complete and most recent analysis can be found by tapping the "More Data" button. Tapping the "Charts" icon at the bottom of your screen will display progress, weight gain and weight loss of all measurements which are analysed over time.

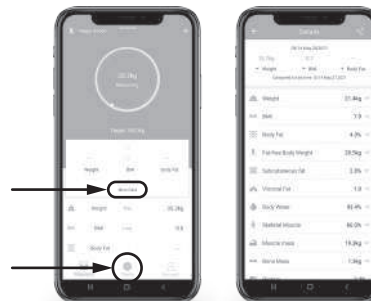TAKING FUTURE MEASUREMENTS

Whether in your personal account menu, or in the main menu, tap the "Measure" icon at the bottom right of your screen, and step on the scale. Please note: always use the scale with bare feet, failure to do so will disallow the scale to measure and analyse all the smart results. Do not step away from the scale until all measurements are calculated and analyzed. The app will send you a notification by sound, the progress-circle will stop moving AND the app will display a message which indicates "Done".

Please note: always wait until "done" appears on your screen in order to fully synchronize the data between the scale and the app. Do not click on any other screens or in-app buttons while the scale is analysing as this will disturb the data transfer.

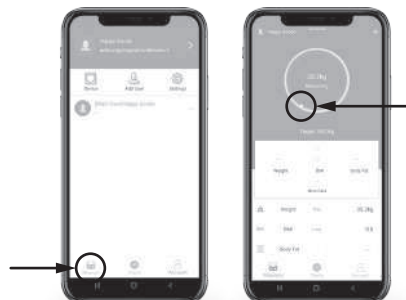
EXTRA FUNCTIONALITIES

Adding multiple users / family members to the app

Tap the "Account" icon and choose "Add user". Enter the personal data and tap "Confirm". All users will now appear in the user list, where you can switch between users. Or simply step on the scale and let the app suggest the already registered users it can recognise.

Alternative method: other users may also install the app on their own mobile device for a full setup and use the app separately.

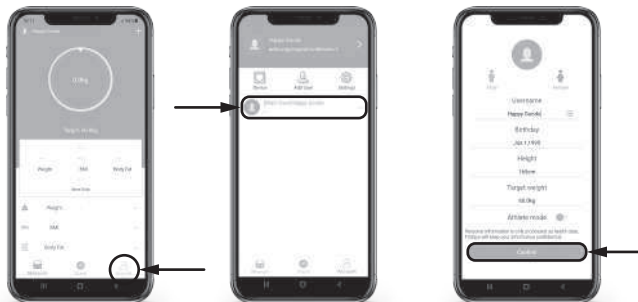
Baby Mode

Tap the "Measure" icon at the bottom of the screen followed by tapping on the "+" icon at the top of the screen. In the dropdown menu which now appears, select "Baby Mode" and follow the provided instructions.

Pair with Fitbit, Google Fit, Samsung Health or Apple Health

Tap the "Account" icon, followed by the "settings" icon. All pairable devices are listed in the settings menu. Choose the preferred device, follow the instructions on the Fitdays app and instructions which appear on your device.

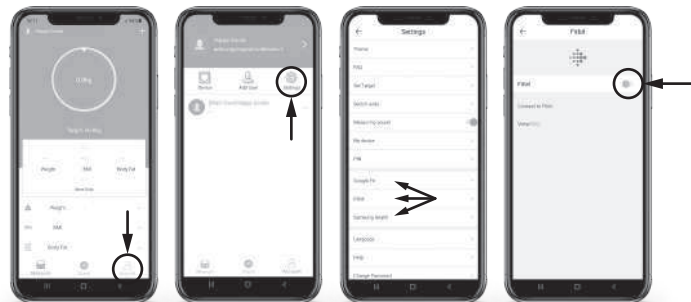FREQUENTLY ASKED QUESTIONS & SOLUTIONS

Bluetooth does not seem to be connecting

- Check if the bluetooth functionality on your mobile device is switched on and working
- Check if the bluetooth icon is displayed on the scale's LCD-screen
- Check your mobile device's operating system and update which must be:
 - iOS 8.0 or higher
 - Android 4.3 or higher
- Open the Fitdays app, tap the "account" icon, followed by the "device" icon to double check that the scale is listed in the app
- Install / pair the scale again if needed

Body fat percentage and other smart data is not displayed

- All additional measurements can only be analysed by combined usage of the scale and the Fitdays app
- All additional measurements can only be analysed by using the scale by stepping on it with bare feet
- Ensure both scale and feet are dry, do not wear socks or other footwear
- Ensure scale and app are paired and used together
- Ensure your mobile device's bluetooth is switched on
- Start the analyses by tapping the "Measure" icon on the app, step on the scale and remain standing on the scale until the app displays "done"
- Do not click on any icons in the app while the analysis is taking place and remain in the main menu

The scale is showing inconsistent results

- Ensure the scale is placed on a flat and hard surface and make sure the scale is not tilted, as this could influence results
- Check the 4 flexible support holders at the bottom of the scale, ensure these are intact and not clogged up by dust or dirt
- Ensure the scale is not placed on carpets or rugs
- Remove and replace the batteries

The app displays incorrect fat-muscle ratios and other smart measurements (BMI / body water % etc.)

- The personal data is entered incorrectly
- Ensure all personal data is accurately entered (gender, age, height/length etc.)

How can I re-calibrate the scale after moving it or change the batteries?

- Touch the sensors with your hand or foot, only to activate the scale (do not yet step on it)
- Wait until the scale automatically drops to 0.0kg

How can I change the displayed weight unit (Kg, Lb, St)

- Press the small button located in the battery compartment 2 -3 times
- Check the changes appearing on the LCD-display

The scale does not turn on

- Ensure the battery protector has been removed
- Ensure the batteries are installed correctly with the polarities pointing in the right direction
- Replace the batteries, they might be faulty or empty
- In case of a battery leakage: erosion could have damaged the battery connectors, try to carefully clean these with a cotton-bud before inserting new batteries

BATTERY ADVICE

- Remove the batteries when the scale won't be used for prolonged periods of time to prevent battery leakage.
- Keep and store batteries out of the reach of children.
- Do not combine old and new batteries, or batteries from different brands to avoid.
- Do not place batteries near / in any type of heat source or open fire.
- Batteries should not be disposed with regular household waste.
- Please check the regulations and solutions offered by your local authority.

Electrical appliances should not be disposed of with regular household waste. Please recycle where facilities exist and contact your local authority for advice and regulations.

PRODUCT DISPOSAL

Parts of our packaging are recyclable. Please dispose of it in an environmentally friendly manner by separating the different materials such as cardboard and plastic. Would you stop using the product while it is still in good working condition, we kindly ask you to consider donation to your local thrift store or charity.

Electrical appliances should not be disposed of with regular household waste. Please contact your local authority for advice on nearby recycling centres for used electronic devices.

Batteries in this product should be removed prior to the disposal of the device and disposed of separately in accordance with the regulations of your local authority. To protect the environment, non-rechargeable and rechargeable batteries may not be disposed of with normal household waste, but should be taken to suitable collection points.

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