

Troubleshooting

State	Possible causes	Treatment measures
Unable to connect power	The host battery is discharged. Use in low or high temperature environment.	Please charge the host. Use at the recommended ambient temperature of (15-35 °C).
No electric stimulation	The power supply is not connected	Briefly press the OK key. Until LED light comes on.
Current stimulation is too weak	The setting intensity is low	Press to increase the strength
Too strong current beam excitation	The setting intensity is too high	Press to weaken and reduce the strength
Sudden interruption of stimulation	The remaining power of the host battery is exhausted	Please charge the host.

⚠ Danger

Do not use this product together with the following medical electronic equipment. It may cause accidents or lead to failure of medical electronic instruments, causing serious damage to the body.


⚠ Warning

Do not use it for the following personnel.
Otherwise, accidents or skin and body problems may be caused.

All electronic rejection products on pregnant women, women just after birth, and patients.

Host usage

- The stimulation intensity is divided into 19 grades (intensity 0-19 grades).
- Please do not set a lower strength when using it for the first time. Please adjust other strength after getting used to low strength.
- While the remaining power of the battery is gradually reduced, the current intensity will also be gradually weakened. If you continue to use it, please charge the host.

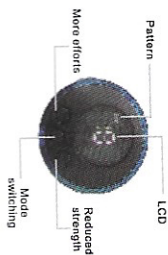


(B-mode)
When the machine is switched on to gear 1-19, it will automatically shut down for 15 minutes normally. If gear 1 is switched to gear 1-19 and is not used, it will automatically switch to gear 0 in 5 seconds, and the switch will automatically shut down in 20 seconds.

Technical parameters

Please set a lower intensity when you use it for the first time. Please adjust other strength after getting used to low strength. For details, please refer to the following message mode guide.

Product name	Magic Paste
Voltage	3.7V
How to use	A-F (8 modes) 19 gears
Frequency	1~100HZ
Working hours	15 minutes
Output	35mA maximum
Size	62x62x16
Weight	31
Battery capacity	200mAh

MAGIC PASTE

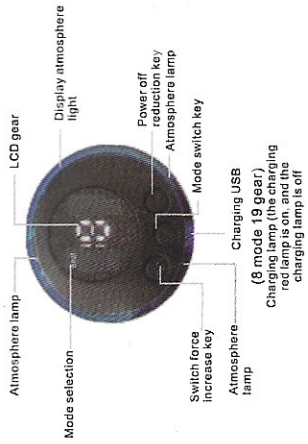

Matters needing attention
 1. When charging, do not unplug the host, just plug in the charging wire for charging.
 2. Two clips must be installed properly when using.

LCD model (8 modes and 19 gears)
 1. It can be used for about 160 minutes after one hour of charging. The greater the strength, the greater the power consumption. If it cannot be used during charging, it can be used after charging.
 2. The button of the host is buttoned on the mat like a button, and pressed on the skin+press twice for startup/increase the force - decrease the force/shut down m switch mode.

What is EMS technology?

Electrical muscle stimulation (EMS), also known as neuromuscular electrical stimulation, is the use of electrical pulses to induce multiple muscle contractions. The central nervous system makes muscles contract, EMS has proven to be a positive muscle strengthening tool.

What is the purpose of EMS?

- Muscle strengthening
- Prevent wasting muscle atrophy
- Muscle relaxation
- Relieve pain
- Relax your muscles
- Cosmetic muscle conditioning

Product warranty card

This card is for the buyer's reference. Please keep it safe. The information on this card is for reference only and does not constitute a warranty. The warranty is subject to the terms and conditions of the warranty.

Model Name	Year, Month, Day	Serial Number
Color/Style	Center	Manufacturer
Manufacturer	Manufacturer	Manufacturer
Product Name	Product Name	Product Name
Product Name	Product Name	Product Name

Warranty Regulations

This card is for the buyer's reference. Please keep it safe. The information on this card is for reference only and does not constitute a warranty. The warranty is subject to the terms and conditions of the warranty.

1. The product is guaranteed for 12 months from the date of purchase. If the product is found to be defective within this period, the manufacturer will provide a replacement or repair at no cost to the buyer.

2. The warranty does not cover damage caused by misuse, accident, or unauthorized repair.

3. The warranty is void if the product is used for purposes other than those specified in the user manual.

4. The manufacturer is not responsible for any consequential or incidental damages.

5. The manufacturer reserves the right to modify the warranty terms without notice.

Certificate

Product Name: **Massage Gun**

Number: **1818181818**

Date of production: **2023-10-27**

Examination date: **2023-10-27**

The product has been inspected and found to meet the quality requirements. It is suitable for use.