MVOLO

Safety and Usage Instructions

TELESCOPIC TUBE

**INFRARED HEAT
LAMP**

Product Introduction

This product has superior performance and is designed and produced in strict accordance with the safety standards of household appliances, safe and reliable, in line with the ideal energy value range of the human body, it is in home life My good helper! Easy to use, start and use.

Usage Method

1. Open the package, take out the product, install the bulb and lamp net cover, and adjust the lamp arm and lighting the required position and angle of the spotlight head to the body.
2. Do not face the human body. After the power is turned on, preheat for about one minute before using, the irradiation time (the same part) is about 30 minutes for adults and 30 minutes for children, about 20 minutes.
3. During auxiliary use, the skin of the irradiated part shall be exposed at a distance of 16 inch, The skin temperature is about 40° F Celsius or varies from person to person. i.e. I feel comfortableMake available.


Matters needing attention

Please read the instructions carefully before use

2. Do not touch the irradiating head, the heating plate in the net cover and other parts during the use, children and people with intellectual disabilities must be supervised by adults to avoid burns or injuryElectric shock accident.
3. When irradiating the face, wear colored glasses or eye masks. Do not irradiating the eyes directly, avoid eye injury.
4. White smoke may appear when the lamp cap is used for the first time or after a long timePhenomenon, which is caused by moisture absorption of thermal insulation material of irradiation lamp cap, which is a normal phenomenon, It will disappear after preheating for a period of time. In case of loss or failure of the product, it is not allowed to disassemble by yourself, contact our company or dealer.
5. Skin allergic reaction or other symptoms during irradiation due to personal physical differences in case of discomfort, stop using it immediately.
6. This product should not be use in strong electromagnetic environment and should avoid electromagnetic interface.
7. Never put water or liquid on the lighted bulb, otherwise the bulb will heap upCold shrinkage, explosion.
8. Avoid scalding when turning on the lamp. Do not touch the bulb.

How to Install:

Note: Do not install in the wrong direction!

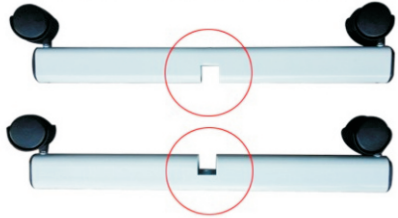- 1 Insert the wheel into the hole

- 2 Fix the expansion pipe and screw it up

- 3 Pull out the stainless steel telescopic rod

- 4 Collect and clamp the lamp cap

Note: It is not too long to pull 3/2 long

- 5 Pull out the stainless steel telescopic rod

Anti scalding Excellent Lamp Holder

Safer, more stable, freely adjustable and convenient to use.

If you have any questions about your product, please contact us at:

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