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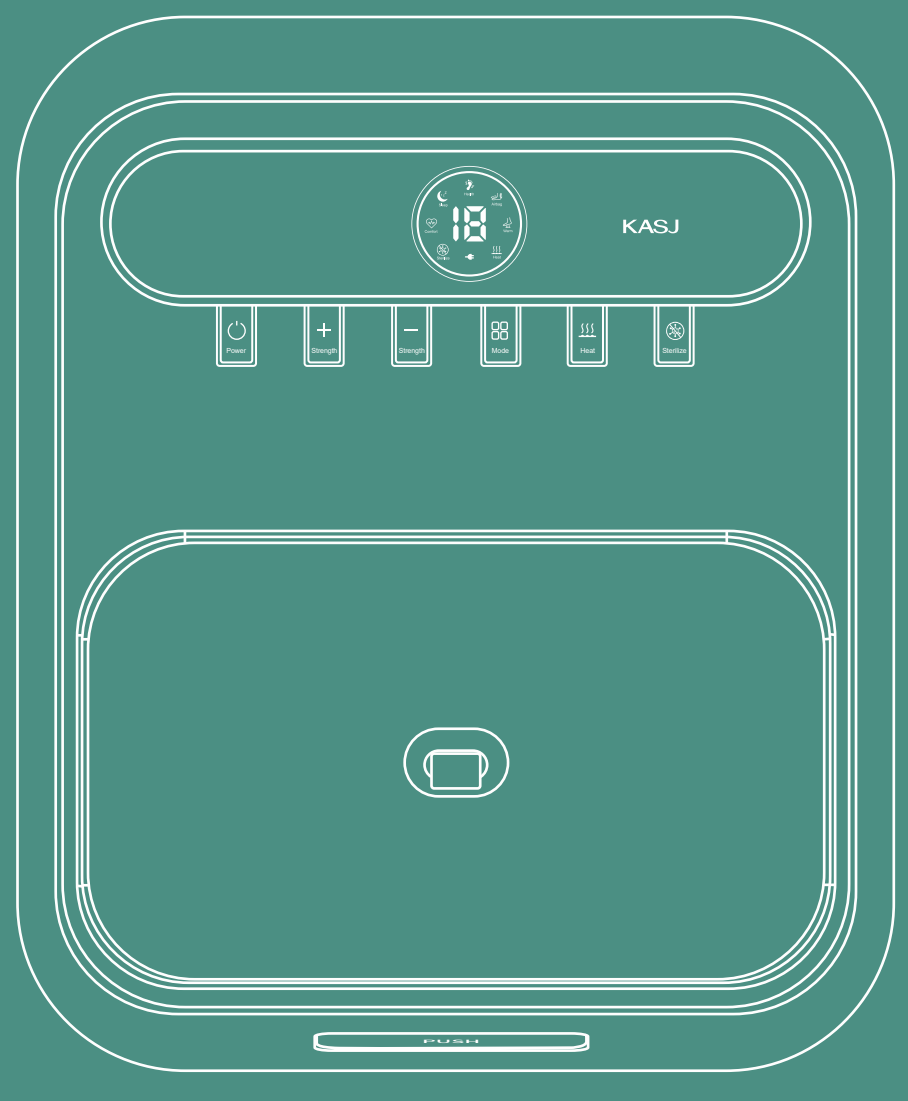
**ENG**

LUNASEA

Thank you for your purchase! Please read the user manual carefully to ensure that you use the foot massager correctly. Keep this manual for future reference if needed.

**PRODUCT PARTS**

**UNPACKING PARTS**

**ENG**

**CONTROL PANEL**

- 1. Control the opening and closing of the whole machine;
- 2. Increase control strength;
- 3. Decrease control strength;
- 4. Control the switching of comfort, sleep, health, air bag and warm modes;
- 5. Control the opening and closing of machine heating;
- 6. Ozon sterilization;

**Attention**  
Please don't press and hold the key too hard, so as to avoid program error.

**CHARGING INSTRUCTIONS**

1. When the power of the machine is low, "K" Flicker, while charging "KK" Flashing from left to right, when charged "KKK"

**NL**

**OPERATING INSTRUCTIONS**

1. After power-on, press "one-key start" once, and the buzzer "beeps", and the display screen shows that the "1S" soothing indicator lights up. The massage intensity defaults to 3 gears, and the working time defaults to 15 minutes.
2. After starting up, press the "Strength Key +" once, and the buzzer "beeps", and the setting intensity increases by n+1. Up to 6th gear. After setting, the remaining time is displayed on the display screen.
3. After starting up, press the "Strength Key -" once, and the buzzer "beeps", and the setting intensity increases by n-1. The lowest is 1st gear. After the setting is completed, the remaining time is displayed on the display screen.
4. After starting up, press the "Mode key" once, the sleep indicator lights up, the buzzer "beeps", and the foot massage and airbag start to work; Push the "mode key" repeatedly, the health indicator lights up and the buzzer beeps. Push the "mode key" repeatedly, the airbag indicator lights up and the buzzer beeps, and the foot massage stops working. Push the "mode key" repeatedly, the warm indicator lights up and the buzzer beeps, and the foot massage stops working and the heating starts working. (Intensity + and heating functions cannot be adjusted in this mode).
5. After starting up, press the "heat key" once in the non-warm foot mode, the buzzer "beeps", the heating indicator on the display screen lights up, and the machine starts heating; Press the "heat key" again, and the buzzer beeps, and the machine stops heating.
6. After starting up, press the "Sterilize Key" once, and the buzzer beeps, and the sterilization indicator on the display screen lights up, and the machine starts sterilization; Press the "Sterilization/Disinfection Key" again, and the buzzer beeps, and the machine stops sterilization.

**ENG**

**ATTENTION**

- It is forbidden to use this machine standing, so as not to damage the machine and cause danger!
- Do not operate the machine with wet hands, and do not use the machine near an open flame.
- Do not place this machine in wet places such as bathroom or swimming pool.
- If you feel uncomfortable (such as dizziness, headache, etc.), please stop using it.
- If the machine is found to be malfunctioning or abnormal, stop using it immediately, and consult or repair with our company or its maintenance centers, and do not disassemble and repair it by yourself.
- If the power cord is damaged, in order to avoid danger, it must be replaced by the manufacturer or its maintenance department or full-time personnel.
- The surface of the machine is hot, and people who are insensitive to heat are forbidden to use it.

**MAINTENANCE AND STORAGE**

- Before cleaning, please make sure that the power switch is turned off and the power plug is unplugged from the socket.
- Please wipe the surface of the product with a soft cloth.
- Dirty places should be cleaned with soft cloth dipped in water or neutral solvent, and then dried with a dry cloth.
- Keep the product clean and away from high humidity or high temperature and environment with volatile gases.
- If you don't use the product for a long time, please put it in the packaging box to prevent dust from falling into it.
- Do not unplug the power plug with wet hands.

**ENG**

- Do not rinse the product with water, or spill water or other liquids into the inside of the product.
- Do not use corrosive detergents such as benzene or diluents to clean products, otherwise the products may fade, and the product structure may be corroded and cracked.
- Do not place the product in a place with direct sunlight or high temperature, so as not to fade the product.

**WARNING TIPS**

- Do not pull the zipper hard when removing the foot cover, so as to avoid damage.
- Please do not use this machine in the bathroom or other places with humid air. If it is soaked in water or damp, please do not use it. It should be repaired by our company or distributor.
- Please use this product in the environment of 0-40 °C.
- Please do not store and use this machine near the heater or electric furnace to avoid sunlight.
- This product is for domestic use only, not for commercial use.

When using it all the year round, please check whether there are the following symptoms:

- Have a burnt smell;
- After touching the wire, sometimes it is powered on and sometimes it is powered off;
- Heating of power cord;
- Other exceptions.

**WARNING**

If the above phenomenon occurs, please stop using the product in order to prevent the occurrence of fault and accident. Turn off the power switch, unplug the power plug, and contact the dealer or our company where you purchased this product. Do not disassemble and repair by yourself.

**ENG**

**COMMON FAULT HANDLING**

Common faults	Reasons	Exclusion method
No power	Low battery	Charging the machine with an adapter
No power	The power adapter is not plugged in properly	Plug in the power plug according to the correct method. Please contact the seller if none of the above methods can solve the problem.

**WARRANTY CARD**

Dear Users,  
Thank you for using the products. Our company will provide you with the following services according to the relevant provisions of the "Consumer Protection Law" with this warranty certificate.

1. According to the regulations, the warranty period of our products is one year (calculated from the date of purchase). During the warranty period, the company will be responsible for free maintenance for any fault caused by the quality problems of the products under normal use.
2. During the warranty period, if the fault is caused by the quality problem of the product itself, please bring the filled warranty card to the authorized maintenance center of our company for free maintenance.
3. Once the warranty card is altered, the warranty will expire immediately.
4. Keep the warranty card as the warranty certificate of this machine. Please keep it properly and don't make up for it if it is lost.
- Repair is not free in the following situations:
  - A. No warranty certificate;
  - B. Fault caused by fault to operate the machine according to the requirements of the manual;
  - C. Damage caused by dismantling by non-special maintenance personnel of the company;
  - D. Fault, scratch or breakage caused by movement or falling;
  - E. Damage caused by improper storage, maintenance and use by users;
  - F. Damage caused by natural disasters.

**ENG**