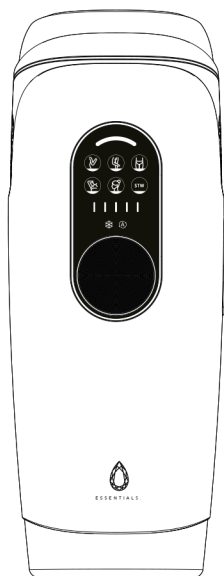




ESSENTIALS.

Instruction Manual

Please read this manual carefully.



IPL Ice Cooling Hair Removal Device

MODEL: LX1000

Contents

Safety Tips	1
Product Introduction	5
Preparation Before Use	9
Operation	12
Maintenance After Use	16
Cleaning and Maintenance	16
Common Problems and Solutions	17

Safety Tips

Before use, do not apply any gel or skincare products (aloe vera gel, etc.) on your skin!
Please assemble and use the product according to the Instruction Manual.
Please read the Instruction Manual carefully and keep it safe for future reference.

- ▲ **Note:** Contents concerning any personal safety issues or property damage that may be brought about by misoperation
- ❗ **Warning:** Contents concerning any danger that may be incurred by misoperation
- ☠ **Danger:** Contents concerning any physical injuries that may be caused by misoperation

⚠ Notes

- 1.Never use the device in a dark environment
- 2.To prevent malfunction, do not use the device in direct sunlight
- 3.Wear goggles when you use the device
- 4.Do not use this device if your skin is naturally dark
- 5.When you first use the device, start from setting 1
- 6.Avoid exposure to the sunlight before and after treatment
- 7.Shave your hair and clean your skin before treatment (do not use chemical methods)
- 8.During treatment, you must press the light exit window against the skin
- 9.Never subject the device to heavy shocks and do not shake it
- 10.To avoid damage, never unplug the device without turning it off

Notice: When using the hair removal device for a break, please make sure that the inlet is facing up and the LED is facing down, which is beneficial to thermal dissipation and cooling.



WARNING

- Before use, clean your skin and make sure that the skin is free of residual hair, oily matters and dirt
- If the skin of the treatment area is not dry or clean, the treatment result of pulsed light will be affected, and non-flashing or unintentional flashing of the device may occur
- Please select the light intensities based on your skin conditions. You must start from setting 1 when you use it for the first time
- Do not use the antiperspirant immediately after treating your armpit
- When you need to treat and cool your back, it is recommended that you find an assistant to help you, as it may be dangerous when you do it alone
- Slight redness appearing on the skin after use, if any, is a normal phenomenon. You may use the cooling function to relieve the discomfort
- Keep the treatment area from direct UV exposure within 24 hours after treatment
- Do not touch the light exit window immediately after use, for the high temperature of it may lead to burns
- Please scrub off the hair and skincare cream on the device after use to prevent device malfunction and skin damage during further use.
- Do not use it around children, for the device will give off bright light that may hurt their eyes
- To prevent eye injuries, do not use the product on the areas around the eyes
- Do not use it near flammable and explosive materials (including sprayers and gas stations), which may cause fire and electric shock injuries, etc.
- To prevent injury, do not use the device on the pets
- To prevent malfunction and damage, do not place the light exit window near anything that is easily demagnetized
- To prevent malfunction and fire, do not use the product when air ventilation slots are obstructed
- To prevent the device from the abnormal operation, do not direct the light exit window at any black objects
- To prevent being burned black, discolored or broken, do not use the device on any clothes, hair or other objects
- Do not leave the product unattended when you turn it on. Turn it off in time after use.
- To prevent burns and malfunction, please check and clean the dirt on the light exit window regularly
- The ambient temperature for the operation of this device is 5°C~40°C. Please check the room temperature before use
- If the product is to be used under a temperature difference between a new environment and a previous one of over 20°C (from a low-temperature environment to a high-temperature environment, or otherwise), please turn off the device, leave it still and use it after 2 hours
- Do not expose the product to sunlight or to UV light for a long time, which will seriously damage the product
- Do not use the device if the built-in fan does not work when the power is turned on (no buzz is heard), so as not to cause smoke, fire and other faults

DANGER

- Wear goggles when you are using the product. While the light won't act directly on your eyes, the high-frequency flash may be harmful to your eyes
- If your hair is longer than 1mm, you should use the razor to shave your hair first. The heat absorbed by excessively long hair on the skin surface can lead to redness and swelling. Unpleasant odor from scorched hair may have a negative effect on your user experience and the treatment result, and cause damage to the product
- Skin tone varies across your body, so different light intensity settings should be applied for actual conditions. Please begin with setting 1, and increase the setting level by level (according to your conditions and the acceptance level of the skin)
- If the skin of the treatment area is not dry or clean, the treatment result may be affected and non-flashing or unintentional flashing of the device may occur
- We suggest that you ask an assistant or partner for help when you are treating or cooling your back as we think it may be dangerous if you do it alone
- Never use the product on the areas around your eyes. Wearing goggles is a common practice that can protect users. When you are using the product on the areas around your eyes, goggles may not be able to fully shield you from the energy. The light will be absorbed by the melanin in the iris at a close range
- Do not use the product close to a bathtub, shower tray, washbasin or other containers filled with water. Please keep the product and the power adapter dry. Do not use the product when you are taking a bath or shower
- Do not insert anything into the product or leave it about
- Do not touch any of the product's internal parts if it's damaged to avoid an electric shock
- If the power line is damaged, in order to avoid danger, contact the customer service for the replacement of it by professionals of maintenance department or like departments of the manufacturer
- Since we lack the data on the use of the product on children under the age of 15, you should never allow any children under the age of 15 to have access to the product. You should also keep a safe distance from any children under the age of 15 when using the product.



THE PRODUCT IS NOT SUITABLE FOR THE FOLLOWING GROUPS OF PEOPLE OR UNDER THE FOLLOWING SITUATIONS

- Expectant mother, women in lactation and menstrual period, and people with superficial scar, diabetes, fragile skin and skin sensitive to lights must not use the product
- People with severe heart disease, moderate or severe hypertension, malignant tumor, epifolliculitis, infectious skin diseases, burns and coagulation disorders must not use the product

- People with abnormal collagen metabolism, including a history of scar and poor wound healing must not use the product
- People with vascular disease, such as varicosity and vascular ectasias at the treatment area must not use the product
- People with skin sensitive to light and easily develops a rash or an allergic reaction must not use the product
- Do not use the product if the treatment area is infected or accompanied by herpes simplex, skin lesion or haematoma, etc
- People with diseases suppressing immunoreaction (including HIV infection or AIDS, lupus erythematosus, porphyria, etc.) must not use the product
- Anyone with hemorrhagic diseases and anyone externally uses isotretinoin, AHA, BHA or azelaic acid must not use the product
- Do not use this product if you are currently on aspirin, immunosuppressive drugs, isotretinoin, etc. Do not use this product if you are on sun-sensitizing drugs recently (such as tetracycline, sulfanilamide, etc.)
- Taking a painkiller will reduce the sensitivity of the skin to the heat, please use the product with caution before its effect disappear
- Anyone who has used chemical hair remover in the past 30 days, and who has used skin whitening products or cosmetics containing vitamin A (shortly) before treatment must use the product with caution
- If you have a sunburn, do not use the product before the skin fully heals. Do not use the product on sensitive areas such as the mucous membrane
- Do not use this product for any types of unhealed wounds. People with hirsutism and hypertrichosis induced by other diseases are not recommended to use the product
- Do not use the product together with essential oil. Please make sure that the skin of the treatment area is clean and dry, and no lotion, lubricant or other cosmetics are applied to the skin. If the skin surface is not clean and dry, the user experience may be affected and the product may be damaged
- If your skin tone or hair color is not proper, do not use the product (if necessary, consult your doctor before using it)

⊗ DO NOT USE THE PRODUCT ON THE FOLLOWING AREAS OR UNDER THE FOLLOWING CIRCUMSTANCES

- | | |
|---|-----------------------------|
| 1、 The inside of the mouth, private parts, mucus membrane, etc. | 8、 Eyebrows and hair |
| 2、 Surgical site of plastic surgery and operative site | 9、 Ears |
| 3、 Laryngeal prominence | 10、 Nipples and areolas |
| 4、 Eyeballs, eyelids, and around the eyes | 11、 The inside of the navel |
| 5、 Head | 12、 Large vein |
| 6、 The areas with acne | 13、 Tattoo, body painting |
| 7、 The areas with spots and moles | |

(Notes: This list is not complete. If you are not sure whether you can use this product, it is recommended that you consult your doctor)

PRODUCT INTRODUCTION

- This is a miniaturized product of IPL freezing point hair removal technology in cosmetic medicine to be used at home. It uses sapphire freezing point skincare technology that can cool down the skin, reduce the heat on the surface and keep the treatment area from burns, so as to better protect the skin. During treatment, press the light exit window against the skin, press the hair follicles to make their roots closer to the surface area with higher energy, which increases the energy absorptivity of hair follicles. Comparing with the traditional hot light hair removal device with a concave surface, the sapphire light exit window can be cleaned easily with a scrub.
- This product has 6 body modes, 5 settings of intensity and an automatic skin tone identification system that satisfy different skin sensitivity.
- The product has an integrated light exit window with a long service life, making it needless to be replaced, saving you the cost for replacement. This means that you don't need to buy accessories frequently. Accordingly, there will be no hidden additional cost for it. In addition, the product has a complete filtering system that can filter unnecessary wavelengths and harmful light

THANK YOU FOR CHOOSING IAM HAIR REMOVAL DEVICE. PLEASE READ THE INSTRUCTIONS CAREFULLY BEFORE USE

Model	Specifications	
LX 1000	Product name	IPL Ice Cooling Hair Removal Device
	Size	170mm×65mm×43mm
	Input	24V ---
	Power	60W
	Optical exposure	2-5J/cm ²
	Light exit window area	3.0cm ²
	Wavelengths	600~1200nm
	Standard	EN IEC 55014-1: 2021 EN IEC 55014-2: 2021 EN IEC 61000-3-2: 2019+A1:2021 EN61000-3-3: 2013+ A1:2019+A2:2021

PACKING LIST


Name	QTY	Name	QTY
Host	1Set	Care mode attachment	1Piece
Power Cord (whit Adapte)	1Piece	Hand mode attachment	1Piece
Instruction Manual	1Piece	Leg mode attachment	1Piece
Face mode attachment	1Piece	Private mode attachment	1Piece
		Armpit mode attachment	1Piece




PRODUCT OVERVIEW

Name of parts


1.Screen	2.Flash button	3.On/off/setting adjustment button	4.Mode adjustment button
5.Air outlet 1	6.Skin tone sensor	7.Light exit window	8.Air outlet 2
9.Air inlet	10.Flash indicator	11.Mode indicator	12.Setting indicator
13.Auto flash icon	14.Cooling icon	15.Power socket	16.Power line of adapter

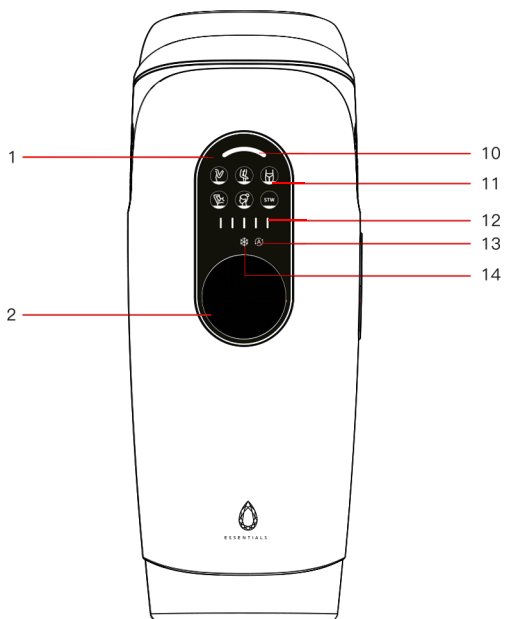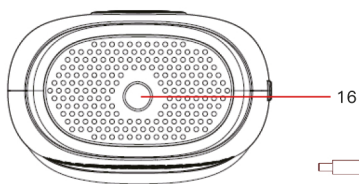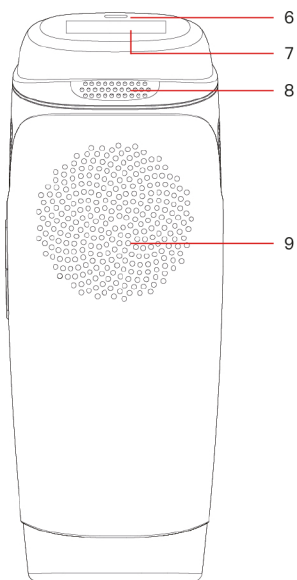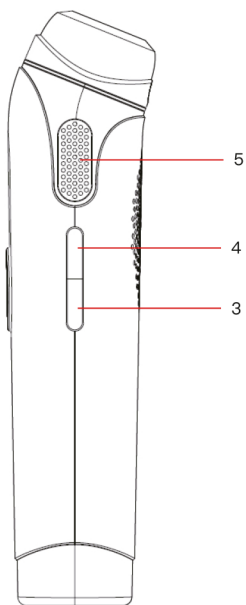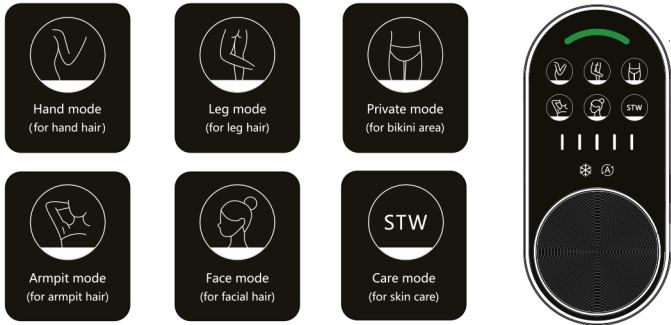(Do not obstruct the air inlet/outlet with anything, otherwise the device will be seriously damaged)






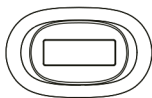
***We have developed 6 treatment modes for different areas of skin and thickness of hair**




***There are 6 different replaceable attachments tailored for different areas**

Note:


- ① Model LX 1000 have 6 attachments: face mode, hand mode, leg mode, private mode, armpit mode and care mode attachment.
- ② Please confirm the machine has exited the automatic light emitting mode, and the symbol  is off before installing the replacement head.
 1. Changing attachment can further pick up the wavelength needed for hair removal, and filter the wavelength unneeded, so as to remove the hair while reducing light irritation
 2. Skin with problems are not recommended to be treated with this product (refer to safety tips on page 1 for detailed information). It is recommended that you place an attachment before using the device on sensitive skin.
 3. To avoid skin injury, do not use care mode on any treatment area within 48 hours after treatment.



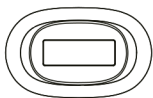
Hand mode attachment
(for hand hair)




Leg mode attachment
(for leg hair)

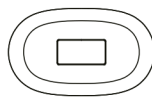

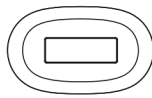
Private mode attachment
(for bikini area)



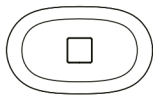Armpit mode attachment
(for armpit hair)



Face mode attachment
(for facial hair)


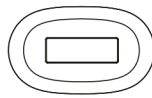


Care mode attachment
(for skin care)

PREPARATION BEFORE USE

1. Clean the skin completely, make sure that the skin is completely dry and free of residual hair and oily matters (residues including perfume, lotion, deodorant, sun cream, etc.)

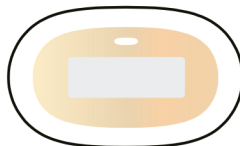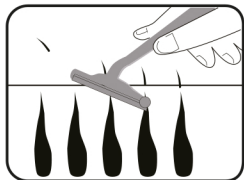2. Make sure that the light exit window is free of residual dust and hair and breakages or abnormalities. If there is dirt, please be sure to clean it to avoid malfunction, deterioration and skin problems.



3. Before use, do not apply any gel or skincare products (aloe vera gel, etc.) on your skin. The product can cool the skin during treatment with the cooling mode, making the treatment process more comfortable and painless.

4. Perform a skin test before use. Check the skin condition after 2 hours. If the treatment area is normal (without redness, blister, etc.), you can use the product without worry.

● LIGHT INTENSITY SETTING














· This product is developed for at-home hair removal, thus its light energy is more tender. The device has 5 different settings of light energy. You may select a proper light energy setting for hair removal based on your skin tone and hair color, which can make the treatment process comfortable. The following recommended setting table can help you in selecting a proper light energy

1. Please refer to the following table to decide which light intensity is suitable for your skin tone and hair color

Notes: ① Place skin tone table on your skin and compare it with your skin tone to confirm your skin tone.

② We recommend you to keep this skin tone table for future reference

2. If you experience any pain or discomfort during treatment, reduce the setting. People with darker skin are prone to skin allergy as the skin will absorb more light energy, so we recommend them to use a lower setting. You may feel warm or hot when the skin is treated by different light intensity, but it should not be painful

Hair color		Skin tone							
		White	Beige	Light brown	Medium brown	Dark brown	Brownish black and darker colors		
White/Grey 	Without attachment	x	x	x	x	x	x	x	
	With attachment	x	x	x	x	x	x	x	
Red 	Without attachment	x	x	x	x	x	x	x	
	With attachment	x	x	x	x	x	x	x	
Pale gold 	Without attachment	x	x	x	x	x	x	x	
	With attachment	x	x	x	x	x	x	x	
Dark gold / light brown 	Without attachment	1/2/3/4/5	1/2/3/4/5	1/2/3/4/5	1/2/3	x	x		
	With attachment	1/2/3/4/5	1/2/3/4/5	1/2/3/4/5	1/2/3/4	x	x		
Brown 	Without attachment	1/2/3/4/5	1/2/3/4/5	1/2/3/4/5	1/2/3	x	x		
	With attachment	1/2/3/4/5	1/2/3/4/5	1/2/3/4/5	1/2/3/4	x	x		
Dark brown 	Without attachment	1/2/3/4/5	1/2/3/4/5	1/2/3/4/5	1/2/3	x	x		
	With attachment	1/2/3/4/5	1/2/3/4/5	1/2/3/4/5	1/2/3/4	x	x		
Black 	Without attachment	1/2/3/4/5	1/2/3/4/5	1/2/3/4/5	1/2/3	x	x		
	With attachment	1/2/3/4/5	1/2/3/4/5	1/2/3/4/5	1/2/3/4	x	x		

Notes: ① Those marked with "x" on the above table are not suitable for the product

② If you want to confirm whether your skin tone is suitable for the product, please consult the skin tone table

③ Your skin may react differently on different days/occasions for a number of reasons

● SKIN TEST

•When you use the product for the first time, or if you haven't used this product for a long time, or if you have accepted tanning recently, perform a skin test on each area to be treated. The skin test is necessary to check your skin's reaction to the light intensity and to determine the correct light intensity setting for each body area

1.Select an area near the area to be treated

(Note: do not try the device on difficult or sensitive areas (corners or bony areas))

2.Press and hold the on/off button for 3 seconds to turn on the device and make sure you select setting 1 and the right mode for the area to be treated

3.Place the light exit window on the skin and make sure it is pressed against the skin

4.When the flash indicator appears on the screen, press the flash button to release a flash

5.Slide the device over the skin to the next area to be treated

6.Increase the setting by one level, apply a flash and slide the device to the next area. Repeat this for all levels within the recommended range of light intensity for your skin

7.Check your skin condition 2 hours after the skin test. If the treatment area is normal (without redness, blister, etc.), you can use the product without worry

● KEYS TO TREATMENT ON EACH AREA

1.Hand and leg treatment

① Before treatment, clean the skin of the hands and legs thoroughly, make sure that the skin is completely dry and free of residual hair and oily matters (residues including perfume, lotion, deodorant, sun cream, etc.)

② You may feel hot on bony and bumpy areas more frequently than on other areas

③ For areas with pigment accumulation and dark color, use the product on a lower setting. Use setting 1 before you are tolerable to it

2.Private treatment

① Before treatment, clean the skin of the bikini line thoroughly, make sure that the skin is completely dry and free of residual hair and oily matters (residues including perfume, lotion, deodorant, sun cream, etc.)

② For areas with loose skin, stretch the skin with your fingers before treatment. It is normal if you feel slightly hot when the skin is under light exit window

③ For areas with pigment accumulation and dark color, use the product on a lower setting. Use setting 1 before you are tolerable to it

3.Ampit treatment

① Before treatment, clean the skin of the ampit thoroughly, make sure that the skin is completely dry and free of residual hair and oily matters (residues including perfume, lotion, deodorant, sun cream, etc.)

② For areas with loose skin, stretch the skin with your fingers before treatment. It is normal if you feel slightly hot when the skin is under light exit window

4.Face treatment

- ① Before treatment, clean the skin of the face thoroughly, make sure that the skin is completely dry and free of residual hair and oily matters (residues including perfume, lotion, deodorant, sun cream, etc.)
- ② Keep it away from areas on the eyes, eyebrows, lips, hair and around the eyes during treatment.
- ③ For areas with pigment accumulation and dark color, use the product on a lower setting. Use setting 1 before you are tolerable to it.

5.Do not use the product on the following areas

Areas with burns, dark color, chilblain, black mole and tattoo, areas around the eyes, and eyebrow, hair, sunburned areas

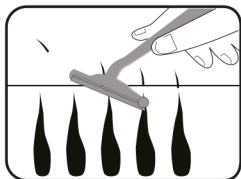
- After treatment, keep away from UV rays 24 hours a day to avoid the skin of the treatment area being sunburned
- Do not look straight at the flash

OPERATION

- If you don't adjust the light intensity setting manually when the device is plugged in and turned on, the device will automatically identify the skin tone and recommend the treatment setting and treatment mode. The setting will be adjusted by the device automatically. You only need to adjust the mode and press the flash button to start treating. You don't need to adjust the setting
- If you adjust the light intensity setting manually when the device is plugged in and turned on, the device will be in the manual mode that requires the manual adjustment of treatment mode. You should adjust the light intensity setting based on your skin tone by consulting the skin tone table, adjust the device to the right mode and press the flash button to start treating

● Steps of Single Flash Release Mode

1. Shave the hair of the area to be treated using a razor, clean the skin and make sure it is completely dry



- 2.Plug the device in, the icons on the screen will be on for 1 second before they go out



3. Press the on/off button for more than 3 seconds to turn on the device. Information will be shown on the display



Information on the screen when the device starts in the mode that automatically adjusts light intensity setting

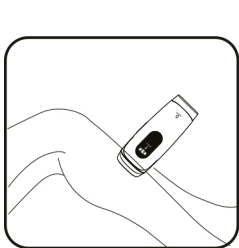


Information on the screen when the device starts in the mode that requires manual adjustment of light intensity setting

4. Compare your skin with the skin tone table. If your skin tone belongs to the part that is marked with x, do not use the device. The skin tone table shows the recommended setting for different skin tones. However, it is suggested that you start from setting 1 when you first use the device
5. Adjust the setting according to your skin tone by pressing the setting adjustment button . You may select among setting 1~5; Press the mode adjustment button  to select the right treatment mode based on the area to be treated. When you first use it, it is suggested that you start from setting 1 to the extent that the skin is not injured. The higher the setting is, the better the treatment result will be
6. It is recommended that you use the corresponding attachment if your skin is sensitive, and select the right modes for different body areas to reduce the irritation of the light to the skin
7. During operation, wear the specialized goggles supplied with the product to avoid the eyes being hurt by the flash



8. Place the light exit window at a 90° angle on the area to be treated. Make sure that the light exit window is pressed against the skin, and the flash indicator on the screen is on before you press the flash button to start treatment (Note: if the light exit window is not pressed against the skin, the flash indicator will be off, and no flash will be released when you press the flash button)




The flash indicator will be "green" when the device is in the mode that automatically adjusts light intensity setting



The flash indicator will be "white" when the device is in the mode that requires manual adjustment of light intensity setting

● STEPS OF AUTOMATIC FLASH MODE

1. When the device is turned on, press and hold the flash button for over 3 seconds, then the symbol  will appear on the screen, indicating that the device is in the automatic flash mode



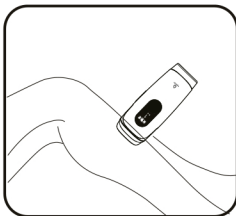
Information on the screen when the device starts in the mode that automatically adjusts light intensity setting



Information on the screen when the device starts in the mode that requires manual adjustment of light intensity setting

2. Press the light exit window against the skin, and slide slowly to the next area. The device will flash automatically without pressing the flash button. Press the flash button again to exit automatic flash mode and enter into manual flash mode



(Note: You may adjust the speed of movement on the skin according to the flash rate to make the coverage more even. The flash rate is higher in setting 1 and slightly lower in setting 2 and so on, the higher the setting is, the lower the flash rate will be)




The flash indicator will be "green" when the device is in the mode that automatically adjusts light intensity setting



The flash indicator will be "white" when the device is in the mode that requires manual adjustment of light intensity setting

- Note: ①Please confirm the machine has exited the automatic light emitting mode and the symbol  is off when testing the skin tone.
- ②Please confirm the machine has exited the automatic light emitting mode and the symbol  is off before installing the replacement head.
- ③Please make sure to wear goggles all the way when it is in service.

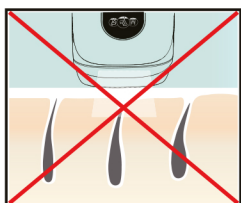
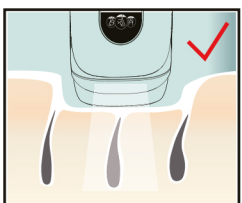
● STEPS OF COOLING CARE MODE

1. After treatment, exit automatic flash mode and place the cooling surface of the device on the treatment area, and slide back and forth to cool and soothe the skin
2. After the skin care, press and hold the on/off button  for 3 seconds to turn off the device

(Without any operation, the device will switch itself off automatically in 5 minutes. After being turned off, the fan will continue to work for 1 minute to cool down the device.)

TIPS FOR USE

1. When you use the device, make sure that the light exit window is pressed against the skin, which allows the full contact between it and the roots of hair follicles, ensuring that the hair follicles absorb higher energy so as to achieve a better result. The cooling surface is in full contact with the skin, which cools it and reduces its heat so as to better protect the skin.



2. Since the continuous flashing of the device will cause the heat to accumulate, in order to make the cooling continuously effective, after finishing treating an area, we recommend that you cool the treatment area for about 30s before treating the next area. It allows you to cool and soothe the skin again and at the same time makes the device's cooling effect better.

● TREATMENT SCHEDULE

Time	Week 1~2	Week 3~4	After week 5
Treatment frequency	3 times a week	Twice a week	Once every month or half a month

(Note: Please treat the skin various times according to the actual growth of the hair. For the hair that is too thick and stiff, increase the setting and frequency of treatment according to the actual conditions)

MAINTENANCE AFTER USE

- It is suggested that you keep away from water for 6 hours after treatment. After that, you can moisturize the treatment area with moisturizer.
- Do not apply perfume, mask and cosmetics containing buffing beads to the skin after treatment.
- Sunbathing may cause your skin to be tanned or sunburned. Treating the skin under such conditions may lead to redness, inflammation, etc. Therefore, avoid exposure to sunlight during the use of this product. The period for sunburn repair should be no less than 1 week
- Never expose the treatment areas to sunlight within 48 hours after treatment. Even if you need to expose the treatment areas to the sunlight after 48 hours, you should pay attention to the skin in case of any undesirable situations. Apply SPE30+ sun cream to the skin after 72 hours to protect it from the sun.
- Artificial tanning during the treatment period is not allowed.
- Consult your doctor if plastic surgery is needed after treatment (only for treatment areas).
- Consult a surgeon if you experience any discomfort.

CLEANING AND MAINTENANCE

In order to ensure the best effect and longest service life of the product, clean the product before and after use to avoid hair and dirt accumulation. If you don't clean it regularly, the product will not be as effective as expected.

- Before cleaning, you must switch it off, unplug it and let it cool
- Scrub and clean the light exit window with a soft, clean and dry cloth
- Please store the equipment in a dry and dust-free environment with a temperature of -10°C - 60°C
- Never clean the product and its accessories using a faucet or dishwasher
- Do not clean the product with steel wool, abrasive cleansers (toothpaste, etc.)
or corrosive liquids (gasoline, machine oil, acetone, etc.)



(Attention: If water does not work well, a few drops of 75% medical alcohol can be used to clean the above accessories)

ENVIRONMENTAL PROTECTION



- Do not discard it together with common household waste. Send it to a designated recycling center, which is a good practice for environmental protection

JUDGEMENT AND SOLUTION ON COMMON PROBLEMS

· This chapter outlines the possible problems when using this product. If you cannot solve the problem based on the following solutions, please contact customer service

- In case of any abnormal situation or failure to troubleshoot, please contact customer service in time
- Please keep your purchase invoice and warranty card. When your product breaks down, you need to take such card.

Problem	Possible Causes	Solutions
Unable to start	The power is not properly connected	Check whether the device is plugged in
Unable to flash	The device is defective	① Check whether the power supply works normally ② Repair or return the product immediately
Smell unpleasant when releasing a flash	① The filter glass is dirty ② Residual hair on body surface	① Wipe the filter glass ② Shave off the hair and reduce the light intensity
The indicator can't work	① Whether the temperature of the device head and body is too high ② Whether the device head is pressed against skin	① Wait a moment, and use the device after it cools down ② Press the device against the skin until the flash indicator is on
The pain is unbearable during treatment	The hair on the part to be removed is not shaved off	Before using the product, please shave off the hair on the area to be treated. If the shaving causes skin allergy, please cut the hair as much as possible and use the device with a comfortable setting
	The product is not applicable to the area to be treated	Do not use the product on inner labia, vagina, anus, papilla and areola. Men should not use it on face or scrotum
	The setting is too high	Lower the setting to a comfortable one
	The light exit window is damaged	Repair or return the product immediately
The treatment result is not satisfactory	The setting you used is too low for you	Please select a higher setting
	You did not release a flash on the area adjacent to the previously treated area	Please release a flash on the adjacent areas and ensure there is some overlap between the areas treated.
	This device is not applicable to your hair color	If your hair is light gold, grey, red or white, the treatment will be ineffective
	The frequency of using this device is insufficient	To successfully remove all hairs, we suggest that you follow the recommended treatment schedule
	Your skin's response to light exposure is slow	You should use the product continuously for at least 6 months, because it can inhibit hair regrowth during this period
After treatment the treatment area turns red	Slight redness is normal, harmless to the body, and will disappear soon	Use aloe essence or the cooling surface to soothe skin
The product becomes hot and damp during use	It may be due to long operation of the device	Turn off the device, and use it after sufficient heat dissipation

STATEMENT FOR HAZARDOUS SUBSTANCES

- As provided by Regulation on the Recovery and Disposal of Waste Electrical and Electronic Products, (hereinafter referred to as Regulation), names and contents of possible hazardous substances in the electrical products of the company are listed as follows:

Name of Part	Hazardous substance					
	Pb	Hg	Cd	Cr6+	PBB	PBDE
use frequency						
Outer Shell	O	O	O	O	O	O
Circuit Board Assembly	X	O	O	O	O	O
Power Line	X	O	O	O	O	O
Adapter	X	O	O	O	O	O

The table is prepared according to the regulations of SJ/T 11364–2014.

O: Indicates that the content of the hazardous substance contained in all of the homogeneous materials used for this part is below the maximum of GB/T 26572.

X: Indicates that the content of the hazardous substance contained in at least one of the homogeneous materials used for this part is above the maximum of GB/T 26572.



The mark "Environment-friendly Use Period" of the product, as shown in the left picture, meets the requirements of China. Restriction of Hazardous Substances (China RoHS) released on July 1, 2016 and China Household Electrical Appliances. Association's Guidance on the Environment-friendly Use Period of Household Electrical Appliances. According to the requirements of China RoHS, the product manufacturer

and importer shall disclose the hazardous substances in the product based on the requirements of and in the form stipulated by relevant applicable standards and mark the environment-friendly use period of the product.

If you have followed the notes on safety and use of the product recorded in the Instruction Manual, and there are no general defenses provided in other laws and regulations, the toxic and hazardous substances or elements in the product will not leak or mutate within the above prescribed years from the date of manufacture. Using the product will not cause severe pollution to the environment or serious damage to the user and his property.

The "Environment-friendly Use Period" is not a period for safe use, particularly differing from the period limited by some factors, including electrical performance safety and electromagnetism safety.

If you want to discard the product after proper usage, please get rid of it as prescribed by the law on recovery and recycling of electronic information product.

Note: The period is "Environment-friendly Use Period", not as quality guarantee period.

Special Instruction

This product is developed according to the skin tones of Asians.

We are solely liable for the revision and explanation of the Instruction Manual, and ensures it's correctness to the greatest extent. The information has been carefully checked. In case of any printing error or omission, the company reserves the right of final interpretation to the extent permitted by law. Meanwhile, the company reserves the right to improve and revise the product and the Instruction Manual. We won't further notify any change in the product structure. Some pictures in the Instruction Manual are schematic diagrams and are for reference only. If the picture does not match the actual product, the latter shall prevail.