bloomfold

Bloomfold Electric Heating Pad
Model Num: BF210

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EN

6. Warnings
7. Product details & instructions
8. Cleaning & storage

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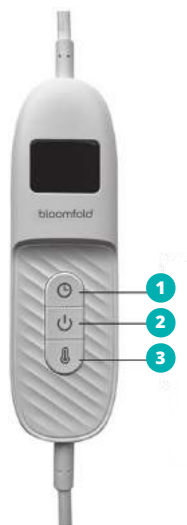
Warnings

Safety precautions before use

1. Read the instructions before use.
2. This product is intended for the intended use only.
3. Always observe precautions when using the electric blanket hereinafter referred to as product.
4. Make sure that packaging material does not get into the hands of children. They may suffocate in it.
5. People who are sensitive or hypersensitive to heat are not advised to use an electric heating pad.
6. Not suitable for children under 3 years of age.
7. Children between the ages of 3 and 8 to be used only under supervision.
8. This product is intended for domestic use.
9. Do not use the product folded or creased.
10. Use the product dry only.
11. Using an electric heating pad may interfere with the operation of a pacemaker. Although with 5,000 V/m, 80 A/m and magnetic induction of 0.1 mT are far below the values, it is advisable to contact your doctor and pacemaker manufacturer.
12. Do not place remote control on the product during use.
13. If the product or cord is damaged or not working properly, do not use it and contact customer service.
14. Remove plug only when product is off.
15. If product is damaged do not use product again, contact customer service to see if repair is possible.
16. Allow product to cool before storing.
17. Keep product out of reach of children.
18. Caution: improper use may result in injury.
19. Do not place objects on the product.

Product details

- 1 Timer
- 2 On/off
- 3 Temperature

Timer setting	Time in minutes
1	20
2	40
3	60
4	80
5	100
6	120
7	140
8	160
9	180

Instructions before use

1. Turn on the product by pressing the power button.
2. Choose one of the 9 heat settings. By default, the product starts at setting 1. After 9, setting 1 starts again.
3. Choose a desired time by clicking on the timer, choose from one of the 9 possible timer options. By default, the product starts at 20 minutes.
4. Stages 2 and 3 can be used to allow the product to heat up faster.
5. If using during sleep, it is advised not to go higher than setting 1.
6. Press the on/off button again to turn the product off manually.

Cleaning & storage

1. Turn off the product and remove the plug from the socket.
2. Disconnect the remote control from the product.
3. Let the product cool down completely before cleaning.
4. Small stains can be washed with a slightly damp cloth and a small amount of detergent.
5. If the product is very dirty, it can be washed by hand at a temperature of up to 30 degrees Celsius.
6. Hand wash the product a maximum of twice a year to prevent damage.
7. After washing let the product dry completely out of the sun, let it dry naturally. Do not place it near a radiator.
8. Never dry the product through use.

Storage

1. If the product will not be used for a longer period of time, it is advised to store the product in its original packaging in a cool dry place.
2. Do not place objects on the product, do not fold the product in half.

Specifications

Measurements	40 x 30 cm	Max. temperature	43 °C
Model number	BF210	Max. timer	180 minutes
Material	Fleece	Voltage	220 - 240V
Cord length	230 cm	Power	100W

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