105mm

148mm

gatson
WEIGHING SCALE
INSTRUCTION MANUAL

Thanks for choosing our product, please read the instruction manual carefully, and keep it for future reference.

1.Features


High precision sensors	Hold function (manual hold & auto hold mode)
High HD LCD display	Tare function
Auto power off	Capacity: 30kg, division: 2g
Overload indication	Display: LCD screen
Low battery indication	Screen size: 79.6x16mm
Unit conversion: kg,g,lb,st:lb	Power: 3x1.5V AAA batteries

2.Operation

Switch button 	1.Press this button to turn on scale 2.Long press this button to turn off scale
Unit button	Press this button to convert units between g/kg/lb/st: lb
Hold button	Press this button to convert modes between <auto hold> and <manual hold >
Tare button	Press this button to clear the weight, back to 0.0kg


2.1 Step for use:

- Step 1. Put the scale on the hard & flat surface.Do not use it on carpet or any soft surface as it may affect accuracy.
- Step 2. Place the tray on the scale.
- Step 3. Press the button to turn on scale, display "HI" about 2 seconds then back to zero, now the scale is ready for weighing.
- Step 4. Put the baby on the scale, the LCD will show the weight of baby.

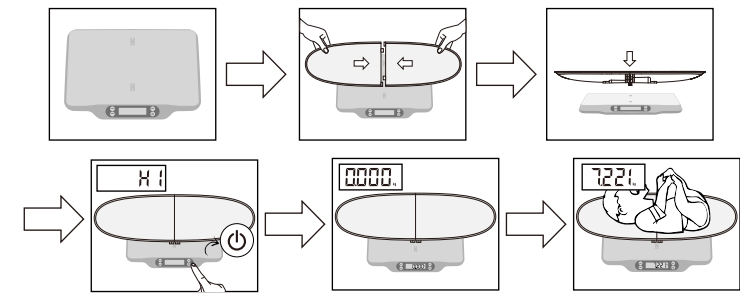
2.2 Reset and Tare Function

A. Press the "TARE" button if you need to put a towel or cloth on the tray before weighing the baby. This will reset the display to "0.00".

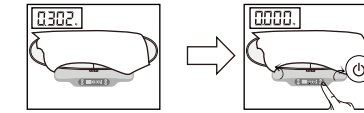B. Press the "TARE" button if the data is locked and you need to weigh another baby. This will reset the display to "0.00", allowing you to weigh the new baby without interference from the previous data. This function can be repeated as needed.

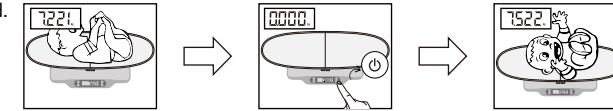- C.Removing the baby on the tray, the LCD will display the minus value; press the <tare> button again, after 2 seconds the LCD will show the weight in zero.
- D.The scale will turn off automatically after approximately 60 seconds while in idle mode.

2.3 Hold mode:

The default mode is auto-hold. Press "Hold" button for about 2 seconds to convert to manual-hold mode. When weighing, press "Hold" button to maintain the weight on the screen.
Scale will display "Hold" only on Auto-Hold Mode

2.4 Unit conversion (g/kg/lb/st:lb)

Under the weight mode, press <Unit> button to convert the unit between "g" or "kg" etc.

2.5 Overload Indication & Low power Indication

When the weight is over the "max capacity", LCD screen will show " 0-Ld", please remove the object immediately, or will damage the device.
When the LCD screen displays "LO", it indicates that the battery is low,please charge the scale.

3.Safety precautions:

- Ensure to place the scale on flat,dry and hard surface when weighing.
- Use this device for its intended purpose at home.
- Do not load object that is over the scale's maximum capacity of 30kg
- Avoid shaking or hitting the scale.
- Do not submerge it into water.
- Do not place anything on the scale when it is not in use
- The optimal temperature for operating and storing the scale is from 5 to 30 degree celsius.

This Scale is meant for domestic use, Not to be used for any transaction or protection.

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