
TORQ™ mini

Percussion Massager

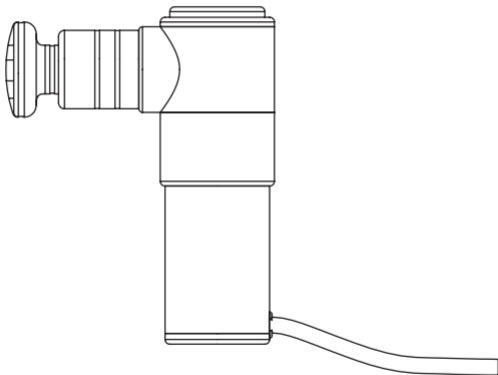mytorq.co

USER MANUAL

TORQ™ mini

USER MANUAL

Thank you for choosing the TORQ mini Percussion Massager.

Before using it, please read this manual carefully to ensure proper operation and keep this manual for future reference.

Contents

Main information.....	2
Product introduction.....	4
Instructions for use.....	6
Product specifications.....	7
Maintenance.....	8
Warranty.....	8
Contact.....	8

PRECAUTIONS

To avoid harm, fault, electric shock, fire, death and equipment damage, please read the instructions and warnings carefully.

CAUTION

1. Use only as instructed. Use TORQ mini as described in the manual. Use only TORQ recommended accessories.
2. The TORQ mini is not intended to be used by small children unless they have been given supervision by a responsible adult. TORQ mini is not a toy.
3. Do not submerge TORQ mini in water.
4. Do not bring TORQ mini into close contact with fire.
5. Do not use TORQ mini on the head, neck, heart or genitals.
6. Use TORQ mini on soft tissue only. When feeling pain or discomfort, stop use immediately.
7. Only use the provided charger to charge your TORQ mini.
8. Do not drop TORQ mini.
9. Do not remove screws or try to disassemble the product.
10. TORQ mini is not an alternative or complement to medical treatment. If you are feeling pain, undergoing medical treatment or using any other medical devices, please consult your doctor before using this product.

11. Keep your fingers, hair or other body parts away from the back of the massage heads to prevent injury.

Do not touch this part when the machine is running without the massage heads being installed.



Do not touch this part when the machine is running with the massage heads installed.

DISCLAIMER

If you have any of the following situations, or any health-related problems, please consult your doctor before using this product.

1. Pregnancy, diabetes, mental disorder, retinal damage and/or other complications.
2. Pacemaker usage, recent surgery, epilepsy, migraine, lumbar disc protrusion, and spondylolisthesis.
3. Recent joint replacement surgery or placement of intrauterine contraceptive devices (IUD), metal implants or other concerns about your own body.

1. PRODUCT INTRODUCTION

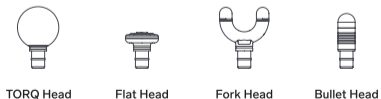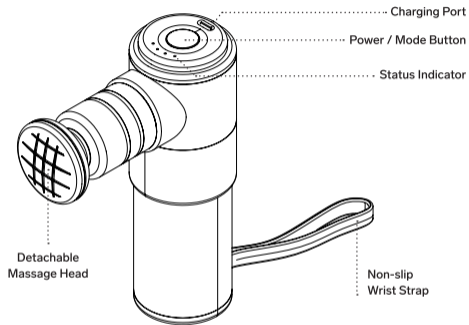Massage Heads

TORQ Head

Your go-to head for your every day massage. Made from soft yet firm foam, the TORQ Head is designed to be used comfortably on all areas of the body. Ideal for breaking down lactic acid from your muscles and perfect for that after work massage.

Intensity Level: 3/10

Flat Head

When you want that deeper tissue massage without any discomfort, TORQ's Flat Head is the accessory for the job. Use anywhere on your body to relax your muscles and tighten your skin.

Intensity Level: 5/10

Fork Head

Designed to pinpoint sore areas, TORQ's Fork Head alleviates areas of tension in no time. Covered in soft, non-porous silicone for that ultimate comfort feeling. Ideal for use on your forearms and feet.

Intensity Level: 8/10

*Do not use longer than 15 seconds in one spot

Bullet Head


When precision matters, the TORQ Bullet Head delivers pinpoint accuracy to those hard to reach spots. Covered in soft, non-porous silicone for that ultimate comfort feeling. Perfect for lower back pain relief and anywhere you would have massaged with your thumb.

Intensity Level: 10/10

*Do not use longer than 15 seconds in one spot

2. INSTRUCTIONS FOR USE

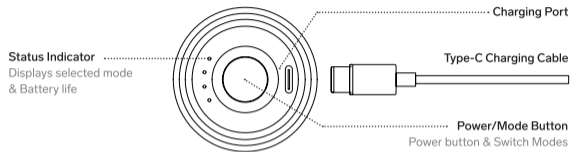
1. Press and hold the Power/Mode button for 2 seconds to turn TORQ mini on. The status indicator lights display the current battery life.
2. Short press the Power/Mode button to switch modes. The status indicator light shows the current mode. Each press will increase the percussion intensity.
3. Press and hold the Power/Mode button for 3 second to turn TORQ mini off.

Charging

1. Fully charge the battery before first use.
2. Turn TORQ mini off when charging.
3. Connect the power supply. Status indicators will flash alternatively.
4. When the four indicator lights are on, the battery is fully charged.

Replacing the Massage Heads

- With TORQ mini turned off, slowly insert the massage head into the barrel.
- A slight twisting movement is recommended.
- Make sure to insert the head in all the way.
- Always insert the massage head perpendicular to the device.
- Do not use excessive force to insert or remove the massage heads.

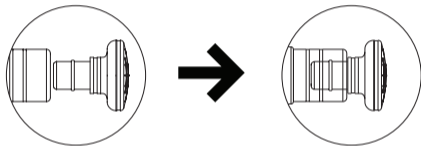3. PRODUCT SPECIFICATIONS

Rated speed: 900-3200 PPM	Stroke: 6 mm
Rated input: 100-240 V~50/60 Hz	Maximum power: 70 W
Rated output: Type-C, 5V 2A	Net weight: 0.5 kg
Battery capacity: 2500 mAh	Size: 142*90*44 mm

4. MAINTENANCE

1. When not using TORQ mini for a prolonged time, store it in a safe, cool and dry place.
2. Charge your TORQ mini at least every 6 months in order to prevent the battery from deteriorating.
3. Only use neutral detergent and a soft cloth for daily maintenance. Do not use engine oil, solvents or other chemicals for cleaning.
4. Turn off the device after use and before cleaning the product.

5. WARRANTY

Visit our website www.mytorq.co for information about warranty.

6. CONTACT

If you have any further questions or concerns and wish to find out more information about TORQ products, please contact us at care@mytorq.co

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