

HANASCO

Quick Guide & User Manual
www.hanasco.com

Portable Electric Water Flosser

Important Notice

⚠ Statement

It is the user's responsibility to read the user manual carefully before using the product. Please read the user manual carefully before using the product.

⚠ Warning

- Do not immerse the product in water when it is not in use.
- Do not disassemble or change the parts of the product to avoid damage or injury.
- Do not get the power adapter wet. Use the power cord and plug to connect to the power outlet.
- Do not use the water when charging, please keep the power plug from getting wet.
- Do not use the product if you have a damaged power cord or if the cord is frayed.
- Do not get the product into the water.

⚠ Precaution

- Please do not use the product for other purposes than oral cleaning.
- Do not use the product for the teeth or gums of children and infants.
- Please do not immerse the product into the water when it is not in use.
- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.


⚠ Important

- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.

⚠ Attention

- Do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.
- Please do not use the product for other purposes than oral cleaning.

Product description

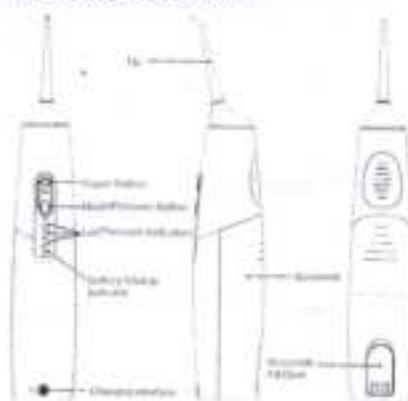

Accessories



Using steps

Please confirm whether the electricity is sufficient before using it. The new machine needs to be recharged before it can be used. When in use or after use, the charging indicator lights up the red light to indicate that the battery power is insufficient and needs to be recharged.

1. Nozzle installation: Insert nozzle with reference to sketch, turn the nozzle 45 degrees in the direction of the arrow and lock it in place as shown in the figure.



2. Open the reservoir cover and hold the flosser horizontally to fill the water tank. Cover the tank after filling the appropriate amount of water.



3. By pressing the mode button to select the desired mode, the product has three modes. The specific functions are as follows:



Standard mode: The first gear indicator is on, and the standard water pressure is efficient and clean for most people.



Soft mode: The second gear indicator is on, and the gentle water pressure is suitable for the first use of water flosser or sensitive mouth.



Pulse mode: The first gear indicator light is on, Pulse water can deeply clean and massage gums.

4. When in use, place the nozzle deep into the mouth, aim at the teeth, gums and gingiva, and ensure that the water flow is perpendicular to the eruption site, while the mouth is slightly open.

* Tip: If you notice itching or slight bleeding when first used, because the gums are not used to the efficient cleaning. The phenomenon usually stops within 1-2 weeks. If it is more than two weeks, it may indicate that your gums are abnormal. Please stop using and consult a dentist.

5. Press the switch button to start the experience. This product has a two-minute timed shutdown function. If you need to continue using it, please press the key.



Product Parameters

Product model	H100
Working Modes	Normal/Soft/Pulse
Pulse	1400-1800 times/minute
Reservoir capacity	230 ml
Nozzle Numbers	4 pcs
Water Pressure	30-110 PSI
Life time	About 80 minutes
Charging time	4 hours
Adapter	DC 5V 1A
Battery parameter	Li-ion battery 2000 mAh
Electricity reminding	Low electricity indicator
Waterproof level	IPX7
Standard	GB4706.58-2008

Battery indicator

Start-up power indicator	
100% - 90%	green light
90% - 20%	blue light
<20%	red light
Charging power indicator	
100%	green light
20% - 100%	blue light
<20%	red light

Environmental protection

The rechargeable lithium battery inside the product may contain harmful substances that pollute the environment. Discard the product or put the product into the designated recycle bin. Be sure to take up the battery and send the battery to the designated recycling.

Cleaning and Maintenance

Clean the product with clean water or neutral detergent instead of corrosive or abrasive detergent. Do not use hot water over 50 degrees to avoid damage to the product.

- Body Cleaning:** Dry cloth or wet cloth can be used to clean the body. Do not immerse the body in water to clean.
- Cleaning of water tank:** Can be wiped clean with a dry cloth or a damp cloth. If you don't plan to use it for a long time, please spray the water inside.
- Nozzle cleaning:** can be wiped with dry cloth or wet cloth, or cleaned with clean water rinse and turn water with cloth, do not touch, pull out from the nozzle edge.
- Cleaning of charging plug:** Please wipe the metal part of plug regularly with cloth to avoid heat caused by dust accumulation.

* Dentist recommend replace the nozzle every six months.

Troubleshooting

Troubleshooting	Reason of issue	Solution
Not working	Not enough charge	Charge
Not water pressure	Blockage of nozzle or too little water	Replace nozzle or clean
Charging indicator abnormal or not on	Block plug power plug not connected properly	Return to basic Correctly connect power plug
	Insufficient charge	Charge