
# AQUACARE WATERFLOSSER

powered by **BRAUN**

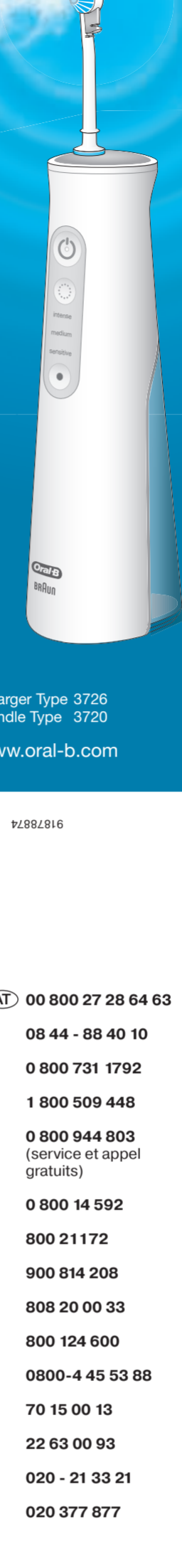Charger Type 3726

Handle Type 3720

[www.oral-b.com](http://www.oral-b.com)

91878874

**(DE) (AT) 00 800 27 28 64 63**

**(CH) 08 44 - 88 40 10**

**(UK) 0 800 731 1792**

**(IE) 1 800 509 448**

**(FR) 0 800 944 803**

(service et appel gratuits)

**(BE) 0 800 14 592**

**(LU) 800 21172**

**(ES) 900 814 208**

**(PT) 808 20 00 33**

**(IT) 800 124 600**

**(NL) 0800-4 45 53 88**

**(DK) 70 15 00 13**

**(NO) 22 63 00 93**

**(SE) 020 - 21 33 21**

**(FI) 020 377 877**

#### **Internet:**

[www.oralb.com](http://www.oralb.com)

[www.braun.com](http://www.braun.com)

[www.service.braun.com](http://www.service.braun.com)

[www.az-oralb.it](http://www.az-oralb.it)

Braun GmbH

Frankfurter Straße 145













61476 Kronberg/Germany

**Charger Type 3726**

**Handle Type 3720**

91878874/X-18

DE/EN/FR/ES/PT/IT/NL/DK/NO/SE/FI
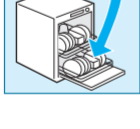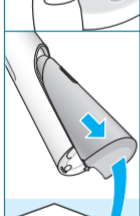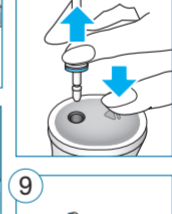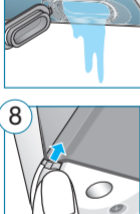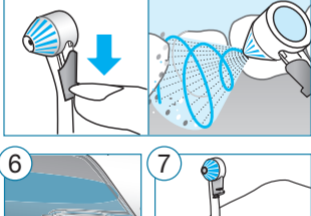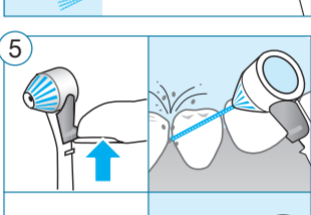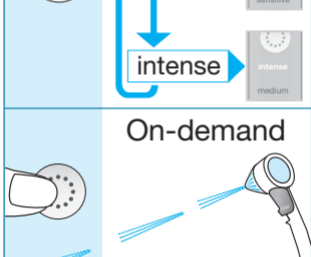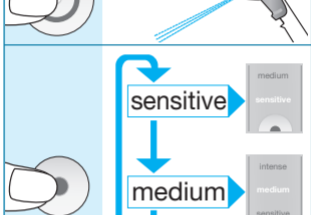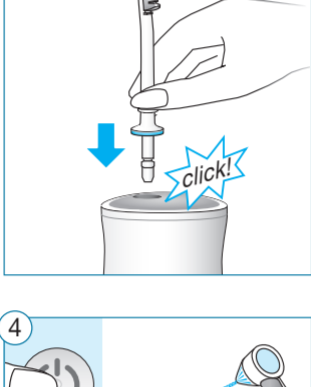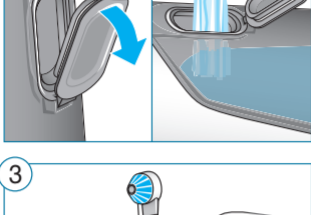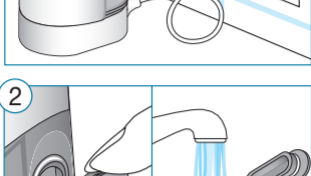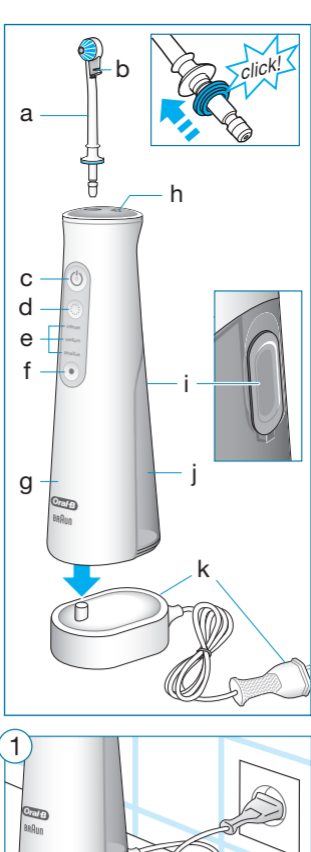
## Deutsch

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## English

Welcome to Oral-B!

Please read the usage instructions carefully before first use and save this manual for future reference.

### WARNING

- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- This appliance contains batteries that are non-replaceable. Do not disassemble the product except when disposing of the battery. When taking out the battery from the handle for disposal, use caution not to short the positive (+) and negative (-) terminals.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This appliance is a personal care device and is not intended for use on multiple persons.
- Small parts may come off, keep out of reach of children.

To avoid irrigator nozzle breakage which may generate a small parts choking hazard or damage teeth:

- If the product is dropped, the irrigator nozzle should be replaced before the next use even if no damage is visible.
- Clean the appliance properly after each use (see «Cleaning recommendations» paragraph). Proper cleaning ensures the safe use and functional life of the appliance.

### IMPORTANT

- Periodically check the entire product/cord for damage. If the product/cord is damaged, take it to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- The oral irrigator is not intended for use by children under age 8. Children aged from 8 years and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge can use the product, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be performed by children.
- Oral irrigators shall be used only with liquids defined in the manufacturer's instructions.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in the appliance. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- Only use the charger provided with your appliance. Do not charge any other device with this charger, as this may cause damage.

### Description

- a Nozzle
- b Switch
- c On/off button (continuous jet)
- d On-demand button

e Intensity level

f Intensity selection button

g Handle

h Release button

i Water reservoir cap

j Water reservoir (150 ml)

k Charger

## Specifications

For voltage specifications please refer to the bottom of the charging unit.

Noise level: 68 dB (A)

## Operating and charging

Plug the charger (k) into an electrical outlet and place the handle (g) on the charger (picture 1).

While charging the on/off button (c) light is fading in and out.

Once it is fully charged the light turns off. A full charge typically takes 14 hours.

If the battery is running low the on/off button is flashing in red.

## First Use

Attach the nozzle (a) to the handle until it locks into place indicated by a click sound (picture 3).

Fill the water reservoir (j) with some fresh and clean water by

lifting the reservoir cap (i) on the

back of the handle (picture 2). Hold the handle over a washbasin,

press the on/off button (c) and

allow the entire water volume to

run through the nozzle. When the

device has not been used for

more than one week repeat this

step to ensure flushing.

## Jet functions:

You can choose between different

jet stream types by sliding the jet

switch (b) down or up (before

doing so, turn off the water

stream with the on/off button (c)

on the handle (g)):

Rotating stream type for general cleaning and gum massage (jet switch in lower position): The

rotation of the water stream with

its micro-bubbles are designed

for easy penetration and gentle

gum massage (picture 5).

Focused stream type (jet switch

in upper position): The focused

water stream with its micro-bub-

bles are designed for targeted

cleaning (picture 5).

Press the on/off button (c) for a

continuous jet stream.

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to trigger a single jet stream. Up-

on releasing the on-demand but-

ton the single jet stream stops.

Press the intensity selection but-

ton (f) successively to adjust the

intensity level (sensitive – regu-

lar/medium – intense) (depend-

ing on model). You may pre-

select the intensity level before

you turn on the handle. The in-

tensity level selected last is

memorized (picture 4).

Note: In the beginning, we rec-

ommend to use the sensitive lev-

el for your teeth and gums to get

accustomed to the usage.

## Using the irrigator

Fill the water reservoir (j) with

lukewarm water by lifting the

reservoir cap (i) on the back of

the handle (g).

Lean over the sink and direct the

nozzle in your mouth. Close your

lips slightly and turn on the

handle by pressing the on/off

button (c).

Guide the jet along the gumline

and pause briefly between teeth.

Alternatively, you may also point

the nozzle tip at a desired area

and use the on-demand button

(d) for a single jet stream.

Keep your mouth slightly open to

allow the water to run into the sink.

Refill water reservoir as many

times as needed to complete

your cleaning routine. It may take

multiple refills to clean your entire

mouth.

## After use

Make sure the water reservoir is

completely empty. Press the on/

off button (c) and leave it turned

on until no more water is running

through the nozzle, to prevent

any deposit build-up or bacteria

from collecting in the standing

water. Turn off the handle by

pressing the on/off button again

once the water is empty. A longer

operation of the handle without

water may damage the motor.

To drain any water out of the han-

dle hold it over the washbasin and

tilt it from side to side (picture 6).

Release the nozzle (a) by press-

ing the release button (h) and

pulling the nozzle at the same

time (picture 7).

Lift up the cap (i) and let the

water reservoir air-dry. Always

wipe the handle and nozzle dry.

Note: Make sure you always

empty the water reservoir before

storing or travelling with the

irrigator.

In general we recommend to re-

place the nozzle every 3 months.

## Cleaning recommendations

Clean and dry handle after each

use and leave the water reservoir

cap (i) open.

The water reservoir (j) is dish-

washer-safe and can be de-

tached via the snap-fit on the

bottom of the handle. Move the

lever and slide the reservoir down

(picture 8). Regularly clean the

charger and handle with a damp

cloth. Store the appliance in a

cool and dry place.

Subject to change without notice.

## Environmental Notice

Product contains batteries

and/or recyclable electric

waste. For environment

protection do not dispose

in household waste, but for

recycling take to electric waste

collection points provided in

your country.

## Guarantee

We grant 2 years guarantee on

the product commencing on the

date of purchase. Within the

guarantee period we will elimi-

nate, free of charge, any defects

in the appliance resulting from

faults in materials or workman-

ship, either by repairing or re-

placing the complete appliance

as we may choose.

This guarantee extends to every

country where this appliance is

supplied by Braun or its appoint-

ed distributor. This guarantee

does not cover: damage due to

improper use, normal wear or

use, especially regarding the

brush heads/nozzles as well as

defects that have a negligible ef-

fect on the value or operation of

the appliance.

The guarantee becomes void if

repairs are undertaken by un-

authorised persons and if original

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