LANAFORM

WELLNESS

2-in-1 Shiatsu Comfort HEATING AND SHIATSU MASSAGER

SHIATSU MASSAGER
WITH 4 MASSAGE HEADS

HEATING FUNCTION
WITH 2 INTENSITIES

ZONES: FEET, SHOULDERS,
NECK, BACK & LEGS

BREATHABLE AND
COMFORTABLE MESH LINING

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VERWARMINGS- EN
SHIATSUMASSEGEKUSSEN

HEIZ- UND
SHIATSU-MASSAGE

CALEFACTOR Y
MASAJEADOR SHIATSU

MASUJACA I
OGRZEWAJACA SHIATSU

ENGLISH

CONTENTS OF THE PACKAGING

- 2-in-1 Shiatsu Comfort massaging and heating pad
- 1 power cable
- 1 instruction handbook

PRECAUTIONARY STATEMENT

▲ Please read all the instructions before using your product, particularly these basic safety instructions.

📷 The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.

▲ Please note

This device is not intended for use by people, including children, with reduced physical, sensory or mental capabilities, or people with a lack of experience or knowledge, unless someone responsible for their safety can supervise or instruct them beforehand on how to use the device. Children should be supervised to make sure they do not play with the device.

Non-compliance with the following instructions may cause physical injury or damage to property (electric shock, skin burns, fire). The aim of the following safety and hazard instructions is not only to protect your health and that of others, but also to protect the product.

Do not use the product on young children or people who do not feel the heat (for example, diabetics, people with skin symptoms caused by an illness or areas of skin with scars over the area of application or who have consumed medicines or alcohol).

This product must not be used by those who are not heat sensitive or by vulnerable people unable to react when it overheats. This product must not be used by very young children (between 0 and 3 years old) because they might not react if it overheats.

This product must not be used by young children (between 3 and 8 years old) unless a parent or guardian has set the remote control and the child has been given enough information to use the device safely.

Use this device only for the purposes for which it was designed and according to the instructions provided for this type of use. It may be dangerous if used inappropriately. The manufacturer cannot be held responsible for damage caused by inappropriate or improper use.

▲ Please note

Read these notes carefully! Non-compliance with the instructions below is likely to cause personal injury or damage to property.

Keep the packaging out of the reach of children. They could suffocate.

Do not use accessories that are not recommended by LANAFORM or that are not supplied with this equipment.

▲ Electric shock

As with any electrical device, this device must be used carefully and with caution to avoid risks due to electric shocks.

Conditions of use for the device:

- use it only at the voltage shown on the device (the rating plate is on the back of the device),
- do not use the device if it is damaged. Do not use this device if the power plug is damaged, if the device is not working properly, has been dropped on the floor or if it is damaged or has been dropped into water. In this event, have the device examined and repaired by the supplier or its after-sales department.

If the device is faulty or breaks down when in use, switch it off immediately and unplug it. Do not pull on the power cable or on the device in order to remove the plug from the wall socket. Do not hold and never carry the device by its power cable. Keep cables away from hot surfaces.

Never pull the mains plug out with damp or wet hands.

Do not expose your device to sunlight or to any hot surface. This will prevent damage to the insulation of the electrical circuit.

The remote control must remain within reach; do not leave the switch hidden under a pillow or any other object.

Do not use this device under a cover or under a cushion. Excessive heat can cause fire, electrocution or injury.

Keep the device away from sources of heat to prevent the casing melting and causing a fire.

Do not use the device close to highly inflammable gas or explosive mixtures.

The device must be switched off and unplugged after each use and before it is cleaned.

Do not immerse the device in water.

Do not attempt to retrieve a device that has fallen into water. Disconnect the device immediately.

Do not use the device if it or its accessories show any signs of damage.

▲ Health risk

Do not use the device on parts of the body that are inflamed, injured or swollen. If in doubt, consult a doctor before using it.

Excessively prolonged use may cause burns.

Do not fall asleep while the device is running.

Electrical and magnetic fields emitted by this device may interfere with the operation of pacemakers, even if the power of the electrical and magnetic fields is clearly below the authorised limits. Consult your doctor and the manufacturer of your pacemaker before using this device.

Please consult your doctor before using it if:

- you are pregnant,
- you have a pacemaker, prosthetic limbs or electronic implants,
- you are suffering from one or more of the following illnesses: circulation problems, varicose veins, open wounds, contusions, scratches, inflammation of the veins.

This product is a non-professional heating device. Do not use this device as a substitute for medical treatment. It is intended only for family use and in a home environment.

If you experience any pain when using this device, stop using it immediately and consult your doctor.

▲ Repairs

This device does not contain any part that can be repaired by the user; never try to repair the device yourself.

▲ Use

Do not power up the device when it is folded or compressed. Do not jam the device between objects and do not bend it too sharply.

Only use the device in a dry environment.

Do not pull, twist or fold the cables by bending them sharply.

Do not carry this device by its power cable or remote control or use the cable as a handle.

Do not insert pins, needles or other pointed objects into this device.

This device is not a toy. Keep out of the reach of children.

Use only the power supply included.

This device is not intended to warm animals up.

Do not use this device with footwear.

Please check this device frequently for any signs of wear and tear. If you detect signs of

wear and tear or damage, or if this device has not been used correctly, it must be returned to the manufacturer or reseller before it is used again.

If the power cable is damaged, it must be replaced by the manufacturer; contact its after-sales department to have it repaired or replaced.

Always take note of the notes on its use, cleaning and maintenance.

Do not place any item on the device.

If you have any questions about using our device, please contact customer service.

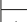
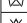
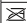
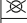
Before using the device, remember to remove the silica gel desiccant.

USE

The 2-in-1 Shiatsu Comfort massaging and heating pad is designed solely for massaging and/or reheating the feet, legs, back, arms or neck of the human body. The device is designed for domestic use only.

Important information

The following symbols are used in these instructions for use, on the packaging and on the data plate of the device and its accessories:

	Read the instructions
	Do not insert needles
	Only use the device in a dry environment
IPx1	Weather proof
	Warning instruction showing a risk of injury or health risks
	Safety note showing a risk of damaging the unit/accessory
	Do not wash
	Do not use bleach
	Do not tumble dry
	Do not iron
	Do not clean
	This product must not be used by very young children (between 0 and 3 years old)

INTRODUCTION

Thank you for buying the 2-in-1 Shiatsu Comfort massaging and heating pad by LANAFORM. The 2-in-1 Shiatsu Comfort massaging and heating pad provides you with a circular kneading massage for releasing tense muscles. This pad's massage nodes move inwards or outwards for a deeper and more authentic massage experience. The versatile design of the massaging and shiatsu pad allows you to target your neck, back, legs, shoulders and feet.

The Shiatsu massage is based on ideas embedded in traditional Chinese medicine, specifically the meridian system, i.e. the channels of the human body in which energy circulates.

The aim of a Shiatsu massage is to promote an individual's physical, emotional and mental wellness by working on energy blockages in the meridians and by stimulating the body's self-regulating capacity with pressure movements exerted along the meridians using fluent movements.

In addition to the Shiatsu massage with two intensities and two rotation speeds, the 2in1 Shiatsu Comfort shiatsu pad has a heating function with two intensities. The two functions can be used at the same time or separately.

DESCRIPTION OF THE DEVICE 1

- 1 Detachable foot warmer
- 2 Detachable shiatsu pad
- 3 Zip fastening
- 4 Remote control
- 5 Power cable
- 6 Adapter

REMOTE CONTROL 2

- 1 On /Off
- 2 Shiatsu massage/massage intensity
- 3 Shiatsu massage rotation direction
- 4 Heating function/heating intensity

AREA OF USE

- Open the cardboard packaging.
- Remove the device upwards, leaving the protective film behind.
- Now remove all the protective film.
- Check that the device, remote control, power cable and mains adapter are not damaged.

INSTRUCTIONS FOR USING THE 2-IN-1 SHIATSU COMFORT MASSAGING AND HEATING PAD FOR THE FEET

- Installing the 2-in-1 Shiatsu Comfort massaging and heating pad on a hard, flat surface.
- Make sure that the shiatsu pad and foot warmer are properly connected using the zip fastening.
- Connect the power cable to the mains adapter and connect this to a power socket, checking that the voltage of the adaptor corresponds to the voltage of your home network.
- Put your feet inside the massaging foot warmer.
- Press the On/Off button to switch the device on. A red LED light is displayed above the corresponding button.

Shiatsu Massage Function

- To activate the Shiatsu massage function, press the "Shiatsu Massage" button. The massage heads will start to turn in an an-

ti-clockwise direction and a green LED light will be displayed.

- The 2-in-1 Shiatsu Comfort massaging and heating pad has 2 Shiatsu massage intensities. To select one of them, keep pressing the "Shiatsu massage" button. A light green LED light corresponds to the lowest speed whereas the dark green light corresponds to the highest intensity.
- To change the direction of the Shiatsu massage heads, press the "Shiatsu Massage rotation direction" button.
- The Shiatsu massage function stops automatically after 15 minutes. If you want to use the massaging and heating pad again, press the "On/Off" button again.
- If you press the massage heads too hard, a system protecting the device will interrupt the massage. To reactivate it, unplug the device, then plug it back in.

Heat Function

- To activate the heating function, press the "Heating Function" button. A red LED light is displayed above the corresponding button.
- The 2-in-1 Shiatsu Comfort massaging and heating pad also has two heating intensities. Keep pressing the corresponding button to select the low intensity (light red LED light) or the high heating intensity (dark red LED light). To stop the heating function, press the "Heating function" button once again. The LED light disappears.
- To reactivate it, press the "Heating Function" button once again.
- The heating function is programmed to stop automatically after one hour. If you want to use the heating function again, press the "Heating Function" button once again.

 Notes:

- If you want to use the function for a lengthy period, we recommend that you select the lowest heat setting on the control mechanism to prevent your feet overheating as this may cause burns;
- Unplug the power cable from the wall socket after each use;
- When first using the device, it may emit a plastic odour. This odour will not persist.

INSTRUCTIONS FOR USING THE 2-IN-1 SHIATSU COMFORT MASSAGING AND HEATING PAD FOR YOUR BODY

- Remove the shiatsu pad from the foot-warmer using the zip fastening.
- Place the shiatsu pad on the area of the body you wish to massage using the instructions in this handbook.
- Refer to the instructions for use in the previous paragraph when using the shiatsu pad.

CLEANING AND MAINTENANCE

⚠ Warning

Before cleaning the device, unplug the power cable from the wall socket.

Do not wash the device and do not submerge it in water or any other liquid. To clean it by hand, used a damp cloth or sponge and a small quantity of liquid detergent, if necessary. Do not put the device in a clothes dryer. The foot warmer must not be wrung out, dry-cleaned, twisted, squeezed or ironed. The removable lining is machine washable at 86°F (30°C). Do not use the foot warmer until it is completely dry.

For long-term storage, put the device back in its original packaging and store it in a dry, uncovered place. Leave the device to dry completely before storing it.

✖ BREAKDOWNS AND SOLUTIONS

⊗ The foot warmer does not heat up or the massage function cannot be activated when:

- the power supply is connected to a plug that is live;
- the power cable is plugged into the foot warmer's power source;

- the two functions (heat +massage) have been activated.

- ⊗ The 2-in-1 Shiatsu Comfort massaging and heating pad is faulty.
- ⊗ Contact customer service.

TECHNICAL DATA

Model	LA110105
Voltage	DC 12 V
Power	30 W
Dimensions	310 × 110 × 310 mm

ADVICE REGARDING THE DISPOSAL OF WASTE

All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.

When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

LIMITED WARRANTY

LANAFORM guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM's control.

LANAFORM may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM's Service Centre. Any attempt to maintain this product by a person other than LANAFORM's Service Centre will render this guarantee void.

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