MEDISANA®

home of wellness

- DE** Massagesitz MC 810
GB Massage seat MC 810
FR Siège de massage MC 810
IT Sedile massaggiante MC 810
ES Asiento para el masaje MC 810
PT Asiento de massagem MC 810
NL Massagekussen MC 810
FI Hierova tuoli MC 810
SE Massagesäte MC 810
GR Κάθισμα μασάζ MC 810

Art. 88937

Gebrauchsanweisung
Instruction Manual
Mode d'emploi
Istruzioni per l'uso
Instrucciones de manejo

Bitte sorgfältig lesen!
Please read carefully!
A lire attentivement s.v.p.!
Da leggere con attenzione!
¡Por favor lea con cuidado!

Manual de instruções
Gebruiksaanwijzing
Καyttöohje
Bruksanvisning
Οδηγίες Χρήσης

Por favor, ler cuidadosamente!
A.u.b. zorgvuldig lezen!
Lue huolellisesti läpi!
Läses noga!
Παρακαλούμε διαβάστε προσεκτικά!

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GB

- 1 Power supply
 - 2 Massage seat
 - 3 Motor arrangement (6 massage motors for shoulders (1), back (2), waist (3/4) and thighs (5/6))
 - 4 Heating
 - 5 Connection for power unit and adapter cable
 - 6 Control unit
 - 7 Illuminated display
 - 8 Function button
- ON/OFF** switches the appliance on/off
-  switches the heating on/off
-  controls the intensity of the massage (high - medium - low)
- MODE** selects the massage programme (1 - 9)
- 9 Remote control pocket
 - 10 Car adaptor cable

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IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use.

Should you give this device to another person, it is vital that you also pass on these instructions for use.

Explanation of symbols

This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user.

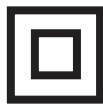
CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

Protection category II

LOT

Lot number

Manufacturer

Safety Instructions

about the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only connect the mains plug to the wall socket when the device is switched off.
- Only use the mains adapter that is provided with the massage seat.
- Keep the power supply, the cable, the control unit and the massage seat away from heat, hot surfaces, moisture and liquids. Never get hold of the power supply with wet hands.
- Do not take hold of a device if it falls into water. Pull out the power plug and the power supply unit from the socket immediately.
- Always unplug the equipment immediately after use.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull or twist the massage seat on the cable, the power supply unit or control unit.
- You must not continue to use the unit if the lead or the mains adaptor are damaged. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the unit in to the service centre.
- Ensure that the cable will not cause tripping. The cable must not be kinked, trapped or twisted.

specific groups of people

- This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless they have been instructed in the use of the device.
- Children must be supervised to ensure that they do not play with the device.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- You should refrain from using the massage seat cover **MC 810**, or should consult your doctor before using it, if
 - you are pregnant,
 - you have a cardiac pacemaker, artificial joints or electronic implants,
 - you suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Never leave the unit unattended if it is connected to the power supply.
- Take care when using the heater. Do not use if you have circulatory problems or on areas not sensitive to heat. **Risk of burning!**
- Consult your doctor if you have any questions concerning the therapeutic use of massage seats.
- If you have unexplained pain, are receiving treatment from a doctor and/or use medical equipment, consult your doctor before using the massage seat.

before using the unit

- Always carefully inspect the power supply, the cable, the control unit and the massage seat for damage before use. Do not attempt to use a defective unit.
- Do not use the unit if there are visible signs of damage on the unit or the cable components, the unit does not work properly or the seat or the control unit have been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.

operating the unit

- The massage seat must only be used for its intended purpose as described in the instruction manual.
- Using the mat for any other purpose invalidates the warranty.
- Never leave the unit unattended if it is connected to the power supply.
- The maximum operating time for one session is 15 minutes.
- Do not allow the massage cushion to come into contact with pointed or sharp objects.
- Never put the unit down next to an electric stove or other sources of heat.
- Never use the massage seat if it is folded up or crumpled.
- Do not stand on the unit.
- Do not touch the power supply unit when you are standing in water and always take hold of the plug with dry hands.

The device is not to be used for commercial use or medical purposes.

If you have health concerns, consult your doctor before using the massage cushion.

SAFETY IN ROAD TRAFFIC

Safety has the utmost priority in road traffic. The driver may not use the massage cushion when driving.

Only use the massage cushion as a passenger or in a parked vehicle.

Only use the massage cushion indoors or in the vehicle.

Do not use the massage seat in damp rooms (e.g. when bathing or showering).

maintenance and cleaning

- The unit is maintenance free. If a fault does occur, just check whether the mains power adapter is properly connected.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all guarantee rights. Ask your specialist dealer and have repairs carried out by the service.

- Do not use the unit if there are visible signs of damage on the unit or the cable components, the unit does not work properly or the massage seat or the control unit have been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.
- Do not immerse the device in water or other fluids.
- If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.

Do not wash!

Do not dry clean!

2 Useful Information

Thank you very much

Thank you for your confidence in us and congratulations!

You have acquired a **MEDISANA** quality product with this massage seat cover **MC 810**. This unit is designed for massaging the back and upper thigh areas. We recommend that you carefully read the following information on use and care of the product so that you experience its benefits and enjoy using your **MEDISANA** massage seat cover **MC 810** for many years.

2.1 Items supplied and packaging

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use it and contact your dealer or your service centre.

The following parts are included:

- 1 **MEDISANA** massage seat **MC 810**
- 1 Control unit
- 1 Power supply
- 1 Car adaptor cable
- 1 Instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required.

If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

3.1 Preparation

- Place the massage seat **2** on a high chair or armchair (**Car seat**) and attach it to this by means of the fasteners on the back.
- Connect the power supply unit cable **1** or the adaptor cable **10** to the remote control **5**.
- Plug the power supply unit **1** into the socket or the adaptor cable **10** into the cigarette lighter.
- Now sit on the seat **2** and familiarise yourself with the controls.

3.2 Operating the unit

The massage seat has four massage zones for specific treatment. These are in the main areas of tension (shoulders, back, waist and thighs). For each massage zone, you can select nine different types of massage with three levels of intensity. You can apply the soothing heat treatment in the waist area separately or use it together with the massage function. The built-in timer switches the appliance off automatically after 15 minutes.

3.3 Function button on the control unit **6**

- Function button **ON/OFF** **8** (Voltage supply):
Switch on: Press **ON/OFF**.
Switching off: Press **ON/OFF** again.
 To reset the timer, press **ON/OFF**.

CAUTION

Make sure that the maximum operating time of 15 minutes is not exceeded.

- Function button **MODE** **8** (Massage programme):
 Press **MODE** to go to the selection of massage programmes (**1 - 9**). The individual massage programmes vary in the type of massage and can be used individually or in combination.

 **Waves:** With this message type the motors are switched in such a way that a wave-like massage (Rolling) feeling is generated.

 **Kneading:** This message type is like a Shiatsu massage.

 **Tapping:** For this type of massage, the motors are switched to generate a tapping massage.

 **Pulsate:** For this type of massage, the massage motors pulsate in a certain sequence.

The nine massage programmes also differ in the sequence of motor use (**1 - 6**) **3**.

Select one of the following programmes **1 to 9**. The background lighting in the display **7** lights up green. Everytime you press **MODE**, you go to the next massage programme.

- Programme 1: **Tapping**
Motor use (1 - 6) ③
6, 5 _ 4, 3 _ 2 _ 1 _ ... (in a cycle)
- Programme 2: **Kneading**
Motor use (1 - 6) ③
2 _ 3, 4 _ 5, 6 _ 1 _ ... (in a cycle)
- Programme 3: **Tapping**
Motor use (1 - 6) ③
1) 5, 6 _ 4, 3 _ 2 _ 1 (5 times)
2) 1, 2, 3, 4 _ 5, 6 (5 times)
3) 1, 2, 3, 4, 5, 6 (5 times) _ ...
- Programme 4: **Tapping**
Motor use (1 - 6) ③
1) 1 _ 2 _ 2 _ 3, 4 (5 times)
2) 3, 4 _ 5, 6 (5 times)
3) 5 _ 6 (5 times) _ ...
- Programme 5: **Pulsate**
Motor use (1 - 6) ③
1 _ 2, 3, 4 _ 5, 6 _ ... (in a cycle)
- Programme 6: **Waves / Kneading / Tapping / Pulsate**
Motor use (1 - 6) ③
1) Waves 1, 2, 3, 4, 5, 6 (5 times)
2) Kneading 1, 2, 3, 4, 5, 6 (5 times)
3) Tapping 1, 2, 3, 4, 5, 6 (5 times)
4) Pulsate 1, 2, 3, 4, 5, 6 (5 times) _ ...
- Programme 7: **Waves**
Motor use (1 - 6) ③
1) 1 _ 2 (3 times)
2) 2 _ 3, 4 (3 times)
3) 5 _ 6 (3 times) _ ...
- Programme 8: **Waves / Kneading**
Motor use (1 - 4) ③
1) Waves 1 _ 2 _ 3, 4 (3 times)
2) Kneading 1 _ 2 _ 3, 4 (3 times) _ ...
- Programme 9: **Tapping / Pulsate**
Motor use (5 - 6) ③
1) Tapping 5 _ 6 _ 5, 6 (3 times)
2) Pulsate 5 _ 6 _ 5, 6 (3 times) _ ...

The display ⑦ shows the number of the set massage programme (1 - 9), the remaining operating time in minutes (15 - 1), as well as the massage intensity in

a bar diagram and the pulse symbol for the type of massage (*Waves / Kneading / Tapping / Pulsate*) being currently used at the bottom left.

The set programme is repeated until the 15 minute operating time has elapsed. Then the built-in timer switches the device off automatically.

- Function button  8 (Massage intensity):
You can choose between high, medium and low for the level of intensity. The medium intensity level is always set when the device is switched on. Press  to select the massage intensity (high - medium - low).
- Function button  8 (Heating on/off):
Press  to switch on the heating 4. The heat function can operate separately or combined with massage function. When the heating is used separately, the background lighting in the display 7 lights up red, when the massage and heating are used together, it lights up orange. The lower back and waist area are pleasantly warmed up. Press  to switch the heating off again.

4.1 Cleaning and maintenance

- Before cleaning the device, make sure that it is turned off and the mains plug is removed from the socket. Allow the unit to cool down.
- The massage seat mat must only be cleaned with a damp sponge. Never use brushes, aggressive cleaning agents, petrol, thinners or alcohol. Rub the seat dry with a clean dry cloth.
- When cleaning, never immerse the unit in water and make sure that no water gets into the unit.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.

4.2 Disposal

This product must not be disposed of together with domestic waste.

All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Consult your municipal authority or your dealer for information about disposal.

4.3 Technical Data

Name and model	:	MEDISANA Massage seat MC 810
Voltage supply	:	100 - 240 V~ 50/60 Hz 1.0 A
Power supply unit	:	12 V= 2500 mA
Car adaptor	:	12 V= 3 A  *2005/83*0011*00
Power consumption	:	approx. 30 W
Operating time	:	max. 15 minutes
Operating conditions	:	only use in dry rooms
Storage conditions	:	in a cool, dry place
Dimensions	:	approx. 106 x 49 cm
Weight	:	approx. 0.86 kg (without power supply)
Power cable length	:	approx. 1.8 m
Adapter cable length	:	approx. 1.8 m
Article No.	:	88937
EAN Code	:	40 15588 88937 0

In accordance with our policy of continual product improvement, we reserve the right to make technical and optical changes without notice.

**Warranty and
repair terms**

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

 **MEDISANA AG**
Jagenbergstraße 19
41468 NEUSS
GERMANY

E-Mail: info@medisana.de
Internet: www.medisana.de

The service centre address is shown on the attached leaflet.

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