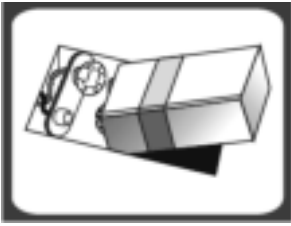
AD 8116

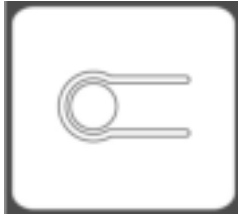
(GB) user manual
(F) mode d'emploi
(P) manual de serviço
(LV) lietošanas instrukcija
(BIH) upute za uporadu

(D) bedienungsanweisung
(E) manual de uso
(LT) naudojimo instrukcija
(EST) kasutusjuhend
(PL) instrukcja obsługi

1

2

3

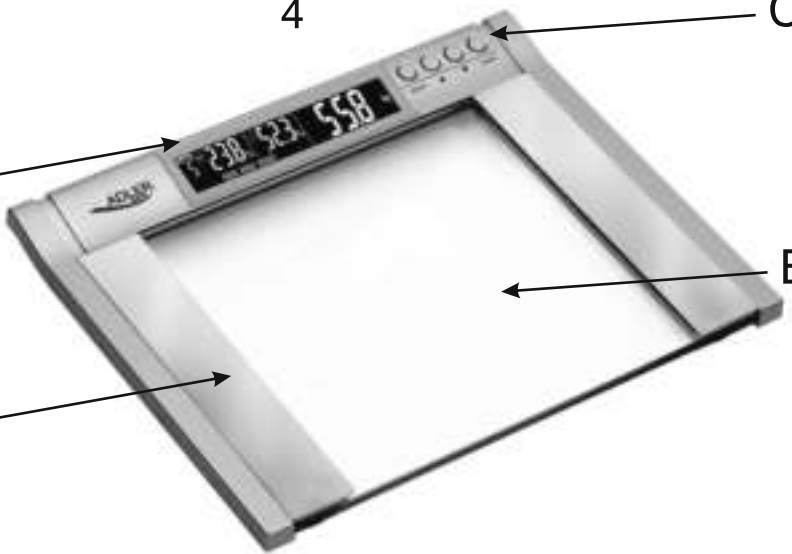
4

C

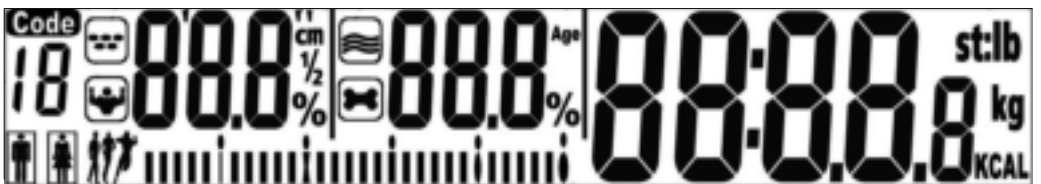
A

B

D

a

FAT/ HYDRATION CONTENT FITNESS ASSESSMENT CHART

After weighting you can use the following chart:

Age	Fat(F)	Hydration(F)	Fat(M)	Hydration(M)	Status	Key Line Indication
≤ 30	4.0-16.0%	66.0-57.8%	4.0-11.0%	66.0-61.2%	Too lean	iii
	16.1-20.5%	57.7-54.7%	11.1-15.5%	61.1-58.1%	Lean	iii iii
	20.6-25.0%	54.6-51.6%	15.6-20.0%	58.0-55.0%	Normal	iii iii iii
	25.1-30.5%	51.5-47.8%	20.1-24.5%	54.9-51.9%	Fat	iii iii iii iii
	30.6-60.0%	47.7-27.5%	24.6-60.0%	51.8-27.5%	Too fat	iii iii iii iii iii
> 30	4.0-20.0%	66.0-55.0%	4.0-15.0%	66.0-58.4%	Too lean	iii
	20.1-25.0%	54.9-51.6%	15.1-19.5%	58.3-55.3%	Lean	iii iii
	25.1-30.0%	51.5-48.1%	19.6-24.0%	55.2-52.3%	Normal	iii iii iii
	30.1-35.0%	48.0-44.7%	24.1-28.5%	52.2-49.2%	Fat	iii iii iii iii
	35.1-60.0%	44.6-27.5%	28.6-60.0%	49.1-27.5%	Too fat	iii iii iii iii iii

M.- Male

F- Female

SAFETY CONDITIONS.

IMPORTANT INSTRUCTIONS ON SAFETY OF USE

PLEASE READ CAREFULLY AND KEEP FOR FUTURE REFERENCE

1. Before using the product please read carefully and always comply with the following instructions. The manufacturer is not responsible for any damages due to any misuse.
2. The product is only to be used indoors. Do not use the product for any purpose that is not compatible with its application.
3. Please be cautious when using around children. Do not let the children to play with the product. Do not let children or people who do not know the device to use it without supervision.
4. This appliance is not intended use by persons (including children) with reduced physical, sensory or mental capabilities, or instruction concerning use of the appliance by a Person responsible for their safety.
5. **WARNING:** This device may be used by children over 8 years of age and persons with reduced physical, sensory or mental capabilities, or persons without experience or knowledge of the device, only under the supervision of a person responsible for their safety, or if they were instructed on the safe use of the device and are aware of the dangers associated with its operation. Children should not play with the device. Cleaning and maintenance of the device should not be carried out by children, unless they are over 8 years of age and these activities are carried out under supervision.
6. Never use the product if it was dropped or damaged in any other way or if it does not work properly. Do not try to repair the defected product yourself because it can lead to electric shock. Damaged device always turn to a professional service location in order to repair it. All the repairs can be done only by authorized service professionals. The repair that was done incorrectly can cause hazardous situations for the user.
7. Never put the product on or close to the hot or warm surfaces or the kitchen appliances like the electric oven or gas burner.
8. Never use the product close to combustibles.
9. Never use this device close to water, e.g.: Under shower, in bathtub or above sink filled with water.
10. Do not touch the device with wet hands.
11. The batteries may leak if they are flat or were not used for an extended period of time. To protect the device and your health, replace them regularly and avoid skin contact with leaking batteries.
12. Do not use the device near water e.g.: in the shower, bath or basin with water.
13. Do not hold the device with wet hands.
14. The device should be turned off each time it is put aside.
15. Please use the device on a stable and flat surface.
16. Please keep the device away from the influence of strong electromagnetic fields.
17. Please ensure that the surface on which the device is used is not damp or slippery.

18. Please ensure that heavy objects do not fall on the display and glass part of the scales.
19. If the scales are not used for an extended period, please remove the battery.
20. DO NOT strike, DO NOT shake and DO NOT drop the scales.
21. Stand still on the scales while the measurement is being taken.
22. Do not use for commercial or medical purposes.
23. It is recommended that you weigh yourself at the same time each day.
24. The measurements given by the scales can be misleading if you weigh yourself after intense exercise, during a strict diet or if you are dehydrated.
25. Children under 10 years of age and the elderly may use only the weight measurement option. The results of measurement options (amount of water, fat or muscle) can be misleading for people over 70 years of age and people who practise bodybuilding.
- 26. The device is not intended for pregnant women, people undergoing dialysis, people with heart pacemakers and other medical implants.**
27. If your measurements differ from the norm, consult a physician.

DESCRIPTION OF THE DEVICE (Fig. 4)

- | | |
|-------------------|------------------------|
| 4A. LCD display | 4B. Glass platform |
| 4C. Function keys | 4D. Measurement plates |

FUNCTION KEYS (Fig. 3)

- | | |
|---|---|
| SET - on button/parameter change button | START- measurement control |
|  - up button |  - down button |

LCD DISPLAY (Fig. a):

USER PARAMETERS:

- | | |
|---|-------------|
|  | - user code |
|  | - gender |
| cm | - height |
| AGE | - age |

MEASUREMENT:

- | | | | |
|--|------------|---|---------------|
|  | - fat % |  | - fat level |
|  | - water % | KCAL | - calories |
|  | - muscle % | st:lb kg | - weight unit |
|  | - bone % | | |

ERROR MESSAGES

- | | |
|--------------------------------------|--|
| Lo- low battery, replace the battery | Err - maximum load exceeded |
| ErrL- low fat level | Errh - high fat level |
| C - error |  - during measurement |

CHANGING WEIGHT UNITS (Fig. 2)

You can change the units that the weight is displayed in using the switch found under the scales, next to the battery container or under the cover of the container. Available units: kilograms (kg) and pounds (lb).

REPLACING THE BATTERY (Fig. 1)

Remove the lid of the battery compartment at the bottom of the scales. Remove the old battery. Insert a new 9V 6F22 battery. Close the lid of the battery compartment.

WEIGHING

- STEP 1. Place the scales on a hard and flat surface (avoid carpets and soft surfaces).
- STEP 2. Gently step on the scales. Stand on the scales evenly distributing your weight and do not move until your weight reading shown on the display stops changing and gets blocked.

SETTING MEASUREMENT PARAMETERS

- STEP 1. Place the scales on a hard and flat surface (avoid carpets and soft surfaces).
- STEP 2. Press the SET button. An option to set the user code will appear.

STEP 3. Use the ▲ and ▼ buttons to choose a user code.

STEP 4. Press SET to confirm.

Repeat steps 3 to 4 to set the gender, height and age. After approximately 7 seconds from the last press of the SET button the scales will switch to the weighing and measuring mode.

WEIGHING WITH MEASUREMENT

STEP 1. Press the SET button. The last user's parameters will appear. If you want to change the user or parameters use the ▲ and ▼ buttons. Wait until the display shows "0".

STEP 2. Gently step on the scales with bare feet. Both feet should touch the measurement plates (4d). Stand on the scales evenly distributing your weight and do not move until your weight reading shown on the display stops changing and the measurement sign: 0000 appears.

Next a sequence of fat, water, muscle, bone and calorie measurements results will appear. They will be repeated three times and saved to the scales memory.

MEASUREMENT CONTROL

The results of the previous measurement can be retrieved using the SET button. Next, select the user code and press START.

AUTOMATIC OFF FUNCTION

When you step off the scales, the device will automatically switch off. The scales automatically switch off when the display shows '0.0' or the same weight reading for approximately 8 seconds.

CLEANING AND MAINTENANCE

1. Clean the scales with a damp cloth. Wipe dry.
DO NOT use chemical/caustic cleaning agents.
2. Store the scales in a cool and dry place.
3. The scales should always be in a horizontal position.

Device is made in class III of insulation.
Device is compliant with EU directives:
Electromagnetic compatibility (EMC)
Device marked CE mark on rating label

TECHNICAL DATA

Power source: battery type: 1x 6LR61/9V
Maximum load: 200 kg

To protect your environment: please separate carton boxes and plastic bags and dispose them in corresponding waste bins. Used appliance should be delivered to the dedicated collecting points due to hazardous components, which may effect the environment. **Do not dispose this appliance in the common waste bin.**

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