


IPL Hair Removal Instrument


Instructions for Use

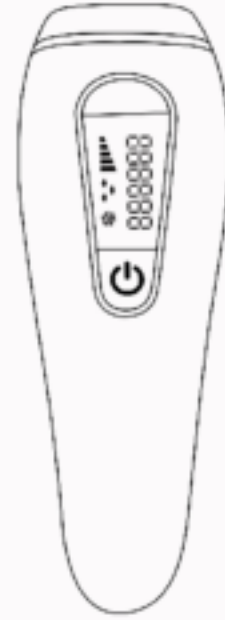

Product Overview


1.1 Instrument Description

- ①Hair Removal Instrument
- ②Contact Surface
- ③Switch button
- ④Screen
- ⑤Light Button

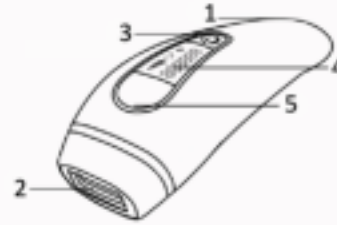

Display Screen

- 1. Flash Residual Number Display
- 2. Indication of Gear Position
- 3. Automatic Continuous Light Output Mode
- 4. Fan

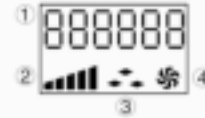

Please read it carefully Before using

It is not anyone that is suitable to use this hair removal instrument. Please refer to the data before use

-Please refer to page 1 [color control diagram]. For skin type VI, do not use this device, otherwise it is very possible to cause adverse reactions to the skin.



1

HOW TO USE

1. BEFORE USE

- 1. Please remove the hair on the skin surface firstly. If the hair does not exceed 2 mm, the instrument can be operated directly without shaving.



2. Steps

- 1. Please clean the skin surface and filter with soft cloth before use.
- 2. Please connect the power adapter and intensity indicator. The flashing "red light" indicates that the circuit is normal.
- 3. Please press the power switch until the screen lights up.
- 4. Put the device at a 90° angle on the skin and make sure the flash window is full contact with your skin.
- 5. If the screen changes from "flashing" to "constantly bright", it indicates the epilator is electrified.



2



- 6. After the machine is pressed against the skin, please gently press the flash button until it emits a flash.
- 7. After using the instrument, please replace the depilated part.
- 8. Please shine as close as possible to the site of previous treatment to ensure that all sites are adequately treated.

Two Mode

Two treatment modes including [press and flash] and [slide and flash] are convenient for you to use on different body parts:

The **press and flash mode** is suitable for small or curved areas, such as knees and armpits. You just need to press the flash button, then the flash will appear.

The **slide and flash mode** is convenient for large areas such as legs. Please hold down the flash and slide the instrument along the skin at one side to emit multiple flashes continuously.

Therapeutic Cycle Reference Table

1-4 week	5-6 week	7-8 week	Above 8 week
3 Times per week	2 Times per week	1 Time per week	1 TIME PER 2 WEEK

3

One depilatory course of treatment refers to 3 times of irradiation. Please do not give light continuously on the same place for three times.

The higher gear position, the better effect. Recommendation is as follows: please use gear 2-4 for the first time use. Then, please choose gear 5. During the interval of treatment, please examine your hair growth speed. The effect differs on various body parts.

After use

①The skin may be slightly inflamed or prickly or warm. These reactions are completely harmless and will soon disappear.

②You can safely apply lotion, cream, deodorant, moisturizer or cosmetics to the treated area. If any tingling or redness appear after the treatment, please wait for these conditions to disappear before applying any other product to your skin.

Clean Preservation

- 1. After use, please turn off the power and pull out the plug to cool down the instrument.
- 2. Please use wet cloth on the following parts: -Filter - Skin Surface
- 3. Please keep the instrument in a dust-free dry place with temperature ranging from 0 to 40°C.



4

The Curative Effect Is Remarkable

The clinically proven results show up to 90% hair reduction after 4-7 times of treatment and 90% customers are satisfied with the treatment result after long-time use. The wavelength over 610nm can close hair follicles and inhibit hair regeneration. So you can farewell repeated hair removal.

When will you see the results?

With just 4-5 times self-treatment sessions, most users can experience visually diminished hair that is comparable to results done by professional administered laser procedures. Optimal results require at least continued use for 6 months.

Where can you use it?

- 1. You can remove hairs on lips, chin, sideburns, armpit, arms, legs, toe-wax, bikini and private parts.
- 2. Level 1 is suggested for the first time use. After that, you can adjust to intensity that you can stand.

Who can you use it?

Please do not use this instrument on naturally dark or recently tanned skin. This beauty tool is ideal for use with lighter skin tones and brown-to-black hair. White, grey, red and blond hairs are not recommended.

5

What are the precautions at the time of using?

Please keep the power cord plugged in; do not pull the plug out during using. The device doesn't need to be charged in advance. Just use it directly after connecting the power cord.

After-treatment Tips

- 1. Please use moisturizer to soothe your skin after the treatment.
- 2. Please avoid ultraviolet rays within 24 hours of treatment.

Contraindications:

- Not for pregnant or breastfeeding women, or during menstrual period;
- Not for individuals with scar diathesis;
- Not for individuals with photosensitization history or diabetes;
- Not for individuals who have depilated in other methods in the past 6 weeks and used skin whitening products or products with vitamin A in the past 1 week;
- Not for individuals who have taken essential oil treatment at the same time;
- Not for individuals who have taken Aspirin and immunosuppressant;
- Not for individuals who have been exposed to strong sunlight in the past 2 weeks;
- Not for sensitive areas, like mucosa, etc and private positions, like upper eyelid, lip, nipple and vagina etc.;
- Not for individuals with hirsutism & other diseases that may trigger complications.

6

7

8

9

10

11