




4D PULSE MASSAGER | INSTRUCTION MANUAL

- 1.Button ON/OFF: total power switch;
- 2.VM key: the vibration switch and the gear key, a total of 5 files, continuous button enhancement; LCD screen display: SPEED said.
- 3.MC key: the pulse current button, a total of 5 files, continuous button enhancement; LCD screen display: M-CURRENT said.
- 4.LED key: LED germicidal lamp switch in and fluorescent detection mode.
- 5.MODE key: mode key, short press a button in the 1/2/3/4/5 file to be replaced, corresponding to the pulse width of 0.5s, 0.25s, continuous discharge and five kinds of vibration mode to switch. LCD screen display: MODE:
- 6.Press the power button after the display lights, according to the need to select a single function or mixed use function, and select the appropriate vibration frequency, pulse width, intensity, moderate rolling back and forth in different parts.
- 7.Operation is completed, please turn off the power supply, 20 minutes without any operation, the automatic shutdown of the product

- ① Pulse ball
- ② LCD display
- ③ Vibration switch/Gear key
- ④ Pulse current button
- ⑤ Bactericidal switch/Fluorescence detection mode
- ⑥ Mode key
- ⑦ Power main switch

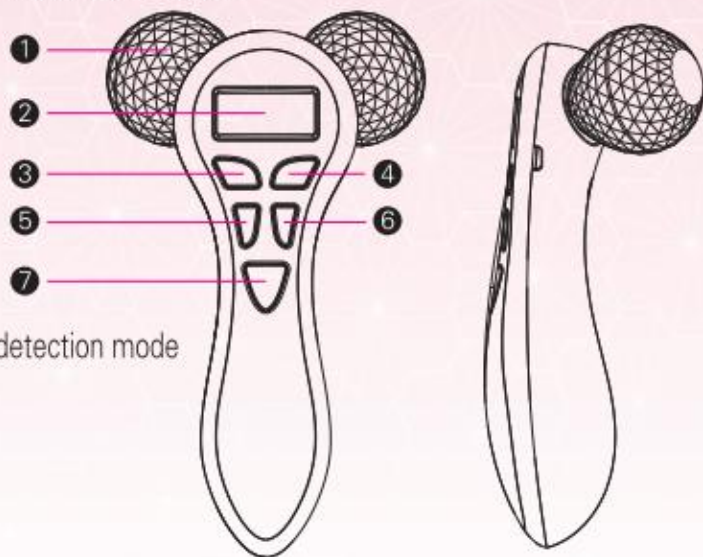



THE MAIN FEATURES OF PRODUCTS

The main features of products: pulse roller massager is high-tech electronic products according to TCM Meridian Theory and modern bio engineering technology combined research made, with five strength of gear and five different function modes, it uses the rules of sphere reasonable ergonomic design, can be in the surface of human skin to skin and move flexibly. The body form a squeeze, massage, with various modes of vibration motor, the purple line light irradiation and electric pulse stimulation, muscle activation of muscle tissue, effectively prevent muscle relaxation, but also can dredge the meridians, relieve pain, is the best product of daily health care, and other clinical treatment has good auxiliary effect.



V face Plastic



Pulse Current



Vibration Massage



Fluorescence
Detection

-
- 1.Name: Pulse roller massager
 - 2.Model:KM-828
 - 3.Product appearance material: ABS
 - 4.Weight: 185g
 - 5.Product specifications: 160x87x69mm
 - 6.Polymer lithium battery: 2200mA voltage: 3.7V output power: 1W
 - 7.Charging time: 5-6 hours



1. Using high frequency vibration mode, promote muscle movement, rapid change of fat granule, achieve partial weight loss, saving time and effort
2. High efficiency current pulse modulation technology, the implementation of the role of the skin, meridians, acupuncture points, clear the meridians to lymphatic drainage
3. The current generation of microwave stimulation of deep subcutaneous tissue to repair muscle, skin elastic fiber and plastic layer, stretch, wrinkles (firming wrinkle)
4. The mechanical vibration and voltage stimulation, accelerate micro vascular blood circulation, exercise muscle tissue, promote formation of different aerobic exercise, enhanced cell permeability, the nutrient supply effectively muscle tissue and skin; the human skin tightening, increase their flexibility, enhance skin cells The new supersedes the old. to slow down aging. Effect
5. The use of **Bluelight** on the skin of irradiated tissue, promote collagen regeneration, remove skin bacteria, promote calm and comfort effect (except sedation)
6. The use of Y type platinum massage wheel design, so as to compact the face, hip and leg, a thin waist can all use and achieve special effect
7. Massage to relax / remove bags under the eyes



Cheek
Massage



Chin
Massage



Neck
Massage



Arm
Massage



Chest
Massage



Abdominal
Massage



Waist
Massage



Massage
For Legs



MATTERS NEEDING ATTENTION / GUARANTEE

- 1 with heart problems, high blood pressure, diabetes, pregnant women, children and other users, please follow the instructions to use.
- 2 if the use of the process of skin pain or other discomfort, please stop using.
- 3 the product can be used together with the mask or essence, the effect is better, can not be washed, so as not to cause short circuit.
- 4 please use the genuine USB charging, full of after the timely removal of the plug.

-
- 1.Free maintain for the breakdown caused by product quality problems;
 - 2.One year guarantee from the date of purchase ;
 - 3.Free to change broken partsm,freight is collected by customers ;
 - 4.If the customer is damaged, the customer shall pay the service parts and the freight;
 - 5.The warranty card to block letters, shall not be altered, please safekeeping, if lost, will not be retroactive;

Product Name		Item No.	
Customer Name		Phone	
Purchase Date			
Address	Means Of Purchase		
Records			
Date	Content	Results	Service Provider