


# 尺寸大小：272mm\*195mm

## 四折纸

Unique and innovative design with 6 led bulbs

68 mm

195 mm



**woosmile**  
Blue Light  
Sonic LED Electric Toothbrush

**INSTRUCTION MANUAL**

The WOOSMILE toothbrush offers you and your home a different, healthy and efficient brushing experience through its innovative design. Please read the product instruction carefully before use and keep it in a safe place.

**BRUSHING MODES**

To select a mode, simply push the arrow button until the desired mode indicator lights up. Once lit, that mode will become the default mode until another mode is selected. Each mode is designed to operate for two minutes, with brief pauses between each 30 second burst.

**NO.1 Premier LED Technology**  
Unique and innovative design with 6 led bulbs, cold blue led light, cleaning & whitening tooth with better effective.

**Clean:**  
Clean is the toothbrush's workhorse mode. At brush strokes per minute, it removes plaque with superior efficiency in a two-minute program. The Clean mode is the base setting on the toothbrush. If you're going to stick to one mode, make it this one.

**Sensitive:**  
If you've got sensitive teeth or gums, or you just find the sonic vibrations a little overwhelming at first, switch to the Sensitive mode. In this mode, the toothbrush vibrates with lower intensity, going easier on your teeth and gums.

**Polish:**  
Improve the gloss of the tooth surface. Swing changes rapidly. The vibration intensity changes rapidly within 0.1 second to stimulate the polishing effect. It can be used for extra treatment of the incisor area to whiten the teeth.

**White:**  
In white mode, the toothbrush works a bit harder to remove surface stains caused by things like coffee and tea.

**TIPS:**  
Use with teeth whitening gel will have better whitening effect.

**WHAT IS INCLUDED**

- 1 x Whitening Sonic Electric Toothbrush
- 2 x LED Electric Toothbrush Heads
- 1 x Charger Base
- 1 x Charging Cable
- 1 x User Manual

**HOW TO USE**

- According to the arrow direction, insert the brush head into the main body of the toothbrush.
- Take a small amount of toothpaste on the brush head and wet it.
- After wetting, choosing different function, put the toothbrush in the mouth and start the toothbrush, then move slowly and uniformly.
- When brushing teeth, the brush and the gum are 45 degree.
- 30 Seconds reminder to replace tooth area, 2 minutes scientific brushing.
- After brushing teeth, rinse brush head and fuselage directly with clean water, then smile.


272 mm

1                      2                      3                      4

正面


**MEMORY FUNCTION**

The toothbrush remembers the mode during the last shutdown, and will repeat it on start.




**INTELLIGENT BRUSHING**

30 second interval pause reminding, 2 minutes smart timer. Quad pacer and Smart timer helps ensure dental professional recommended 2 minutes brushing time and encourages brushing in each quadrant of the mouth

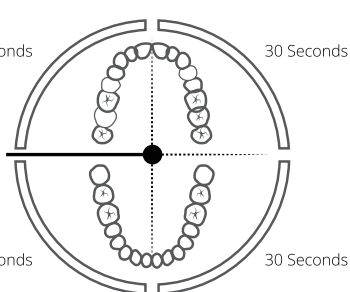


2 minutes total ( 2 x daily )

**LONG-LASTING POWER**  
Fully charged in 10 hours, 4 weeks of usage.



**CHARGING & CLEANING**



LED charging indicator will turn RED during charging. After fully charged, the indicator will turn GREEN.

**Charging The Toothbrush**

Be sure to turn the toothbrush off before charging. Charge by plugging the base into an outlet and placing the handle on the charger. The charger indicator light will turn RED, indicating the handle is charging. After the battery is fully charged, the charge indicator light will turn GREEN, indicating the toothbrush is fully charged.

- The battery will be fully charged after 10 hours. At full charge, the toothbrush can be used for approximately 4 weeks, pre-timed brush cycles.
- The toothbrush must be charged for 10 hours if:
  - this is the first time using the toothbrush
  - the toothbrush has not been used for more than three months
  - the charge is drained from use
- To keep the battery fully charged, we recommend that you keep the toothbrush on the charger when not in use.
- It is normal for the toothbrush to be warm when charging.

**Cleaning & Storing**

After brushing, rinse the toothbrush head under running water for a few seconds to remove any leftover toothpaste or debris. To maintain smooth operation, occasionally remove the toothbrush head and rinse the handle shaft under running water.

Avoid covering the toothbrush head after cleaning. This can cause bacteria growth on head and bristles.

Store toothbrush in an upright position to allow proper moisture drainage and air drying.

If not using for an extended period of time, unplug from outlet, clean and store in a cool, dry area, away from direct sunlight.

**PARAMETERS**

Rated Voltage: DC 3.7V  
Input Voltage: DC 5V/1A  
Standby Power Consumption: <0.5W  
Battery Type: Polymer 800mA Li-Ion Battery  
Battery Life: Over 4 weeks after fully charged  
Vibration Frequency: 38400 rpm  
Waterproof Level: IPX7

**DISCLAIMERS / WARNINGS**

- Discontinue use if this product appears damaged in any way (brush head, handle, charger.) This product contains no serviceable parts.
- Never use the charger if cord/plug is damaged.
- Keep the cord away from heated surfaces.
- This product is designed to clean your teeth, gums and tongue only. Discontinue use of this product and contact a physician/dentist if discomfort or pain occurs.
- Keep out of reach of children

**CAUTIONS REGARDING USE**

- To prevent electric shocks and fire, do not charge the product in the bathroom or put it into water or other liquids, and do not wash the charging port with water.
- Do not use the product, when the cable or the plug is damaged, or the socket is not tightly inserted.
- Do not put the product into hot water.
- Keep the appliance away from children.
- Always ensure the appliance is operated on an electric power source matched to its rated voltage. Failure to do so may result in fire or electric shock.
- It is suggested to change a new brush head when using for 3 months.
- When the cable is damaged, change a new one so that you can charge it.

**DANGER**

In order to reduce the risk of electrocution:

- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not wash the charger.
- Do not use while bathing.

**IMPORTANT**

- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Do not modify the products. Do not disassemble or repair. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact a service center. Do not disassemble the product except when disposing of the battery.
- This product is not intended for use by children.

**WARRANTY**

WOOSMILE guarantees its products for one years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at WOOSMILE expense provided that convincing proof of purchase in the qualifying period is provided. Under following cases WOOSMILE won't be responsible for warranty.

- Brush heads
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

**IMPLIED WARRANTIES**

Any implied warranties, including implied warranties of merchant ability and fitness for a particular purpose, shall be limited to the duration of the express warranties set forth above. In some states limitations on duration of implied warranties do not apply.

1                      2                      3                      4

反面