

DOUBLE-HEAD STETHOSCOPE

(Form & Structure)

1. The contents of Packing

Double type Chestpiece, Y-Type Tube, Binaural

2. Assembly completed chart. Fig 1

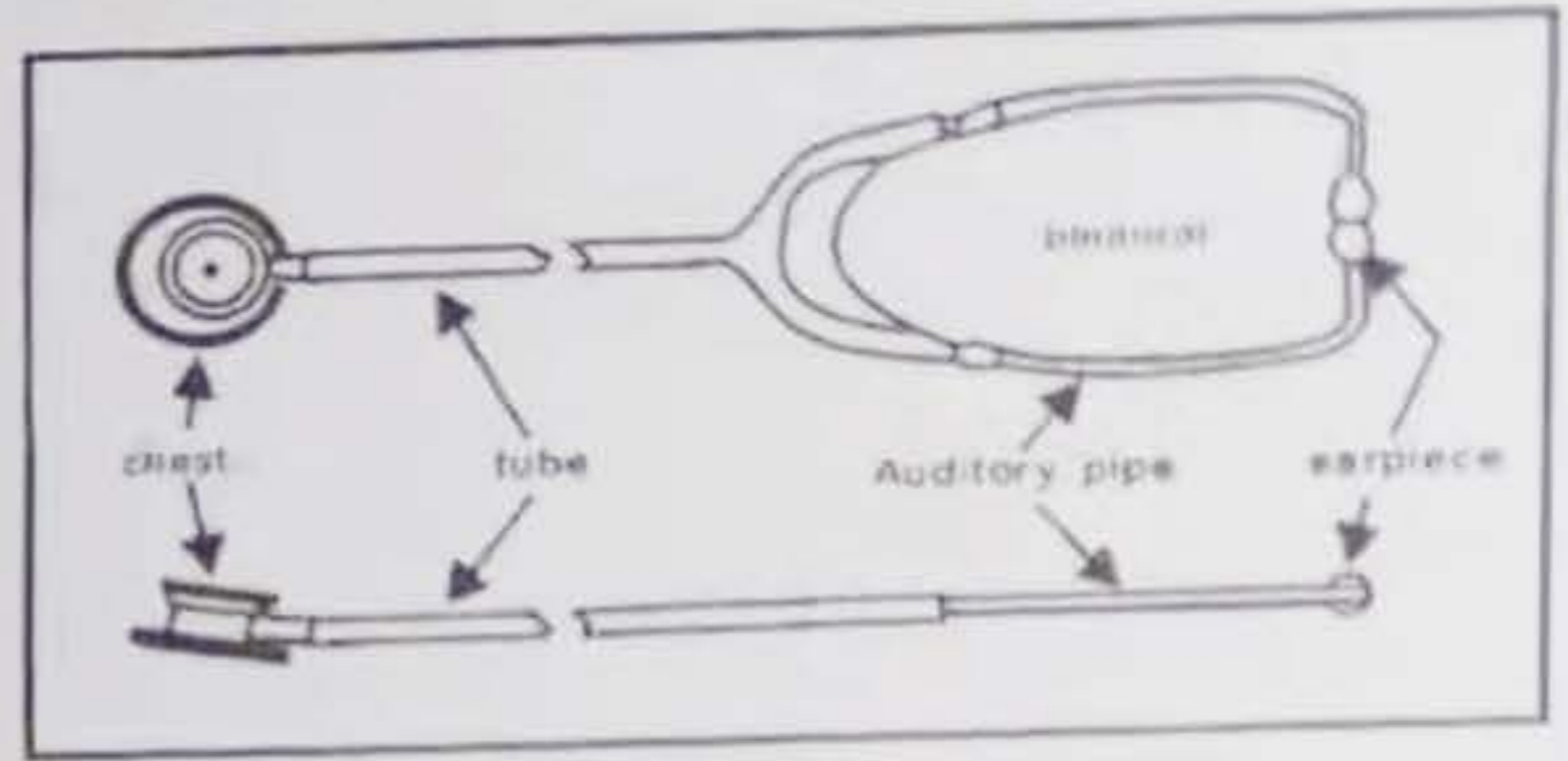


fig 1

(How to use)

1. How to fix the Auditory pipe

Auditory pipe are rotated as shown in the following figure. Fig 2.3

Please adjust the holes of earpieces to align with direction of auditory pipe.

For reference, the standard angle is 18 degree.

As the auditory pipe are assembled a little bit harder

there may be a friction noise "creak", but it is not a sign of defect.

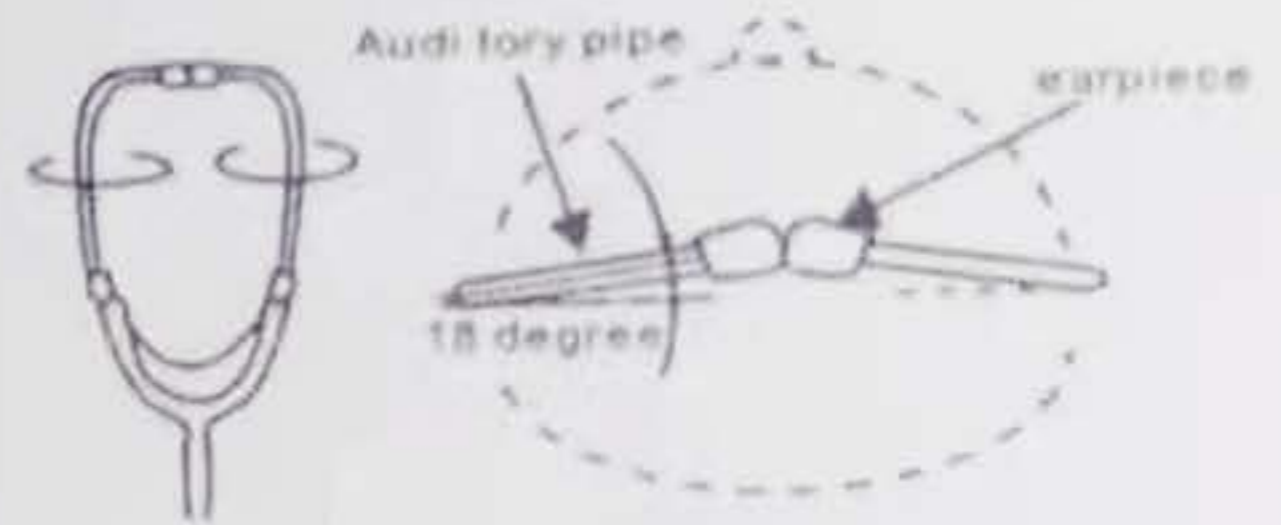


fig 2

fig 3

2. How to changeover between Diaphragm plane and Bell plane

It is possible to changeover the diaphragm plane to the bell plane by rotating the shaft as shown in following figure. Fig.4

Please turn until you feel a click. When the changeover mark is visible, it apply the equipment to the surface of auscultation.

In case of the status shown in the following figure

the auscultation is possible from the diaphragm plane

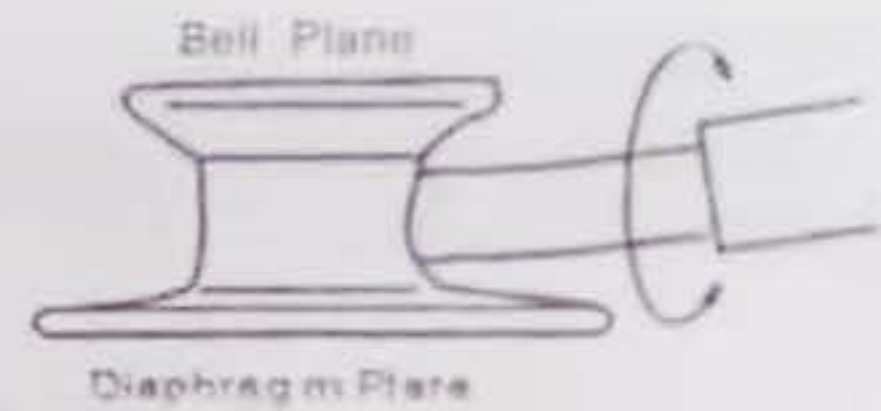


fig 4

3. How to loose the Binaural

Auditory pipe are removed as shown in the following figure. Fig 5

When you remove an auditory pipe, please extend the both sides of an auditory pipe simultaneously

(CAUTIONS: Be careful of the obstacle to the face by the return of spring part.)

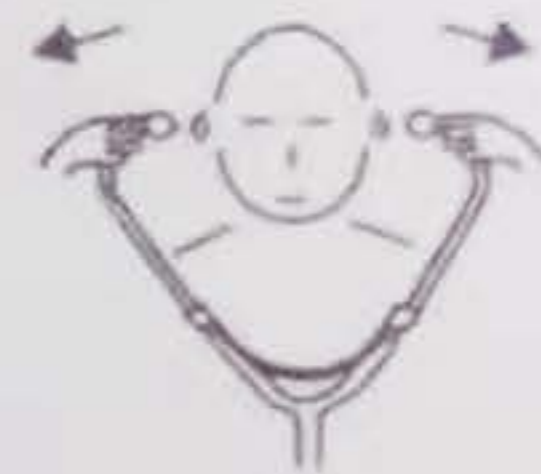


fig 5

4. How to exchange the parts.

I. The attachment & detachment of earpiece

The earpiece is taken off as shown in the following figure. Fig. 6

II. The attachment & detachment of diaphragm

The diaphragm is taken off as shown the following figure. Fig.7

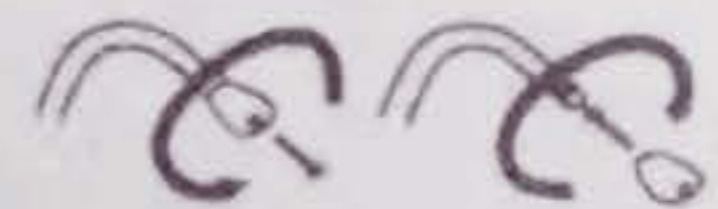


fig 6



fig 7