FITAGE

USER MANUAL

Requires mobile app to experience all the features.
Compatible with iOS and Android devices.
iOS 8.0 or later / Android 4.3 or later.

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WELCOME TO **FITAGE** GET IN SHAPE

Thank you for choosing FITAGE!

With this smart scale we will be your personal health assistant. The device uses bioelectrical impedance analysis (BIA) to provide you the data you need to keep track of your personal health metrics such as weight, BMI (Body Mass Index), body fat and muscle mass.

With the bathroom scale you can use our exclusive FITAGE app, which will be your own personal trainer in your back pocket. Begin the beautiful journey of becoming healthier with FITAGE.

SPECIFICATIONS

Model:	FITAGE-21
Dimensions:	30x30x2.4 cm (11.8"x11.8"x0.94")
Weight:	1.65 kg (3.6 lb)
Maximum weight capacity:	225 kg (500 lb)
Minimum weight measurement:	0.2 kg (0.4 lb)
Power source:	3x AAA Batteries

WHAT'S INCLUDED?

1. FITAGE smart scale
2. 3x 1.5V AAA batteries
3. User manual

APP INSTALLATION & FIRST MEASUREMENT

The smart scale connects with your telephone or tablet via Bluetooth. To do so, your mobile device needs to have the FITAGE app installed. The easiest way to download the app is to scan the QR code below. Follow these steps to experience the full functionality of your FITAGE smart scale!

STEP 1: DOWNLOAD THE APP

1.


2.

3.

1. Download the FITAGE app by scanning this QR code
2. Create an account or continue as guest
3. Install the batteries in the scale (at the backside)


STEP 2: PAIR THE SMART SCALE WITH OUR APP

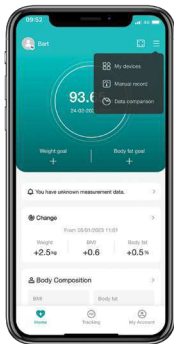1.

2.

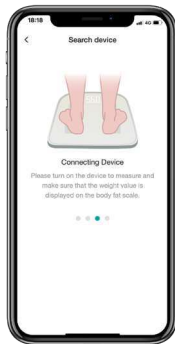3.

4.

1. Turn on Bluetooth on your mobile device and add a new scale in the app
2. Select the FITAGE scale which is turned on
3. Step with dry and bare feet on the scale
4. Wait till the binding is successful

Caution: using older mobile operation systems (like Android 6.0) will require you to enable Location Permission in order for the Bluetooth connection to work. Please check this in the general mobile device settings when pairing does not work like instructed.

Need more information about data processing?
Check the privacy policy on our website.

STEP 3: MEASUREMENTS & BODY ANALYSIS REPORTS

1.

2.

3.

1. Quickly step on and off the scale to calibrate the device (can be done anytime)
2. Step on the scale and stay on it until the displayed digits are stabilized and the numbers flash twice
3. Check the mobile app for your personal body analysis report

When the displayed number flashes twice, it indicates that the weight measurement is complete. Once the light bar stops moving back and forth, the body fat and other body composition measurements have also finished and the data is sent to the phone.


MEASUREMENT TIPS

Read these tips to get the most accurate measurements from your FITAGE smart scale.

- This device uses bioelectrical impedance analysis (BIA) technology. In order to have accurate and consistent measurements, please weigh every day on the same moment. For example, every morning when you get up.
- Please avoid taking measurements:
 - During or immediately after heavy exercise
 - After a heavy meal
 - Immediately after a shower, bath or sauna
 - During or after alcohol consumption
 - During a fever
- Please be barefooted and stand straight while standing on the smart scale. Wearing socks or boots may


prevent certain metrics from being measured. Dirty feet can also affect the accuracy of the measurement.

- Weight units can be changed from kilogram (kg) to pound (lb) or vice versa via the mobile app under Settings. The height indicator will change accordingly with it (meter or foot/inch).
- Place the smart scale on a flat and hard floor before taking any measurements. Calibrate the device if needed, then the weight should display 0.0 kg/lb before stepping on the smart scale.

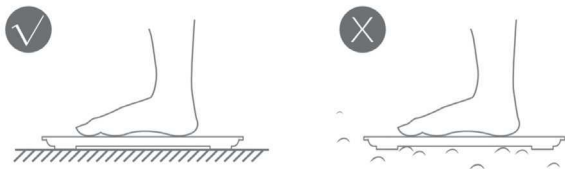Always use the device on a hard surface

- Measure body fat, muscle mass and other metrics correctly. Have both feet uniformly touch the electrodes as showed below.

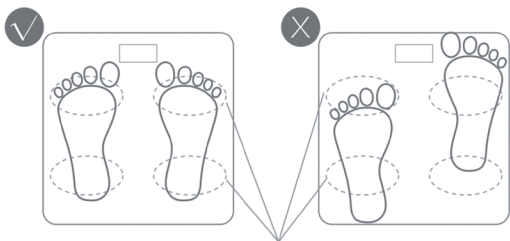Conductive areas


POWERING UP THE SCALE

Your scale shipment includes 3 AAA batteries. Insert these batteries into the designated compartment. Once powered on, your scale should display 0.0, indicating that it is ready for use.

With the batteries installed you can measure your body weight by stepping on it. However, to measure all the body metrics (such as your body fat percentage) you need to use the FITAGE app via a mobile device.

Please note that the included batteries are demo versions. Replace them within two months in order to avoid damage and leakage to the smart scale.


ADD MULTIPLE USERS

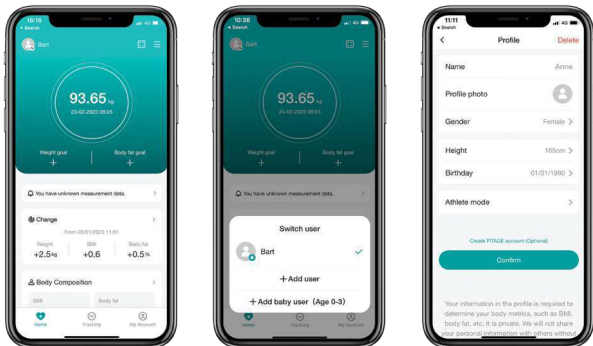1. Click on the avatar in your app home screen to add or switch users
2. You can delete a user by selecting the profile and go to "My Account"

You can also add an extra user by installing the app on another phone or tablet. The bathroom scale can connect to multiple mobile devices. This way you can achieve your health goals together!

DISPLAY DESCRIPTION

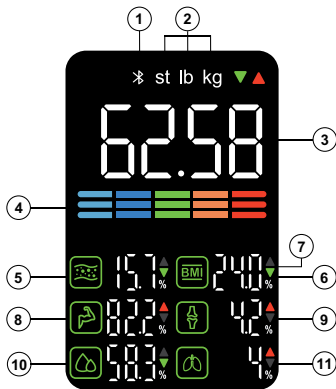1. Bluetooth
2. Weight Units
3. Weight
4. BMI Colour Indication
5. Body Fat
6. BMI
7. Lower or Higher Arrow indication
8. Muscle Mass
9. Bone Mass
10. Body Water
11. Visceral Fat

1. The BMI color indication, shown in five colors from left to right, correspond to different BMI levels: very low (light blue), low (dark blue), normal (green), high (orange), and very high (red).
2. The green and red arrows indicate the current measurement value compared to the last measurements:
 - If the result has increased, the red arrow lights up.
 - If the result is the same, no arrow lights up.
 - If the result has decreased, the green arrow lights up downward.

3. Different indicators are displayed based on the age of the user.

Metric	Required Age	Metric	Required Age
Weight	Effective at any age	Visceral fat grade	≥ 18 years old
BMI (Body Mass Index)	≥ 4 years old	Skeletal muscle percentage	≥ 18 years old
Body Fat %	≥ 10 years old	Muscle mass	≥ 18 years old
Body Water %	≥ 10 years old	Bone mass	≥ 18 years old
Subcutaneous Fat	≥ 10 years old	Protein	≥ 18 years old
BMR (Basal metabolic rate)	≥ 10 years old	Metabolic age	≥ 18 years old
Fat-Free Body Weight	≥ 10 years old		

Note: Some of the metrics vary according to different ages in our system.

MEMORY FUNCTION

Your LCD bathroom scale offers a memory function that can store up to 20 measurements for 8 users. It automatically saves measurements when a measurement is made with the app closed, overwriting the oldest entry when the memory is full. When you open the app, all stored data transfers to your phone. The smart scale recognizes users and may send notifications when new data is detected, streamlining your tracking experience. When the weight difference between two users is less than 2kg, the scale can't identify the user. When connecting the app, you can choose which user the data belongs to.

BLUETOOTH TROUBLESHOOTING

If you experience frequent Bluetooth connection failures, please try these steps:

1. Reset your FITAGE scale by taking out one battery for 5 seconds and then put it back in.
2. Reset the Bluetooth connection on your smartphone or tablet.
3. For older operating systems (like Android 6.0): ensure that Bluetooth and Location Permissions are enabled on your phone. This can be done in the general settings of the mobile device.


Then:

- Open the FITAGE app
- Select the menu on the home screen
- Select My Devices in the top right "☐●".
- Find your smart scale and slide it left to delete the device
- Tap "☐●" on the top right to pair it again.
- Step with one foot on the scale to turn it on.
- After a few seconds the scale will appear.
Tap on your device.
- Step on the scale with both feet to connect the device with the app via Bluetooth.

If this does not work, please go to “My Account” and select “Feedback”. Send your feedback. Then mail our customer support team with your problem and account mail address. Together we will try to solve your issue.

VIEW REPORTS

Click “Tracking” to view your progress history. Tab the clock icon in the to see the calendar and individual measurements. Here you can delete measurements and export your data. Share data with your health buddies via Facebook, Instagram or other social media channels.

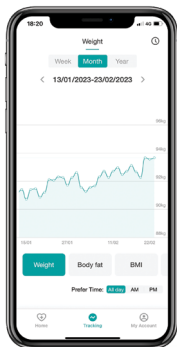SYNC TO FITNESS APPS

Apple Health



1. Enter Settings and select "Apple Health" (iOS only)
2. Turn on the Apple Health button
3. Turn on all the desired metrics that you want to share

Health Connect



1. Enter Settings and select “Health Connect” (Android only)
2. Turn on the Health Connect button, login to your Health Connect account
3. You will be redirected to the Health Connect page
4. After the connection has been made, the data will be synched to Health Connect
5. On Android 13 and below, Health Connect is an app without a desktop icon. On Android 14 and above, Health Connect is a system component.

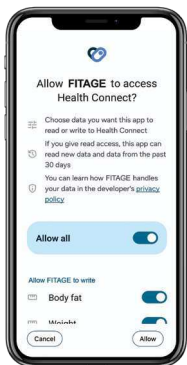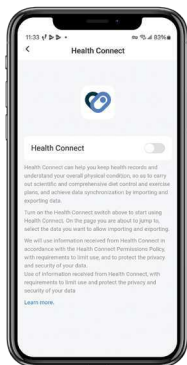Fitbit App

1. Enter Settings and select "Fitbit"
2. Turn on the Fitbit button, login to your Fitbit account
3. Select "Allow" to access the Fitbit app
4. After connection, the data Weight, Body Fat and BMI will synch to the Fitbit app. Check your data via their app or website.

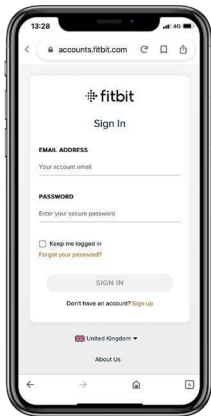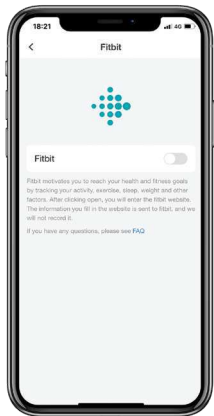FREQUENTLY ASKED QUESTIONS

1. How do I get the most accurate results?

- Stand firmly on the scale, do not move while measuring
- Place your smart scale on a hard and flat surface
- Fill in correct information on your profile setup, accurate body composition matters when measuring specific metrics
- Stand uniformly with both feet on the scale and make sure your feet touch the electrodes
- Measure your weight each day at the same time under the same conditions (for example, directly after waking up)
- Measure with bare and dry feet

2. Why is my smart scale not working?

- Please check if the batteries are installed properly, replace the batteries if necessary
- When data disappears with a flash; restart the smart scale by taking out one battery for five seconds and put it back in

3. The smart scale only measures my weight?

- Please step on the scale with bare and dry feet. Ensure that the Bluetooth connection is working, and ensure that your personal data like gender and height is entered. Double check if your feet are positioned correctly on the smart scale.

4. Are my body fat or muscle mass measurements accurate?

- When measuring body fat or muscle mass there is always an error factor. Even when you undergo a professional DEXA scan, the results will not be completely accurate. It is extremely difficult to measure accurate body fat percentage due to many biological factors. Use the smart scale as an indicator and do not rely completely on its body fat measurement.
- When the measured body fat is very off, please turn on/off the Athlete mode in the app under your profile. This will correct the measurement accordingly.

5. I have an 8 in my display and the scale is not working. What can I do?

Please reinstall or replace the batteries with new ones. If that does not work, contact our support team.

6. How can I delete one of the measurements?


When you are in the FITAGE app, you can click on tracking (bottom menu bar middle option). Then click on the history icon (top right). Then you can delete the specific data.

7. Can I manually add weight to the app?

Yes, this is possible in the home screen. Select the menu and tab on "Manual Record".

8. I have 3 stripes on the display and it does not turn off?

Please contact our support team for further assistance.

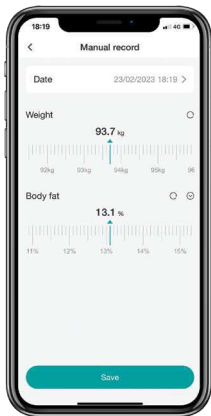WARNING

Please read carefully before use.

- This Bluetooth scale should not be used to diagnose or treat any medical issues. You should always consult a doctor when undertaking any diet, exercise program or other medical treatments.
- If the scale has a malfunction, please check its batteries first and replace them if needed.
- Do not use this product if you wear a pacemaker, artificial lung, ECG or any other medical device.
- Do not use the bioelectrical impedance analysis (BIA) while pregnant.
- Check the device before use. Do not use the scale if it is damaged. Continuous use of a damaged scale may cause improper results or injuries.
- Do not use the scale on a slippery floor and do not use it when your body is wet. Keep your body and the scale dry when using the device.

- When the FITAGE scale is not in use, please remove the batteries and place the device in a dry room that protects it against moisture, heat, lint, dust and direct sunlight for optimal durability. Never place any objects heavier than 180 kg (396 lbs) on the weight scale.
- The FITAGE scale is a high-precision device. Do not jump or stomp on the device. Please handle the scale with care to avoid damage.
- This Bluetooth scale is not fully waterproof. Never put the scale into water. If the scale accidentally gets wet after usage, clean it. Do not use soap or other cleaning chemicals.
- The data measurements may not be accurate when using the scale immediately after heavy physical exercise.
- This scale is only suited to full extend for users over the age of 16. There are restrictions in usage for people that are younger.

WE WANT TO LISTEN TO OUR CUSTOMERS

Your opinion is everything for FITAGE.

Your experience drives us to help contributing to improve people's health and fuels our community.

If you like your FITAGE product, please share your experience with us!

1. Go to Amazon / Bol.com / fitage.com
2. Click "Your Orders" and search for the FITAGE product.
3. Click on "Write A Review" and tell us about your experience!

ANY PROBLEMS?

If you encounter any problems,
please contact us via Email.

We are available on
Monday - Friday 09.00-17.00 (+1 GTM)
in English or Dutch.

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