FITAGE

USER MANUAL

ENGLISH	<u>2</u>
NEDERLANDS	<u>11</u>
FRANÇAIS	<u>20</u>
DEUTSCH	<u>29</u>
ESPAÑOL	<u>38</u>
ITALIANO	<u>47</u>

WELCOME TO

FITAGE

GET IN SHAPE

Introducing the FITAGE massage gun, a revolutionary device designed to help you achieve optimal health and wellness. This massage gun is engineered to provide deep tissue massage and relaxation, helping you recover faster and perform better.

Whether you're an athlete looking to improve your performance, or simply seeking relief from daily stress and tension, the FITAGE massage gun is the perfect tool to help you reach your goals. With its advanced features and user-friendly design, this massage gun is a must-have for anyone looking to become healthy again. So why wait?

Invest in your well-being today with the FITAGE massage gun.

SPECIFICATIONS AND CERTIFICATIONS

Product Name	FITAGE Massage Gun
Product Model	FITAGE-04
Massage Heads	12 Pieces
Battery Type	18650 Lithium Battery
Motor	Brushless Motor
Exterior shell	ABS
Power Supply	USB-C
Input Voltage	12V
Certifications	CE, RoHS

WHAT'S INCLUDED

- 1x Percussion Massage Gun
- 1x Travel Case
- 12x Massage Heads
- 1x USB-C Cable
- 1x Massage Ball
- 1x Instruction Manual

INSTRUCTIONS FOR USAGE

Install and remove massage heads:

- Choose the massage head you want to use and insert it directly into the device. Don't twist the massage head while inserting it.
- Before using the massage gun, make sure that the massage head is firmly inserted and secured in place.
- To remove the massage head, pull it straight out of the massage gun. Avoid twisting while removing.

How to turn the unit on and off:

- To turn the unit on, move the power switch located at the bottom of the battery to the ON position.
- Once the unit is turned on, press the fingerprint button at the center of the LCD screen to activate the first gear. The display will show "1".
- To increase the speed to the next gear, press the fingerprint button at the center of the LCD screen.
- To turn the unit off, press and hold the fingerprint button at the center of the LCD screen. By turning the unit off, the gear resets to 1 when turned on.
- To completely turn off the power of the device, move the power switch located at the bottom of the battery to the OFF position.

1. Pressure indicators.

The more pressure there is been given, the more indicators will turn on.

The 4th and 5th dots are red and represent the highest pressure range. Be careful with applying high pressure for a long amount of time.

2. Speed level indicator.


3. Decreasing speed button.

4. Increasing speed button.


**5. Start/Stop button,
pause the massage gun
and reset it to '00'**

USB-C
Charging port

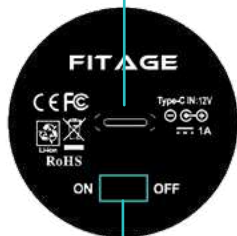Power switch

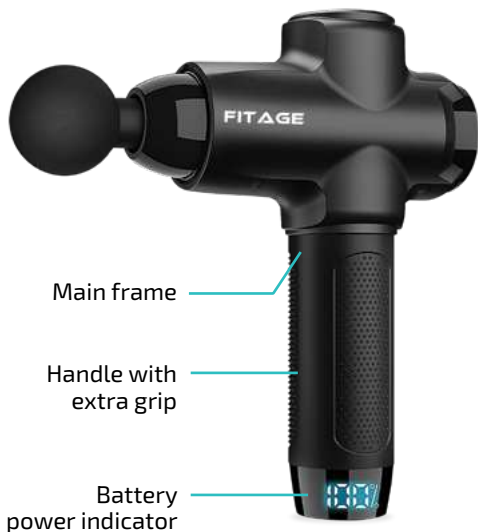Main frame

Handle with
extra grip

Battery
power indicator

HOW TO USE THE MASSAGE GUN

1. Start massaging the muscles by setting the massage gun to the lowest setting.
2. Place the massage gun on the muscle and hover over it with the massage gun. Gradually increase the speed if you want the massage to be stronger. You don't have to apply a lot of pressure; the massage gun does most of the work.
3. Massage back and forth across the muscle and avoid keeping the massage gun in the same spot for longer than five seconds. If you are massage a more tense area, you can hold a little longer. Be sure not to get too close to harder bone area's.

BATTERY INSTRUCTIONS

The device offers approximately 5 to 7 hours of battery life on the lowest setting. When the battery is low, it takes about 2-3 hours to fully charge it. The total working time of the battery depends on the intensity level on which the massage gun is being used.

To start using the massage gun, charge the battery fully before the first use. Connect the USB-C charging cable to the charging port of the unit. While charging, the battery power indicator on the bottom of the device will flash to indicate the current battery level. The indicator will show 100% when the battery is fully charged. You can recharge the battery at any time and at any battery level, but it is recommended to recharge it before it is fully depleted. Do not turn on the device while the battery is charging.

MASSAGE HEADS & TARGETED MUSCLE GROUPS

ROUND HEAD  Palm, Back, Buttocks, Legs	MUSHROOM HEAD  Joints, Sensitive Muscles	U-HEAD  Neck, Back, Achilles
SPADE HEAD  Back, Buttocks, Legs, Calves	FLAT HEAD  Back, Buttocks, Legs, Calves	THUMB HEAD  Palm Plantar
SEMI-CIRCLE HEAD  Legs, Buttocks, Arms	T-SHAPED HEAD  Back, Buttocks, Legs, Calves	FAN-SHAPED HEAD  Back, Buttocks, Legs, Calves
DISPERSED HEAD  Shoulders, Buttocks	MUSHROOM HEAD  Joints, Sensitive Muscles	GYRO HEAD  Joints, Sensitive Muscle

BENEFITS

EFFECTIVE AND RELIEVING

IMPROVES
BLOOD CIRCULATION

REJUVENATES

RELIEVES PAIN

MAINTENANCE

- To clean the massage gun and massage heads, use a slightly damp cloth and then gently dry them with a soft cloth.
- Make sure to regularly clean the motor heat sink and keep it free of any debris or dust.
- When storing or traveling with the massage gun, make sure to place all the items in the travel case and securely close it.

FREQUENTLY ASKED QUESTIONS

1. The massage gun automatically turns off after 10 minutes of use. Is there any way to prevent this?

Yes, the massage gun is designed to automatically turn off after running for more than 10 minutes or if placed on a hard surface with too much pressure. To prevent this, simply turn the massage gun back on if you need to continue using it.

2. How loud is the massage gun when set to its highest level?

The FITAGE Massage Gun is designed to provide high power while keeping

noise levels low, ranging from 35dB to 55dB. It is more quiet than standard massage guns on the market.

3. What is the warranty period for this product?

We offer a two-year repair warranty for this product. Please refer to the next page for any exceptions.

4. What should I do if my massage gun won't turn on?

1. First, make sure you have properly charged the device using the included power charger. The charger light will turn green once the battery is fully charged.
2. If that doesn't work, check that the power button is switched to the "ON" position.
3. If the massage gun still won't turn on, please contact our customer support team at **support@fitage.nl**.

NOTE:

SAFETY WARNING:

1. To prevent damage to your property or injury to yourself and to ensure that the product continues to function properly over time, please read and follow these instructions carefully while using the product.
2. This unit is not suitable for individuals who are 12 years old or younger.
3. In order to minimize the risk of injury, do not directly massage the skin with the device. Wear clean and dry clothing over the massage area and use gentle pressure for no more than 60 seconds to prevent abrasions.
4. Use high speed and high pressure settings only on soft tissue areas of the body to prevent pain and discomfort. Do not use high speed and pressure on the head or any bones to avoid injury.

5. Keep your fingers and hair away from the telescopic rod and motor heat dissipation hole to avoid pinching injury.
6. Use the device in a dry environment with normal room temperature to maintain its long-term capability.
7. Handle the device gently and avoid dropping or hitting it, as this can damage the device.
8. This device is for household use only and should not be used commercially. Do not run or charge the device in an unmanned state.
9. Do not use the device to diagnose or treat any medical issues, and consult a doctor before using it if you have any health problems. Do not use the device if you are pregnant, have diabetes, neuropathy, retinal damage, wear a pacemaker, or have undergone surgery.
10. If the device malfunctions, check the battery level first and charge it if necessary. Do not attempt to open or modify the device. Use the original charger and remove it after the battery is fully charged.
11. To avoid damaging the device, do not put it in water or let water enter the openings. If the device becomes wet, use a dry cloth to clean it and do not use soap or any other cleaning agents. Also, remember to clean the motor heat sink frequently to keep it free of dust and debris.
12. Do not use the unit continuously for more than one hour.
13. To prevent any damage to the device, do not use or charge it in high-temperature environments. Also, make sure to remove the charger from the unit before using it and after it has fully charged.

WARRANTY

We offer a two-year repair warranty on our product, except in the following cases:

- The equipment is damaged due to incorrect operation or transportation by the user.
- The equipment has been tampered with or repaired without authorization.
- The customer did not refer to the instructions manual for assistance.
- The product has been overused beyond its recommended capacity.
- Damage is caused due to poor storage or maintenance by the customer.
- The customer fails to provide proof of purchase.

WE WANT TO LISTEN TO OUR CUSTOMERS

At FITAGE, we value your opinion and strive to use your feedback to improve global health and build our community. If you're happy with your FITAGE product, we'd appreciate it if you could share your experience.

Here's how you can do it:

1. Go to Amazon / Bol.com / fitage.com
2. Click "Your Orders" and search for the FITAGE product.
3. Click on "Write A Review" and share your thoughts with us!

ANY PROBLEMS?

If you experience any issues, please don't hesitate to contact us via email.

Our customer service team is available
Monday to Friday from 09.00-17.00 (+1 GTM)
and can assist you in English or Dutch.

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