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V-Goods Premium Ultrasonic tooth cleaner User Manual

• Specifications

Charging Voltage: 5V
Charging Current: ≤1A
Charging Time: 2-3 Hours
Battery Capacity: 3.7V D.C.: 1500 mAh
Power consumption: ≤2.8W
Frequency: 36-40 kHz
Product Size: 1.3 × 7.4 Inch
Net Weight (about): 5.6 oz/159g
Gross Weight (about): 11.6oz/330g

• Package list

1 × Electric dental calculus remover
2 × Pointed teeth cleaning head
1 × Flat teeth cleaning head
1 × Dental mouth mirror
1 × Spanner
1 × User manual
1 × USB cable

• Charging and Store


1. Charge the device with 5V 1A adapters.
2. Charge the device when the power indicator turns yellow during use. Note: The battery will run out if the yellow light is flashing.
3. Please charge the device every 6 months for better battery maintenance.
4. Place the device in the cool place. Avoid direct sunlight or high temperature.

• How to use

1. Power off and fully charged before use.
2. Remove the transparent dust cover.
3. Choose the proper teeth cleaning head.
Pointed teeth cleaning head for the slit of the teeth
Flat teeth cleaning head for the tooth surface


4. Install and TIGHTEN the head with the spanner
5. Press the power button to turn on the device.
6. Press the power button to choose a suitable mode.
The corresponding mode indicator will turn green.
Note: With the smart memory function, the device can retain the previous selection mode.
7. Press the light button to turn on the built-in spotlight.
8. Hold the dental mouth mirror. And apply the teeth cleaning head to the needed area on your teeth to remove the tartar or stain.
Note:

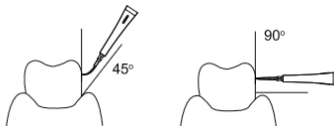A: Keep an angle of 45 degrees between the pointed teeth cleaning head and teeth to remove the tartar.
B: Keep an angle of 90 degrees between the flat teeth cleaning head and teeth to polish or remove the stain.
9. The device will shut down after 10 minutes of operation. Please restart it for continuous use.
10. Press the power button for 2 seconds to turn off the device after use.
11. Disinfect teeth cleaning head with alcohol or boiling water for 5 minutes, and then dry it for the next use.

• Product Description

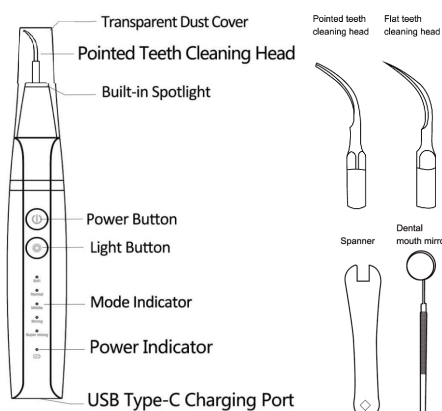• Safety Instructions

1. Do not drop the device or insert any objects into any opening or hose.
2. Use this device only for its intended use as described in this manual.
Do not use attachments not recommended by the manufacturer.
3. This device is not intended for use by children (Less than 16 years old).
4. Keep the power cord away from the heated surface.
5. Do not operate the device with a damaged cord or not work improperly.
6. The teeth cleaning head is sharp, please be aware of injured.
7. Do not use the device while sleeping or drowsy.
8. Do not touch the teeth cleaning head while the device is operating.
9. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
Do not use in the following situations
A. Patients with periodontitis, gingivitis, and other dental and oral diseases
B. Bleeding gums
C. Dental treatments include crowns, bridges, root canals, fillings

• FAQ

1. My gums bleed when I use the device.
A: The device cleans the tartar with the high-frequency vibration, which stimulates the gums. Gums bleeding is normal and not harmful to the gums.
2. The device shuts down during use.
A1: Power indicator light is green: operating more than 10 minutes, the device will shut down. Restart it for continuous use.
A2: Power indicator light is yellow or the device cannot operate: the battery is run out, please recharge the device.
3. The device seems not powerful.
A: Please ensure the connection between the teeth cleaning head and the device is well. If not, tighten them with the spanner.
4. The teeth cleaning head is getting hot when using.
A: It's normal as the ultrasonic energy accumulates in the head for the high-frequency vibration.
5. The device is getting hot when using it.
A: The slight temperature rise is normal during operation.

• Warranty

The product supports a one-year warranty from the date of purchase. This warranty does not apply to products which have been lost, altered, broken, or damaged by misuse or failure to follow the instruction provided with the product.
The consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, and return the entire product to support a warranty claim.
To ensure proper use and trouble-free operation, please read this manual carefully and store it in a safe place for future reference.

Waterproof level: IPX 6