LANAFORM

WELLNESS

Massage Mattress

VIBRATING MASSAGE

4 MASSAGE AREAS

2 HEATING ZONES

3 INTENSITIES

5 PRESET MODES

MATELAS
MASSANT

MASSERENDE
MATRAS

MASSAGE
MATRATZE

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
MATERAC
DO MASAZU

1

- EN Control panel, 4.
- FR Commande de contrôle, 6.
- NL Bediening, 8.
- DE Fernbedienung, 10.
- ES Panel de control, 11.
- IT Comando di controllo, 13.
- PL Urządzenia sterujące, 15.
- CS Dálkové ovládání, 16.
- SK Ovládací príkaz, 18.
- HU Távirányító, 20.
- SL Daljinski upravljalnik, 21.
- HR Upravljač, 23.
- RO Comenzi, 24.
- SV Manöverkontroll, 25.

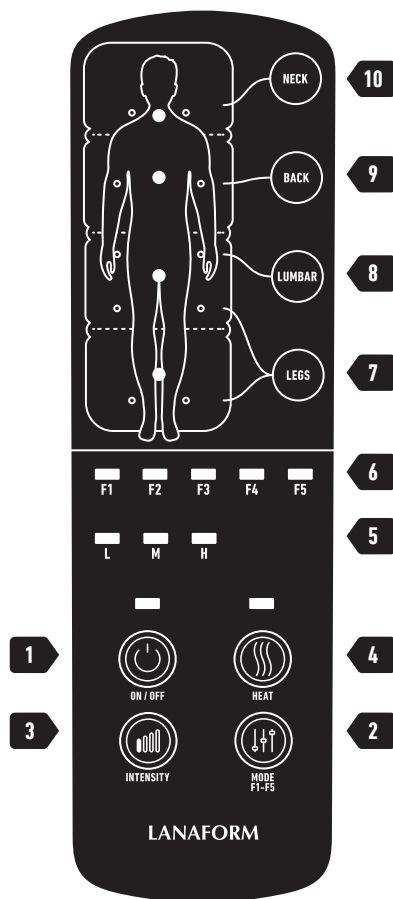EN Instructions

📷 *The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.*

INTRODUCTION

Thank you for choosing the LANAFORM *Massage Mattress*, a quality massage device equipped with ten motors, enabling you to find the massage you want. They can work independently or together over the length of your body, using the programmed options. This massage mattress also has two heating zones. Its relaxing, stimulating and/or invigorating effects will help reduce tension and muscular tiredness.

PLEASE READ ALL THE INSTRUCTIONS BEFORE USING YOUR MASSAGE MATTRESS, ESPECIALLY THIS BASIC SAFETY INFORMATION

- Only use this device according to the instructions given in this manual.
- Before using the appliance, remember to remove the silica gel desiccant.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have

been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

- Check that your mains voltage is the same as that of your device.
- Do not use accessories that are not recommended by LANAFORM or that are not supplied with this device.
- If the power cable is damaged, it must be replaced by a similar one available from the supplier or its after-sales department.
- Do not use this device if the power plug is damaged, if the unit is not working properly, has been dropped on the floor or if it is damaged or has been dropped into water. Have the device examined and repaired by the supplier or its after-sales department.
- If there is an unusual smell when it is being used

as described in these instructions for use, switch the device off, unplug it and have it examined by the supplier or the latter's after-sales department.

- Do not carry this device by its power cable or use the cable as a handle.
- Always unplug the device after use or before cleaning it.
- Unplug a device that has fallen into water immediately before retrieving it.
- To disconnect the device from the mains, ensure that it is in the "ON/OFF" position and then unplug it.
- Do not connect or disconnect the plug of the device with wet hands.
- Do not take apart either the device as a whole or even part of it.
- An electrical device should never remain connected unsupervised. Unplug it when you are not using it.
- Keep the power cable away from hot surfaces.
- Do not use this product before going to sleep. Massage acts as a stimulant

that may hinder sleep.

- Never use this massage device when sleeping.
- Never let anything fall into and never insert any object (tacks or other metal fastening or clasp) into any of the openings of this device.
- Never use this device in a room where aerosols (sprays) are used or in a room or where oxygen is being administered.
- Do not use this device under a cover or under a cushion. Excessive heat can cause fire, electrocution or injury.
- Do not leave the device exposed to direct sunlight.
- Do not use this device when it is folded or rolled up.
- Do not place this device on any white or light-coloured seat or armchair to avoid any colour transfer from this massage mattress.
- If you are concerned about your health, consult a doctor before using this appliance.
- If you experience any pain when using this device, stop

using it immediately and consult your doctor.

- Never use this device on a body part that is swollen or inflamed or if you have a skin rash.
- This product is a non-professional massage appliance designed to relax tired muscles. Never use this massage device as a substitute for proper advice from your doctor or as a substitute for medical treatment. This device is intended only for use in the home.

▲ *Any pain and any persistent symptoms must be regarded as indications of a problem that may be more serious. If any of these appear, please consult your doctor before continuing to use the Massage Mattress.*

- Do not use this equipment without consulting your doctor in the following cases:

If you are pregnant, lactating, asthmatic, if you have a physical, neurological or mental disability, if you are suffering from psychological trauma, if you have a

pacemaker, a prosthesis or hernia, back problems, disc disease, circulation problems, stomach problems, if you are epileptic, diabetic, if you are prone to thrombosis, if you suffer from migraines or if you are under the influence of alcohol or drugs. This list is not exhaustive. If in any doubt, consult your doctor before using the *Massage Mattress* or if doubts or health problems arise. In any event, you must be in good health to use this equipment.

- Do not use this device for longer than 15 minutes at a time.

INSTRUCTIONS FOR USE

Features

- Vibration-based massage
- 10 motors
- 5 preset modes combining different massage areas
- 3 levels of massage intensity
- 4 massage areas (shoulders/upper back/lower back/legs)
- 1 "heating" function made up of 2 heating zones
- Switches off automatically after 15 minutes

Description of the Massage Mattress

System of straps and positioning of the device

The *Massage Mattress* has a system of straps so that you can fix it to any long seat. Just thread the straps across the back of the chair, joining them up and fix them to each other with a plastic buckle. You can also place this massage mattress on your sofa or settee, your bed or even on the ground if you wish.

Control panel, ■ 1

- 1 To select on/off
- 2 To select the preset massage mode
- 3 To select the massage intensity level
- 4 To select the heating function
- 5 Displays massage intensity level selected
- 6 Displays the preset massage mode selected
- 7 To select the massage area legs
- 8 To select the massage area lower back
- 9 To select the massage area upper back
- 10 To select the massage area shoulders

Using the Massage Mattress

The *Massage Mattress* is used with a transformer. First of all, connect the end of the power cable to the female input connector located on the right-hand side of the massage mattress. Then connect the transformer to a plug socket. Once connected to the plug socket, the *Massage Mattress* is ready for use.

- ON/OFF

Turn the appliance on or off by simply pressing the "ON/OFF" button (diagram 1, point 1).

- OPERATION

Very easy to use, enabling you to select

- A preset massage mode from the five modes on offer (diagram 1, point 2).
- One of three different levels of vibration intensity (diagram 1, point 3).
- The heating function (diagram 1, point 4).
- One or more of four massage areas to suit your needs and/or wishes (diagram 1, points 7, 8, 9, 10).

MAINTENANCE

Cleaning

- Unplug the mattress from the mains and wait a few minutes before cleaning it.
- Clean it using a soft, slightly damp sponge.
- Never let water or any other liquid come into contact with the mattress.
- Do not immerse it in liquid to clean it.
- Never use abrasive detergents, brushes or a thinning agent to clean it.

Storage

- Place the massage mattress in its box and in a dry place away from direct sunlight.
- Avoid contact with sharp edges and other pointed objects, which might cut or pierce the fabric surface.
- To prevent the power cable breaking at any point, do not wind it round the mattress.
- Do not hang the mattress up by the manual control lead or power cable.

ADVICE REGARDING THE DISPOSAL OF WASTE

All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.

When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

LIMITED WARRANTY

LANAFORM guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM's control.

LANAFORM may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM's Service Centre. Any attempt to maintain this product by a person other than LANAFORM's Service Centre will render this guarantee void.

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