
beper

IT	Massaggiatore cervicale	pag.2
EN	Neck Massager	pag.7
FR	Masseur cervical	pag.12
DE	Zervikales Massagegerät	pag.17
ES	Masajeador cervical	pag.22
GR	Μασάζ Αυχένα	pag.27
RO	Dispozitiv pentru masarea gâtului	pag.32
CZ	Masážní přístroj na krk s technologií EMS	pag.37
NL	Nek massage toestel	pag.43
LV	Kakla masieris	pag.48

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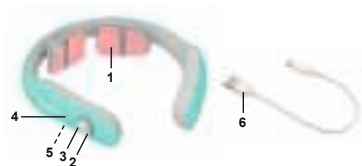Warnings for use

Do not use this device in combination with other electronic medical devices, such as pacemakers, artificial hearts and lungs, and other life functions maintenance and monitor devices.

The use of this device in combination with other high frequency surgical appliances can cause burns and damages to the massager.

Consult a doctor before using this device in case of pregnancy, high skin sensitivity, heart or cardiovascular diseases, tumors, or other medical conditions or chronic diseases that require periodic medical assistance.

Do not use this device under the shower, while bathing and/or sleeping.

This device cannot be used by children and/or by individuals with reduced physical, psychological, or sensory capabilities.

This device cannot be used by individuals who are insensible to heat and/or with skin perception issues.

Do not immerse in water or other liquids.

Tips for use

This massager is designed to ergonomically adapt to the neck area, to ensure correct stimulation of the area adjacent to the cervical vertebrae.

It is advisable to use this device twice a day, for max.15 minutes per session.

Before use, it is better to slightly wet the skin on the neck area (with warm water or water-based gel or cream), to ensure the perfect adhesion of the electrodes to the skin.

The intensity of the electric impulses should be gradually set, starting from the lowest level 1, and slowly increasing to the next levels, always paying attention to your own needs and feelings; if the massage does not feel right or creates discomfort, it is better to lower the intensity level or to change massage mode.

If the skin is not hydrated enough, or if the electrodes are not touching the skin, it is possible to see a small rash on the skin, which will generally disappear in 1-2 days.

Description of product Fig.A

1. Electrodes
2. ON/OFF and mode selection key
3. Intensity selection knob
4. Heating function key
5. USB connector
6. USB cable

Charging

Before first use, it is better to fully charge the device. Insert the USB cable inside the connector, and connect the massager to a charging device (powerbank, laptop, etc.) or to a USB adapter (not included). The LED indicator's blue light will flash during charging, and it will stay on when charging is complete.

Disconnect the USB cable from the massager once charging is completed.

If the LED indicator flashes during operation of the device, this means the battery is running out; proceed to charge it as described above.

Operation

Before use, it is better to slightly wet the skin on the neck area (with warm water or water-based gel or cream), to ensure the perfect adhesion of the electrodes to the skin, and to avoid rashes.

Press the ON/OFF key for a few seconds to activate the device.

The device is equipped with a guide voice that will indicate selected mode and intensity level.

When switched on, the device will start operating in auto mode: the blue LED indicator is on, and the guide voice states "Mode 1" (mode one); the device operates combining all other 5 massage modes available, alternating them in cycles of 1 minute.

By pressing the ON/OFF key, it is possible to select the desired massage mode, the guide voice will state the selected mode (Mode 1 – mode one, Mode 2 – mode two, Mode 3 – mode three, Mode 4 – mode four, Mode 5 – mode five, Mode 6 – mode six)

	Mode 2	Mode 3	Mode 4	Mode 5	Mode 6
Pulse width	100µs	100µs	500µs	500µs	500µs
Pulse frequency	50Hz	40Hz	10Hz	25Hz	100Hz
On-Time	20s	10s/2*	10s	4s	6s/2*
Off-Time	2s	2/10s**	1s	0,5s	5/4s**

On time (s) indicated switch-on time cycle in seconds (contraction).

Off time (s) indicates switch-off time cycle in seconds (relaxation).

* There are 2 individual rounds in cycle 3 and 6.

** The off-time is 10 seconds for every 2nd round in cycles 3 and 6.

To adjust pulse intensity, rotate the intensity selection knob (clockwise to increase intensity, counter clockwise to reduce it). For heating function, press the heating function key, the low heating setting will activate (temperature between 38 and 41° C); by pressing the key again, the max heating level will activate (temperature between 42 and 45°C); press the key once again to deactivate heating function. To switch the device off, press the ON/OFF key for a few seconds, until the blue LED indicator will turn off.

Note: the device automatically shuts off after 15 minutes of continuous operation

Cleaning and maintenance

Make sure that the device is off and disconnected from USB cable before proceeding with cleaning.

Remove any dust and/or stain found on the device using a soft cloth slightly dampened with warm water.

⚠ Warning: do not immerse in water or other liquids; do not use solvents or other detergents

Store the device in a cool, dry place, away from direct sunlight, preferably inside its box.

Technical data

Lithium battery 900mAh

Input: 5V DC USB

For any improvement reasons, Beper reserves the right to modify or improve the product without any notice.

The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment. The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

GUARANTEE CONDITIONS

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- The appliance has not been tampered with, otherwise it is not maintainable.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) the aesthetic parts are excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, carelessness, wrong or improper installation, damage during transport and any other damage not attributable the supplier. For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge.

In any case, if the part to be replaced for defect, breakage or malfunction is an accessory and/or a detachable part of the product, Beper reserves the right to replace only the very part in question and not the entire product For technical assistance and/or repairs after the guarantee period you can directly contact the below address:

Contact your distributor in your country or after sales department Beper.

E-mail **assistenza@beper.com** which will forward your enquiries to your distributor.

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