## **BenjaStore**

Info@benjastore.com

www.benjastore.com

## **Manufacturer:**

**Shenzhen Mison Technology Co., Ltd**

Email: info@mison.com

Web: www.mison.com

MS-306C

5 IN 1 Multifunctional Beauty Device

Use Manual

- English 2
- Español 12
- Français 22
- Deutsch 32
- Italiano 42
- Nederlands 52

# CONTENTS

Product Introduction .....	2
Packing List .....	2
Product Schematic .....	3
Safety Precautions .....	3
Product Use Method .....	5
Product Maintenance .....	9
Common Problems And Solutions .....	9
Product Specification .....	10
Product Warranty .....	10

## ► PRODUCT INTRODUCTION


This multifunctional beauty device is using the 4 most popular beauty technologies: RF(Radio frequency),EMS(Micro current), Acoustic vibration, Led light therapy, so that you can enjoy professional skin care at home like beauty institutions.

At the same time our company carefully study the five major beauty modes: clean, Lift, Lead in, Anti-aging, Remove acne, you can choose to use which mode according to your needs, saying goodbye to the traditional, complex choices of mode beauty device, making your using more easily. The common skin care products and only affect in the skin surface, our device can stimulate the deep skin, so you can effectively improve the skin problems.

## ► PACKING LIST

- User Manual x1
- Main Body x1
- USB Charging Cable x1
- Charging Base x1

## ► PRODUCT SCHEMATIC
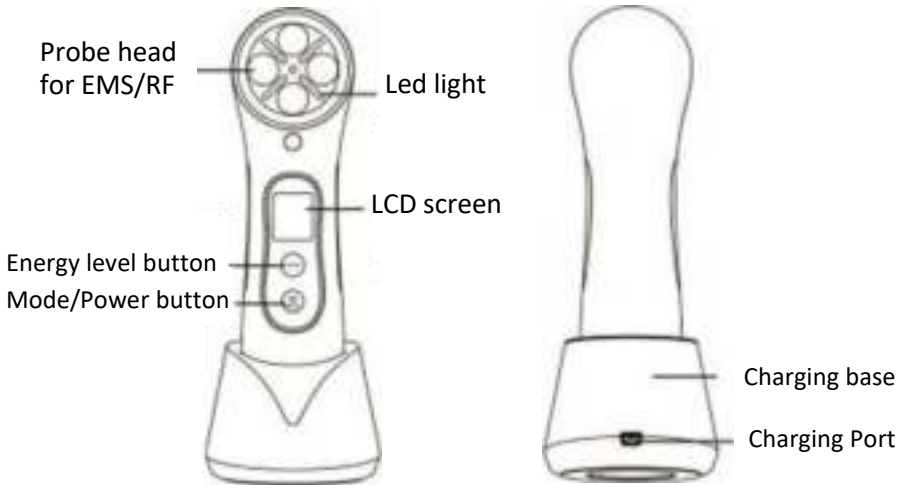
## ► SAFETY PRECAUTIONS

Please read these safety precautions carefully, then use the device according to the manual. In order to make sure you can use the device correctly and avoid the damage to others, please follow these precautions.

### Warning

- Please do not use other ways besides the manual to use the device, otherwise it may break down and cause an accident.
- If users are suffering from physical, sensory and psychiatric disease or undergoing dermatological treatments, then please use our product when you are monitored by some special person or responsible person (such as a skin specialist).
- Please do not twine the charging cable on the power plug.

- Please do not destroy the charging cable, and do not bend, twist, stretch it, or add heavy things on it, or put it on high temperature place, also do not bundle and process it to prevent from.
- Using the device for a long time or in wrong way will cause the charging cable broken, if you find following conditions, please stop using the device:
  1. The cable can not be original shape.
  2. The head of cable is damaged.
- Please do not throw the device in the fire.
- Please do not put the device into the water or wash it.
- Please do not give the device to the baby and child and put it away from them.

### Attention

- This product is only for daily use, please do not disassemble, remould, repair it.
- Please do not several people share this device.
- Please stop using the device and consult the professional doctors when the skin is pain and uncomfortable.
- Please do not use or store the device and its accessories in some damp places like shower room.
- Please do not throw the device strongly to avoid damage it by collision.
- Please stop using the device when it is broken.
- Please do not touch the device when it is in the water, and disconnect power, preventing from electric shock.

## ► PRODUCT USE METHOD

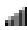
### Before Using

Please make sure the device charged 2 hours and enough energy.

### 1. Steps Of Using

Please press  and hold it for 3 seconds to turn on the device, light the Lcd screen, now the mode is acquiescent, that is dean mode, energy level is 1, time is 5 minutes, and time is 5 minutes countdown state.

### 2. Energy Level Adjustment

In any mode of beauty device, you can press  to choose which level, the order is from 1 to 5 level, which is showed by the  on the LCD screen.

Show 	Energy Level:1
Show 	Energy Level:2
Show 	Energy Level:3
Show 	Energy Level:4
Show 	Energy Level:5

### 3. Description of Modes

Press  to choose the mode, when the mode is selected, on the screen there will be a  after that mode, the mode which is not selected not have .

Using RF(Radio Frequency), EMS(Micro current), Acoustic Vibration, Led Lights technology together to study 5 beauty modes:

### **3.1 RF(Radio Frequency)Function**

RF is a high frequency and a thermal penetration therapy, which is through RF conductive head.it transmits a kind of electromagnetic energy caused by high frequency ac change to the subcutaneous tissue. It can also stimulate the dermis collagen regeneration, increasing the thickness and the density of the dermis. It will lift and fill wrinkles, eliminate the scars and restore skin elasticity and luster.

### **3.2 EMS(Micro current)Function**

It is a special technology of electrical perforation method, and it can in short time create space between the cells so that the essence of liquid components can be directly delivered to the skin surface and dermal layer, compared with Ions, its ability is 10 times on the 1MM from the surface of skin. 4mm is 25 times. So it can make sure the essence can be absorbed effectively.

### **3.3 Acoustic Vibration Function**

It can massage and smooth nerves, promote blood circulation, speed up the metabolism, make the skin more compact.

### **3.4 Led Light Function**

#### **Red light - Wavelength:620nm(±5nm)**

It can promote collagen regeneration, make the skin more elastic, suitable for relaxation or wrinkled skin use.

#### **Green light - Wavelength:520nm (±5nm)**

It can solve pigmentation, discoloration, freckles, sun spots and age spots. It can also dilute and reduce the existing pigmentation, improve the skin color, and protect the skin to from future potential pigmentation crisis. Melanin grows on epidermal surface, green light can inhibit the excessive production of melanin, preventing the pigment growth.

**Blue light - Wavelength: 465nm ( $\pm 5$ nm)**

It can kill the bacteria, promote blood circulation. It also can penetrate the skin below the 0.25MM, destroying bacteria living environment, eliminating acne on the skin, pain, congestion, which can make skin soothing and resting.

**Yellow light - Wavelength: 590nm ( $\pm 5$ nm)**

It will stimulate the lymph and nervous system, which is used to regulate the muscles and stimulate the immune system, and relax and restore the balance of sensitive skin.

**Pink light - Wavelength: 700nm ( $\pm 5$ nm)**

It will promote blood circulation, strengthen the beauty solution absorption, whitening skin, suitable for dark skin color.

**Five modes:**

Clean Mode	Clean up harmful substances in the skin, metal ions, reduce pigmentation, improve skin tone, accelerate blood circulation and skin metabolism
Lift Mode	Restore skin flexibility, pulling wrinkles, adjust the muscles and stimulate the lymphatic system, reshape the face curve
Lead in Mode	Create space between cells in short time, led in the essence of nutrition effectively
Anti-aging Mode	Promote skin collagen regeneration, enhance skin vitality, compact skin, improve skin relaxation
Remove Acne Mode	Promote blood circulation, anti-inflammatory, sterilization, effectively remove the skin acne, shrink pores

## 4. Skin Care Time

After research, the following is the best time for each mode of use:

- Clean mode                                      Default time 5 minutes
- Lift mode                                         Default time of 10 minutes
- Lead in mode                                    Default time 10 minutes
- Anti-aging mode                                Default time is 10 minutes
- Remove Acne mode                            Default time 5 minutes

Once you choose a mode, the LCD screen at the bottom of the time will be based on the default time to show the countdown state, such as 5mins, 4mins... until it turns off,

Long press  3 seconds, can also be turn off.

## 5. Operation Method

After cleaning the skin, take the essence.

Operate from bottom to top, from inside to outside of the skin, pulling on the skin to move.

This product must be used in conjunction with the beauty liquid, the use of normal skin is two or three times a week.

This product must be used in conjunction with the beauty liquid, the use of normal skin is two or three times a week.

## 6. Power Query Function

Power is 3 grid display, power display is as follows:

1 Grid Battery	Battery remaining capacity: 0-30%
2 Grid Battery	Battery remaining capacity: 30%-70%
3 Grid Battery	Battery remaining capacity: 70%-100%

When the battery is low, the upper right corner of the battery icon has been flashing, it is prompt to charge the device.

## 7. Charging

### Charging Method

- Connect the power cable and the adapter to the device base for charging, or connect the power cable and the adapter to the device directly.
- The battery icon flashes during charging.
- When the charge is complete, the battery icon returns to normal.

## ► PRODUCT MAINTENANCE

- After using, please use paper, towels or dry soft doth to clean the head of the device, when the dirt is difficult to remove, please try use soft cloth with a little water.
- When you do not use this product for a long time, it is recommended that you charge the device regularly and store it in a dry, cool place.
- The device is not completely waterproof, please do not use water to clean the instrument to avoid damage.

## ► COMMON PROBLEMS AND SOLUTIONS

- No Feeling on Face: 1 . Please put the instrument head parallel to the face.  
2. Adjust the energy power.
- The device is not work: Lack of battery remaining capacity, please charge it.
- The energy level is not available: 1. Internal parts of the device damaged  
2. Button stuck

- If this product can't work normally, please stop using it then contact the factory after-sales service personnel

## ► PRODUCT SPECIFICATION

Charging Voltage	5V, 500mA
Charging Hours	2 hours
Battery Capacity	3.7V, 1000 mAh
Vibration Frequency	8600 ( ± 10 ) rpm
Net Weight	135g
Charging Base	DC 5V 1A

## ► PRODUCT WARRANTY

1. This device warranty time is 12 months, and it counts from the day you buy.
2. In warranty, following the manual and normal use, then if the device appears materials defect, use troubles and craft problem, we will repair it for free.
3. Due to the improper storage, crash, unauthorized disassemble, the battery normal loss and other causes of product failure, these conditions is beyond the range of warranty.
4. This warranty card can not reissued under any circumstances, please be sure to keep it and present it at the time of warranty. In the after-sales service, if there are problems, please contact us.

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1. This device warranty time is 12 months, and it counts from the day you buy.
2. In warranty, following the manual and normal use, then if the device appears materials defect, use troubles and craft problem, we will repair it for free.
3. Due to the improper storage, crash, unauthorized disassemble, the battery normal loss and other causes of product failure, these conditions is beyond the range of warranty.
4. Assicurati di contattarci prima di intraprendere qualsiasi azione di restituzione: **Info@benjastore.com**  
 Si prega di allegare il numero dell'ordine di Amazon, Ti risponderemo e ti forniremo una soluzione soddisfacente in 24 ore.

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