TUBBLE[®]
AMSTERDAM

TUBBLE MASSAGE GUN USER MANUAL
EN / DE / FR / NL

TUBBLE

Enjoy your special moment of relaxation.

Follow us and share your Tubble moments:

Instagram: @letstubble/ Facebook: letstubble

If you have any questions, complaints or compliments, please e-mail us at **hello@tubble.com** and we'll be happy to help you!

TABLE OF CONTENTS

1. ENGLISH.....4

2. GERMAN9

3. FRENCH.....14


4. DUTCH.....19

To reduce the risk of electric shock, fire and personal injury, please follow the instructions below:

- Read all safety warnings and instructions.
- This massage gun is suitable to be used only by adults.
- Failure to follow the warnings and instructions may result in electric shock, fire and or injury.
- Keep this manual.
- Due to the force of the massage gun we recommend to wear clothing on the massage areas. Do not let the massage gun make direct contact with skin.
- Make sure you have the right massage head before you start using the muscle massage gun.
- Massage time on the same muscle part should not exceed 1 minute.
- Do not use on your head nor any bony parts of your body.

- Frequent high speed and high pressure on the same site can cause abrasion.
- Stop use immediately if you experience any pain or discomfort during massage.
- Keep fingers and hair away from the massager to prevent injury.
- Use only the original charger to charge the massage gun.
- Keep the massage gun away from water and from heat sources.
- Inspect the massage gun for any damage prior to each use. If damaged, do not use the massager.
- Modifications of this massage gun is strictly prohibited.
- Do not use during pregnancy or in case of diabetes, pacemaker, recent surgery, or with any other health issues.
- When unsure, always consult a certified doctor or physician before using the massage gun.

PRODUCT INFORMATION

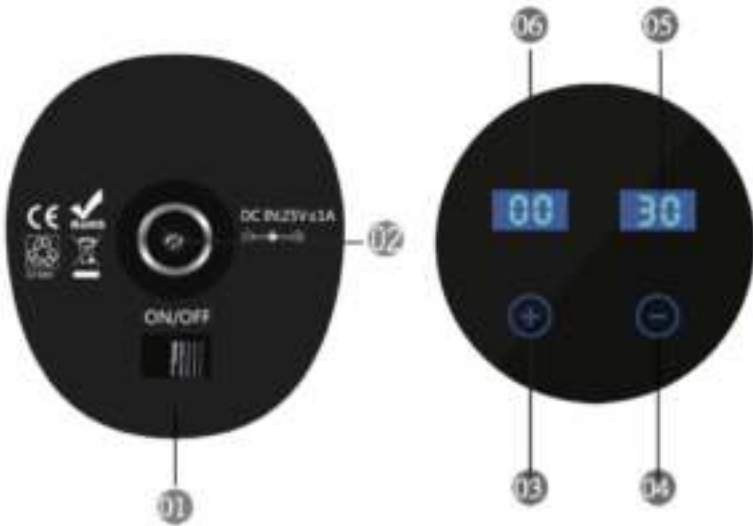No	Name	No	Name
01	Power Switch	02	Charge Port
03	Increase Intensity	04	Decrease Intensity
05	Battery Indicator	06	Level Indicator (1-30)

FUNCTION INFORMATION

1. Power Switch	Turn the power switch left (ON) to power on the message gun Turn the power switch right (OFF) to turn the power off
2. Charge Port	When using for the first time, make sure the battery is fully charged Use the original charger and make sure the machine is not powered on when charging
3. Increase Intensity	Click to increase the massaging intensity
4. Decrease Intensity	Click to decrease the massaging intensity
5. Battery Indicator	Shows the remaining battery life
6. Level Indicator	Shows the intensity level Be careful when you use higher intensity levels

MUSCLE MASSAGE GUN HEADS

Every muscle massage gun head is suitable for massaging a different muscle group

Photo	Name	Function
	60MM round head	Large muscle groups
	50MM round head	Large muscle groups
	Y head	Neck, spine and Achilles tendon
	Small flat head	Muscle relaxation and for all muscle groups
	A head	Joints, trigger points, deep tissue massage
	D head	Relaxation and shaping of various muscle parts such as shoulders, neck, IT band, spine, etc.

HOW TO USE THE BATTERY

- This massage gun can work up to 6 hours at the lowest level.
- For security reasons and battery life optimization the massage gun stops working automatically after 8-9 minutes. You can switch in on again immediately and continue the massage.
- The battery should be charged 4-5 hours when it runs out.
- Never leave the massage gun charging unattended.
- Be sure to turn the power off when the massage gun is not in use.

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